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DOMINIQUE DUTSCHER SAS

## TABLE OF CONTENTS

<b>1.</b>	<b>SYMBOLS, ABBREVIATIONS AND STICKERS .....</b>	<b>5</b>
<b>2.</b>	<b>GENERAL DATA .....</b>	<b>6</b>
2.1.	Aim of this manual .....	6
2.2.	Application .....	6
<b>3.</b>	<b>DESCRIPTION OF THE FREEZE DRYER.....</b>	<b>7</b>
3.1.	Technical features .....	7
<b>4.</b>	<b>INSTALLATION INSTRUCTIONS .....</b>	<b>8</b>
4.1.	Transport and unpacking .....	8
4.2.	Placing on site .....	8
4.3.	Electrical connection .....	9
<b>5.</b>	<b>BEFORE START-UP .....</b>	<b>10</b>
5.1.	Vacuum pump .....	10
5.2.	Refrigerating unit .....	11
<b>6.</b>	<b>DETAILS OF THE FREEZE-DRYER.....</b>	<b>11</b>
6.1.	Description .....	11
<b>7.</b>	<b>RECOMMENDATIONS .....</b>	<b>12</b>
7.1.	Ice condenser .....	12
7.2.	Greasing .....	12
7.3.	Long storage .....	12
<b>8.</b>	<b>WORKING PRINCIPLE.....</b>	<b>13</b>
8.1.	Freeze .....	13
8.2.	Lyophilization .....	13
8.2.1.	<i>In manifold</i> .....	13
8.2.2.	<i>In chamber</i> .....	14
8.2.3.	<i>End of the process</i> .....	15
8.3.	Lyophilization modes .....	15
8.4.	Selectable processes .....	15
<b>9.</b>	<b>OPTIONAL ACCESSORIES.....</b>	<b>16</b>
9.1.	Chambers .....	16
9.2.	Manifolds .....	18
9.3.	Flasks condensations collector (Drip tray) .....	18
9.4.	Lyologger software .....	19
<b>10.</b>	<b>OPERATION AND CONTROL INSTRUCTIONS.....</b>	<b>20</b>
10.1.	Start up .....	20
10.1.1.	<i>Main window</i> .....	20
10.2.	Freeze-drying .....	22
10.2.1.	<i>Probes</i> .....	23
10.2.2.	<i>Recipe</i> .....	23
10.2.3.	<i>Semiautomatic cycle</i> .....	25
10.3.	Auxiliary cycles .....	26
10.3.1.	<i>Defrosting</i> .....	26
10.3.2.	<i>Alarms</i> .....	27
10.3.3.	<i>Alarms list</i> .....	28
10.4.	Synoptic .....	29
10.5.	Settings .....	31
10.5.1.	<i>Language</i> .....	32
10.5.2.	<i>Date/hour</i> .....	32
10.5.3.	<i>IP Settings</i> .....	33
10.5.4.	<i>Technical Service</i> .....	33
10.6.	Recipe edition .....	34
10.6.1.	<i>Load recipe</i> .....	35
10.7.	Users .....	36
10.7.1.	<i>User's levels</i> .....	36
10.7.2.	<i>User login</i> .....	36
10.7.3.	<i>Password change</i> .....	37
10.7.4.	<i>Users management</i> .....	38
10.7.5.	<i>User's log out</i> .....	40

<b>11. MAINTENANCE</b> .....	<b>41</b>
11.1. Cleaning and disinfection of the unit.....	41
11.1.1. <i>Cleaning</i> .....	42
11.1.2. <i>Disinfection</i> .....	42
11.2. Vacuum unit maintenance.....	43
11.3. Refrigerating system maintenance .....	45
11.4. Maintenance table.....	45
<b>12. TROUBLESHOOTING</b> .....	<b>46</b>
12.1. Vacuum unit.....	46
12.2. Refrigerating unit.....	47
12.3. Freeze-dryer .....	47
<b>13. ANNEXES</b> .....	<b>48</b>
13.1. Equipment spare part codes .....	48
13.2. Examples of recipes .....	49
13.3. Maintenance records .....	50
13.4. Decontamination declaration .....	52
13.5. Refrigerating system P&ID (LyoQuest -55) .....	53
13.6. Refrigerating system P&ID (LyoQuest -85) .....	54
13.7. Refrigerating system P&ID (LyoQuest Arctic -85) .....	55
13.8. Vacuum system P&ID .....	56
13.9. Declaration of conformity .....	57
13.10. Accessories configuration .....	58

Translation of the original manual

**1. SYMBOLS, ABBREVIATIONS AND STICKERS**


Abbreviations used in this manual:

V	Volts	s, h	Seconds, hours
A	Amperes	Hz	Hertz
W, kW	Watts, kilowatts	dB	Decibel
mm, cm, m	Millimetres, centimetres, metres	°C	Celsius
kg	Kilogram	GLP	Good Laboratory Practice

The European Union issued the Directive 2012/19/EU on Waste of Electrical and Electronic Equipment (WEEE).

This product is required to comply with the European Union's Waste Electrical & Electronic Equipment (WEEE) Directive. Mark symbol is shown at the right.

Azbil Telstar sells products through distributors throughout Europe. Contact your local distributor for recycling/disposal.

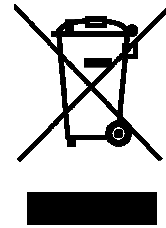The name plate & CE mark sticker is located at one side of the freeze-dryer.

## **2. GENERAL DATA**

### **2.1. Aim of this manual**

The aim of this manual is to give the user guidelines for the handling and operation of the **LyoQuest / LyoQuest Arctic laboratory freeze dryers**. This equipment is designed to be a multipurpose unit in the research and development field, of biological, pharmaceutical and food products. This instruction manual is valid for the following models:

**LyoQuest -55 / LyoQuest -55 Plus**

**LyoQuest -85 / LyoQuest -85 Plus**

**LyoQuest Arctic -85 Plus**

#### **WARNING**

- This manual must be read and fully understood before its installation, maintenance and operation of the equipment.
- All safety instructions contained herein must be respected.
- Any use or application failing to comply with the recommendations expressed in this manual may cause danger and loss of rights to manufacturing guarantee as well as freeing the manufacturer from all responsibility arising.
- Before carrying out any maintenance work on the equipment, it must be disconnected from the energy source.

### **2.2. Application**

This model of freeze-dryer can be used for the freeze-drying in flasks, vials, ampoules or in bulk, by means of the different configurations of the basic unit and the possibility of coupling several accessories.

#### **WARNING**

- This freeze dryer is NOT a sterilizable unit.
- This freeze dryer consists of methacrylate parts. Not following the instructions about the disinfection of the unit (chapter 11.1) or in the event of freeze-drying methacrylate corrosive products can damage the parts, being excluded of the guarantee of the unit.
- Under no circumstances may the unit be used with corrosive, flammable or explosive gases.
- In the event of freeze-drying products containing corrosive or explosive solvents, please contact the Customer Service, to take the specific preventive actions for the specific product.
- Certain freeze-drying cycles increase the required preventive maintenance actions on the vacuum pump. See chapter 11.2

#### **MANUFACTURER: AZBIL TELSTAR TECHNOLOGIES, S.L.U.**

Av. Font i Sagué, 55  
Parc Científic i Tecnològic Orbital 40  
08227 TERRASSA (SPAIN)  
Telf (+34) 93 736 16 00  
Fax (+34) 93 786 13 80  
e-mail: [telstar@telstar.com](mailto:telstar@telstar.com)  
<http://www.telstar.com>

### 3. DESCRIPTION OF THE FREEZE DRYER

#### 3.1. Technical features

<b>VACUUM PUMP (optional)</b>		ULVAC version	ILMVAC version
Dimensions (length x width x height)	mm	487.5x170x249.5	500 x 300 x 340
Weight	kg	27	25.4
Vacuum at the pump inlet	mbar	See user manual of the pump	
Motor power	kW	0.4	0.4

<b>ICE CONDENSER</b>		<b>LQ -55</b>	<b>LQ -85</b>	<b>LQ Arctic -85</b>
Condensation surface	m <sup>2</sup>	0.25		0.25
Minimum temperature (T≤+21°C)	°C	-55	-85	-85
Ice capacity	kg	8 (4 kg/24 h)	8 (6 kg/24 h)	8 (6 kg/24 h)
Material		AISI – 316L		AISI – 316L

<b>REFRIGERATING UNIT</b>		<b>LQ -55</b>	<b>LQ -85</b>	<b>LQ Arctic -85</b>
Number of compressors		1	2	2
hermetic type		One stage		One stage
Refrigeration power (at 35°C)	W	495	800	800
Refrigerant (free of CFCs)		R-507	R-507/ R-23	R-290/ R-170
Refrigeration of the unit		By air circulation		

<b>EQUIPMENT</b>		<b>LQ -55</b>	<b>LQ -85</b>	<b>LQ Arctic -85</b>
Process vacuum. Unloaded, T≤ 21°C, condenser temperature: -55°C (LyoQuest -55), condenser temperature: -85°C (LyoQuest -85), condenser temperature: -85°C (LyoQuest Arctic -85)	mbar	< 5 x 10 <sup>-2</sup>	< 3 x 10 <sup>-2</sup>	< 3 x 10 <sup>-2</sup>
Dimensions (height x width x depth)	mm	480 x 425 x 575		480 x 506 x 575
Weight	kg	65	70	71

## **4. INSTALLATION INSTRUCTIONS**

### **4.1. Transport and unpacking**

When our products leave our warehouse, they have been tested thoroughly and packed with precision. At delivery, examine the packing and the handling indicators meanwhile the forwarder is present.

The freeze-dryer is delivered together with the optional accessories and documentation in a single pack; it is packed in a suitable way to prevent it from being damaged during transportation.

The unit has to be shipped in vertical position according to the crate indications (adhesive labels).

The pack includes a tube of silicone grease, to carry-out the maintenance of the equipment seals. This manual is included in the scope of supply.

When unpacking the freeze-dryer, check if the supply includes accessories, separately packed to ensure suitable protection.

Red ShockWatch or Tiltwatch indicators means impact or unacceptable tilting. In this case mark this on the shipping documentation while the forwarder is present and state an official complaint to the forwarding company.

If there is no visible damage to the unit, unpack inspect (and test) the equipment within 8 days of delivery. Do not return goods without notifying and authorization.

### **4.2. Placing on site**

The unit must be placed on a table that has a smooth and horizontal surface. The basic unit must be firmly settled to prevent vibrations.

#### **WARNING**

- When assembling the chamber to the LyoQuest / LyoQuest Arctic, do not forget to assemble the methacrylate base. It is required for proper working when chamber is installed.
- The optimal working conditions are:
  - Temperature: < 25°C (recommended 23°C±2°C)
  - Relative humidity: 50%±20%

If any of these conditions is not fulfilled, optimum equipment performance cannot be guaranteed.

### 4.3. Electrical connection


- a) Ensure that the power supply is the correct one for the unit.
- b) Ensure that the connection to the power supply suits the freeze-dryer's power.
- c) Connect the electrical cable of the unit to the power supply.

**WARNING** 

- Before switching the Freeze-Dryer on, check that the supply voltage is the same as the one mentioned on the specification plate. It should be connected to a suitable socket base with a grounded plug.
- Whenever there is a power cut or the supply fails for a long period of time, stop the freeze-dryer.
- Installation in the place where the Freeze-Dryer is connected to must meet the safety regulations in force and be according to the local regulations of the country or customer facilities.

Connect the cable of the vacuum pump to the unit.

Connect the electrical cable of the basic unit to the power supply (example for plug wire supplied with a European plug).

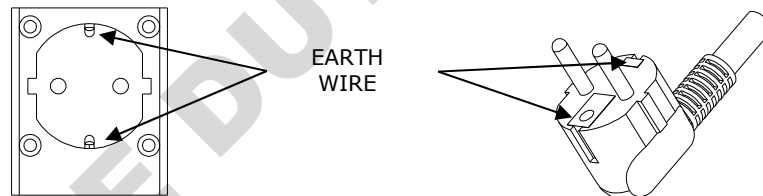Electrical supply: 230V monophase + PE, 50 Hz  
230V monophase + PE, 60 Hz

Installed power (\*): LyoQuest -55: 2.5 kW  
LyoQuest -85: 3.5 kW  
LyoQuest Arctic -85: 3.5 kW

Consumption (\*): LyoQuest -55: 8.6A (11A max)  
LyoQuest -85: 10.6A (15A max)  
LyoQuest Arctic -85: 10.6A (15A max)

\* 230V 50 Hz

## 5. **BEFORE START-UP**

Before proceeding with the start-up of the freeze-dryer the following instructions must be followed step by step.

### **WARNING**


- **Prior to begin start up, electrical connection must meet all requirements according to chapter 4.3.**

### **NOTE**

Before beginning a new operation:

- Open the drain valve in order to check that all water from the last operation has been totally eliminated. Close the valve again.
- Check that the condenser is covered by its transparent cover (manifold or chamber).
- Check that all valves of the chamber are closed correctly.

Drain valve (open)
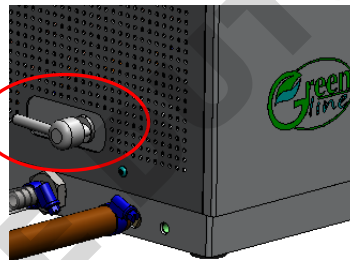
### 5.1. **Vacuum pump**

The pump is supplied with the oil loaded. Verify the oil level (see section Vacuum unit maintenance).


In order to prevent the non-condensable vapors of the product from contaminating the oil of the pump, the system has Gas-Ballast function (see section Vacuum unit maintenance). When the pump is switched off, a closing valve isolates the pump from the freeze-dryer, preserving the vacuum.

The vacuum pump is equipped with an exhaust filter, whose function is to retain oil vapors.

The electrical connection of the vacuum pump has to be connected to the base placed at the side of the unit. Furthermore, the vacuum pump admission has to be connected to the basic unit through the supplied pipe.

### Connection of the vacuum pump expulsion:

To direct the vapors outside of the pump, connect the pipe to vacuum pump exhaust, removing the installed filter.

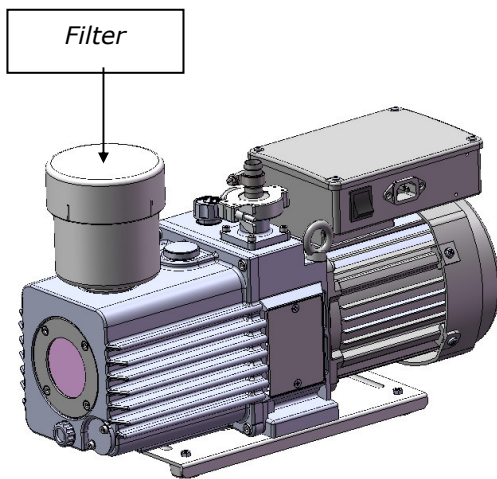**ULVAC pump**

**ILMVAC pump**

### 5.2. Refrigerating unit

In order to assure a suitable cooling of the refrigerating unit (via forced air circulation by means of a fan), all vents must be separated by at least 20cm from other equipment or the walls.

## 6. DETAILS OF THE FREEZE-DRYER

### 6.1. Description

The LyoQuest / LyoQuest Arctic freeze-dryers are constructed as a compact desktop version. The basic unit consists of a metallic box with lateral access doors to perform maintenance tasks. The vacuum pump is supplied separately and is optional. The control panel is located at the front.

A refrigerating system cools the condenser coil to the working temperature. The condensation surface has been dimensioned ensuring the condensation of the water vapor released by the products that are being freeze dried. This prevents vapors from entering the vacuum pump. The refrigerating system is formed by:

- Compressor:

LyoQuest -55: one-stage hermetic compressor, condensed by air.

LyoQuest -85 / LyoQuest Arctic -85: two one-stage hermetic compressors, condensed by air.

- Condenser:

Material: AISI 316L stainless steel.

#### **Models:**

- LyoQuest: Prepared to work with all accessories, but not with heated plates.
- LyoQuest **Plus** / LyoQuest Arctic **Plus**: Prepared to work with all accessories and heated plates.

## **7. RECOMMENDATIONS**

### **7.1. Ice condenser**

Before each operation it is necessary to open the water drain valve of the condenser in order to assure that the ice condensed during the previous operation has been totally eliminated. If water residues remain in the condenser or in the drainpipe, they will be frozen as a result of the vacuum, thus possibly compromising work conditions.

### **7.2. Greasing**

Slightly grease the airtight seals with silicone grease. The contact surfaces should always be kept clean and greased. Avoid dirt accumulation.

### **7.3. Long storage**

#### **WARNING**

- In case of extended shutdown of the freeze-dryer, the vacuum pump oil must be changed to avoid any corrosive waste that could damage it. Disconnect the unit from the electricity supply. Clean all the components of the equipment (manifold, chamber, ...). Clean the inside of the chamber, check that the condenser has been completely emptied and lubricate the airtight gaskets of the chamber, valves and condenser.

## 8. WORKING PRINCIPLE

### 8.1. Freeze

The product inserted into the equipment can either be previously frozen or can be frozen in the ice condenser. To freeze it in the condenser, the product is placed into the plate stack (optional accessory), and the condenser is cooled selecting the option "freezing" on the touch screen, both in semiautomatic mode and in automatic mode.

The freeze temperature is proportional to the solidification point of product. Freezing speed influences the size of the crystals and, consequently, sublimation speed. Therefore, the most suitable freeze parameters for each product must be determined. It is advisable to freeze products at a temperature below  $-20^{\circ}\text{C}$  for a period of 3-4 hours.

### 8.2. Lyophilization


#### 8.2.1. In manifold

It is used to lyophilize bulk products in glass flasks with a grinding mouth ( $29/32\text{mm } \varnothing$ ). The products contained in the flasks (filling a quarter up to maximum half the flask) have to be previously congealed in a rotating bath of alcohol, freezing the product to the walls of the flask.

When the condenser has reached the programmed temperature, the following steps must be performed:

1. Connect the vacuum unit.
2. Connect the flasks or bottles (with previously frozen product) to the valves connected to the condenser.
3. Open the valves.

Then the vacuum is done inside flasks, and the sublimation of the product starts. A pressure of less than 10mbar is reached within 10 minutes. If this does not happen, the pump will stop, and the alarm of the vacuum system will appear. In this case the valve-flask connections must be reviewed in order to assure that there is no leak.

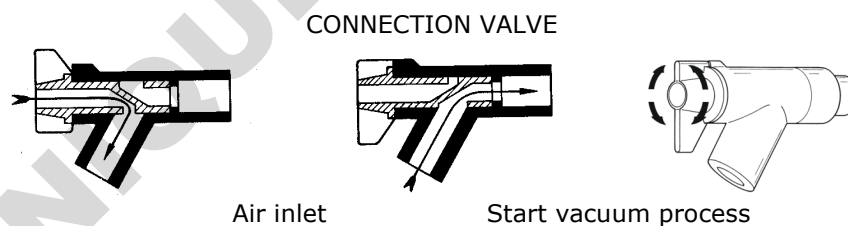If it is necessary to connect additional flasks when the freeze-drying process has already begun, the following steps must be performed:

1. Close all valves.
2. Connect the new flask and open its corresponding valve.
3. Wait 10-15 seconds in order for the pressure to be restored in the unit.
4. Open the valves of the other flasks.

It is possible to provide external heat by means of an electric heat source, in order to speed up the drying process. The heat source should be at the suitable distance to apply the radiation heat that the product requires. A hot air heater directed directly to the flasks may also be used.

At the beginning of the process a layer of ice will form on the outside wall of the flask; this will gradually disappear as drying out continues. When the exterior of the flask is at ambient temperature, no ice remains in the product (the process has finished). Once the process has finished, the corresponding valves of the manifold are closed, and air enters in the flasks. So, they can be removed.

### 8.2.2. In chamber

The product may be freeze-dried inside the condenser or in the chamber:

- **Inside the condenser** (see Fig. A): It is advisable to load the product only on the lower and middle plates. The higher plates do not receive cold as do the other plates because the vacuum supply is placed at top of the condenser. The equipment allows you to program recipes (see chapter 10.6) to run a complete cycle without making any manual actions.
- **Outside the condenser:** This is recommended when the quantity of product allows you to load the three plates. It is advisable to freeze the product inside the condenser because the heated plates will be cold and heating rates from negative temperatures may be reached. Once the product is frozen, remove the plates of the condenser and place them in the chamber (see Fig. B). Also, it is possible to load a product that has been frozen in other equipment. The equipment allows you to program recipes (see chapter 10.6) to run a complete cycle without making any manual actions once the product is frozen and the chamber has been placed.

#### NOTE

- The product layer should be as thin as possible (recommended height <20mm). The time of the freeze-drying process depends for an important part from this height. It is very advisable to load the freeze-dryer uniformly and that the height of product is the same in all zones.
- The time of the process also depends on:
  - The heat provided to product (external if the manifold is used or by means of an electrical resistor placed inside the plates).
  - The maximum temperature accepted by the product.
  - The final degree of humidity desired.
  - The percentage of solid material in the product.
  - The vacuum degree obtained.
  - The total amount of product to be dried each time.
  - The temperature of the condenser.

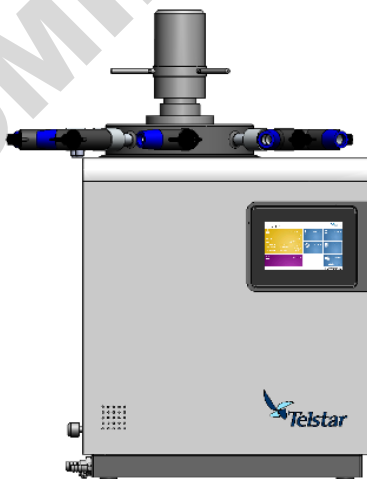
Fig. A

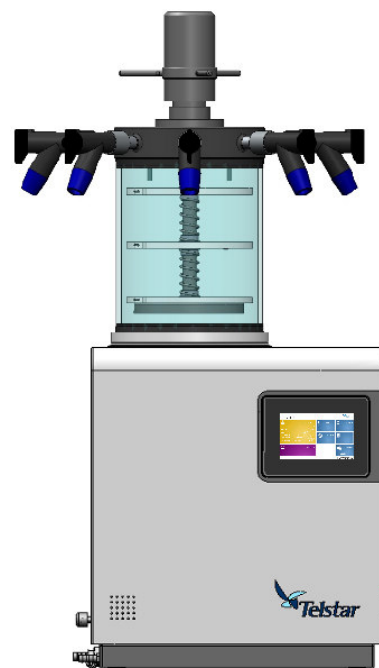
Fig. B

### 8.2.3. End of the process

Stop the pump and the refrigerating circuit. Open the air inlet of the valve on the side of the chamber. The air inlet has to be opened very slowly to avoid the formation of turbulences that will spread about dust in the flasks or plates. Then, the defrosting function is activated.

At the end of the process, the pump and the refrigerating unit will stop (the compressor stop is delayed in order to extract all gases from the circuit). Next, a connection valve of the manifold or the venting valve placed at the top cover has to be opened in order to ventilate the inside of the equipment.

It is necessary to connect a container to the outlet of the outflow pipe. Once this is done the manual outflow valve placed at the left side can be opened. Start defrosting mode (see chapter 10.3.1) in order to heat the surface of the condenser which allows you to extract the ice block in a few minutes. Also, it is possible to heat the condenser until the entire ice block has been defrosted and all the water has been drained through the outflow valve. In both cases it is advisable to remove the chamber or the manifold to accelerate the defrosting.

Before the start of a new operation, check that the condenser has been totally defrosted, that there are no water traces in the overflow pipe, and that the overflow valve is closed again.

### 8.3. Lyophilization modes

There are two possible lyophilization modes: the AUTOMATIC mode (see chapter 10.2) and the SEMIAUTOMATIC mode (see chapter 10.2.3). The automatic mode allows you to establish the processes, the parameters and the time by means of a RECIPE (see chapter 10.6).

### 8.4. Selectable processes

**Freezing:** It starts up the refrigerating system and cools the condenser. The cooling system requires a preparation time (5 minutes approx.) in which the temperature may increase slightly. After this the temperature will quickly decrease. This step must be longer than 30 min to make sure that the system can reach the minimum temperature.

**Cool+Vacuum:** It starts up the refrigerating system, cools the condenser and starts up the vacuum unit. The vacuum begins when the temperature of the condenser is lower than the established parameter (normally -40°C). The control of the vacuum level will depend on the programmed Set Point. The system will need a time longer than 20 min to reach the minimum vacuum level.

**Heat plates:** It starts up the electrical resistors of the plates. They warm up until the programmed temperature has been reached. Also, it starts up the vacuum unit. The vacuum begins when the temperature of the condenser is lower than the established parameter (normally -40°C). The control of the vacuum level will depend on the programmed Set Point.

Each one of the processes may have parameters to set. In the automatic cycle it is also necessary to set the time for each step. The parameters are the following ones:

	Semiautomatic	Automatic
<b>Freezing</b>	-	Time (HH:mm)
<b>Cool + vacuum</b>	Vacuum (mBar)	Vacuum (mBar)
		Time (HH:mm)
<b>Heat plates</b>	Vacuum (mBar) Plates temperature (°C)	Vacuum (mBar)
		Plates temperature (°C)
		Time (HH:mm)

Depending on the temperature of the condenser probe, the electrically operated valve of the refrigerating circuit may take a certain time to start working. As a result, the temperature may increase some degrees initially.

The pump starts to work when the Set Point has been reached. The vacuum valve will open few seconds later.

## 9. OPTIONAL ACCESSORIES

(Compilation at the end of the manual)

### 9.1. Chambers

Model	Description	Diameter	Distance between plates
C220 SP3	Cylindrical chamber equipped with 3 standard plates	220mm	70 mm standard
C220 HP3	Cylindrical chamber equipped with 3 heated plates	220mm	70 mm standard / 140 mm large spacer (LS)
C220 HP3 S	Cylindrical chamber equipped with 3 heated plates and vial stoppering device	220mm	
C220 HP3 M8	Cylindrical chamber equipped with 3 heated plates and 8 inlets manifold	220mm	
C220 HP3 S M8	Cylindrical chamber equipped with 3 heated plates, vial stoppering device and 8 inlets manifold	220mm	
C350 HP3	Big size cylindrical chamber equipped with 3 heated plates	350mm	71 mm standard / 142 mm large spacer (LS)
C350 HP3 S	Big size cylindrical chamber equipped with 3 heated plates and vial stoppering device	350mm	
C350 HP4	Big size cylindrical chamber equipped with 4 heated plates	350mm	71 mm standard
C350 HP4 S	Big size cylindrical chamber equipped with 4 heated plates and vial stoppering device	350mm	


Material: Methacrylate (PMMA).

220mm versions optionally available in Pyrex®.

The (under vacuum) vial stoppering devices are operated manually, and they require vials provided with closure.

C220 SP3

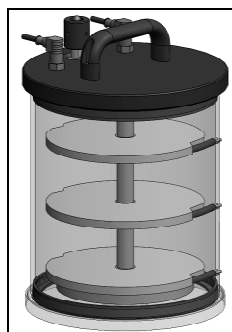C220 HP3


C220 HP3 S


C220 HP3 M8


C220 HP3 S M8

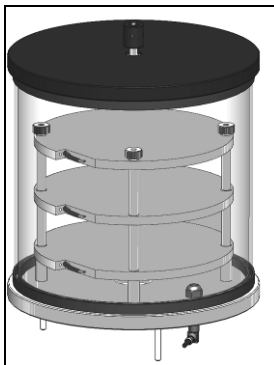
C350 HP3

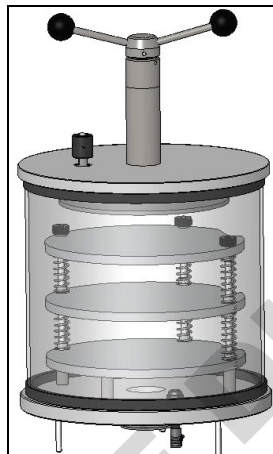C350 HP3 S

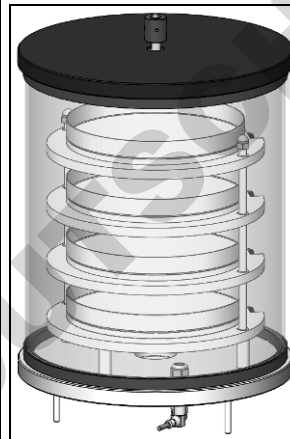C350 HP4

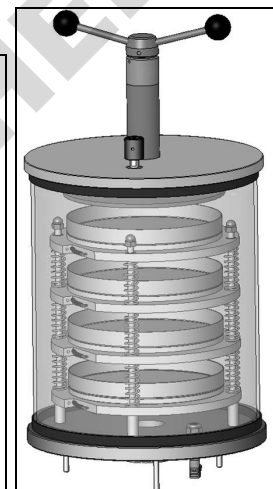C350 HP4 S

Model	Description	Standard plates	Distance between plates
CC SP	Cubic chamber	2	145 mm
		4	70 mm
		6	56 mm

Material: Stainless steel AISI 316  
 Material: Methacrylate (PMMA).  
 Internal dimensions: 484x340x376 mm  
 External dimensions: 542X460X617 mm  
 (height x width x depth)

CC SP

## 9.2. Manifolds

### Manifolds with 8-12-16 inlets

Horizontal manifold. Each one of the manifold inlets is provided with a 3-way valve which allows isolating and restoring the pressure in each one of the flasks independently. The valves allow the junction to grinding flasks by means of the conical adapter that they incorporate (external diameter: 29-32 mm; internal diameter 18 mm). By removing this adapter, tubes or flasks can be inserted directly in the valve.

Circular manifold with 8 inlets (fig C) and straight one with 12 inlets (fig D).

By joining together two 8 inlets manifolds it is possible to get 16 inlets (fig E).

Fig. C

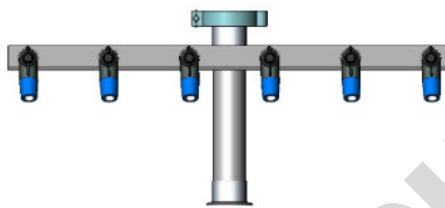
Fig. D

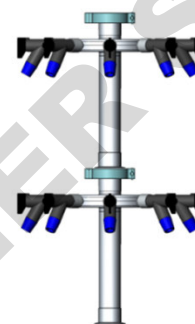
Fig. E

### Manifolds with 40-80 inlets

Manifold with 40 inlets for 13mm diameter pipes (fig F). By means of a distributor it is possible to get 80 inlets (fig G).

Fig. F

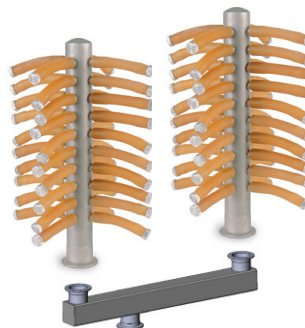
Fig. G

## 9.3. Flasks condensations collector (Drip tray)

Cover to collect the condensations generated on the outer part of the flasks of the manifold (Fig. H).

Fig. H

#### 9.4. Lyologger software

Software for monitoring on PC that provides data for the vacuum and the temperatures of the condenser, plates and product.

The equipment is connected by a RS-232 serial port on PC and it sends information in real time.

A data packet is sent every second in order to know the freeze dryer state.


The serial port is automatically configured by the software.

The dataflow is unidirectional from the freeze dryer to the PC. It means that it is only possible to monitor, and not to control, the equipment connected to the PC.

Minimum computer requirements:

- Processor Intel Pentium II or AMD Athlon 32 bits.
- 128 megabytes RAM.
- Windows 98SE.
- Hard disk IDE fast.
- Graphic board SVGA (color 16 bits at 800 x 600 Pixels).

See instructions in the User's guide of LyoLogger by clicking the software help button.
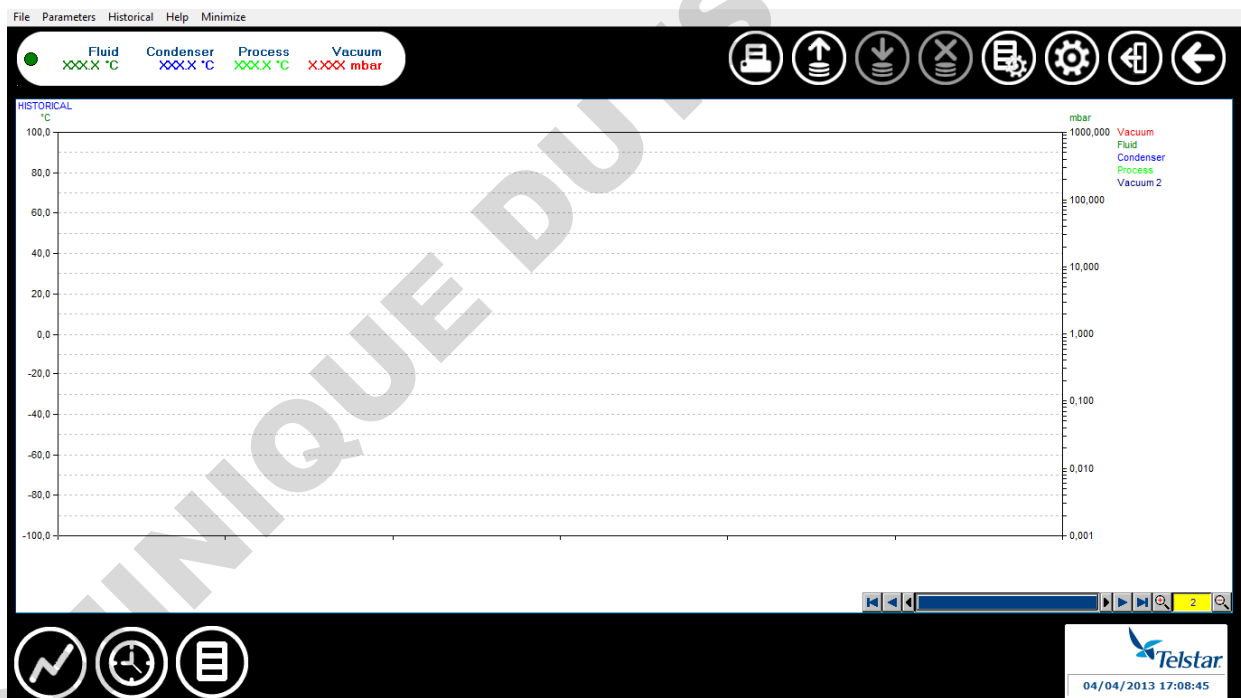
## 10. OPERATION AND CONTROL INSTRUCTIONS

### 10.1. Start up



To start up the equipment, press the switch located at the right of the equipment, moving it to the position I.

The equipment will start up and the introduction screen will appear. Press the screen to access to main screen of the application.
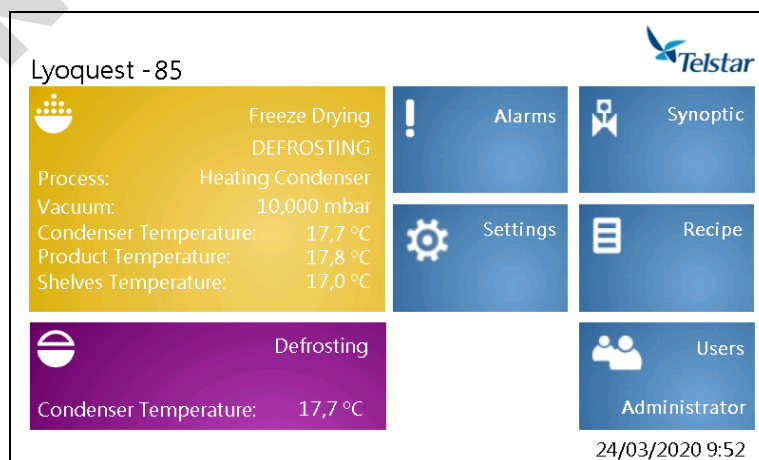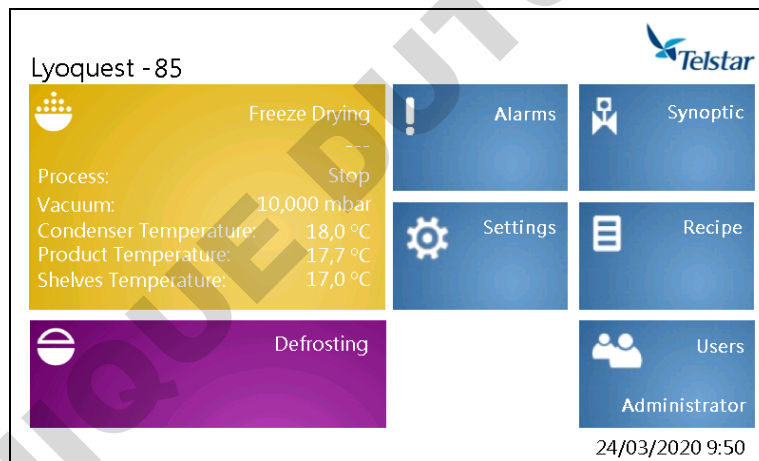
#### 10.1.1. Main window

Once in this menu, the user can navigate through their screens by clicking over each icon.

-  **Freeze drying.** This box also shows the process in execution and real time values of the equipment sensors.
-  **Defrosting.** When this process is active, the condenser temperature is displayed in this button.
-  **Alarms.** Exclamation symbol will turn into red when any alarm has been triggered in the equipment.
-  **Synoptic.** Access to see the status of the equipment components.
-  **Settings.** They can be seen and modify them, when applicable.
-  **Recipe.** Access to see and edit them.
-  **Users.** Access to manage them. This box also shows the current logged user.

DOMINIQUE DUTSCHER SAS

## 10.2. Freeze-drying

When pressing this button, placed at the main window, user can access the *Freeze-drying* window.

### Freeze drying

Recipe	Receta1		    
Process	Stop	Step 0	
Step time	SP - : -	PV - : - : -	
Vacuum	- mbar a	10,000 mbar a	
Shelves Temperature	- °C	17,0 °C	
Condenser Temperature	17,8 °C		

It consists of:

- info of the recipe in execution in automatic mode, shown in the central part:

**Recipe:** Name of the loaded recipe.

**Process:** Process under execution.

**Step:** Step number of the executing recipe.

**Step time (SP):** Time Set point for recipe step.

**Step time (PV):** Current elapsed time.

**Vacuum (SP):** Vacuum Set point of the recipe.

**Vacuum (PV):** Current vacuum value.

**Shelves temperature (SP):** Recipe temperature setpoint.

**Shelves temperature (PV):** Current shelves temperature value.

**Condenser temperature:** Current condenser temperature value.

- Buttons located at the lower and right part of the screen. Depending on level Access, user will visualize some of them:

Automatic cycle activation.

Semiautomatic cycle activation.

Window to visualize the values of the probes.

Window to select and visualize the required recipe to be executed.

Return to the previous window.

The following buttons appear when Automatic cycle is on:

By pressing this button, the PLC loaded recipe will start.

By pressing this button, current recipe stops and Automatic mode turns off.

The following button appears when Semiautomatic cycle is on:

By pressing this button, the window to execute processes in semiautomatic mode will appear.

### 10.2.1. Probes

 When pressing this button, placed at the *Freeze-drying* window, user can access the *probes* window.

The value of each probe of the equipment are displayed.  
Shelves and product temperature only appear when shelf option is active.

**Sensors**

Condenser Temperature		17,7	°C
Vacuum		10,000	mbar
Shelf Temperature		17,0	°C
Product Temperature		17,7	°C

### 10.2.2. Recipe

#### 10.2.2.1. LOADED RECIPE

 When pressing this button, placed at the *Freeze-drying* window, user can access the *Loaded recipe* (active) window.

The recipe currently loaded in the PLC is displayed. It will be executed when Automatic mode is on.

**Loaded Recipe (Steps 1 to 5)**

Number	Recipe			
1	Receta1			
Step	Process	°C	mbar	hh:mm
1	Cool + Vacuum		0,500	1 : 0
2	Heat Shelves	10,0	0,001	0 : 30
3	Heat Shelves	30,0	0,001	1 : 0
4	Freezing			0 : 30
5	Stop			

Data of °C, mbar and hh:mm are only shown when applicable.

By pressing this button, recipe steps 6 to 10 will be displayed.

Open the window for loading the recipe. The recipes saved in the equipment will appear (4 possible recipes).


Return to the previous window.

### 10.2.2.2. LOAD RECIPE

Press this button to access the *Load Recipe* window:

The recipe to be loaded in the PLC can be selected here. That chosen recipe will be executed when the automatic mode is on.

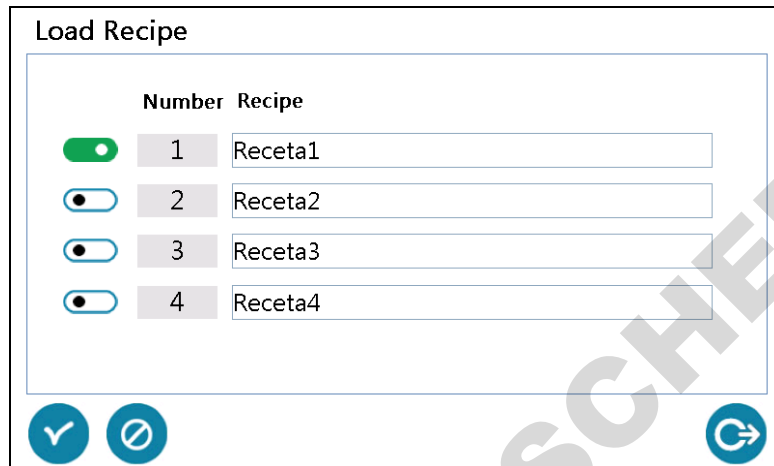Number	Recipe
1	Receta1
2	Receta2
3	Receta3
4	Receta4

Recipe selection button. White if not selected. Green if selected.

Selected recipe confirmation.

Cancel and go to the previous window.

Return to the previous window.

### 10.2.3. Semiautomatic cycle

When pressing this button, placed at the *Freeze-drying* window, user can access the *Semiautomatic cycle* window.

The user can select the chosen process by pressing on the corresponding process box.

It consists of:

- Set point and real time values of the shelf temperature and the condenser vacuum, shown in the upper part:


**Condenser Vac. SP:** Condenser vacuum set point.

**Condenser Vac:** Current condenser vacuum value.

**Shelves temperature SP:** Temperature set point.

**Shelves temperature:** Shelves temperature value in real time.

- The processes selectable by the user (see chapter 8.4), shown in the lower part.

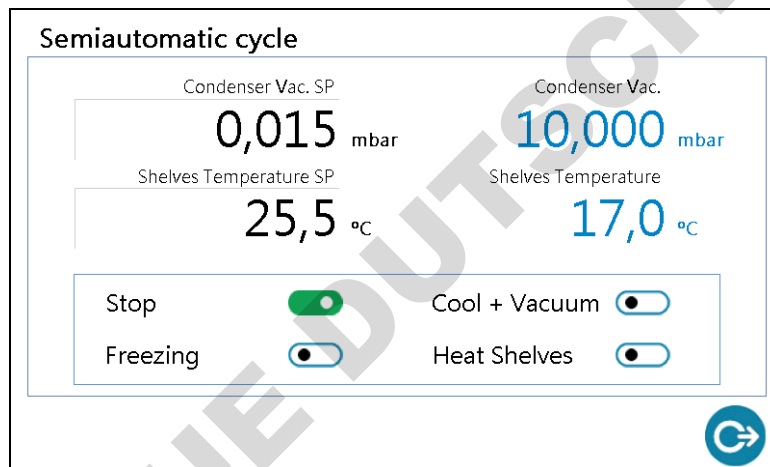Process selection button. White if not selected. Green if selected or under execution.



### 10.3. Auxiliary cycles

#### 10.3.1. Defrosting

When pressing this button, placed at the main window, user can access the *Defrosting* window.

Status of this process (ON/OFF) and condenser temperature are displayed in this screen.

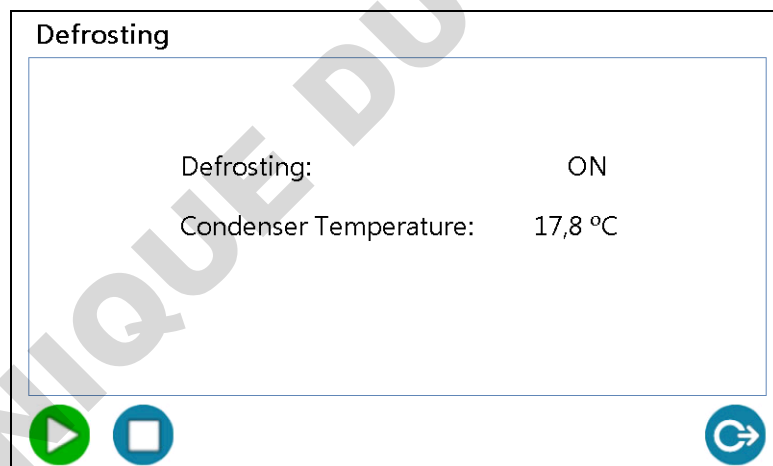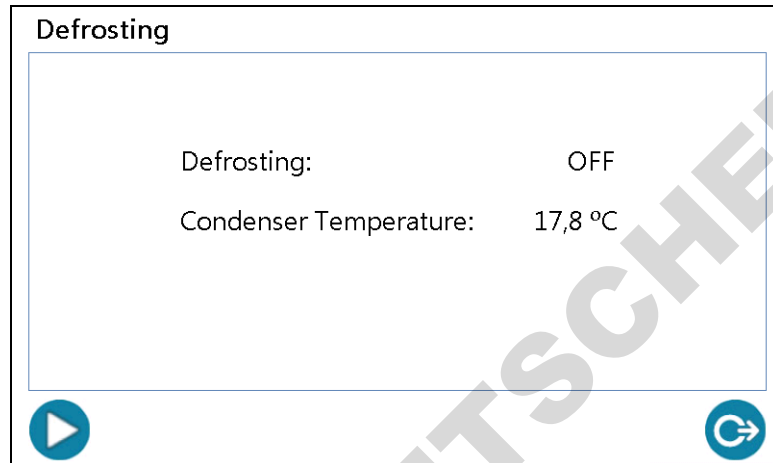By pressing this button, defrosting starts.

By pressing this button, defrosting stops.


Return to the previous window.

### 10.3.2. Alarms

 When pressing this button, placed at the main window, user can access the *Alarms* window.


All alarms will be displayed indicating the status, date and hour of occurrence and a descriptive text.

Active alarms:

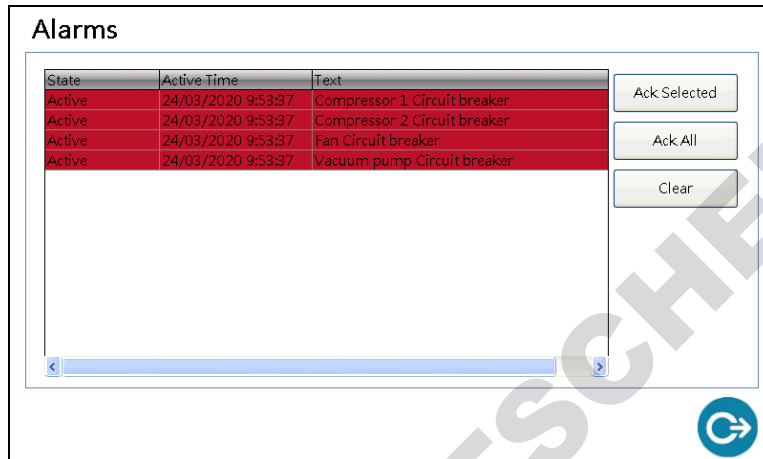
State	Active Time	Text
Active	24/03/2020 9:53:37	Compressor 1 Circuit breaker
Active	24/03/2020 9:53:37	Compressor 2 Circuit breaker
Active	24/03/2020 9:53:37	Fan Circuit breaker
Active	24/03/2020 9:53:37	Vacuum pump Circuit breaker

Muted alarms:

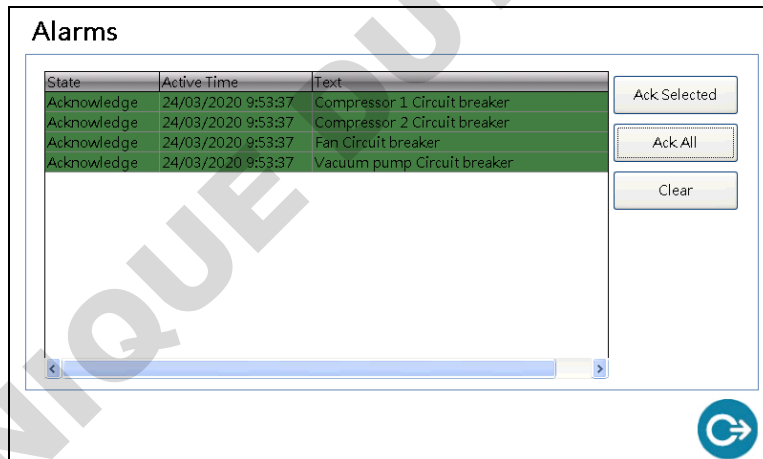
State	Active Time	Text
Acknowledge	24/03/2020 9:53:37	Compressor 1 Circuit breaker
Acknowledge	24/03/2020 9:53:37	Compressor 2 Circuit breaker
Acknowledge	24/03/2020 9:53:37	Fan Circuit breaker
Acknowledge	24/03/2020 9:53:37	Vacuum pump Circuit breaker

Already corrected alarms (problem was solved before Acknowledgement):

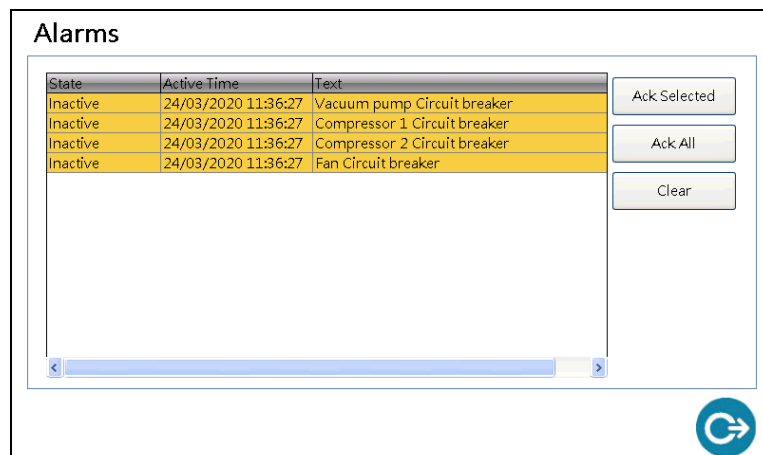State	Active Time	Text
Inactive	24/03/2020 11:36:27	Vacuum pump Circuit breaker
Inactive	24/03/2020 11:36:27	Compressor 1 Circuit breaker
Inactive	24/03/2020 11:36:27	Compressor 2 Circuit breaker
Inactive	24/03/2020 11:36:27	Fan Circuit breaker


Solved alarms (not removed ones):

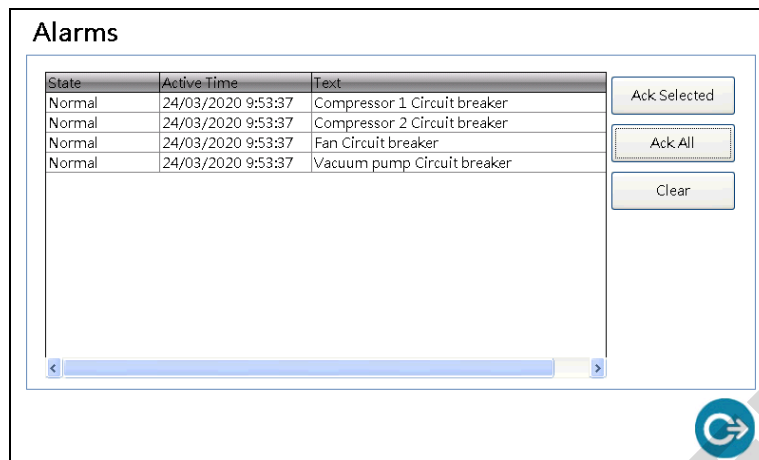
The following buttons appear at the window:

It acknowledges and mutes the selected alarm.

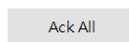It acknowledges and mutes all the alarms at once.

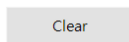It erases all the alarms. But active alarms will keep appearing.

Return to main menu.

### 10.3.3. Alarms list

Nº	MESSAGE	Description
1	Compressor 1 High Pressure	High pressure in compressor 1
2	Condenser Temperature probe fail	Condenser temperature probe
3	Vacuum probe fail	Vacuum probe
4	Compressor 2 MCB	Compressor 2 circuit breaker
5	Fans MCB	Fans circuit breaker
6	Compressor 1 MCB	Compressor 1 circuit breaker
7	Vacuum pump MCB	Vacuum pump circuit breaker
8	Vacuum Fail	Vacuum not reached
9	Shelf Heater CB	Shelf heater circuit breaker
10	Shelf Temperature Probe	Shelf temperature probe
11	Compressor 2 High Pressure	High pressure in compressor 2

#### 10.4. Synoptic


By pressing this button, placed at the main window, the user can access the *Synoptic* window.


The value and status of the equipment components will be displayed.

Depending on the versions of the equipment (with heated shelves or not), the screen will be displayed as follows:

Without shelves:

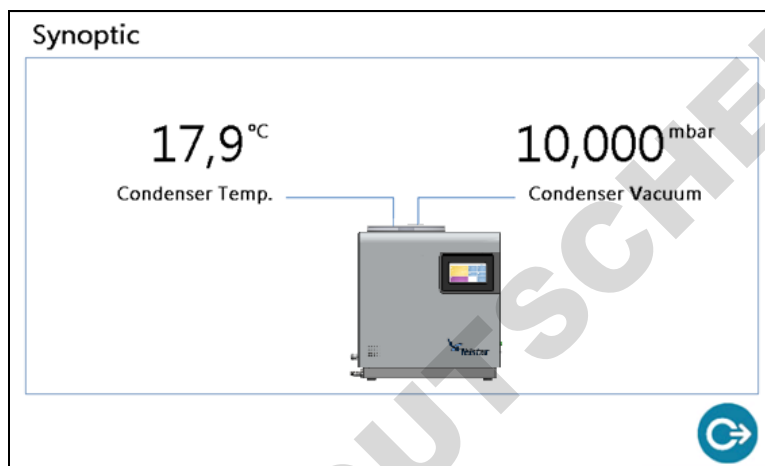With heated shelves:

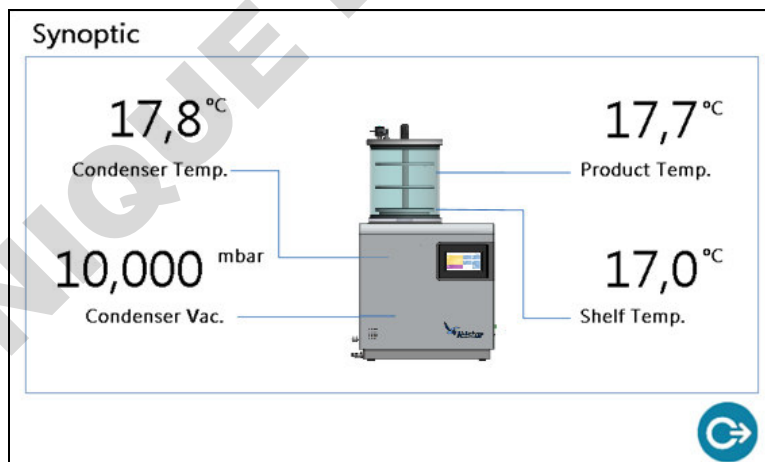The following values will be displayed:

- Condenser temperature
- Condenser vacuum
- Product probe temperature (*only if shelves option is enabled*)
- Shelves probe temperature (*only if shelves option is enabled*)


By pressing over the base unit in the middle of the picture, it will be displayed the equipment cooling system synoptic.

Depending on the base unit version, the corresponding synoptic appears.


Condenser temperature probe and vacuum probe values and status will be displayed.

Status: Green if active. Grey if stopped.


LyoQuest -55:

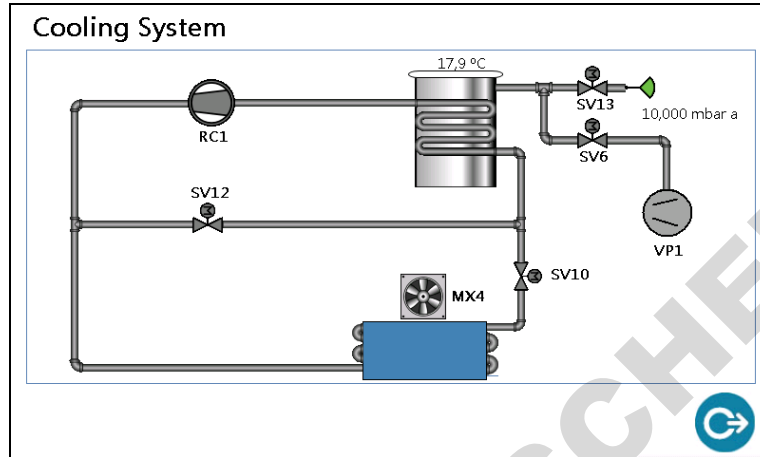LyoQuest -85:

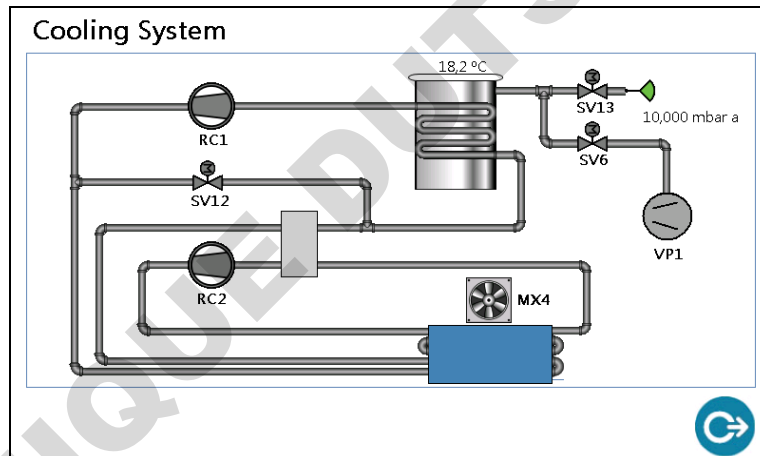LyoQuest -85 Arctic:
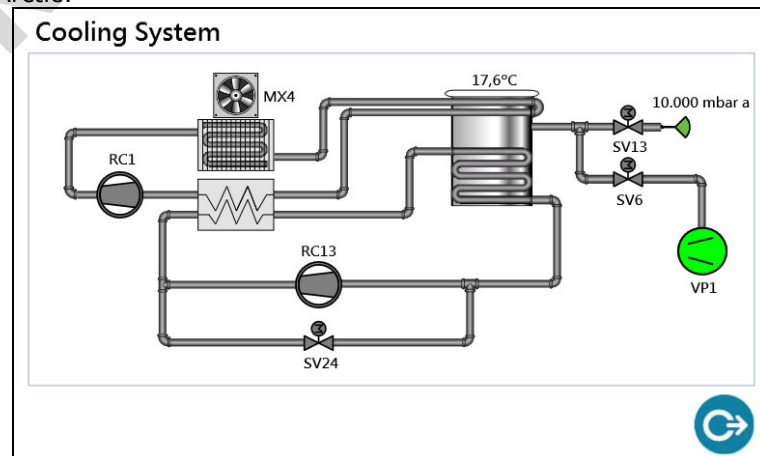
## 10.5. Settings


When pressing this button, placed at the main window, the user can access the *Settings* window:

The equipment and display settings will be shown.

The following actions will be available:

- See the vacuum pump elapsed time (hours).
- Set up the condenser temperature value that determines the vacuum pump start up.
- Set up the manifold/shelves option.
- See the PLC and HMI versions.

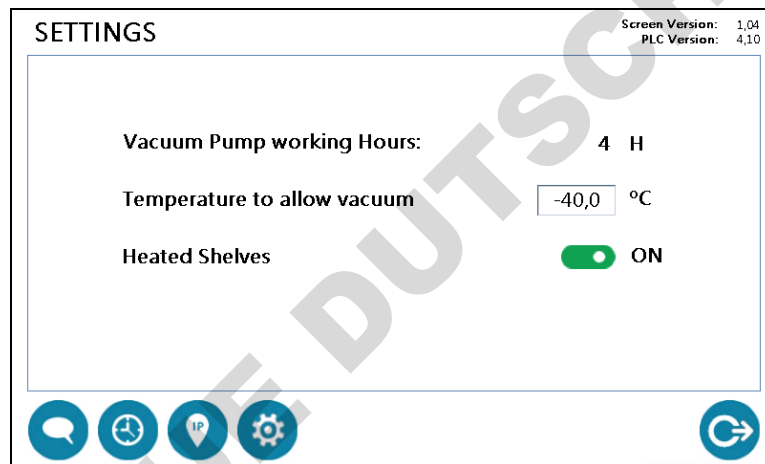Selection button. White if not selected. Green if selected.

- By pressing the icons displayed in the lower area, the user can access additional settings windows:

Languages setting window.

Time/date system setting window.

IP address of the equipment system window (by default 192.168.9.248)

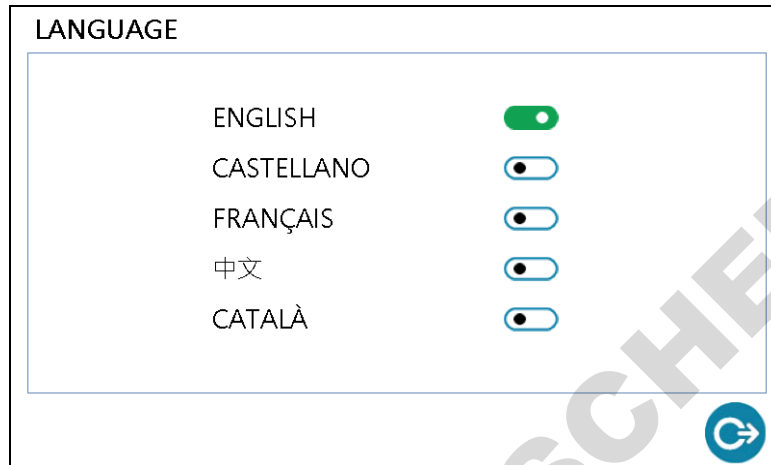
Window for the Technical Service.

Return to main menu.

### 10.5.1. Language

 When pressing this button, placed at the *Settings* window, the user can access the *Language* window:

 Selection button. White if not selected. Green if selected.

### 10.5.2. Date/hour

 When pressing this button, placed at the *Settings* window, the user will be able to set the date/hour of the system.

### 10.5.3. IP Settings

 By pressing this button, placed at the *Settings* windows, the user can configure the local network settings where the equipment will be connected to.


The image shows a dialog box titled "IP Settings" with a "LAN" tab. It contains two radio buttons: "Obtain an IP address via DHCP" (unselected) and "Specify an IP address" (selected). Below these are input fields for "IP Address" (192.168.9.248), "Subnet mask" (255.255.255.0), "Default Gateway", "Primary DNS", and "Secondary DNS". At the bottom are "OK", "Cancel", and "Apply" buttons.

Default IP: 192.168.9.248. Default subnet mask: 255.255.255.0

### 10.5.4. Technical Service

 By pressing this button placed on the *Settings* window, the system opens the window restricted to the Technical Service:

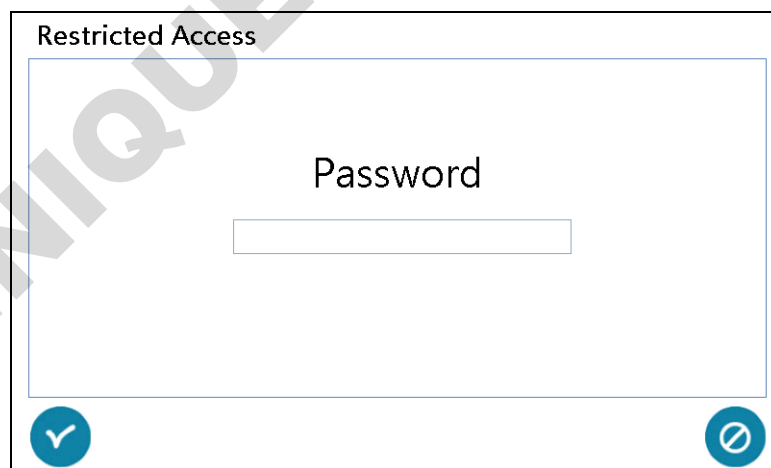Telstar password is required to Access.

The image shows a dialog box titled "Restricted Access". It contains a large text input field with the word "Password" centered above it. At the bottom left is a blue circular button with a checkmark, and at the bottom right is a blue circular button with a slash through it.

-  Confirms the password. If valid, the system will grant the access.
-  Cancel and go to the previous window.

The factory parameters can be configured in this restricted access window. The vacuum pump elapsed time can also be reset.

## 10.6. Recipe edition

 By pressing this button, placed at the main window, the user can access the *Recipe Edition* (steps 01 to 05) window to create and edit the recipes of the equipment.

Recipe can be edited in the central part. Other function icons are found in the lower part of the window. Each recipe consists of 10 steps.

**Recipes Edition (Steps 1 to 5)**

Number	Recipe			
1	Receta1			
Step	Process	°C	mbar	hh:mm
1	Cool + Vacuum		0,500	1 : 0
2	Heat Shelves	10,0	0,001	0 : 30
3	Heat Shelves	30,0	0,001	1 : 0
4	Freezing			0 : 30
5	Stop			

**Recipes Edition (Steps 6 to 10)**

Number	Recipe			
1	Receta1			
Step	Process	°C	mbar	hh:mm
6	Stop			
7	Stop			
8	Stop			
9	Stop			
10	Stop			

**Number:** Recipe number to be edited. (*Non-editable*)

**Recipe:** Recipe name. (*Editable*)

**Step:** Recipe step number.

**Process:** 4 selectable process types (see chapter 8.4): Stop, Freezing, Cool + vacuum and Heat Shelves

**°C:** Shelves temperature set point. (*Editable*)

**mbar:** Condenser vacuum set point. (*Editable*)

**hh:mm:** Hours and minutes that this step selected process will last. (*Editable*)

°C, mbar y hh:mm only appear when the selected process requires such parameters.

Access to the steps 6 to 10 of the recipe edition.

Return to the steps 1 to 5 of the recipe edition.

Access the load recipe edition in order to select the recipe to be edited.

Return to main menu.

Save the edited recipe. The following confirmation window will appear.

Save Recipe

Are you sure?

### 10.6.1. Load recipe

By pressing this button, placed at the *Recipe Edition* window, the user can access the *Load Recipe* window in order to choose the recipe to be edited (4 possible recipes).

Load Recipe

Number		Recipe
<input checked="" type="checkbox"/>	1	Receta1
<input type="checkbox"/>	2	Receta2
<input type="checkbox"/>	3	Receta3
<input type="checkbox"/>	4	Receta4

Recipe selection button. White if not selected. Green if selected.

Selected recipe confirmation and load.

Cancel and go to the previous window.

Return to the previous window.


## 10.7. Users

In order to log in, the user must introduce a password by pressing the *user's* button from the main menu.

It allows:

- User log in
- Change user passwords
- Create and edit users
- User log out
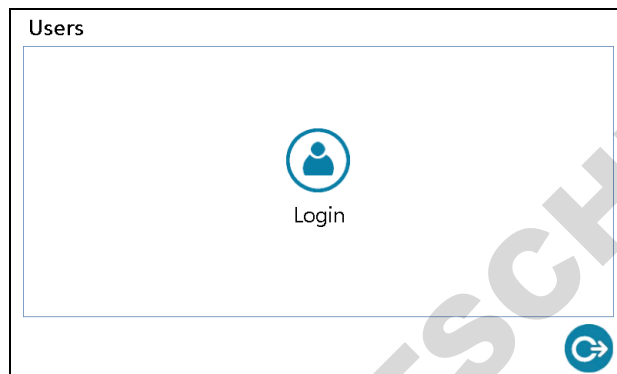
### 10.7.1. User's levels


Two user's levels are defined:

- **Administrators:** they can access all windows of the equipment except the advance settings one.
- **Operators:** They can only access to Freeze-drying, defrosting, alarms, synoptic and users. They are not allowed to create neither change equipment settings.

### 10.7.2. User login

Press this button to access the user's management window in order to select the user and input the password

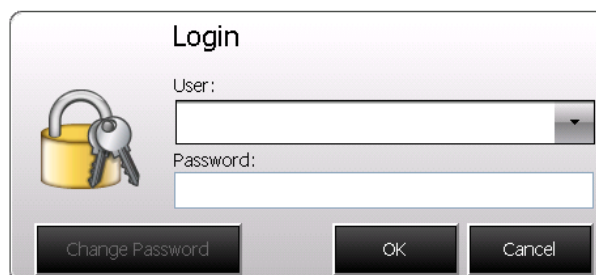**User:** Selection of the user. By default it is defined as Administrator.

**Password:** Introduce the selected user's password. By default the password is "a" for Administrator and "o" for Operator.

Change the password of the selected user.

User and password confirmation. If correct, user is logged in.

Cancel and return to previous window.


### 10.7.3. Password change

Press this button to change the password of the selected user.

To change the password, select the user and then press

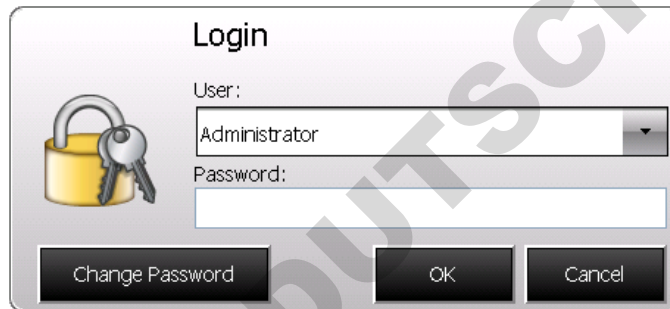The Login dialog box features a yellow padlock icon on the left. It contains a 'User:' label above a dropdown menu with 'Administrator' selected. Below that is a 'Password:' label above an empty text input field. At the bottom, there are three buttons: 'Change Password', 'OK', and 'Cancel'.


A new pop-up window will appear in order to input the old password and create the new password.

The Change Password dialog box features a yellow padlock icon on the left. It contains three text input fields labeled 'Old Password:', 'New Password:', and 'Confirm Password:'. At the bottom, there are three buttons: 'Change', 'OK', and 'Cancel'.

### 10.7.4. Users management

To manage the users of the equipment it is required to log in as administrator.

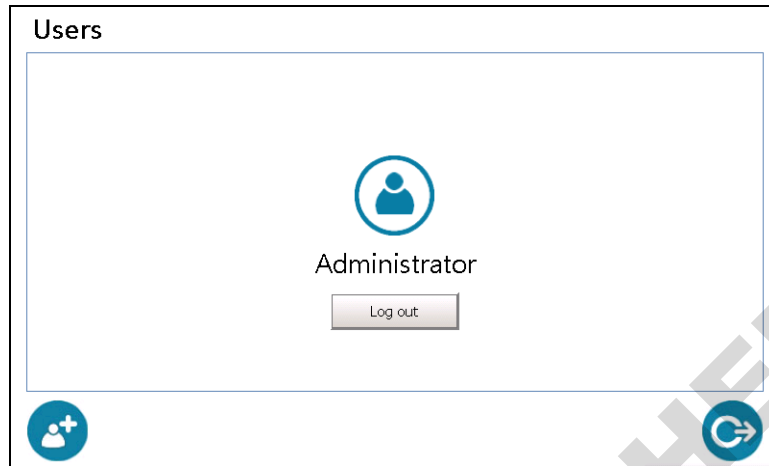A new icon will appear at the left lower part of the screen:  
By pressing that icon, a new pop-up window will allow to manage the users.


The full list of the users will appear. Several icons will be showed at the lower part:


Create a new user.

Edit the selected user.

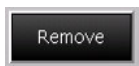
Remove the selected user.

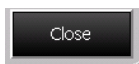
Exit and return to previous window.

### 10.7.4.1. CREATE NEW USER

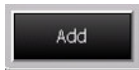Press this button to create a new user.

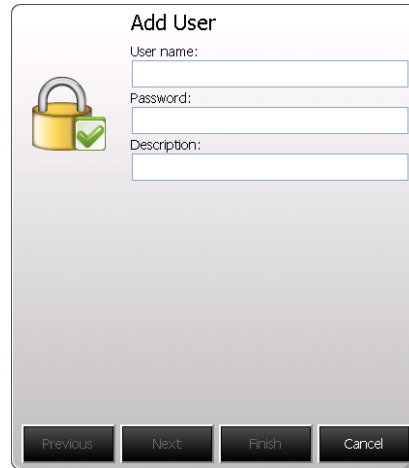The 'Add User' dialog box contains the following fields and controls:

- User name:** Text input field
- Password:** Password input field
- Description:** Text input field
- Navigation buttons:** Previous, Next, Finish, Cancel

**User name:** New user name.  
**Password:** New user password.  
**Description:** User description (optional).


Return to previous window.

Go to next window.


Save new user and exit.

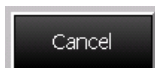Cancel process and exit.

Once data edition is finished, press "next" to go to the next window. Then the desired user level is to be defined.

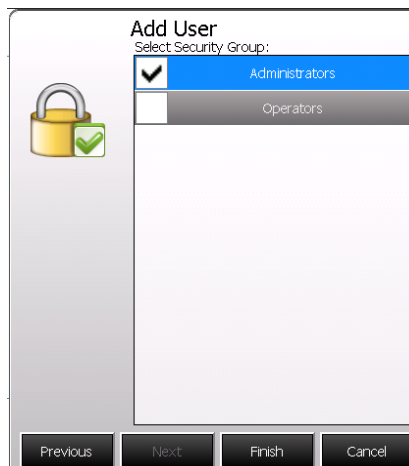The 'Add User' dialog box now shows the 'Select Security Group' step:


- Select Security Group:** A list with two options:
  - Administrators
  - Operators
- Navigation buttons:** Previous, Next, Finish, Cancel

After selecting the desired group, press Finish to save and exit.

### 10.7.5. User's log out

Log out is made in two ways:

- Automatically: After 5 minutes without activity in the equipment.
- Manually: A window will be displayed in the user menu to log out.
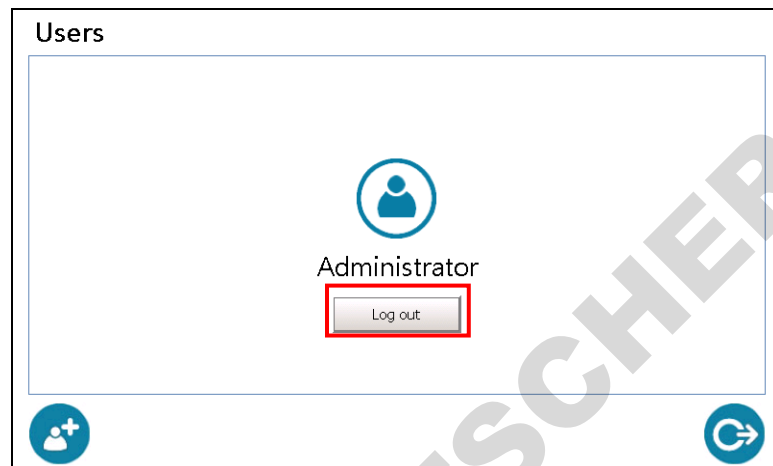
## **11. MAINTENANCE**

The aim of this chapter is to provide a guide for operating on all installed components of the freeze-dryer that require maintenance.

The maintenance plan is based in a series of inspections and operations that must be carried out at intervals depending on the specific task that range from a daily basis up to once every year. The corresponding information should be compiled and noted in detail in order to be used for the following maintenance operations.

### **WARNING**

- All maintenance actions on this equipment must be performed by a qualified technician who is familiar with the proper maintenance procedure required for this equipment. In case that these actions are not carefully carried out, the equipment could work in abnormal conditions.
- Before beginning any maintenance operation or work on the freeze-dryer, disconnect it from the power supply.
- In general, a regular inspection of the correct operation of all the components of the unit should be performed once a year.

### **11.1. Cleaning and disinfection of the unit**

#### **WARNING**

- In case of using any solvent (isopropanol, ethanol, etc.), it is a must to follow these rules:
  - Empty the condenser.
  - Disinfect all the surfaces with a cloth impregnated with a disinfectant product.
  - Dry immediately the surfaces with a clean and dry cloth.
- It is recommended NOT to use disinfectant products containing chlorine or halogens because they can damage the stainless steel and painted surfaces.

#### **WARNING**

- After using solvents in the chamber and condenser, it is a must to ventilate during 24 hours in order to evaporate any remainder. Otherwise the methacrylate surfaces will be damaged due to solvent action.

The normal procedure consists of a previous cleaning prior to disinfection.

For both, unit and accessories, use sterile or non-use disables cloths, which do not shed particles or fibers.

To disinfect use cloths dampened with a disinfecting solution which do not damage or affect the paintwork, stainless steel or chamber glass. A previous clean to disinfecting can also be done with water.

To access to corners or difficult access zones, another type of tools will have to be used.

### 11.1.1. Cleaning

- To be performed at two levels:

- Routine cleaning: To be performed always between two consecutive runs with same product

- Radical cleaning: To be performed always when changing of product or when unit has been out of use for a long time

- According to the level of dirtiness of the unit, cleaning products more or less aggressive will be used, i.e.:

Cool or hot water, mild detergent solutions, organic solvents. In all cases, use of non-aggressive tools is recommended: particle free clothing or paper, sponge, brushes...

• The target of the cleaning is to perform a mechanical work able to remove the residues and organic material stuck to the surface of the equipment. This dragging effect is also useful to remove a big part of microbial load. After cleaning, a complete drying is required.

- Nonaggressive tools recommended:
  - particle free clothing or paper
  - sponge
  - brushes, etc.
- Cleaning products more or less aggressive:
  - Cool or hot water
  - mild detergent solutions
  - organic solvents

### 11.1.2. Disinfection

• After cleaning, disinfecting of surface to remove microbiological contamination is necessary. This is not a sterilization but sanitization, by means of suitable disinfectants.

• It is recommended to use several types of disinfectants with different mechanism, and use them alternatively, changing from one to another every 1-2 months to avoid development of bacterial resistance.

- Some recommended disinfectants are:

- Alcohol: Isopropanol 70%, or ethanol 70%
- Aldehydes: Formaline, Glutaraldehyde
- Cationic: Benzalconium chloride, quaternary ammonium
- Peroxide and oxides: Hydrogen peroxide and peracetic acid

• Once cleaning and disinfecting is done, a test is required to check that procedure has been properly performed (with a small broom or with a piece of sterile cotton). A visual inspection and a microbiological control of collected samples are recommended.

## **11.2. Vacuum unit maintenance**

### **WARNING**

- Respect all service alarms and all the warnings included in the user manual of the vacuum pump, specially about the required oil level and its quality.
- Certain freeze-drying cycles\*, basically those in which the solvent is not water, require increasing the preventive service operations on the vacuum pump.

**\* Not performing the required inspections in such cases may result in severe damage to the unit and/or the product to be freeze-dried.**

Check level and condition of the oil. Add oil if necessary, to reach the correct level.

In order to remove non-condensable vapors from the product and avoid oil contamination of the pump, it is suggested to leave the Gas-Ballast open  $\frac{1}{4}$  during the freeze-dryer cycle.


**Inspection of oil condition.** The usual state of the oil is clear and transparent. The visual inspection is performed through the sight glass at the front of the pump. Oil draining and changing is recommended when it looks dirty or milky. Waste disposal of used oils should be done only following the current prohibitions and safety standards about protection of the environment and according the local regulations of the country or customer facilities.

**Inspection of oil level.** The correct level for oil is halfway up the viewing sight glass, with the pump switched off. If the oil does not reach this level, add the required quantity of oil, through the filling hole, with the pump switched off.

**Visual inspection of oil leaks.** Make sure that there is no oil near the vacuum pump. The presence of drops or puddles of oil suggests possible oil leaks via gaskets or caps or also a low oil level. The area must be cleaned, and the origin of the oil leak must be identified.

**Oil change.** See manual of the corresponding pump version

**ULVAC**

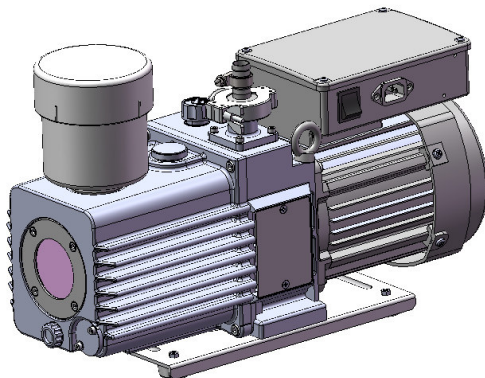
**ILMVAC**

**Vacuum pump internal cleaning.** Each time a change of oil is done, it is advisable to clean the inside of the pump body. It is harmful for the pump to be stored for a long period of time filled with dirty oil; therefore, is advisable to change the oil beforehand.

It is advised to use a funnel when emptying the oil of the pump, as follows:

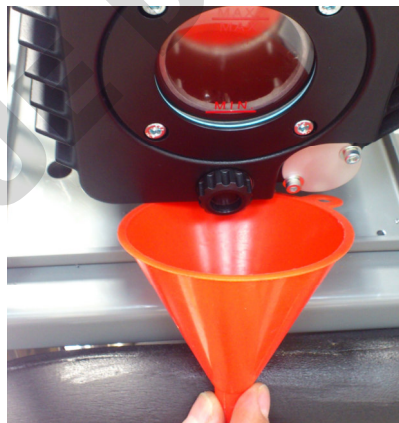See the manual of the pump for the manufacturer recommended oils.

Recommended Azbil Telstar Technologies oils are:

- **AV-30**, mineral oil for general applications in industry and laboratory.
- **CHEMOIL-500**, oil to be used at high temperatures.
- **WETOIL-8r**, oil to be used in aspiration of large quantity of water steam.

About the **ILMVAC pumps**, LABOVAC 14 is recommended for standard applications. See the manual of the pump for special applications recommended oils.

**NOTE:**

- Follow the current environmental regulations about the waste and disposal of used oil.

### 11.3. Refrigerating system maintenance

The refrigerating system requires minimal maintenance. If due to accidents or breakages oil or refrigerant needs to be added, it will be necessary to get in contact with Azbil Telstar Technologies Customer Service.

To assure the proper refrigeration of the system, all the grids must be at a minimum distance of 20cm away from any other equipment and from the wall.

It is recommended to perform a periodic cleaning of the dust accumulated in the condenser batteries of the refrigerating system. To proceed, the lateral doors of the freeze dryer must be opened to clean the battery by means of a small paintbrush, and if it is possible by means of compressed air. This should be done particularly before summer.

When the switch is off, the compressor continues working for the time set by the collecting timer (2 minutes) in order to inhale the gas of the evaporator. If the time has passed and the compressor has not stopped, the equipment will have to be checked, for it's not advisable that the compressor works without gas circulation.

#### Refrigerant loads

	LyoQuest -55	LyoQuest -85		LyoQuest Arctic -85
<b>R-23</b>	-----	95 g	<b>R-170</b>	70 g
<b>R-507</b>	420 g	390 g	<b>R-290</b>	125 g

### 11.4. Maintenance table

The frequencies for each operation, particularly the oil change in the vacuum pump, that are shown in the following Maintenance Table, are suggested values for a normal operation of the freeze-dryer. Extremely severe working conditions or the use of aggressive fluids can reduce the intervals mentioned in the Table.

OPERATION	PERIODICITY
<b>VACUUM UNIT</b>	
- Oil condition control	Daily
- Change of oil filter	1 year
- Oil level control	Daily
- Visual inspection of oil leaks	Daily
- First oil change	After the first 150 hours of service
- Oil change	Minimum every 500 hours of service
- Pump body cleaning	1 year
- Check the performance of the vacuum unit, particularly the standard time for reaching vacuum throughout the unit with cold condenser $P_{atm}$ to 1mbar < 5min	6 months
<b>REFRIGERATING UNIT</b>	
- Condenser panel cleaning	6 months
- Check the performance of the refrigerating system, verifying the normal times to cool the condenser. LyoQuest -55: +20 to -30 <10 min LyoQuest -85: +20 to -40 <5 min LyoQuest Arctic -85: +20 to -40 <5 min	6 months
<b>FREEZE DRYER</b>	
- Joint cleaning and greasing	6 months

## 12. **TROUBLESHOOTING**

### 12.1. **Vacuum unit**

**Problem: The vacuum pump does not start.**

<b>CAUSE</b>	<b>SOLUTION</b>
The pump motor is not receiving voltage.	Check the connection of the on/off switch located upon the motor of the pump.
	Check the electrical connection terminals of the pump.
	Check the corresponding circuit-breaker.
	Check that the environmental temperature is not too high (<35°C), this might have caused the circuit breakage.
The pump's starter protector is blocked because the condenser does not reach the correct temperature. Minimum temperature to switch on the pump: LQ-55 <-30°C / LQ-85 <-40°C LQ Arctic-85 <-40°C	Check the temperature of the condenser.
Pump wheel is seized up. The motor axle does not turn.	Contact Azbil Telstar Technologies Customer Service.

**Problem: Vacuum cannot be obtained.**

<b>CAUSE</b>	<b>SOLUTION</b>
The motor turns but does not create vacuum.	Contact Azbil Telstar Technologies Customer Service.

**Problem: The pump does not reach the vacuum pressure limit.**

<b>CAUSE</b>	<b>SOLUTION</b>
Lack of oil.	Check the level of oil of the pump, and add oil up to the suitable quantity
The oil is contaminated.	Empty the oil from the pump and fill it with new oil.

**Problem: Abnormal heating of the pump.**

<b>CAUSE</b>	<b>SOLUTION</b>
Lack of ventilation or environmental temperature too high (>35°C).	Place the pump in an area with better ventilation.

**Problem: Mist leak from the emission channel.**

<b>CAUSE</b>	<b>SOLUTION</b>
Gas-Ballast is opened too much.	Close the Gas-Ballast.
There is a leak.	Locate the leak.

### 12.2. Refrigerating unit

**Problem: The refrigerating unit does not start up.**

CAUSE	SOLUTION
The compressors do not receive voltage.	Check electrical connection terminals of the compressors.
	Check the circuit breaker corresponding to each compressor.

**Problem: The condenser does not reach the minimum temperature.**

CAUSE	SOLUTION
Refrigerant leak.	Contact Azbil Telstar Technologies Customer Service.
The refrigerating fans of the refrigerating unit do not work.	Check the circuit breaker corresponding to each fan.

**Problem: Intermittent stoppage of the compressor.**

CAUSE	SOLUTION
Insufficient refrigeration of the condenser component.	Place the freeze-dryer in an area with better ventilation.
	Check that the temperature in the room is lower than 25°C.
	Clean the area formed by the condenser component with compressed air and a paintbrush.
	Check that there is a suitable circulation of the forced air in the fan.

**Problem: The compressor stops after alternative stoppages.**

CAUSE	SOLUTION
The thermal protector of the motor of the compressor is disconnected.	The same solutions than the previous problem.
	Wait about 30 minutes to cool down the motor and unblock the protector, allowing again its operation.

### 12.3. Freeze-dryer

**Problem: Working vacuum cannot be reached.**

CAUSE	SOLUTION
There is a leak in the equipment.	Check that the unit is watertight.
	Lubricate the watertight joints with silicone.
	Check that all the manifold keys are closed.
	Check that the condenser discharge valve is closed.
	Check there is no defrosting water of the condenser remaining from the last operation in the discharge tube.

**Problem: The vacuum gauge indicates "gauge short-circuit" on the display (-100).**

CAUSE	SOLUTION
Vacuum gauge short-circuit	Replace the gauge.

### 13. ANNEXES

#### 13.1. Equipment spare part codes

Only original parts and accessories will guarantee correct operation of the freeze-dryer.

Please do not hesitate to contact the Customers Service pointing out the model, operating frequency (Hz) and serial number of the unit.

DESCRIPTION	CODE
-------------	------

#### VACUUM UNIT

##### ULVAC version (standard)

Oil rotary vacuum pump ULVAC GLD137CC	636214
Oil mist filter OMT-200A ULVAC	622106
Major kit ULVAC GLD137CC	633910
Oil AV-30 for ULVAC vacuum pump (1 litre)	15263

##### ILMVAC version (optional)

Vacuum pump ILMVAC P6Z-101	72324
Cartridge for oil mist 800160 for P6Z-101	630249
Service kit ILMVAC P6Z-101 302076-01 (Rotary vane kit)	630248
Service kit ILMVAC P6Z-101 402008-01 (Diaphragm kit)	644895
Oil Labovac-14 for ILMVAC Vacuum Pump (1 litre)	644897

Vacuum transducer THYRACONT	57403
Silicone grease	4716

#### REFRIGERATING SYSTEM

##### LyoQuest -55 / LyoQuest -85

Compressor R-507 SC15CL Danfoss (-55°C & -85°C, 50 Hz)	58411
Compressor R-507 SC15CLX Danfoss (-55°C & -85°C, 60 Hz)	62773
Compressor R-23 SC10CL Danfoss (-85°C, 50 Hz)	58410
Compressor R-23 SC10CLX Danfoss (-85°C, 60 Hz)	62771

##### LyoQuest Arctic -85

Compressor R-290 NLE12.6CNL SECOP (-85°C, 50 Hz)	639869
Compressor R-170 DLE7.5CN SECOP (-85°C, 50 Hz)	639880

#### GENERAL

Condenser tightness gasket	1728821
----------------------------	---------

### 13.2. Examples of recipes

The following 2 exemplary recipes are intended to facilitate the use of the equipment and provide the user with an initial starting point to start lyophilization.

**Attention:** These recipes are merely intended as a tutorial and therefore they cannot be applied to any other product. The recipes were developed for use with the intended products only (i.e. mushrooms and skimmed milk) and regarding pressure and temperature conditions. Times may also vary depending on the container used, product height, etc. It is important that the products are well frozen before placing them in the equipment. Evidently every other product will have its own optimal recipe which will depend on its intrinsic nature. By consequence, these recipes can only be used in order to get a general idea of how to take a start with the lyophilization process. Further experimentation based on these recipes can nevertheless lead to the optimal recipe for each specific product.

#### ANNEX A.1: RECIPE FOR MUSHROOMS LYOPHILIZATION

Step N°	Process	Parameters		
		Time (hh:mm)	Pressure (mBar)	Temperature(°C)
1	Cool + vacuum	00:30	0,2	-
2	Cool + vacuum	01:00	0,2	-
3	Heat plates	01:00	0,2	30
4	Heat plates	99:00	0,2	30

**Fig. A1: lyophilized mushrooms.** As can be seen the lyophilized mushrooms maintain the original volume and appearance. Noticeably this product has practically lost all its humidity, which results in a strongly reduced density compared to the original product. Drying products by lyophilization involves two major advantages: on the one hand it is possible to conserve the product during a much longer period of time, while on the other hand it entails a considerable reduction in transportation costs (due to its reduced weight). In order to maximize the conservation of the mushrooms it is recommended that once lyophilized, the mushrooms are inserted in a bag that previously has been vented with nitrogen. Once the mushrooms have been inserted, the bag needs to be closed immediately to avoid contamination by atmospheric humidity.

Fig. A1

Fig. A2

#### ANNEX A.2: RECIPE FOR THE SKIMMED MILK LYOPHILIZATION

Step N°	Process	Parameters		
		Time (hh:mm)	Pressure (mBar)	Temperature (°C)
1	Cool + vacuum	00:30	0,2	-
2	Heat plates	01:00	0,2	20
3	Heat plates	24:00	0,2	20
4	Heat plates	00:15	0,001	35
5	Heat plates	06:00	0,001	35

**Fig. A2: lyophilized skimmed milk.** The lyophilized skimmed milk has the appearance of a well-formed, solid, consistent tablet with minor shrinkage. In order to obtain optimal conservation, it is recommended to lyophilize the milk in vials, which are to be closed under vacuum atmosphere once the lyophilization process has come to an end. In this way contamination of the product by atmospheric humidity is avoided.

**13.3. Maintenance records**

Write down in the attached tables the information about the maintenance operations performed in the vacuum pump and in the refrigerating compressor. It is necessary to indicate the number of hours operated before resetting the time meters.

<b>VACUUM PUMP</b>			
<b>DATE</b>	<b>HOURS</b>	<b>OPERATION</b>	<b>DONE BY:</b>

<b>COMPRESSORS</b>			
<b>DATE</b>	<b>HOURS</b>	<b>OPERATION</b>	<b>DONE BY:</b>
<b>OTHER MAINTENANCE RECORDS</b>			
<b>DATE</b>	<b>HOURS</b>	<b>OPERATION</b>	<b>DONE BY:</b>

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*By means of this document the customer declares that the equipment has been decontaminated and therefore it is not a hazard for Telstar's personnel that have to proceed with its manipulation.*

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
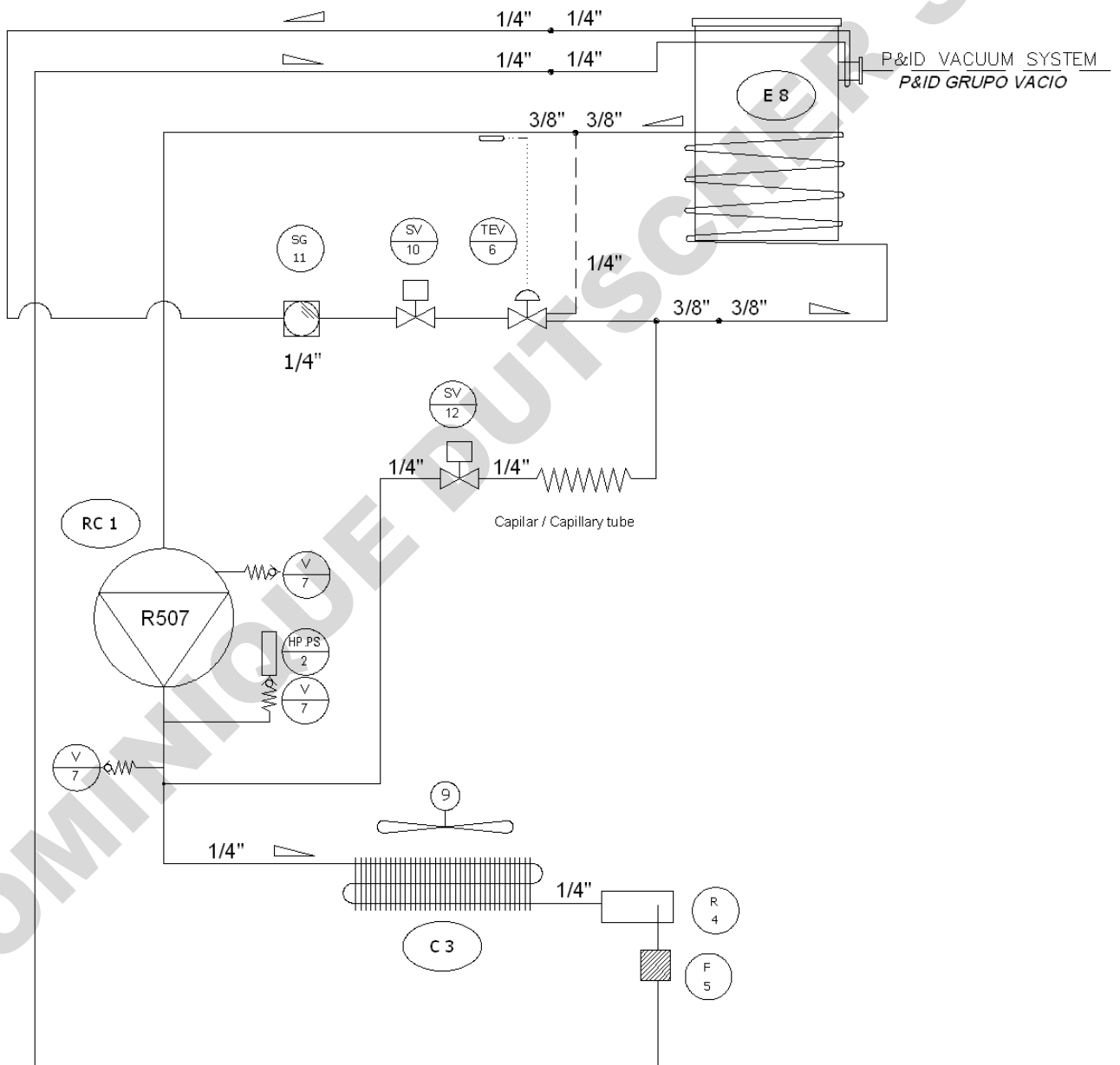
**13.6. Refrigerating system P&ID (LyoQuest -85)**


No.	DESCRIPTION	CODE	No.	DESCRIPTION	CODE
RC1	Compressor R-507	58411	V9	Access valve	3023
HPPS2	High-pressure controller	57311	E10	Evaporator	55190
C3	Air condenser	37223	TVE11	Expansion valve R-23	19657
R4	Liquid receiver	1740404	SV12	Solenoid valve	3510
F5	Dryer filter	19813	HE13	Heat exchanger	19534
F6	Filter	60625	SC14	Subcooler	---
TEV7	Expansion valve R-507	3485	15	Fan motor	57208
RC8	Compressor R-23	58410	SG16	Sight glass	20661

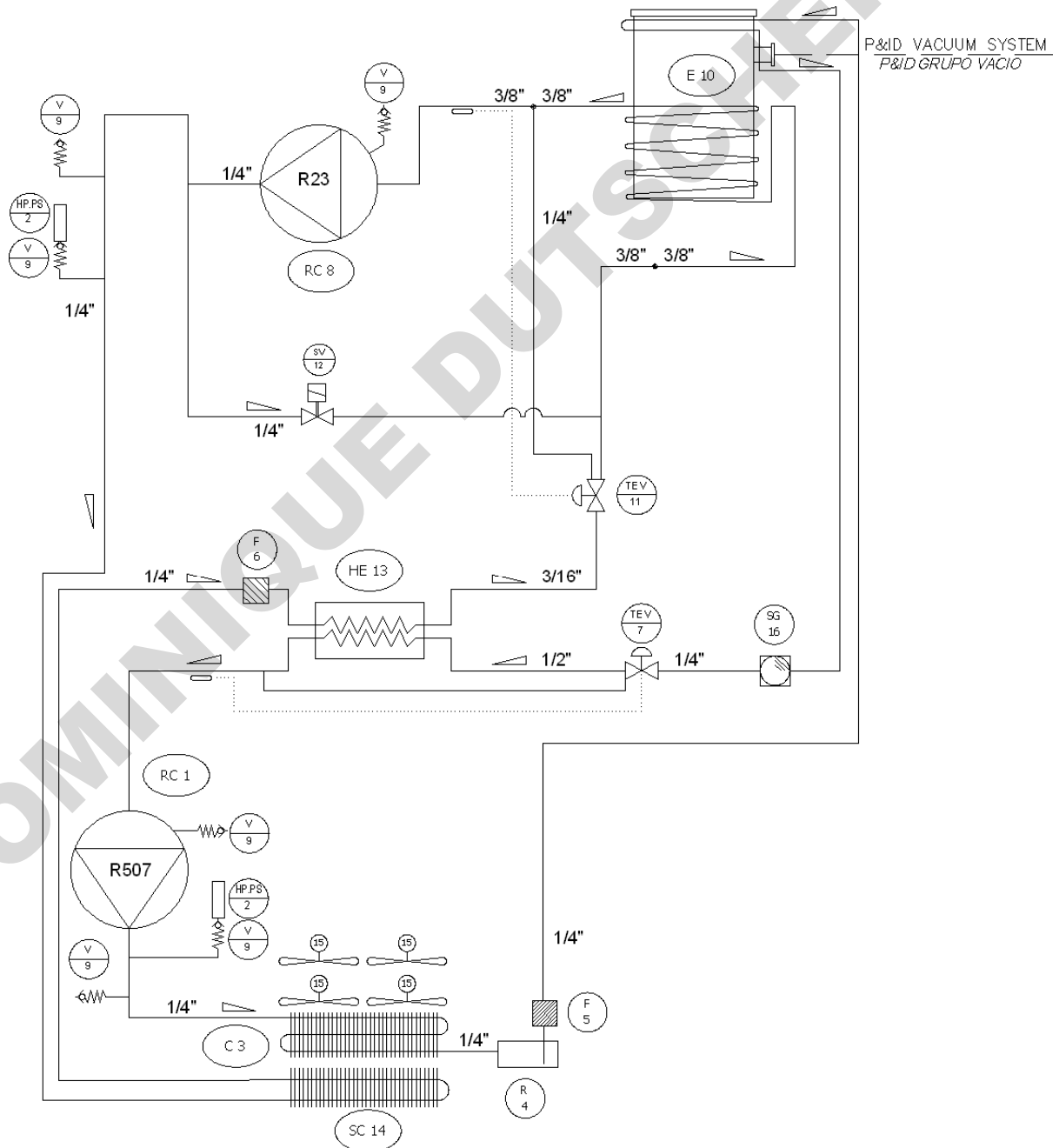
**13.7. Refrigerating system P&ID (LyoQuest Arctic -85)**

No.	DESCRIPTION	CODE	No.	DESCRIPTION	CODE
<b>RC1</b>	Compressor R-290	639869	<b>HV11</b>	Access valve	629121
<b>PSH2</b>	High-pressure controller	57311	<b>RC13</b>	Compressor R-170	639880
<b>HE3</b>	Air condenser	639885	<b>HV14</b>	Access valve	629121
<b>MX4</b>	Fan motor	639899	<b>PSH16</b>	High-pressure controller	57311
<b>HV5</b>	Access valve	629121	<b>F17</b>	Dryer filter	639883
<b>SG8</b>	Sight glass	645114	<b>EMV19</b>	Expansion valve R-170	639887
<b>F7</b>	Dryer filter	639883	<b>SV24</b>	Solenoid valve	3510
<b>EMV9</b>	Expansion valve R-290	645772	<b>HV22</b>	Access valve	629121
<b>HE10</b>	Heat exchanger	616008	<b>HE20</b>	Evaporator	1699738

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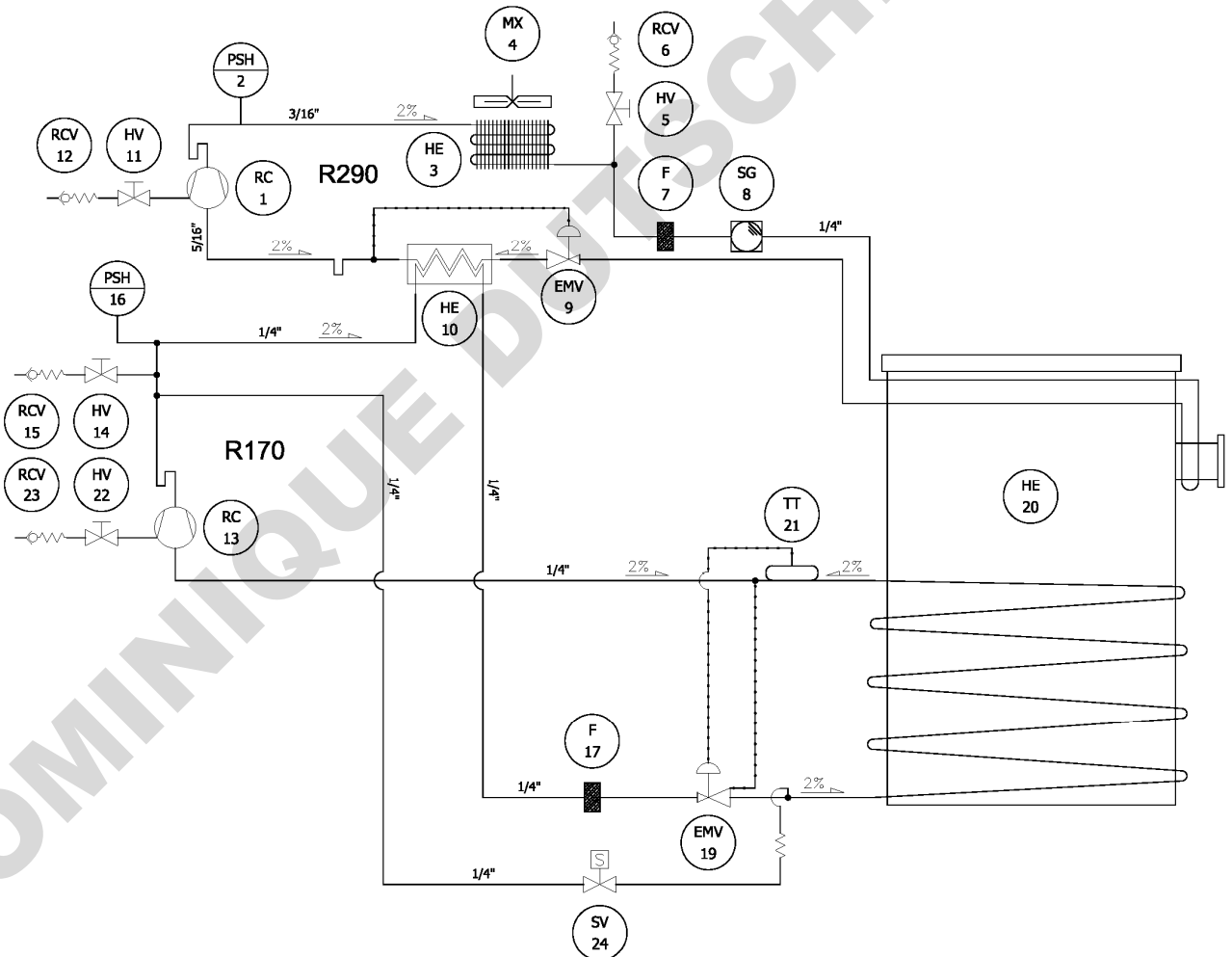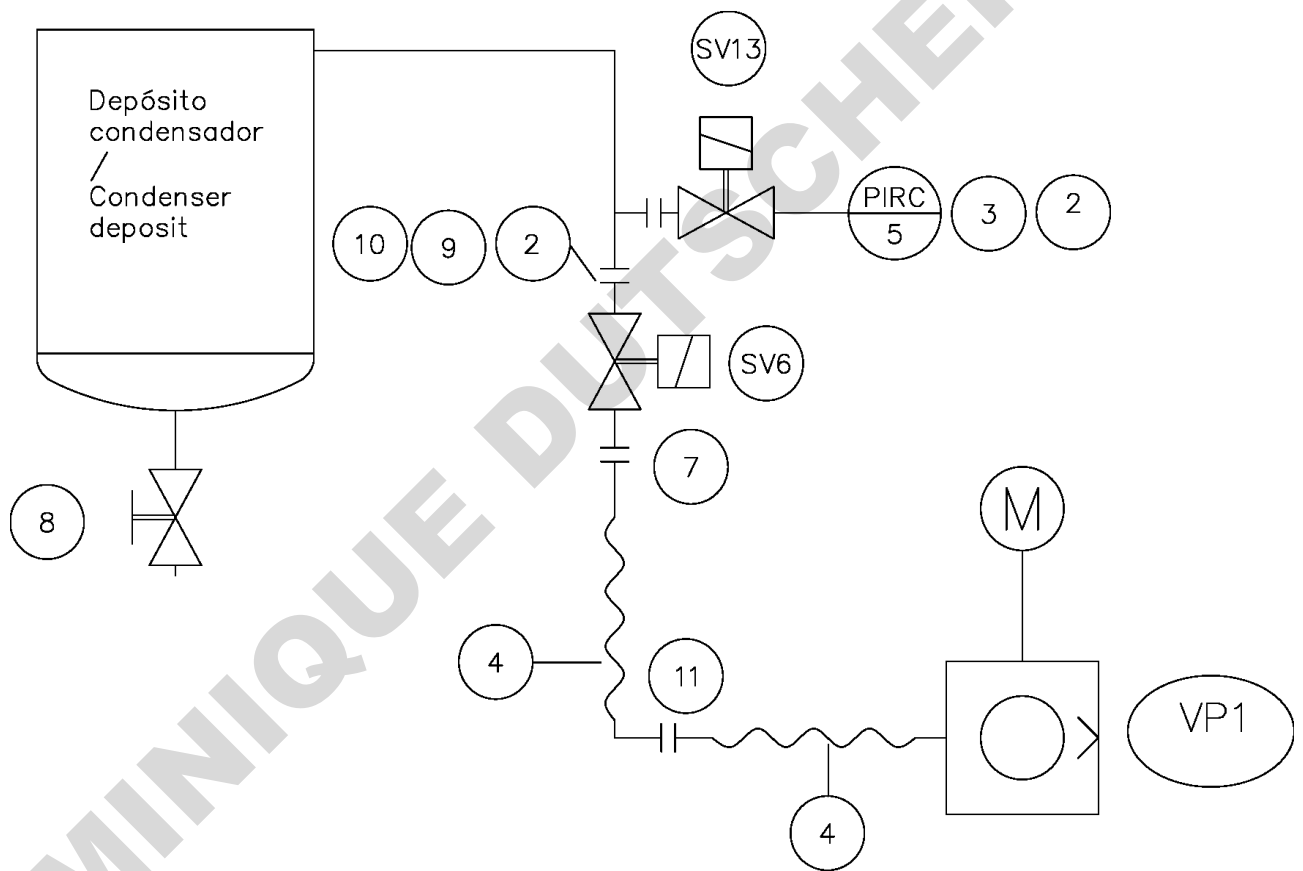
### 13.10. Accessories configuration

## LYOQUEST ACCESSORIES CONFIGURATION

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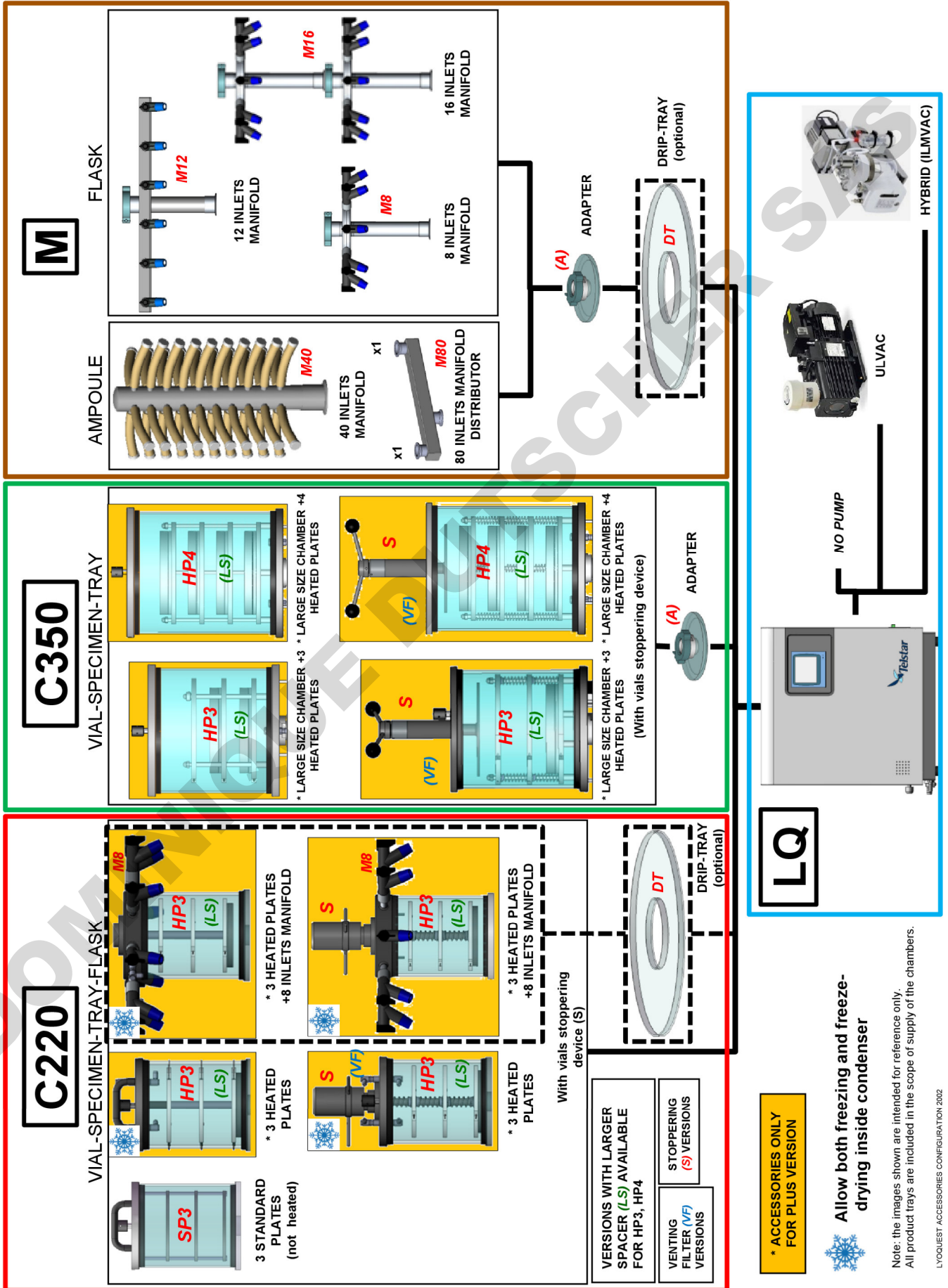The image shows a dialog box titled "IP Settings" with a "LAN" tab. It contains two radio buttons: "Obtain an IP address via DHCP" (unselected) and "Specify an IP address" (selected). Below are input fields for "IP Address" (192.168.9.248), "Subnet mask" (255.255.255.0), "Default Gateway", "Primary DNS", and "Secondary DNS". At the bottom are "OK", "Cancel", and "Apply" buttons.

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The image shows a screen titled "Acceso restringido" (Restricted Access). It features a large text label "Password" above a single-line text input field. At the bottom left is a blue circular button with a checkmark, and at the bottom right is a blue circular button with a crossed-out circle.

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The image shows a 'Login' dialog box with a yellow padlock icon on the left. It contains a 'User:' dropdown menu with 'Administrator' selected, and a 'Password:' text input field. At the bottom, there are three buttons: 'Change Password', 'OK', and 'Cancel'.

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
The image shows a 'Change Password' dialog box with a yellow padlock icon on the left. It contains three text input fields: 'Old Password:', 'New Password:', and 'Confirm Password:'. At the bottom, there are three buttons: 'Change', 'OK', and 'Cancel'.

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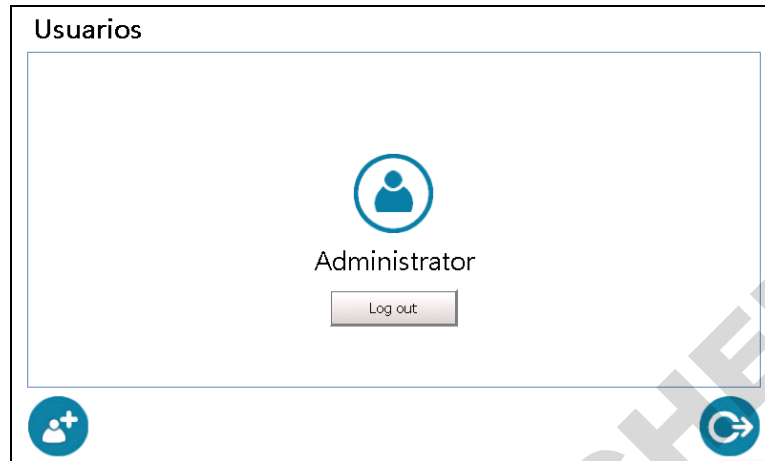
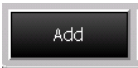
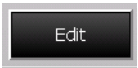
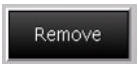
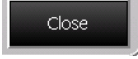
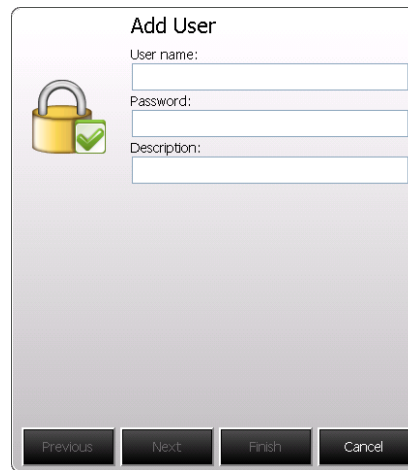
- User name:** Text input field.
- Password:** Text input field.
- Description:** Text input field.
- Navigation buttons:** Previous, Next, Finish, and Cancel.

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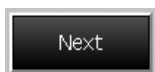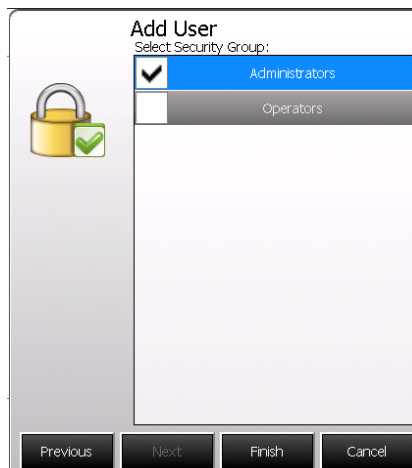
The 'Add User' dialog box now shows a list of security groups:

- Select Security Group:**
  - Administrators
  - Operators
- Navigation buttons:** Previous, Next, Finish, and Cancel.

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