
Cole-Parmer®

CB-200 Series Block Cooling Unit

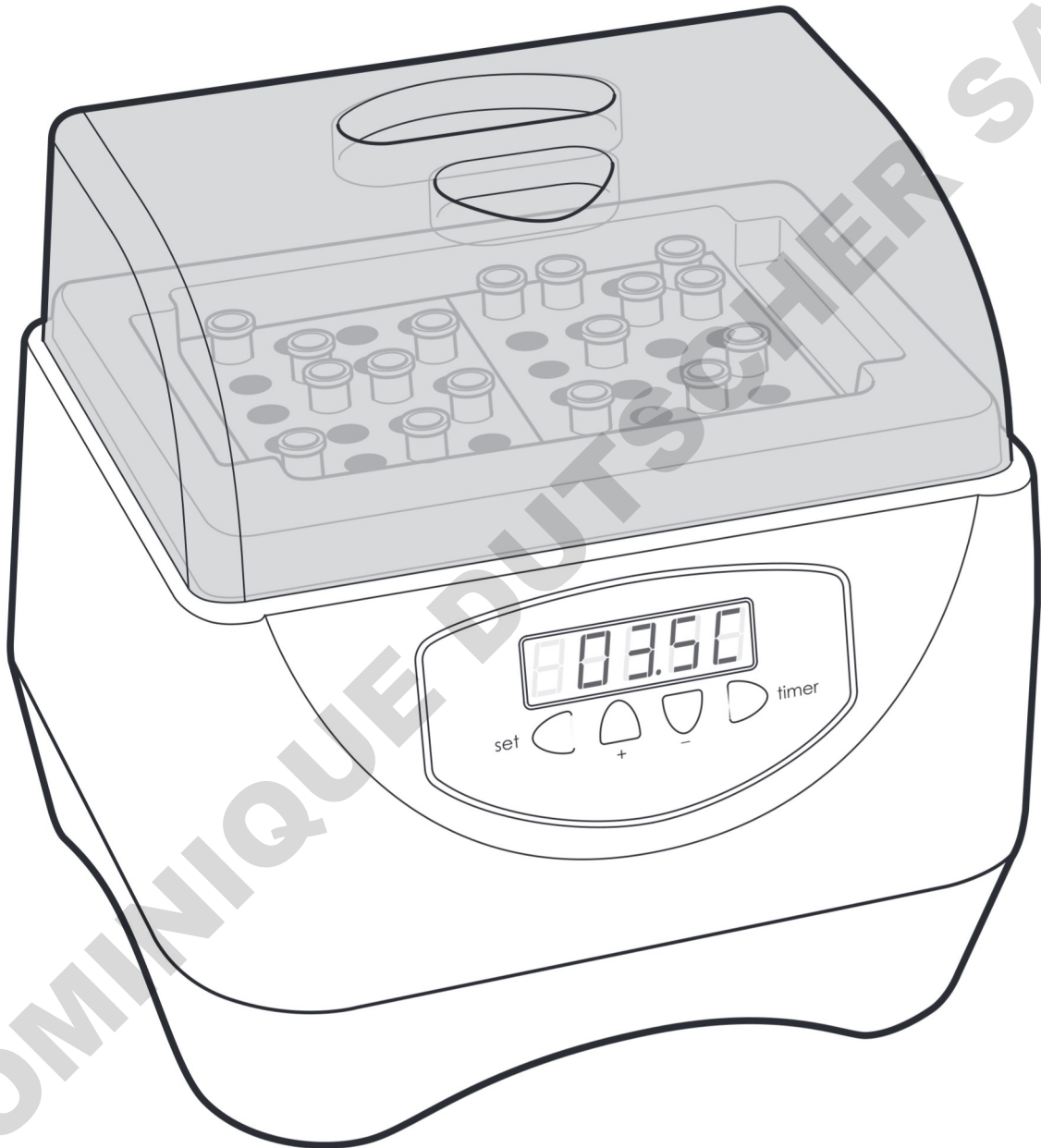Instruction Manual
6106334 Version 4.2

Cole-Parmer®
essentials

CONTENTS


	PAGE
INTRODUCTION	2
BEFORE USE	3
Unpacking	3
SAFETY AND INSTALLATION	4
SÉCURITÉ ET CONSIGNES D'INSTALLATION	6
SICHERHEITS- UND INSTALLATIONSINFORMATIONEN	8
INFORMAZIONI SULLA SICUREZZA E	10
INFORMACIÓN DE SEGURIDAD E INSTALACIÓN	12
INSTALACIÓN CONTACT INFORMATION	14
SPECIFICATION	15
OPERATION	16
Installation	16
Power up Screen	16
Front Panel Controls	16
Setting the Operating Temperature	17
Changing the Timer Mode	17
Setting the Count Down Time	18
Starting the Timer	18
Changing the Temperature Units	19
After Use	19
User Maintenance	19
FAQS	20
ADDITIONAL INFORMATION	21
Fault Finding	21
Fuses	21
Insulation Tests	21
ACCESSORIES	22
REPLACEMENT PARTS	23

INTRODUCTION

The CB-200D is a Peltier controlled cooling/warming device which is ideal for incubating sample tubes at temperatures in the range of 0 to 40°C. The unit has an integral count down/count up timer with audible alarm to set and monitor reaction times.

The CB-200D can hold two interchangeable aluminium insert blocks or one 96-well plate block for cooling a variety of tube formats. The blocks help to maintain a temperature stability of $\pm 0.2^{\circ}\text{C}$ for precise control of reaction temperature.

The unit is supplied with a dark plastic lid which maintains the cooled environment and helps to protect light-sensitive reagents.

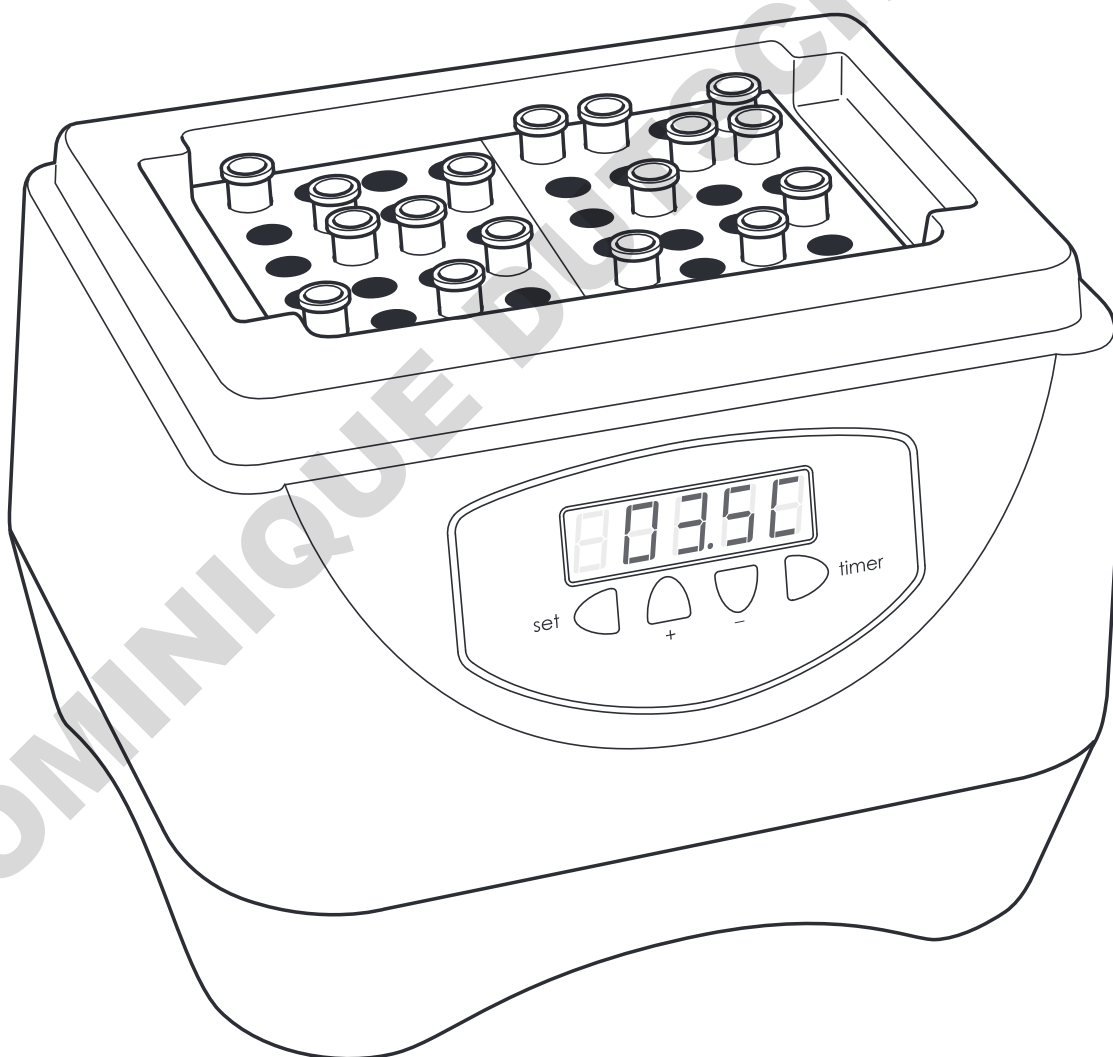BEFORE USE

Before using the CB-200D please make sure you have read this manual carefully. If there is any doubt relating to the proper use of this equipment, the staff at Antylia Scientific or your supplier will be happy to assist you.

UNPACKING

When unpacking the unit please ensure that the following have been removed from the packaging:

- CB-200D
- Mains cables (UK, EU and US)
- Operator's manual

Note: Block inserts should be ordered separately.

The user is advised to keep the original packaging in case the instrument ever needs to be returned for service or repair. Antylia Scientific accepts no responsibility for damage incurred unless the unit is correctly packed and transported in its original packaging.

SAFETY AND INSTALLATION

Please read all the information in this manual before using the CB-200D.

WARNING

Contact with low-temperature objects may cause cold burns. Users should be aware of the following potential hazards:

- USE CARE AND WEAR PROTECTIVE GLOVES TO PROTECT HANDS
- DO NOT place any liquid directly into the instrument.
- DO NOT position the product so that it is difficult to access the ON/OFF switch or disconnect it from the mains supply using the mains plug.

OPERATOR SAFETY

All operators of Cole-Parmer equipment must have available the relevant literature needed to ensure their safety. It is important that only suitably trained personnel operate this equipment, in accordance with the instructions contained in this manual and with general safety standards and procedures. If the equipment is used in a manner not specified by Antylia Scientific the protection provided by the equipment to the operator may be impaired.

All Cole-Parmer instruments are designed to conform to international safety requirements and are fitted with an over-temperature cut-out. If a safety problem should be encountered, switch off the unit at the mains socket and remove the plug from the electricity supply.

If the equipment is not used in the manner described in this manual and with accessories other than those recommended by Antylia Scientific Ltd, the protection provided may be impaired.

INSTALLATION

The instrument should be carried using both hands. Never move or carry the instrument when in use or connected to the mains electricity supply.

1. All Cole-Parmer instruments are supplied with a power cable; this may be integral or plug-in.
2.  Before connecting the instrument to the mains electricity supply, check the voltage against the rating plate (located on the back of the unit). **Please note that the unit must be earthed to ensure proper electrical safety.** Connect the mains cable to a suitable plug according to the table below.

Connection	220/240V, 50/60Hz Supply	110V/120V Supply, 50/60Hz
Live	Brown	Black
Neutral	Blue	White
Earth	Green/yellow	Green

3. UK ONLY: The fused plug supplied with the mains cable is fitted with a 10 amp fuse to protect the instrument and the operator.
4. The units are rated to operate at 100-230V, 50/60Hz.
5. Place the unit on a suitable flat bench or in a fume cupboard if required, ensuring that the air vents on the rear and underside are free from obstruction.

6. Plug the mains cable into the socket on the back of the instrument.

7. Switch on the instrument. Symbols on or near the power switch of the unit have the following meanings:

ⓘ Mains Switch On

○ Mains Switch Off

REPLACEMENT CABLE

Should the mains lead need replacement, a cable of 1mm² of harmonized code H05VV-F connected to an IEC320 plug should be used. **IF IN DOUBT CONSULT A QUALIFIED ELECTRICIAN.**

ENVIRONMENTAL CONDITIONS

The CB-200D is designed to operate under the following conditions:

- Indoor use
- Ambient temperature range +5°C to +40°C
- Altitude to 2000m
- Relative humidity not exceeding 80%
- Mains supply fluctuations not exceeding 10%
- Over voltage category II IEC 60364-4-443
- Pollution degree 2

Note: The control specifications are quoted at an ambient temperature of 20°C. The specification may deteriorate outside an ambient temperature of between 10°C and 30°C.

The instrument has been tested for radio frequency interference and is certified under EN61326.

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SPECIFICATION

Temperature

Temperature range	0.0 to 40.0°C*
Temperature accuracy	±1°C
Temperature units	°C or °F
Cooling technology	Peltier
Temperature display resolution	0.1°C
Temperature display	Orange LED, 5 digits
Timer	m/s, h/m, d/h up to 9 days

Dimensions

Height	225mm
Width	240mm
Length	190mm
Weight*	3.0kg

*Shipping weight excludes block inserts which must be ordered separately

Power

Power supply	100 to 230V, 50/60Hz, 75W
Fuse	Glass 5x20mm, T5A, 250V

* Please note the primary design of BL°CKICE is to maintain biological samples safely between 4°C and 37°C. CB-200D will achieve any temperature in this range within 30 minutes in an ambient temperature of 20°C. CB-200D will also reach 0°C or 40°C, but this may take several hours. Reaching 0°C can be speeded-up by pre-cooling the aluminium blocks before use.

OPERATION

INSTALLATION

See also the sections on Safety and Installation.

1. Place the unit on a suitable flat bench or in a fume cupboard if required, ensuring that the air vents on the rear and underside are free from obstruction.
2. Plug the mains cable into the socket on the back of the instrument.
3. Place the required block inserts inside the unit. Up to two standard block inserts or one 96-well block insert can be used.
4. Connect to the mains electricity supply with the plug provided or one wired correctly for the supply. Switch the power ON from the switch located next to the mains input lead. The front display will then light up.

POWER UP SCREEN

When first switched on, the display will show the name of the product and then show the current temperature of the unit. The unit will then proceed to heat or cool to the last set temperature.

FRONT PANEL CONTROLS

The front panel consists of four buttons for setting the unit temperature and timer and a 5 digit LED display.

The SET Button

Use this button to view and set the temperature, set the timer to count up or down, set the count down time and set the temperature units.

The PLUS Button


Use this button to increase a temperature or time value and to start the timer.

The MINUS Button

Use this button to decrease a temperature or time value and to pause the timer.

The TIMER button


Use this to view the set count down time and to start the timer.

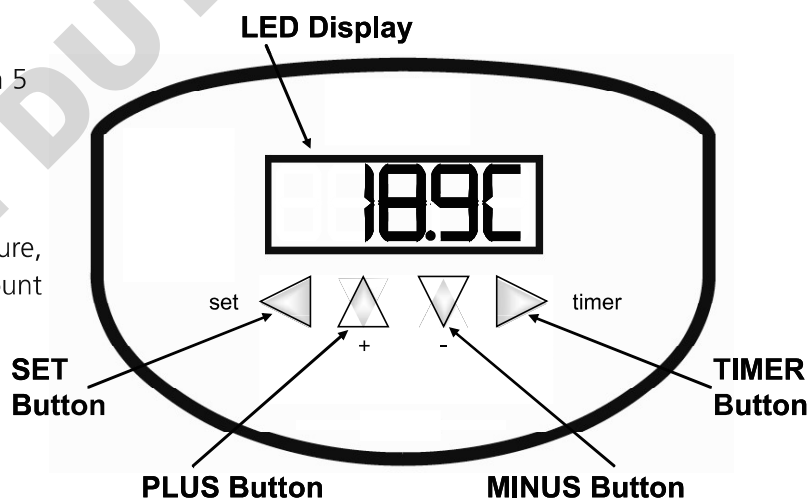
SETTING THE OPERATING TEMPERATURE

1. To view the current set temperature, press the SET button. The set temperature will be displayed and is indicated by the set point indicator in the top left hand corner of the display.
2. To edit the set temperature, press and hold the SET button until the set point indicator moves one position to the right. This is the edit position. Release the SET button.
3. To **increase** the set temperature, whilst the indicator is in the edit position, press the PLUS button. Holding down the PLUS button will cause the value to increase rapidly. A single press of the PLUS or MINUS buttons will increase or decrease the set temperature by 0.1°C.
4. To **decrease** the set temperature, whilst the indicator is in the edit position, press the MINUS button. Holding down the MINUS button will cause the value to decrease rapidly. A single press of the PLUS or MINUS buttons will increase or decrease the set temperature by 0.1°C.
5. Once the set temperature has been changed, the unit will heat or cool to the new set temperature. A buzzer will sound once a stable temperature has been reached.

Set Point Indicator

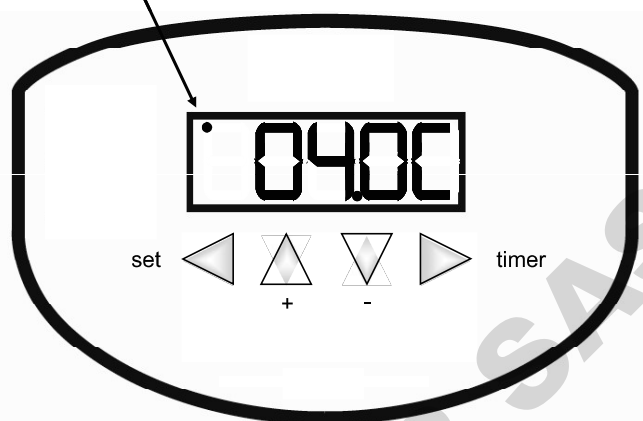Edit Set Point Indicator

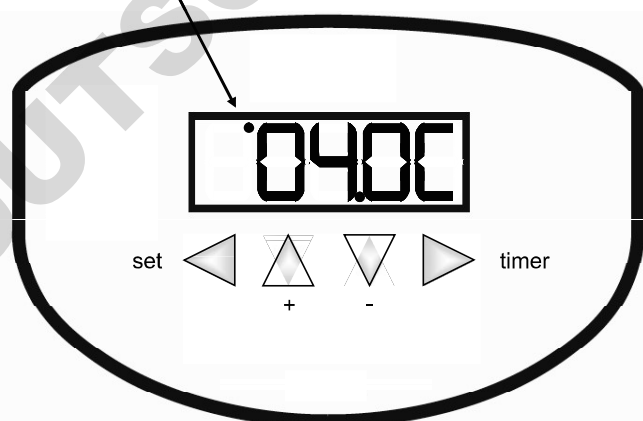
CHANGING THE TIMER MODE

The timer can be set to count up, count down or be disabled.

6. To change the timer mode, press the SET button twice, holding the button down on the second press until the indicator point appears.
7. Press the PLUS button to toggle between **d** for count down, **u** for count up and - - - for timer disabled.
8. A few seconds after the buttons have been released the unit will automatically revert to the temperature display.


Indicator point

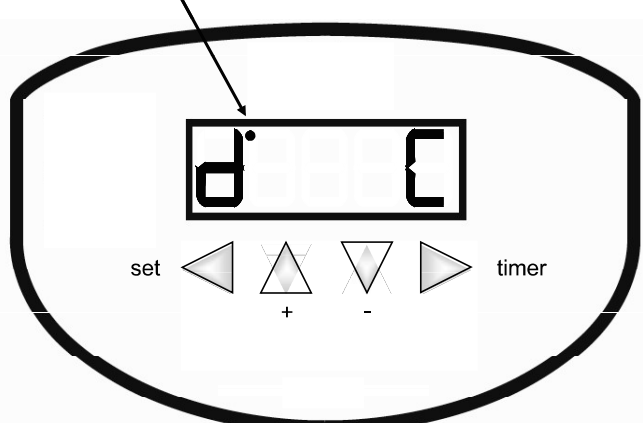SETTING THE COUNT DOWN TIME

The timer can be set to count down time in minutes and seconds, hours and minutes or days and hours up to a maximum of 9 days.


9. First set the timer to count down mode as described in steps 6 to 8.
10. To set the count down time, press the SET button 3 times, holding the button down on the third press until the indicator point appears next to the letter t.
11. To **increase** the count down time, press the PLUS button. Holding down the PLUS button will cause the value to increase rapidly. A single press of the PLUS or MINUS buttons will increase or decrease the time by 1 unit (s, m, h or d).
12. To **decrease** the count down time, press the MINUS button. Holding down the MINUS button will cause the value to decrease rapidly. A single press of the PLUS or MINUS buttons will increase or decrease the time by 1 unit (s, m, h or d).

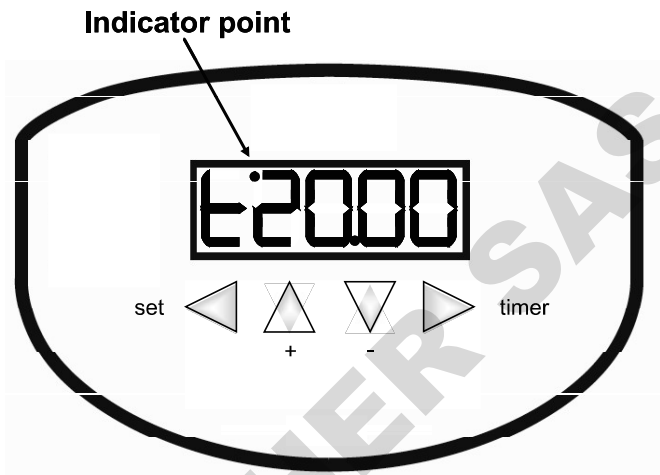STARTING THE TIMER

The timer is activated using the TIMER button. The timer will not begin to count up or down until a stable temperature has been reached.

13. To view the timer setting, press the TIMER button. This will show whether the timer is in count down mode (**d**), count up mode (**u**) or disabled (- - -). If the timer is in count down mode, the set time will also be shown. To change the count down time, follow steps 9-12.
14. To activate the timer, press the TIMER button twice, holding down on the second press until the letter **t** appears.
15. To start the timer, press the PLUS button. The letter **t** will flash while the timer is in operation (either counting up or counting down).
16. To pause the timer, press the MINUS button.

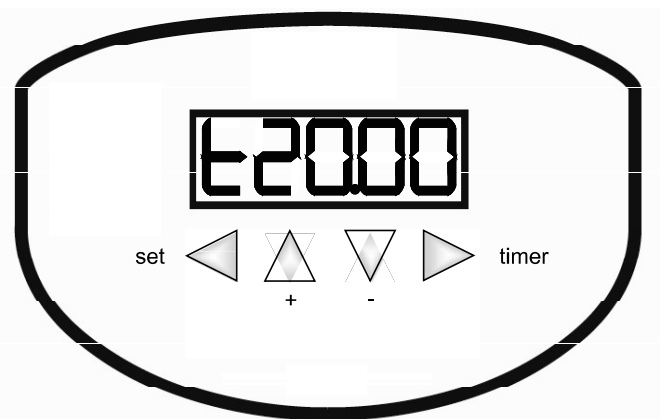
17. To restart the timer after a pause, press the PLUS button.
18. To stop the timer and reset, press the TIMER button.
19. When the count down timer reaches zero, a buzzer will sound and the timer will begin to count up to indicate the over-run time. To stop the buzzer, press the TIMER button.

CHANGING THE TEMPERATURE UNITS

The temperature of the display can be shown in either °C or °F.

20. To change the temperature units, press the SET button 4 times, holding the button down on the fourth press until the indicator point appears next to the letter **U** (units).
21. Press the PLUS button to toggle between °C and °F.
22. A few seconds after the buttons have been released the unit will automatically display the current temperature in the selected units.

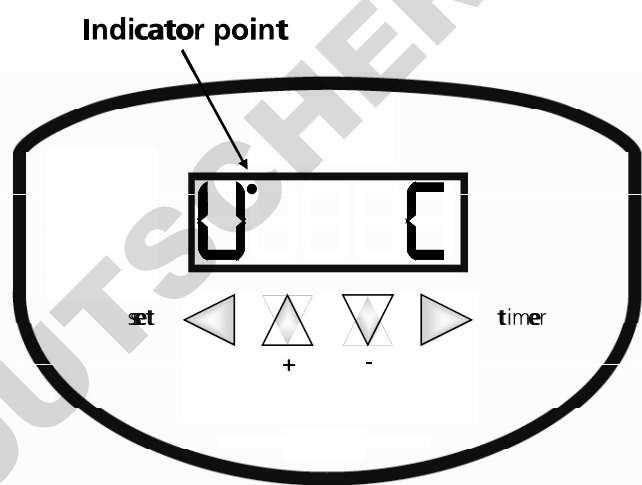AFTER USE

After cooling samples, the bucket and beads may be very cold. Allow to warm to ambient conditions before lifting from the unit or wear protective gloves. Avoid handling cold objects with wet hands.

USER MAINTENANCE

Before cleaning the unit, disconnect it from the power supply. The outer case of the CB-200D may be wiped clean with a damp cloth; ethanol, methanol or formaldehyde can also be used. No part of the case or cover should be immersed in the solvents. Do not use aggressive solvents such as acetone or abrasive cleaners. Allow the unit to dry before using.

Before using any cleaning or decontamination method except those recommended here, the responsible person should check with Antylia Scientific that the proposed method will not damage the equipment.

FAQ'S

Q1 How can I cool the beads more quickly?

A1 **If you intend to use the CB-200D to keep samples cool, then the block inserts can be stored in a refrigerator and then inserted into the unit when required.**

Q2 I have a Cole-Parmer Block Heater; can I use the block inserts I already use with this?

A2 **Yes. The CB-200D will accept any two standard block inserts from the Block Heater range or a single 96-well block insert.**

If you require further technical or application assistance please contact Antylia Scientific at:

Antylia Scientific Ltd.
Service Department,
Beacon Road,
Stone,
Staffordshire,
ST15 0SA,
Email: cpSERVICE@antylia.com
Phone: +44 (0)1785 810475

If you require further technical or application assistance please contact:

Email: cpSERVICE@antylia.com
Phone: +44 (0)1785 810433

We are continually striving to improve our products and software. If you have any comments and suggestions on how we can do things better please send them to us at: cpSERVICE@antylia.com

ADDITIONAL INFORMATION

Note that **this equipment should only be dismantled by properly trained personnel**. Removing the outer cover exposes potentially lethal mains voltages. There are no user serviceable parts within this equipment.

FAULT FINDING

Should you have any problems with your CB-200D which cannot be easily remedied, you should contact your supplier and return the unit if necessary. Please include details of the fault observed and remember to return the unit in its original packing. Antylia Scientific accepts no responsibility for damage to units which are not properly packed for shipping: if in doubt, contact your supplier, giving the full serial number of the unit.

FUSES

If the display on the front panel is not lit, one of the two fuses may have blown. Check that there is no external cause, such as a faulty plug or lead. Check both fuses and replace the faulty fuse with a new one of the correct value (fuse values are given on the label next to the power inlet). Note that **fuses should only be replaced by a qualified electrician.**

The holder for the two fuses is built into the mains input socket. First remove the power cable and then gently prise the fuse drawer open with a flat-bladed screwdriver or similar tool. Each fuse can be removed by using the screwdriver as a lever.

Exchange the faulty fuse in the fuse holder for a working fuse of the correct value. Finally, replace the fuse drawer in the fuse compartment and push the drawer shut. Fuses which blow repeatedly indicate a serious fault and you should return the unit to your supplier for repair.

INSULATION TESTING

This equipment is fitted with RFI suppression circuitry. Any check of the electrical insulation by means of high voltage dielectric testing (for example as in BS EN 61010-1) must be carried out using only a DC voltage.

This unit contains semiconductor components which may be damaged by electric field effects.

REPLACEMENT PARTS

The following replacement parts can be obtained from Antylia Scientific or your Cole-Parmer dealer:

Product Code	Description
HH179(S)	UK 230V mains lead with plug
HH180(S)	European mains lead with plug US
HFCABLEUS	mains lead with plug
700738(s)	Fuse Type T5A, 250V (5A)

ACCESSORIES

The following block inserts are available for use with the CB-200D and can be ordered from Antylia Scientific or your Cole-Parmer dealer:

Please visit www.coleparmer.com for a full list of available accessories.

Product code	Tube size (diam.)	Number of holes	Block dimensions (mm)
F3501	Plain block	None	95 x 76 x 51
F3502	6mm	30	95 x 76 x 51
F3503	10mm	20	95 x 76 x 51
F3504	12mm	20	95 x 76 x 51
F3505	13mm	20	95 x 76 x 51
F3506	15mm	12	95 x 76 x 51
F3507	16mm	12	95 x 76 x 51
F3508	19mm	8	95 x 76 x 51
F3509	25mm	6	95 x 76 x 51
F3510	10mm square cuvettes	2 channels	95 x 76 x 51
F4460	Plain block	Thermometer hole only	95 x 76 x 51
F4461	7 and 9mm	20/10	95 x 76 x 51
F4462	24mm	6	95 x 76 x 51
F4463	26mm	6	95 x 76 x 51
F4464	1.5ml microfuge tubes	20	95 x 76 x 51
F4465	0.5ml micro tubes	30	95 x 76 x 51
F4467	Hi-Temp 96-well plate	96	95 x 151 x 61
F4468	Falcon round bottom plate	96	95 x 151 x 61
F4469	Falcon fl at bottom plate	96	95 x 151 x 61
F4470	2.0ml microfuge tubes	20	95 x 76 x 51
F4471	0.2ml micro tubes	72	95 x 76 x 51
F4473	96 x 0.2ml micro tubes or 96-well	96	95 x 151 x 61
F4474	plate 96 x 1.0ml Porvair plate	96	95 x 151 x 48
F4476	Block for Gelation timer	1 sample cup	95 x 76 x 51

This product meets the applicable CE Directives and UKCA Legislation for radio frequency interference and may be expected not to interfere with, or be affected by, other equipment with similar qualifications. We cannot be sure that other equipment used in its vicinity will meet these standards and so we cannot guarantee

that interference will not occur in practise. Where there is a possibility that injury, damage or loss might occur if equipment malfunctions due to radio frequency interference, or for general advise before use, contact the manufacturer.

Declaration of Conformity is available to view online at www.coleparmer.com

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Ordering Information

Order No.	Series	Model	Legacy SKU
36620-17	CB-200	CB-200D	FBLOCKICE

Warranty Registration

UK

T: +44 (0) 1480 272279
E: uk.sales@antylia.com
W: coleparmer.co.uk

Germany

T: +49 (0) 9377 92030
E: de.sales@antylia.com
W: coleparmer.de

France

T: +33 (0) 1486 37800
E: fr.sales@antylia.com
W: coleparmer.fr

Italy

T: +39 (0) 284349215
E: it.sales@antylia.com
W: coleparmer.it

India

T: +91 22 61394444
E: info@antylia.com
W: coleparmer.in

China

T: +1 847 549 7600
E: sales@antylia.com
W: coleparmer.com

USA

T: +1 847 549 7600
E: sales@antylia.com
W: coleparmer.com

Canada

T: +514 355 6100
E: info@antylia.ca
W: coleparmer.it

Other

T: +1 847 549 7600

Cole-Parmer®
essentials

Antylia Scientific Ltd
Beacon Road,
Stone,
Staffordshire,
ST15 0SA
United Kingdom