# Adventurer™ Balances Instruction Manual

## TABLE OF CONTENTS

|   |           |
|---|-----------|
| <b>1. INTRODUCTION</b> .....  | <b>3</b>  |
| 1.1 Description .....   | 3         |
| 1.2 Features .....  | 3         |
| 1.3 Definition of Signal Warnings and Symbols .....                                     | 3         |
| 1.4 Safety Precautions .....  | 3         |
| <b>2. INSTALLATION</b> .....  | <b>4</b>  |
| 2.1 Unpacking .....   | 4         |
| 2.2 Selecting the Location .....  | 4         |
| 2.3 Leveling the Equipment .....  | 4         |
| 2.4 Connecting Power .....  | 5         |
| 2.5 Connecting the Interface .....  | 5         |
| 2.6 Initial Calibration .....   | 5         |
| <b>3. OPERATION</b> .....   | <b>6</b>  |
| 3.1 Overview of Display, Home Screen .....  | 6         |
| 3.2 Principal Functions and Main Menu .....   | 7         |
| 3.3 Overview of Parts and Features – Draft Shield Models .....                          | 8         |
| 3.4 Overview of Parts and Features – Non Draft Shield Models .....                      | 8         |
| <b>4. APPLICATIONS</b> .....  | <b>9</b>  |
| 4.1 Weighing .....  | 9         |
| 4.1.1 Item Settings .....   | 10        |
| 4.1.2 Sample Name .....   | 11        |
| 4.1.3 Save to USB .....   | 11        |
| 4.2 Parts Counting .....  | 11        |
| 4.2.1 Item Settings .....   | 12        |
| 4.3 Percent Weighing .....  | 14        |
| 4.3.1 Item Settings .....   | 14        |
| 4.4 Dynamic Weighing .....  | 15        |
| 4.4.1 Item Settings .....   | 15        |
| 4.5 Density Determination .....   | 16        |
| 4.5.1 Measuring the Density of a Solid Using Water (default) .....                      | 16        |
| 4.5.2 Measuring the Density of a Buoyant Material Using Water .....                     | 17        |
| 4.5.3 Measuring the Density of a Solid Using an Auxiliary Liquid .....                  | 17        |
| 4.5.4 Measuring the Density of a liquid using a Calibrated Sinkers (not supplied) ..... | 18        |
| 4.5.5 Measuring the Density of Porous Material Using Oil .....                          | 19        |
| 4.6 Check Weighing .....  | 22        |
| 4.6.1 Item Settings .....   | 22        |
| 4.7 Display Hold .....  | 23        |
| 4.7.1 Display Hold .....  | 23        |
| 4.7.2 Peak Hold .....   | 23        |
| 4.7.3 Item Settings .....   | 24        |
| 4.8 Totalization .....  | 24        |
| 4.9 Formulation .....   | 25        |
| 4.9.1 Item Settings .....   | 25        |
| 4.10 Additional Features .....  | 26        |
| 4.10.1 Weigh Below .....  | 26        |
| <b>5. MENU SETTINGS</b> .....   | <b>27</b> |
| 5.1 Menu Navigation .....   | 27        |
| 5.1.1 Changing Settings .....   | 28        |
| 5.2 Calibration .....   | 29        |
| 5.2.1 Calibration sub-menu (InCal models) .....   | 29        |
| 5.2.2 Internal Calibration .....  | 29        |
| 5.2.3 Automatic Calibration .....   | 29        |
| 5.2.4 AutoCal™ Adjustment .....   | 30        |
| 5.2.5 Span Calibration .....  | 30        |
| 5.2.6 Linearity Calibration .....   | 31        |
| 5.2.7 Calibration Test .....  | 31        |
| 5.3 Balance Setup .....   | 31        |
| 5.3.1 Balance Set-up sub-menu .....   | 31        |
| 5.3.2 Language .....  | 31        |
| 5.3.3 User Setting .....  | 32        |

|            |   |           |
|------------|---|-----------|
| 5.3.4      | Filter Level.....   | 32        |
| 5.3.5      | Auto Zero Tracking.....   | 32        |
| 5.3.6      | Auto Tare.....  | 32        |
| 5.3.7      | Graduations.....  | 32        |
| 5.3.8      | Date & Time.....  | 33        |
| 5.3.9      | Approved Mode.....  | 33        |
| 5.4        | Weighing Units.....   | 34        |
| 5.4.1      | Units Sub-menu.....   | 34        |
| 5.5        | Data Maintenance.....   | 34        |
| 5.5.1      | Data Maintenance sub-menu.....  | 34        |
| 5.6        | Communication.....  | 35        |
| 5.6.1      | Communication Sub-menu.....   | 35        |
| 5.6.2      | Baud Rate.....  | 35        |
| 5.6.3      | Transmission.....   | 36        |
| 5.6.4      | Handshake.....  | 36        |
| 5.6.5      | Print Output.....   | 36        |
| 5.6.6      | Auto Print.....   | 36        |
| 5.6.7      | Print Content.....  | 36        |
| 5.6.8      | Feed.....   | 36        |
| 5.6.9      | Format.....   | 36        |
| 5.6.10     | Print Calibration Data.....   | 36        |
| 5.7        | GLP and GMP Data.....   | 38        |
| 5.7.1      | GLP Data Sub-menu.....  | 38        |
| 5.7.2      | Header.....   | 38        |
| 5.7.3      | Balance Name.....   | 38        |
| 5.7.4      | User Name.....  | 38        |
| 5.7.5      | Project Name.....   | 38        |
| 5.8        | Factory Reset.....  | 38        |
| 5.9        | Lockout.....  | 38        |
| <b>6.</b>  | <b>LEGAL FOR TRADE (LFT).....</b>                                       | <b>39</b> |
| 6.1        | Settings.....   | 39        |
| 6.2        | Verification.....   | 39        |
| 6.3        | Sealing.....  | 39        |
| <b>7.</b>  | <b>PRINTING.....</b>  | <b>40</b> |
| 7.1        | Connecting, Configuring and Testing the Printer/Computer Interface..... | 40        |
| 7.2        | Output Format.....  | 40        |
| 7.3        | Printout Examples.....  | 41        |
| <b>8.</b>  | <b>MAINTENANCE.....</b>   | <b>44</b> |
| 8.1        | Calibration.....  | 44        |
| 8.2        | Cleaning.....   | 44        |
| 8.3        | Troubleshooting.....  | 44        |
| 8.4        | Service Information.....  | 44        |
| <b>9.</b>  | <b>TECHNICAL DATA.....</b>  | <b>45</b> |
| 9.1        | Specifications.....   | 45        |
| 9.2        | Drawings and Dimensions.....  | 54        |
| 9.3        | Parts and Accessories.....  | 54        |
| 9.4        | Communication.....  | 55        |
| 9.4.1      | Interface Commands.....   | 55        |
| 9.4.2      | RS232 (DB9) Pin Connections.....  | 56        |
| 9.5        | The USB Interface.....  | 56        |
| <b>10.</b> | <b>SOFTWARE UPDATES.....</b>  | <b>57</b> |
| <b>11.</b> | <b>COMPLIANCE.....</b>  | <b>58</b> |

## 1. INTRODUCTION

### 1.1 Description

The Adventurer balance is a precision weighing instrument that will provide you with years of service if properly cared for. The Ohaus Adventurer balances are available in capacities from 120 grams to 8,200 grams.

### 1.2 Features

**Touch Controls:** Quick, graphical access to all control functions, over a dozen applications and many features.

### 1.3 Definition of Signal Warnings and Symbols

Safety notes are marked with signal words and warning symbols. These show safety issues and warnings. Ignoring the safety notes may lead to personal injury, damage to the instrument, malfunctions and false results.

#### Signal Words

|                       |  |
|-----------------------|--|
| <b>WARNING</b>        | For a hazardous situation with medium risk, possibly resulting in injuries or death if not avoided.  |
| <b>CAUTION</b>        | For a hazardous situation with low risk, resulting in damage to the device or the property or in loss of data, or injuries if not avoided. |
| <b>Attention Note</b> | For important information about the product<br>For useful information about the product  |

#### Warning Symbols

General Hazard

Electrical Shock Hazard

Alternating current

Direct current

### 1.4 Safety Precautions

**CAUTION:** Read all safety warnings before installing, making connections, or servicing this equipment. Failure to comply with these warnings could result in personal injury and/or property damage. Retain all instructions for future reference.

- Verify that the AC adapter's input voltage range and plug type are compatible with the local AC power to be used.
- Make sure that the power cord does not pose a potential obstacle or tripping hazard.
- The balance is for indoor use only.
- Use the balance only in dry locations.
- Do not drop loads on the pan.
- Use only approved accessories and peripherals.
- Operate the equipment only under ambient conditions specified in these instructions.
- Disconnect the equipment from the power supply when cleaning.
- Do not operate the equipment in hazardous or unstable environments.
- Service should only be performed by authorized personnel.
- Do not position the balance such that it is difficult to reach the power connection.

## 2. INSTALLATION

### 2.1 Unpacking

Carefully remove your Adventurer balance and each of its components from the package. The included components vary depending on the balance model (see table below). Save the packaging to ensure safe storage and transport. Please read the manual completely before installing and using the Adventurer balance to avoid incorrect operation.

#### Included Components

- Balance
- Power Adapter
- Wind Ring (only for 0.1 mg and 1 mg models)
- Warranty Card
- Software Compact Disk

### 2.2 Selecting the Location

Avoid excessive vibrations, heat sources, air current, or rapid temperature changes. Allow sufficient space.

### 2.3 Leveling the Equipment

The Adventurer has a level bubble in a small round window beside the display.


To level the balance, adjust the 4 **Leveling Feet** until the bubble is centered in the circle. See the Level Assist section below on information about how and which feet to turn. Be sure the equipment is level each time its location is changed.

#### Level Assist

A level assist function is available to help leveling the Adventurer. There are two ways to access the function:

1. **Weighing Application** -> **Item Settings** -> **Level Assist**. See section 4.1.1 for more information.
2. **Main Menu** -> **Balance Setup** -> **User Settings** -> **Level Assist**. See section 5.3.3 for more information.

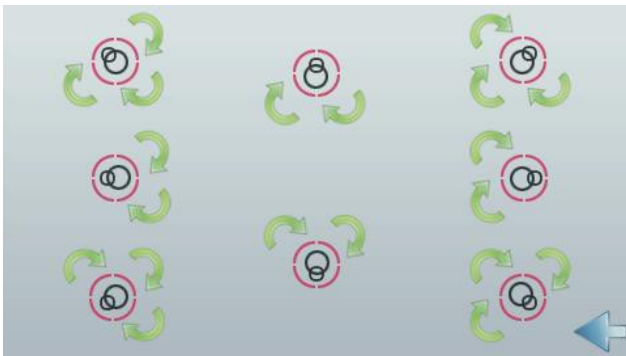Rotate the feet according to the image above depending on the location of the level bubble until the bubble is centered.

## 2.4 Connecting Power

Connect the DC output connector to the power receptacle on the rear of the balance. Then connect the AC power cord to a suitable electrical outlet.

**CAUTION:** For use with CSA certified (or equivalent approved) power source, which must have a limited current output.


**Attention:** For optimal weighing performance, allow the balance to warm up for 60 minutes prior to use.

## 2.5 Connecting the Interface

Use the built-in RS-232 Port to connect either to a computer or a printer with a standard (straight-through) serial cable. Or connect using the scale's USB port.

Interface connections on the rear of the balance:

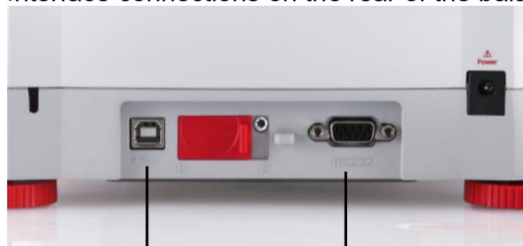USB1

RS232

USB connection on the front of the balance:

USB2

USB1: Used to connect to PC only

USB2: Used to connect a USB flash driver only

RS232: Used to connect to PC or Printer

**Note:** For configuration and interface commands, see the Communication Menu Settings section. For Connecting, Configuring and Testing the Printer/Computer Interface, and for sample Print Output Formats, see the Printing section.

## 2.6 Initial Calibration

When the Balance is first installed, and when it is moved to another location, it must be calibrated to ensure accurate weighing results. Most Adventurer Balances have built in AutoCal which can calibrate the balance automatically and does not require calibration masses. If preferred, the balance can be manually calibrated with external masses. Have the appropriate calibration masses available before beginning calibration. Refer to the Calibration Section for masses and calibration procedure.

### ▣ AutoCal™

Fully automatic internal calibration system that assists with routine maintenance by automatically calibrating the balance daily is available on most models.

Automatically calibrates the system when it senses a temperature change sufficient enough to affect weighing accuracy ( $>1.5^{\circ}\text{C}$ ), or every 11 hours.

### ▣ External Calibration


Select precision models feature traditional external calibration in which external weights (user's choice of calibration weight values) are used to calibrate the balance to ensure accuracy.

### 3. OPERATION

#### 3.1 Overview of Display, Home Screen

This equipment utilizes a touch-sensitive display with *Touch* areas and Buttons to control the equipment's functions.

#### CONTROLS
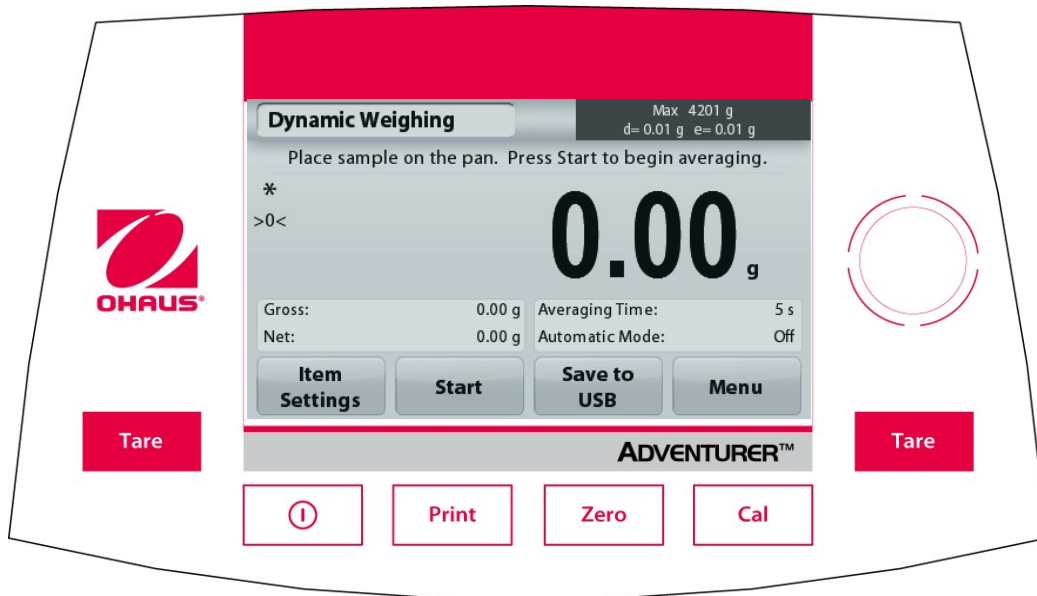
| Button | Action   |
|--------|--|
|        | Short Press (if powered Off): Turns on the scale<br>Long Press (if powered On): Turns off the scale<br><b>Note:</b> The balance will automatically power on when power is connected. |
|        | Short Press: Prints the present data to a printer or a computer.   |
|        | Short Press: Perform Zero operation  |
|        | Short Press: Perform Calibration operation   |
|        | Short Press: Perform Tare operation  |

#### Main Application Screen

|   |  |  |
|---|--|--|
| <p>Application</p>  |  | <p>Capacity and readability</p>  |
| <p>Instructional Messages</p> <p>Stability (*), Net (NET), Gross (G) and/or center of zero (&gt;0&lt;) indicators</p> |  | <p>Result Field: Information varies by application<br/>Touch <b>g</b> to change unit</p> |
| <p>Reference Fields</p>   |  | <p>Application Buttons: Functions vary by application</p>                                |

### 3.2 Principal Functions and Main Menu

**Weighing:** Press **Zero** to set the display to zero. Place an item on the pan. Display indicates gross weight.

**Taring:** With no load on the pan, press **Zero** to set the display to zero. Place an empty container on the pan and press **Tare**. Add material to the container and its net weight is displayed. Remove container and container's weight appears as a negative number. Press **Tare** to clear.

**Zero:** Press **Zero** to zero the balance

#### MENU & SCREEN NAVIGATION

Touch **Menu** to open the menu list.

#### Calibration:

Touch to view calibration options.

#### Balance Setup:

Touch to view and change balance settings.

#### Weighing Units:

Touch to view and change weighing units.

#### Data Maintenance:

Touch to view data maintenance settings.

#### GLP and GMP Data:

Insert user data for traceability.

#### Communication:

Touch to view COM Device Settings and Print Settings.

#### Factory Reset:


Touch to do a Factory reset of menu settings.


#### Lockout:

Touch to view lockout options.

### 3.3 Overview of Parts and Features – Draft Shield Models

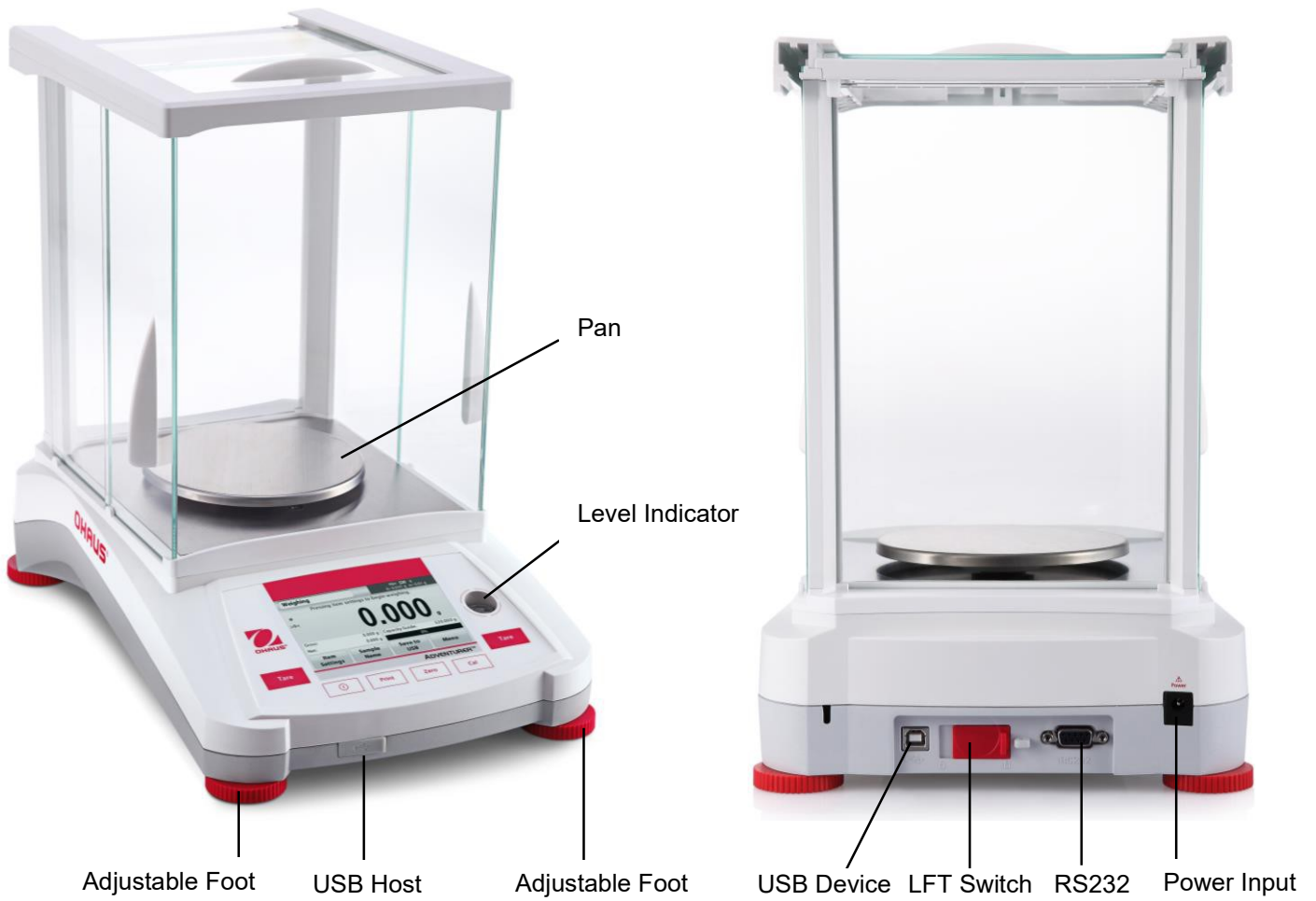
### 3.4 Overview of Parts and Features – Non Draft Shield Models

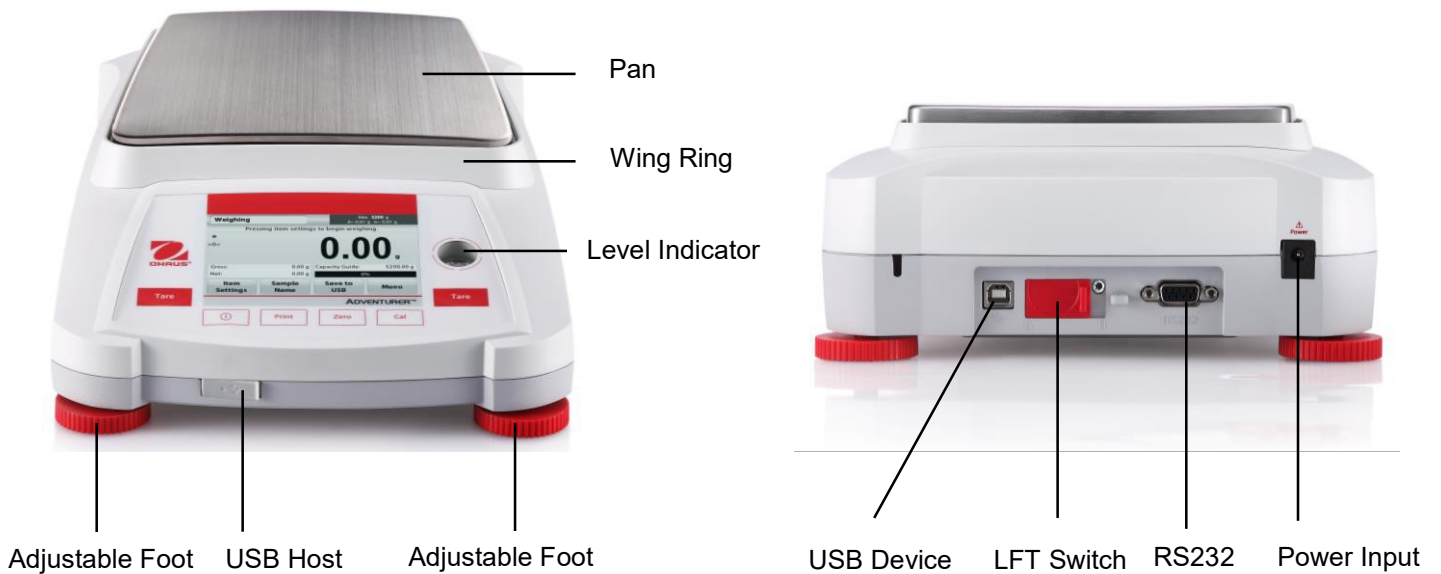
## 4. APPLICATIONS

The balance can be configured to operate in various Application modes. Touch the top left Application field (weighing in the example below):

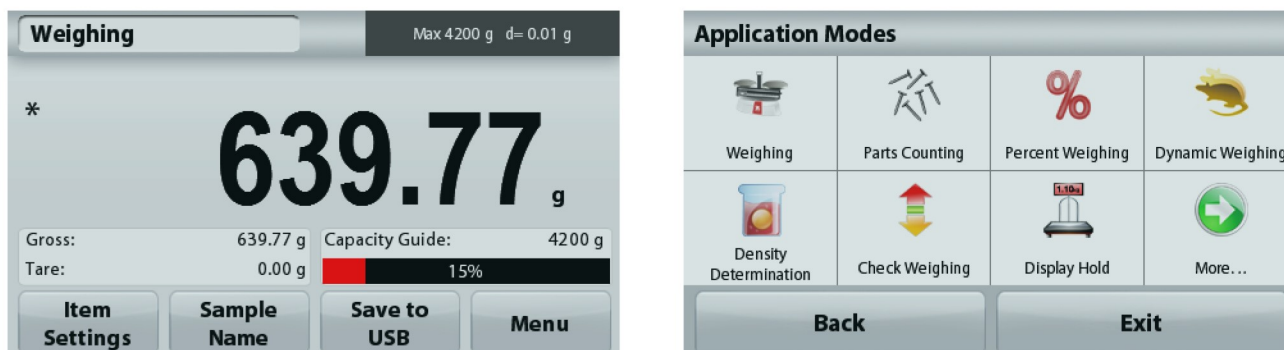The Adventurer Balance has 9 application modes, as follows:
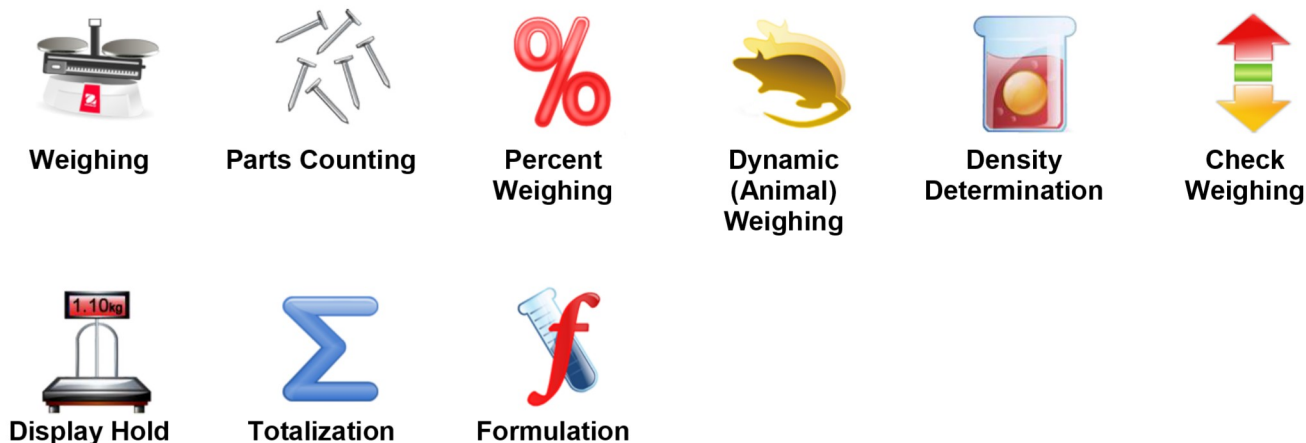
### 4.1 Weighing

**Note:** Before using any application, be sure the balance has been leveled and calibrated.

Use this application to determine the weight of items in the selected unit of measure.

#### Weighing

1. In the upper left portion of the home screen, select Weighing (this application is the default).
2. Press **Tare** or **Zero** if necessary to begin.
3. Place objects on the pan to display the weight. When stable, the \* appears.
4. The resulting value is displayed in the main Display Line in the active unit of measure.

The **WEIGHING** Home screen

Main Display Line

Touch **g** to change unit

Reference Fields

Application Buttons

Application Icon

4.1.1 Item Settings

To view or adjust the current settings

Touch the **Item Settings** button. The Settings screen appears.

**Minimum Weight:** establish a minimum weight value to be used to verify a reading. If an actual weight is below the established Minimum Weight value, it is flagged by a color change: **yellow**.

To adjust the Minimum Weight value, touch the **Minimum Weight** button.

A numeric input window appears.

Use the keys to enter the desired Minimum Weight, then press **Save**.

The display reverts to the previous screen.

To return to the Weighing home screen, touch **Exit** at the bottom of that screen.

**Capacity Bar:** When set to ON, a capacity bar is displayed in the reference field. The capacity guide will show the current weight as a percentage of balance capacity.

If Capacity Bar is set to OFF, the reference field will show Minimum Weight and Sample Name.

**Weighing Units:** Change the displayed unit. See section 5.4 for more information


**Note:** Touching the weighing unit from application home screen will also open the Weighing Units screen.

**Filter Level:** Change Filtering level. See section 5.3.4 for more information

**GLP & GMP Data:** See section 5.7 for more information

**Print settings:** Change printing settings. See section 7 for more information.


**Level Assist:** Instructions on how to move the balance feet to level the balance.

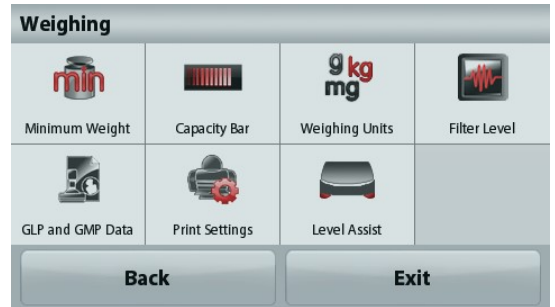
### 4.1.2 Sample Name

Press this button to add a sample name. An alphanumeric input window appears. Press  to alternate between Lower and Upper case characters.

#### Lower Case:
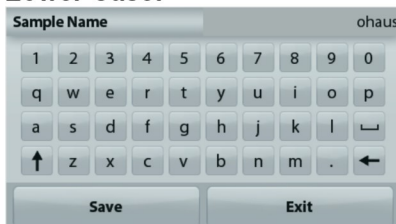
#### Upper Case:

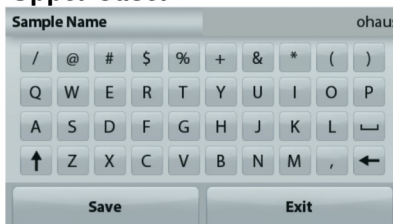Key in the desired sample name and press Save to save the name and return to weighing home screen.

### 4.1.3 Save to USB

Insert the USB flash drive into the USB slot located on the front of the balance. Next, press the Save to USB button to save the data to the USB flash drive. Once saved, the button will momentarily change color to orange.

**Notes:** The first time a USB flash drive is connected to the balance there might be some delay time before the button **Save to USB** works. This is due to that the balance is creating the necessary folders on the USB flash drive where the data will be stored. The Density Determination and Check Weighing applications do not have a Save to USB button.

#### CAUTION:

The weighing data will be saved to USB every day. However, if different weighing modes are used the data will be separately saved to individual files. Depending on the USB drive used, all data might not be transferred from the balance or the display might freeze. If this happens, unplug the USB flash drive and try another USB flash drive. Ohaus takes no responsibility if data on USB flash drive is erased or if the USB flash drive breaks while it is connected to the balance. To minimize the risk of problems arising, Ohaus suggests using a high quality USB flash drive.

## 4.2 Parts Counting

**Note:** Before using any application, be sure the balance has been leveled and calibrated.

Use this application to count samples of uniform weight.

### Parts Counting

1. In the upper left portion of the home screen, select Parts Counting
2. Press **Tare** or **Zero** if necessary to begin.
3. Place objects on the pan to display the weight. When stable, the \* appears.
4. The resulting value is displayed in the main Display Line in pieces (PCS).

The **Parts Counting** Home screen

Main Display Line

Reference Fields

Application Buttons

Application Icon

4.2.1 Item Settings

To view or adjust the current settings

Touch the **Item Settings** button. The Settings screen appears.

**Samples:** The sample size can be 1 to 10 000 pieces. The default sample size is 10. Once a sample size is changed, the balance will immediately open the recalculate APW screen, expecting to establish a new APW.

To adjust the sample size, touch the **Samples** button.

A numeric input window appears.


Use the keys to enter the desired sample size, then press **Save**.

The next screen appears, with the message to place the reference weight on the pan.

Place the reference weight on the pan, then touch **Accept** to capture the value, the screen shows number of pcs.

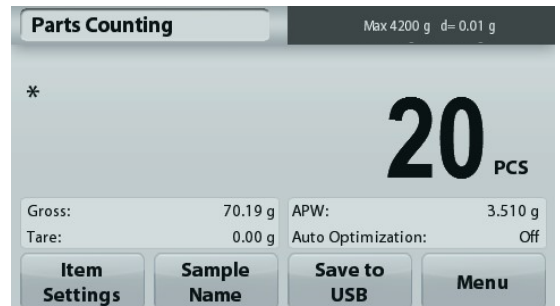**Establish an Average Piece Weight (APW):**

Each time a new type of part is counted, the nominal weight of one piece (Average Piece Weight or APW) must be established using a small quantity of pieces. This APW is stored until replaced by another APW.

There are two methods to establish the APW value:

1. The actual APW is known
2. The APW must be derived by weight. For this case the current sample size will be used

**Set a known Average Piece Weight (APW)**

To adjust the APW value directly, touch the **APW** button. A numeric input window appears.

Key in the desired APW Weight, then press **Save**.

The display returns to the application home screen with the new APW value in the reference field.

**Set a new Average Piece Weight (APW) – Derived**

To establish a new APW, place the reference weight on the pan and press **Accept** to continue.

**Note:** The sample size that is displayed will be used. To use a different sample size, change that first. (See above.)

The home screen shows the number of pieces at the new APW

**Auto Optimization:** When set to **On**, improves counting accuracy by re-calculating the piece weight automatically as parts are added. Default is **Off**.

**Notes:**

APW Optimization occurs only when the number of pieces added to the pan is between one and three times the number already on the pan.

If the APW is entered manually by the numeric keypad, APW auto optimization will not occur.

**Print settings:** Change printing settings. See section 7 for more information.

**Note:** The **Save to USB** button will only appear after a USB flash drive has been connected to the balance. See section 4.1.3 for more information.

### 4.3 Percent Weighing

**Note:** Before using any application, be sure the balance has been leveled and calibrated.

Use Percent Weighing to measure the weight of a sample displayed as a percentage of a pre-established reference weight.

The default (or last) reference weight is displayed.

#### Percent Weighing

1. In the upper left portion of the home screen, select Percent Weighing.
2. Place an object on the pan. The difference between the sample and the reference weight is displayed as a percentage.

The **Percent Weighing** Home screen

Main Display Line

Reference Fields

Application Buttons


Application Icon

#### 4.3.1 Item Settings

To view or adjust the current settings

Touch the **Item Settings** button. The Settings screen appears.

**Recalculate Ref. Wgt:** To establish a new reference weight value, touch the recalculate reference weight button.

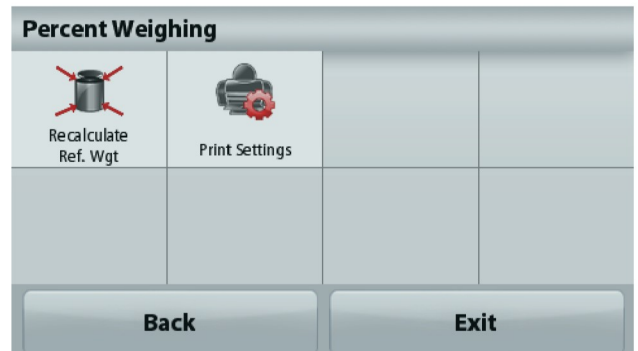Follow the screen instructions to establish a new reference weight.

Alternatively, press the **Ref. Weight** button from the Percent Weighing Recalculate Ref. Wgt screen to establish a new reference weight manually through a numerical keypad.

**Print settings:** Change printing settings. See section 7 for more information.

**Note:** The **Save to USB** button will only appear after a USB flash drive has been connected to the balance. See section 4.1.3 for more information.

## 4.4 Dynamic Weighing

**Note:** Before using any application, be sure the balance has been leveled and calibrated.

Use this application to weigh an unstable load, such as a moving animal. Two different start/reset modes can be selected: **Manual** (start and stop via key press) and **Automatic** (start and stop automatically).

### Dynamic Weighing – Manual (default)

1. In the upper left portion of the home screen, select Dynamic Weighing
2. Place objects on the pan and press the **Start** button.

The **Dynamic Weighing** Home screen

Main Display Line

Reference Fields

Application Buttons

Application Icon

3. The balance begins a countdown (averaging process). During the countdown, the information line displays the time remaining.  
If necessary, press **Stop** to quit.
4. When the countdown ends, the result line is displayed and held. Press **Clear** to clear the held weight and return to the initial screen.

**Note:** Clear the pan before beginning a new Dynamic weighing cycle.

### 4.4.1 Item Settings

To view or adjust the current settings

Touch the **Item Settings** button. The Settings screen appears.


**Averaging Time:** Set the averaging time to a value between 1 and 99 seconds. Default is 5 seconds.

**Automatic Mode:** When set to On, the cycle begins when an object is placed on the pan, and the held value is automatically reset when the object is removed from the pan.

**Sample Name:** Assign a name to the sample.

**Print settings:** Change printing settings. See section 7 for more information.

**Note:** The **Save to USB** button will only appear after a USB flash drive has been connected to the balance. See section 4.1.3 for more information.

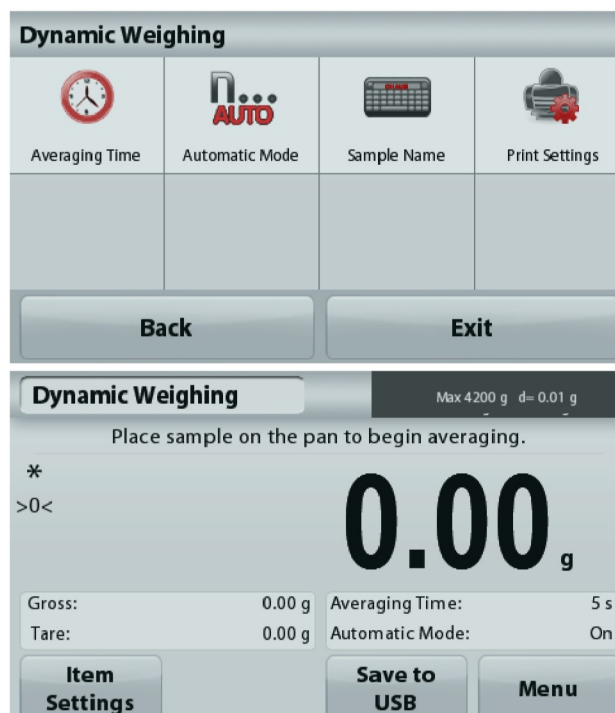
## 4.5 Density Determination

**Note:** Before using any application, be sure the balance has been leveled and calibrated.

Use this application to determine an object's density. Four types of density determination can be made:

1. Solids more dense than water
2. Solid less dense than water
3. Liquid density
4. Porous material (impregnated with oil)

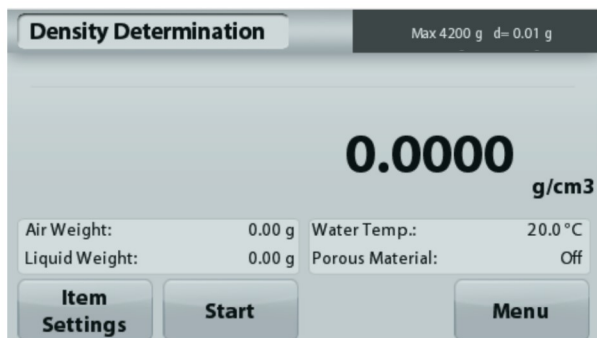The **Density Determination** Home screen

Main Display Line

Reference Fields

Function Buttons

Application Icon

### 4.5.1 Measuring the Density of a Solid Using Water (default)

Press the **Item Settings** button to open the Density Determination **Settings** screen.

Confirm the following **Setups** are selected:

- ✓ **Density Type: Solid**
- ✓ **Auxiliary Liquid: Water**
- ✓ **Porous Material: Off**



To adjust the water temperature value, touch the **Water Temp.** button.

The balance calculates water density based on the water temperature value entered.

Measure the actual water temperature using a precision thermometer.

A numeric input window appears.

Enter the actual water temperature, then press **Save**.

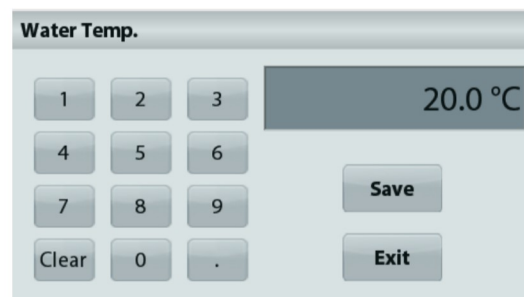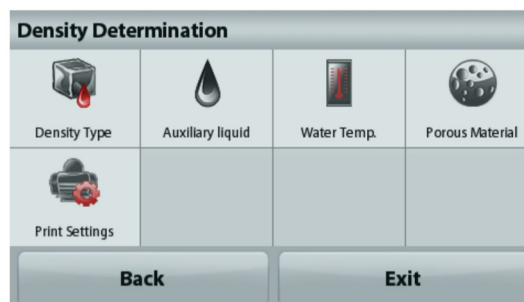To return to the Density Determination home screen, touch **Back**.

**Step 1 of 2 – Weigh the Sample in Air.**  
Press **Start**. Follow screen instructions, then press **Accept** to store the dry sample weight (“in air”).

**Step 2 of 2 – Weigh the Sample Submerged in the Liquid.**  
 Follow the screen instructions, then press **Accept** to store the sample weight (submerged in liquid).

Once the necessary weights have been determined, the density of the sample is displayed in **g/cm<sup>3</sup>** (along with the weight in air, weight in water) on the Application screen.

Press **Start** to reset all the weight values and restart the process.

**4.5.2 Measuring the Density of a Buoyant Material Using Water**

Press the **Item Settings** button to open **Settings** screen.  
 Confirm the following **Setups** are selected:

- ✓ **Density Type: Solid**
- ✓ **Liquid Type: Water**
- ✓ **Porous Material: Off**

Press **Back** to return to the Density Determination home screen.

Follow the same procedure as Solid Material above, except in Density Determination step 2, **push the sample down** into the liquid until it is fully submerged.

**4.5.3 Measuring the Density of a Solid Using an Auxiliary Liquid**

To enable this feature, enter the Density Determination Setup menu and select the following: **Density Type: Solid; Liquid Type: Other; Porous Material: Off.**

Confirm the default values displayed (Liquid Density, etc) are correct.

To adjust the Liquid Density value, touch the **Auxiliary liquid** button and select Other.


A numeric input window appears.

Key in the density in g/cm<sup>3</sup>, then press **Save**.

The display returns to the previous screen.

To return to the Density Determination home screen, touch **Back**.

Begin the Density determination process as per above.

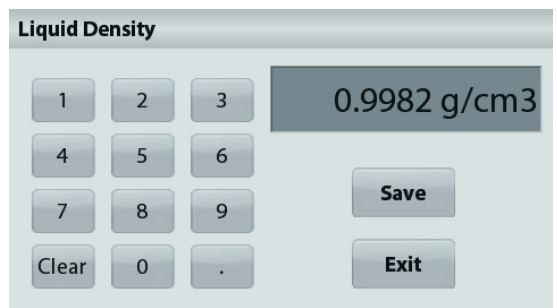**4.5.4 Measuring the Density of a liquid using a Calibrated Sinker (not supplied)**

To enable this feature, enter the Density Setup menu and select the following; **Density Type: Liquid**.

**Note:** when the Density Type is set to Liquid, the Liquid type and Porous material selections are disabled.

The **DENSITY DETERMINATION – LIQUID** Home screen

Main Display Line

Reference Fields


Function Buttons

Application Icon

Confirm that the default value displayed (sinker volume) is correct. To edit the default values, touch **Sinker Volume**.

To adjust the Sinker volume value, touch the **Sinker Volume** button.

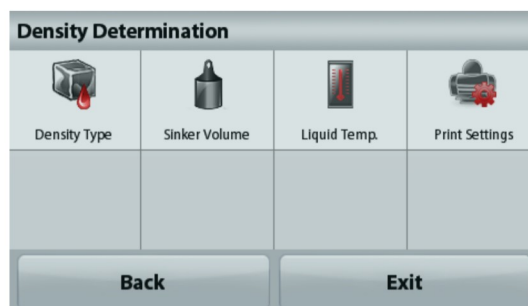
A numeric input window is displayed.

Key in the desired value, then press **Save**.

The display returns to the previous screen with the new value highlighted.


To return to the Density Determination home screen, touch **Back**.

Press **Start** to start the process.

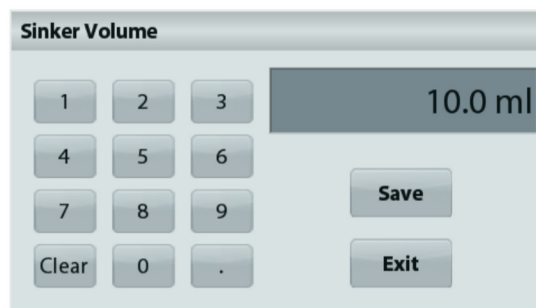
**Step 1 of 2 – Weigh the Sinker in Air.**

Follow the screen instructions, then press **Accept** to store the sinker weight (“in air”).

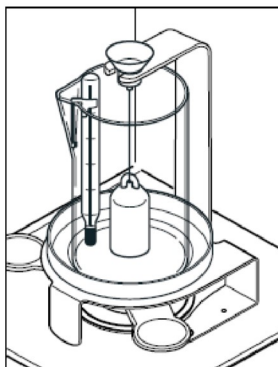**Step 2 of 2 – Weigh the Sinker Submerged in the Test Liquid.**

Follow the screen instructions, then press Accept to store the sinker weight (submerged in liquid).

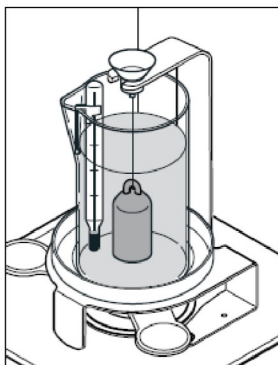Once the necessary weights have been determined, the density of the Liquid sample is displayed in **g/cm<sup>3</sup>** (along with the weight in air, weight in water) on the Application screen.

Press **Start** to reset all the weight values and restart the process.

**4.5.5 Measuring the Density of Porous Material Using Oil**

To enable this feature, enter the Density Determination **Setup** menu, and set the following:

- ✓ **Density Type: Solid**
- ✓ **Liquid Type: Water**
- ✓ **Porous Material: On**

The **DENSITY DETERMINATION – POROUS** Home screen

Main Display Line

Reference Fields

Functions

Application Icon

Confirm the default values displayed (Water Temp) are correct.

To edit the default values, touch **Item Settings**.

The Settings screen appears.

The balance calculates water density based on the water temperature value entered (look-up table).

Measure the actual water temperature using a precision thermometer.

To adjust the Water Temperature or Oil Density values, touch the **Water Temp** or **Oil Density** button.

Numeric input windows appear.

Key in the desired value, then press **Save**.

The display returns to the previous screen with the new value highlighted.

To return to the Density Determination home screen, touch **Exit**.

Press **Begin Density Calculation**.

**Step 1 of 3 – Weigh the Un-Oiled Sample in Air.**

Follow the screen instructions then press **Accept** to store the dry sample weight (in air).

**Step 2 of 3 – Weigh the Oiled Sample in Air.**

Follow the screen instructions then press **Accept** to store the sample weight (oiled).

**Step 3 of 3 – Weigh the Oiled Sample Submerged in Liquid.**

Follow the screen instructions then press **Accept** to store the oiled sample weight (submerged in liquid).

Once the necessary weights have been determined, the density of the sample is displayed in **g/cm<sup>3</sup>** (along with the weight in air, un-oiled and oiled, and weight in water) on the Application screen.

The value stays on the display until **Start** is touched.

Press **Start** to reset all the weight values and restart the process.

### 4.6 Check Weighing

**Note:** Before using any application, be sure the balance has been leveled and calibrated.

Check Weighing is used to compare the weight of a sample against target limits.

#### Check Weighing

1. In the upper left portion of the home screen, select Check Weighing
2. The default (or last) Check weight limits are displayed.
3. Place objects on the pan.
4. The Under/Accept/Over status is shown in the progress bar area while the actual weight of the item is shown on the main Display Line.

The **Check Weighing** Home screen

Main Display Line

Reference Fields

Function Buttons

Application Icon

To set the *Over Limit value*, touch the **Over Limit** button  
 To set the *Under Limit value*, touch the **Under Limit** button  
 A numeric input window appears.

Enter the desired Limit Weight, then press **Save**.

To return to the CHECK WEIGHING home screen, touch **Exit**.


#### 4.6.1 Item Settings

To view or adjust the current settings

Touch the **Item Settings** button. The Settings screen appears.

**Sample Name:** Assign a name to the sample.

**Print settings:** Change printing settings. See section 7 for more information.
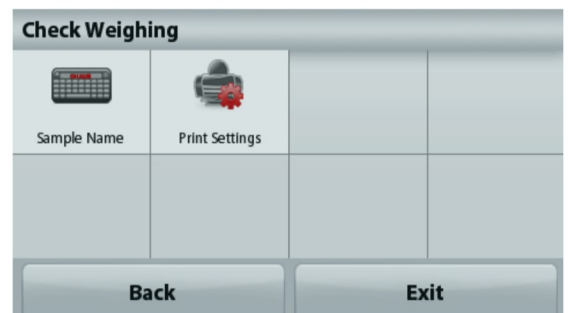
## 4.7 Display Hold

**Note:** Before using any application, be sure the balance has been leveled and calibrated.

Two Modes are available:

**Display Hold** - allows the user to capture and store a stable weight.

**Peak Hold** - allows the user to capture and store the highest stable weight.


The **Display Hold** Home screen

Main Display Line

Reference Fields

Function Buttons

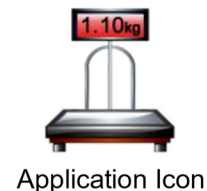Application Icon

## Display Hold

### 4.7.1 Display Hold

1. In the upper left portion of the home screen, select Display Hold
2. Place the sample on the pan and press **Hold** at any time while the weight is being captured.
3. The Main Display Line now shows the first stable weight.
4. Press **Clear** to remove the hold and return to Display Hold Home screen.

## Peak Hold

### 4.7.2 Peak Hold

1. In the upper left portion of the home screen, select Display Hold
2. Choose Peak Hold Mode in Item Settings (see section 4.7.3).
3. Place sample on the pan and press **Start** to begin.
4. Continue to weigh samples. The highest stable weight will be held.
5. To remove the hold and return to normal operation press **Stop**.

### 4.7.3 Item Settings

To view or adjust the current settings

Touch the **Item Settings** button. The Settings screen appears.

**Mode:** Choose between Peak Hold and Display Hold (default).

**Sample Name:** Assign a name to the sample.

**Print settings:** Change printing settings. See section 7 for more information.

**Note:** The **Save to USB** button will only appear after a USB flash drive has been connected to the balance. See section 4.1.3 for more information.

### 4.8 Totalization

**Note:** Before using any application, be sure the balance has been leveled and calibrated.

Totalization measures the cumulative weight of a sequence of items. The cumulative total may exceed the capacity of the Balance. The maximum number of samples (n) is 99.

The **Totalization** Home screen

Main Display Line

Reference Fields

Application Buttons

Application Icon

### Totalization

1. In the upper left portion of the home screen, select Totalization
2. Place item on the pan to begin. The sample weight is shown on the Main Display Line.
3. Press **Accumulate** to add the weight (when stable) of the item to the total.
4. Remove the item from the weighing pan, then add the next item and continue as above.
5. Press **Result** to view the results from the totalization.
6. When finished, press the Clear Total button to reset the accumulated weight to zero.

| Totalization |          |
|--------------|----------|
| Item         | Result   |
| Samples      | 3        |
| Total        | 506.95 g |
| Average      | 168.98 g |
| σ (stdev)    | 38.90 g  |
| Minimum      | 117.00 g |
| Maximum      | 210.57 g |
| Range        | 93.57 g  |

7. Press **Save to USB** to save the result to a USB flash drive or **Exit** to return to the Totalization Home screen.

**Note:** Changing units converts the Accumulation results to the selected unit.

## 4.9 Formulation

**Note:** Before using any application, be sure the balance has been leveled and calibrated.

Use this application for compounding and recipe making. The maximum number of components is 50.

The **Formulation** Home screen

Main Display Line

Reference Fields

Application Buttons

Application Icon

### Formulation

1. In the upper left portion of the home screen, select **Formulation**
2. Press **Start** to begin the compounding process.
3. Place the first ingredient on the pan (or in a tared container) and press **Accept** to store the component.
4. Continue adding components and pressing **Accept** to store the weight of the individual components until the formula is complete. The **Total** line shows the total weight of all the components.
5. Press **Stop** to finish the Formulation. The Formulation results are displayed:

**Note:** If Filler is set to On (see section 4.9.1 below), the balance will ask to add a filler material to complete the formulation. Add the filler material and press **Accept** to complete the formulation and display the results.

| Formulation |          |
|-------------|----------|
| Comp.       | Weight   |
| 1           | 241.76 g |
| 2           | 272.95 g |
| Comp. Total | 514.71 g |

Save to USB      Exit

6. Press **Save to USB** to save the results to a USB flash drive or **Exit** to return to the Formulation Home screen.

**Note:** The formulation results will be cleared when a new formulation starts.

### 4.9.1 Item Settings


To view or adjust the current settings

Touch the **Item Settings** button. The Settings screen appears.

**Filler:** If set to On, a filler material is asked for at the end of the formulation.

**Automatic Mode:** If set to On, the balance will automatically Tare after the component weight has been accepted.

**Print settings:** Change printing settings. See section 7 for more information.
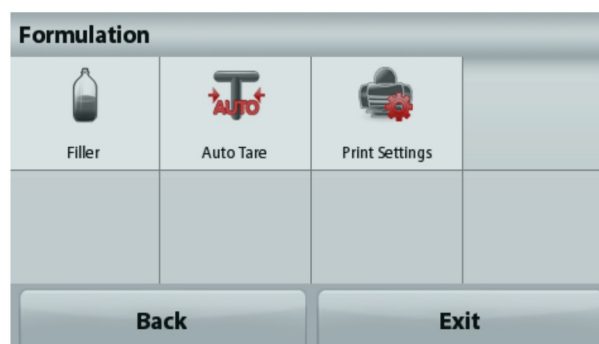
## 4.10 Additional Features

### 4.10.1 Weigh Below


The Adventurer balance is equipped with a weigh below hook for weighing below the balance.

Note: Before turning the balance over, remove the pan and draft shield elements (if present) to prevent damage

**Attention: Do not place the balance on the pan support cone or Load cell Pins**

To use this feature, remove power from the balance, then remove the protective cover for the weigh below opening.

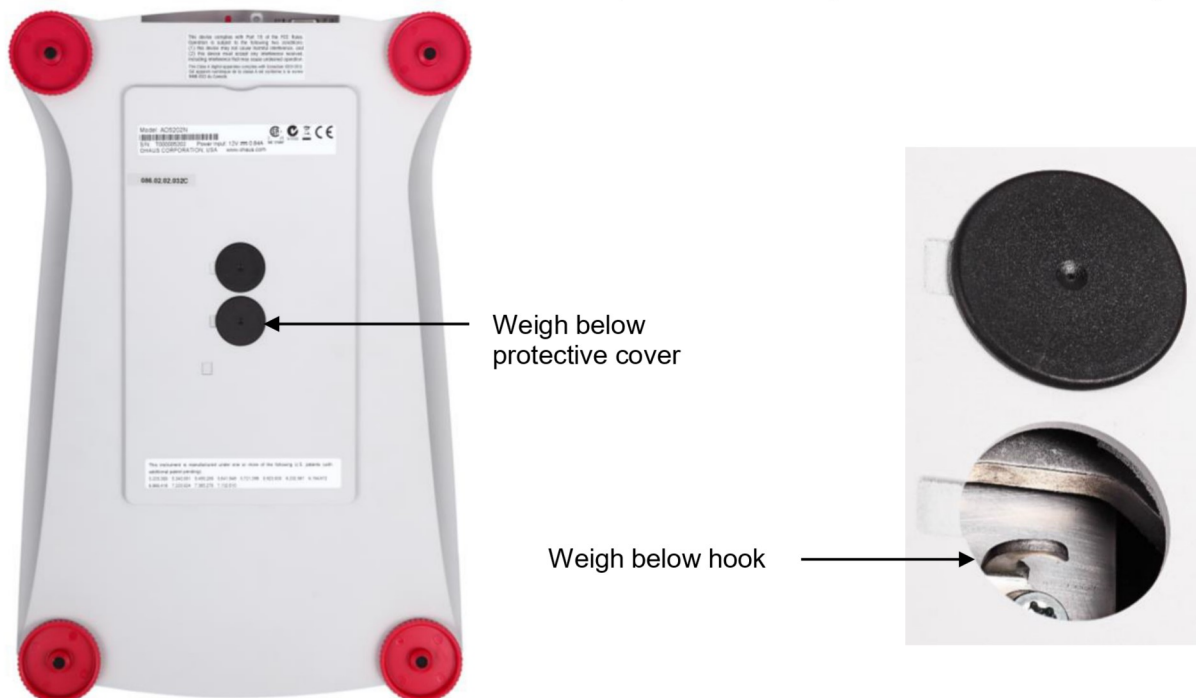The balance can be supported using lab jacks or any other convenient method. Ensure the balance is level and secure. Power on the balance, then use a string or wire to attach items to be weighed.

5. MENU SETTINGS

5.1 Menu Navigation

User menu structure:

| Application Modes   | Main Menu        | Calibration   | Balance Setup      | Weighing Units   | Data Maintenance   | Communication      | GLP and GMP Data | Factory Reset | Lockout          |
|---------------------|------------------|---------------|--------------------|------------------|--------------------|--------------------|------------------|---------------|------------------|
| Weighing App        | Calibration      | Internal Cal  | Language           | Milligram        | Export to USB      | RS232 Standard     | Headers          |               | Calibration      |
| Min Weight          | Balance Setup    | Automatic Cal | User Settings      | Gram             | App. Mode Settings | Baud Rate          | Header 1         |               | Balance Setup    |
| Capacity Guide      | Weighing Units   | AutoCal™      | Touch Calibrate    | Kilogram         | Menu Settings      | 2400               | Header 2         |               | Weighing Units   |
| Units               | Data Maintenance | Span Cal      | Brightness         | Carat            | Import from USB    | 4800               | Header 3         |               | Data Maintenance |
| Filter Level        | Communication    | Linearity Cal | Beep               | Ounce            | App. Mode Settings | 9600               | Header 4         |               | Communication    |
| GLP and GMP Data    | GLP and GMP Data | Cal Test      | Auto Dim           | OunceTroy        | Menu Settings      | 19200              | Header 5         |               | GLP and GMP data |
| Print Settings      | Factory Reset    |               | Level Assist       | Pound            | Balance_Info       | 38400              | Balance Name     |               | Factory Reset    |
| Level Assist        | Lockout          |               | Filter Level       | Pennyweight      |                    | Transmission       | User Name        |               |                  |
| Counting App        |                  |               | Auto Zero Tracking | Grain            |                    | 7 E 1              | Project Name     |               |                  |
| Sample Size         |                  |               | Auto_Tare          | Newton           |                    | 7 E 2              |                  |               |                  |
| APW                 |                  |               | Graduation         | Momme            |                    | 7 N1               |                  |               |                  |
| Auto Opt            |                  |               | Date & Time        | Mesghal          |                    | 7 N2               |                  |               |                  |
| Print Settings      |                  |               | Date               | HKTael           |                    | 7 O 1              |                  |               |                  |
| Percent App         |                  |               | Time               | SGTael           |                    | 7 O2               |                  |               |                  |
| Recalculate Ref Wgt |                  |               | Approved Mode      | TWTael           |                    | 8 N1               |                  |               |                  |
| Print Settings      |                  |               |                    | Tical            |                    | 8 N 2              |                  |               |                  |
| Dynamic App         |                  |               |                    | Tola             |                    | Handshake          |                  |               |                  |
| Average Time        |                  |               |                    | Baht             |                    | Print Settings     |                  |               |                  |
| Auto Mode           |                  |               |                    | Custom1          |                    | Print Output       |                  |               |                  |
| Sample Name         |                  |               |                    | Unit Name        |                    | Stable Weight Only |                  |               |                  |
| Print Settings      |                  |               |                    | Factor           |                    | Numeric Value Only |                  |               |                  |
| Density App         |                  |               |                    | Exponent         |                    | Single Header Only |                  |               |                  |
| Density Type        |                  |               |                    | 10 <sup>-3</sup> |                    | Print Options      |                  |               |                  |
| Auxiliary Liquid    |                  |               |                    | 10 <sup>-2</sup> |                    | Auto Print         |                  |               |                  |
| Water Temp.         |                  |               |                    | 10 <sup>-1</sup> |                    | Auto Print Off     |                  |               |                  |
| Porous Material     |                  |               |                    | 10 <sup>0</sup>  |                    | On Stability       |                  |               |                  |
| Oil Density         |                  |               |                    | 10 <sup>1</sup>  |                    | Interval (seconds) |                  |               |                  |
| Sinker Volume       |                  |               |                    | 10 <sup>2</sup>  |                    | Continuous         |                  |               |                  |
| Liquid Temp.        |                  |               |                    | 10 <sup>3</sup>  |                    | Print Content      |                  |               |                  |
| Print Settings      |                  |               |                    | LSD              |                    | Selection          |                  |               |                  |
| CheckWeighing App   |                  |               |                    | 0.5              |                    | Header             |                  |               |                  |
| Sample Name         |                  |               |                    | 1                |                    | Date & Time        |                  |               |                  |
| Print Settings      |                  |               |                    | 2                |                    | Balance ID         |                  |               |                  |
| DisplayHold App     |                  |               |                    | 5                |                    | Balance Name       |                  |               |                  |
| Display Hold Mode   |                  |               |                    | 10               |                    | User Name          |                  |               |                  |
| Sample Name         |                  |               |                    | 100              |                    | Project Name       |                  |               |                  |
| Print Settings      |                  |               |                    |                  |                    | Application Name   |                  |               |                  |
| Totalization App    |                  |               |                    |                  |                    | Sample Name        |                  |               |                  |
| Formulation App     |                  |               |                    |                  |                    | Result             |                  |               |                  |
| Filler              |                  |               |                    |                  |                    | Gross              |                  |               |                  |
| Automatic Mode      |                  |               |                    |                  |                    | Net                |                  |               |                  |
| Print Settings      |                  |               |                    |                  |                    | Tare               |                  |               |                  |
|                     |                  |               |                    |                  |                    | Information        |                  |               |                  |
|                     |                  |               |                    |                  |                    | Signature Line     |                  |               |                  |
|                     |                  |               |                    |                  |                    | Line Feed          |                  |               |                  |
|                     |                  |               |                    |                  |                    | Save To USB        |                  |               |                  |

All menu navigation is performed by touching the display. To enter the Menu, touch **Menu** from any Application Home screen. The Main menu appears, with buttons for **Back** and **Exit**. Continue touching the appropriate list item to navigate to the Menu items.

**5.1.1 Changing Settings**

To change a menu setting, navigate to that setting using the following steps:

**Enter the Menu**

From any Application screen, Touch **Menu**. The Main Menu List appears on the display.

**Select the Sub-Menu**

Find the item of the Main Menu List and touch it. The Sub-Menu appears.

**Select the Menu Item**

Continue until the desired setting is chosen in the Menu list. Touch the setting to change it. The changed setting will be displayed as highlighted yellow for about 1 second to confirm the changed value.

**Exit the Menu and Return to the Current Application**

After the setting is confirmed, touch **Exit** to return to the Application.

**Note:** at any time the **Back & Exit** buttons can be touched to navigate to the desired area of the menu or return to the current Application. Continue until the desired setting is chose in the menu list.

The Adventurer balance Main menu structure is illustrated below.


## 5.2 Calibration

Adventurer Balances (InCal models) offer a choice of six calibration methods: Internal Calibration, Automatic Calibration, AutoCal™ Adjustment, Span Calibration, Linearity Calibration and Cal Test (Calibration Test).

**Note:** The calibration unit is always in grams.

**Attention:** Do not disturb the balance during any calibration.

### 5.2.1 Calibration sub-menu (InCal models)

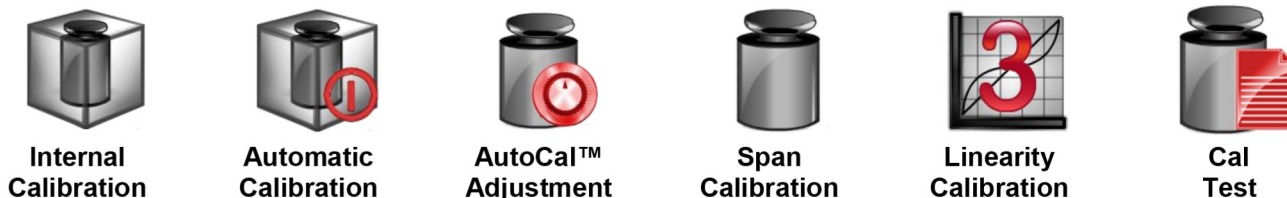**Note:** /E models only have Span Calibration, Linearity Calibration and Cal Test.

### 5.2.2 Internal Calibration

Calibration is accomplished with the internal calibration mass. Internal Calibration can be performed at any time, provided the balance has warmed up to operating temperature and is level.

With the Balance turned ON and no load on the pan, touch **Internal Calibration**. Alternatively, press the Cal key on the balance.

The Balance begins to calibrate. And no load on the pan, touch Internal Calibration.

The display shows the status, then returns to the current application.

To cancel at any time, press **Save**.

### 5.2.3 Automatic Calibration

When **Automatic Calibration** is set ON, the balance performs a self-calibration:

- when it senses a temperature change of 1.5°C
- or every 11 hours

AutoCal will automatically calibrate the Balance (using the internal mass) each time there is a change in temperature significant enough to affect accuracy or every 11 hours.

An information screen will appear when an Automatic Calibration is about to start. Three option buttons will be displayed:

- Now** – Press to perform the calibration directly.
- 5 min** – Press to perform the calibration after 5 minutes.
- Deactivate** – Press to deactivate the Automatic Calibration function.

#### 5.2.4 AutoCal™ Adjustment

Use this calibration method to adjust the span calibration point, without affecting the span or linearity calibration.

Calibration Adjust may be used to adjust the result of the Internal Calibration by  $\pm 100$  divisions.

**Note:** Before making a calibration adjustment, perform an Internal Calibration. To verify whether an adjustment is needed, place a test mass equal to the **span calibration value** on the pan and note the difference (in divisions) between the nominal mass value and the actual Balance reading. If the difference is within  $\pm 1$  division, calibration adjustment is not required. If the difference exceeds  $\pm 1$  division, calibration adjustment is recommended.

*Example:*

|                              |                           |
|------------------------------|---------------------------|
| Actual weight reading:       | 200.014                   |
| Expected weight reading:     | 200.000 (Test mass value) |
| Difference Weight (d):       | 0.014                     |
| Difference weight in digits: | -14 (Adjust value)        |

To perform a Calibration Adjustment, touch AutoCal Adjustment from the Calibration Menu; Enter the value (positive or negative divisions) to match the difference noted earlier in the procedure.

Recalibrate using Internal Calibration. After calibration, place the test mass on the pan and verify that the mass value now matches the displayed value. If not, repeat the procedure until Internal Calibration reading agrees with the test mass.

Once completed, the balance stores the Adjustment value and the display returns to the current application.

#### 5.2.5 Span Calibration

Span calibration uses two calibration points, one at **zero load** and the other at **specified full load** (span). For detailed calibration mass information please see the specification tables in section 9.

With the balance turned ON and no load on the pan, touch Span Calibration to initiate the procedure. Additional calibration values to be used are shown on the display. The best accuracy is achieved using the mass closest to the full span value.

**Note:** To change the span calibration point, touch the alternate weight shown on the display. Follow the screen instructions and place the specified calibration weight on the scale when prompted to do so. When complete, the display shows the Span calibration status and returns to the current application.

### 5.2.6 Linearity Calibration

Linearity calibration uses three calibration points, one at zero load and the others at specified loads. With no load on the scale, press Linearity Calibration to begin the process. The balance captures the zero point, and then prompts for the next weight. Continue to follow the instructions on the display until the calibration is completed.

To cancel at any time, press **Cancel**.

### 5.2.7 Calibration Test

Use Calibration Test to compare a known calibration weight against the stored span calibration data.

With no load on the scale, press **Cal Test** to begin the process. The balance captures the zero point, then prompts for the span value.

The display shows status, followed by the difference between the current calibration weight and the stored calibration data.

## 5.3 Balance Setup

Enter this sub-menu to customize Balance functionality.

### 5.3.1 Balance Set-up sub-menu

Language

User Settings

Filter Level

**AZT**

Auto Zero Tracking

Auto Tare

**I/O**

Graduations

Date & Time

Approved Mode


Factory default settings are shown below in bold.

### 5.3.2 Language

Set the language displayed for menus and displayed messages.

English  
German  
French  
Spanish  
Italian  
Polish  
Turkish  
Czech  
Hungarian
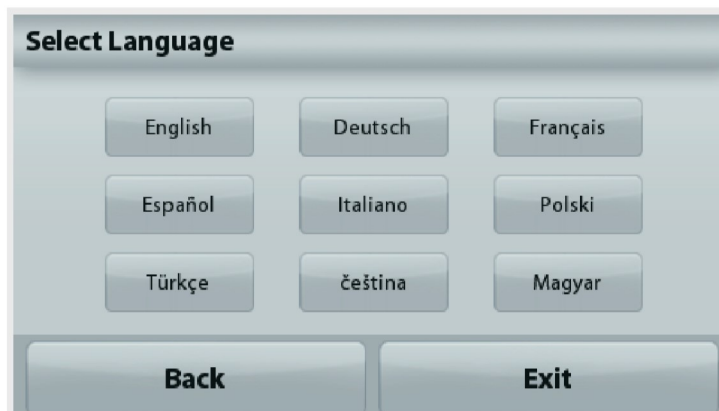
### 5.3.3 User Setting

Use this sub-menu to change the setting for:

#### Touch Calibrate

“Runtime calibration, please touch the screen at the center of the ring”  
(First top-left, then bottom-right.)

#### Screen Brightness:

LOW = low screen brightness.

**MEDIUM** = normal screen brightness.

HIGH = high screen brightness.

#### Beep:

**OFF** = disabled

ON = enabled

#### Auto Dim (Dims the display if no Screen Activity for x minutes)

**OFF** = disabled

10 min

20 min

30 min

**Level Assist:** Instructions on how to move the balance feet to level the balance.

### 5.3.4 Filter Level

Set the amount of signal filtering.

LOW = faster stabilization time with less stability.

**MEDIUM** = normal stabilization time with normal stability.

HIGH = slower stabilization time with more stability.

### 5.3.5 Auto Zero Tracking

Set the automatic zero tracking functionality.

OFF = disabled.

**0.5 D** = display maintains zero up to a drift of 0.5 graduation per second

1 D = display maintains zero up to a drift of 1 graduation per second.

3 D = display maintains zero up to a drift of 3 graduations per second.

### 5.3.6 Auto Tare

Set the automatic Tare.

**OFF** = disabled.

ON = enabled.

‘Place container on the pan’ will be displayed when an Automatic Tare is about to start.

A **Deactivate** button is displayed underneath the text. Press this button to deactivate the Automatic Tare function

### 5.3.7 Graduations

Set the displayed readability of the balance.

**1 Division** = standard readability.

10 Divisions = readability is increased by a factor of 10.

For example, if the standard readability is 0.01g, selecting 10 Divisions will result in a displayed reading of 0.1g.

### 5.3.8 Date & Time

Set the current Date and Time.

Change the format (if desired), then enter the current value.

Press **Save** to confirm the new value.

### 5.3.9 Approved Mode

Use this menu to set the Legal for Trade status.

**OFF** = standard operation.

**ON** = operation complies with Legal Metrology regulations.

**Note:** When Approved Mode is set to ON, the menu settings are affected as follows:

Calibration Menu:

AutoCal internal Calibration is forced to ON and hidden. Internal Calibration and Calibration Test functions are available. All other functions are hidden.

Balance Setup Menu:

Filter Level is locked at the current setting.

Auto Zero Tracking is limited to 0.5 Division and OFF. The selected setting is locked.

Auto Tare is locked at current setting.

Graduations is forced to 1 Division and the menu item is hidden.

Communication Menu (Communication->Print Settings->Print Output):

Stable Weight Only is locked ON.

Numeric Value Only is locked OFF.

Communication Menu (Communication->Print Settings->Auto Print):

Auto print mode selections are limited to OFF, On Stability, and Interval.

Continuous is not available.


Data Maintenance Menu:

Export to USB is hidden

Import from USB is hidden

Lockout Menu:

Menu is hidden

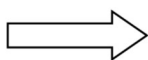**Note:** The security switch located at the rear of the balance must be in the locked position to set Approved Mode to ON. The security switch must be in the unlocked position to set Approved Mode to OFF. See Section 6.

Weighing application main screen with LFT turned ON.

## 5.4 Weighing Units

Enter this sub-menu to activate the desired units of measure. This menu can also be accessed by pressing the unit symbol in an application home screen.

**Note:** Due to national laws, the balance may not include some of the units of measure listed.

### 5.4.1 Units Sub-menu

|                            |                           |                         |                       |                     |                          |
|----------------------------|---------------------------|-------------------------|-----------------------|---------------------|--------------------------|
| <b>mg</b><br>Milligram     | <b>g</b><br>gram          | <b>kg</b><br>Kilogram   | <b>ct</b><br>carat    | <b>oz</b><br>ounce  | <b>ozt</b><br>ounce troy |
| <b>lb</b><br>Pound         | <b>dwt</b><br>pennyweight | <b>Grain</b><br>Grain   | <b>N</b><br>Newton    | <b>mom</b><br>momme | <b>msg</b><br>mesghal    |
| <b>HKt</b><br>Tael (HK)    | <b>SGt</b><br>Tael (SG)   | <b>TWt</b><br>Tael (TW) | <b>tical</b><br>tical | <b>tola</b><br>tola | <b>bht</b><br>baht       |
| <b>C1</b><br>Custom Unit 1 |                           |                         |                       |                     |                          |

**Note:** If Approved Mode is set to **ON**, some units will not be displayed.

Use the Custom Unit to display weight in an alternative unit of measure. The custom unit is defined using a conversion factor, where the conversion factor is the number of custom units per gram expressed in scientific notation (Factor x 10<sup>Exponent</sup>).

*For example:* To display weight in troy ounces (0.03215075 troy ounces per gram) enter a Factor of 0.3215075 and an Exponent of -1.

The Custom Unit's name can be customized up to 3 characters.

## 5.5 Data Maintenance

Enter this sub-menu to customize data transfer settings.

### 5.5.1 Data Maintenance sub-menu

The Export and Import functions

Setting up multiple balances is simple by exporting the profile from a master balance via a USB drive. The data maintenance tool allows you to save user and application settings to a USB, which can be easily transferred to other Adventurer balances. The data can then be used to configure additional Adventurer balances with the data imported from original balance.


### 5.5.2 Export to USB

Export weighing data to a USB flash drive. Two types of data can be exported:

- Application settings (APW, Ref. weight and etc.)
- Menu settings (balance setup function and etc.)

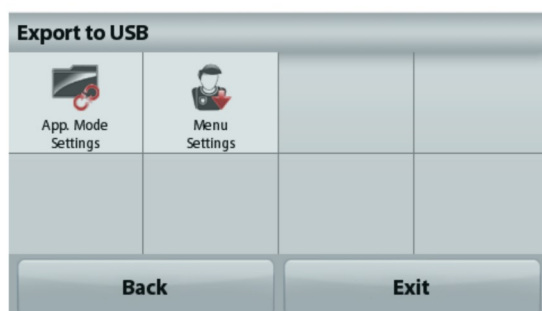**Note:** The function Save to USB needs to be set to ON to enable data transfer to USB. Please see section 5.6 for more information.

### 5.5.3 Import from USB

Import weighing data from a USB flash drive.

### 5.5.4 Balance Info

Enter to view information about the balance. Information displayed includes: Balance Type, Balance ID, Capacity, Readability and Software Version.

## 5.6 Communication

Enter this menu to define external communication methods and to set printing parameters. Data may be output to either a printer or PC.

Factory default settings are shown in bold. Enter to view information about the balance.

### 5.6.1 Communication Sub-menu

**RS-232 Standard**

**Print Settings**

**Save to USB**

### RS-232 Standard:

Enter this sub-menu to customize RS-232 Standard settings.

#### 5.6.2 Baud Rate

Set the baud rate (bits per second).

|             |             |
|-------------|-------------|
| 2400        | = 2400 bps  |
| 4800        | = 4800 bps  |
| <b>9600</b> | = 9600 bps  |
| 19200       | = 19200 bps |
| 38400       | = 38400 bps |

**5.6.3 Transmission**

Set the data bits, stop bit, and parity.

- 7 EVEN 1 = 7 data bits, even parity, stop bit 1
- 7 ODD 1 = 7 data bits, odd parity, stop bit 1
- 7 EVEN 2 = 7 data bits, even parity, stop bit 2
- 7 ODD 2 = 7 data bits, odd parity, stop bit 2
- 7 NONE 1 = 7 data bits, no parity, stop bit 1
- 8 NONE 1** = 8 data bits, no parity, stop bit 1
- 7 NONE 2 = 7 data bits, no parity, stop bit 2
- 8 NONE 2 = 8 data bits, no parity, stop bit 2

**5.6.4 Handshake**

Set the flow control method.

- NONE = no handshaking
- XON-XOFF** = XON/XOFF handshaking
- HARDWARE = hardware handshaking

**Print Settings**

Enter this sub-menu to customize data transfer settings.

**Print Settings sub-menu**

|  |   |   |   |
|--|---|---|---|
| <br>Print Output | <br>Auto Print              | <br>Print Content | <br>Feed |
| <br>Format      | <br>Print Calibration Data |   |   |

**5.6.5 Print Output**

**Stable Weight Only**

Set the printing criteria.

- OFF = values are printed immediately, regardless of stability.
- ON = values are printed only when the stability criteria are met.

**Numeric Value Only**

Set the printing criteria.

- OFF = All Result and G/N/T data values are printed. See section 7.2 for more information
- ON = Only numeric data values are printed

**Single Header Only**

Set the printing criteria.

- OFF = Headers will be printed for every print requirement
- ON = Headers will be printed once a day

**Print Options**

Set the printing criteria.

- Printer = Print data to a printer
- PC = Print data to a PC

### 5.6.6 Auto Print

Set the automatic printing functionality.

- OFF** = disabled
- ON STABILITY**<sup>1</sup> = printing occurs each time the stability criteria are met.
- INTERVAL**<sup>2</sup> = printing occurs at the defined time interval.
- CONTINUOUS** = printing occurs continuously.

<sup>1</sup>When ON STABILITY is selected, set the conditions for printing.

- LOAD** = Prints when the displayed load is stable.
- LOAD ZERO** = Prints when the displayed load and zero reading is stable.

<sup>2</sup>When INTERVAL is selected, set the time interval using the numeric keypad.

Settings of 1 to 3600 seconds are available. Default is 0.

### 5.6.7 Print Content

Define the content of the printed data.

#### Selection

Set the status.

**Deselect All** = all are set to OFF

**Select All** = all are set to ON

#### Header

**Date & Time**

**Balance ID**

**Balance Name**

**User Name**

**Project Name**

**Application Name**

**Sample Name**

**Result**

**Gross**

**Net**

**Tare**

**Information**

**Signature Line**

### 5.6.8 Feed

Set the paper feed.

**1 LINE** = move the paper up one line after printing.

**4 LINE** = move the paper up four lines after printing.

### 5.6.9 Format

Set the printing format.

**Single Line** = prints on a single line.

**Multiple Lines** = prints on multiple lines.

### 5.6.10 Print Calibration Data

Set the printing criteria.

**Off** = disabled

**On** = enabled.

### Save to USB

Set the status.

**OFF** = the data will not be saved to USB

**ON** = the data will be saved to USB

## 5.7 GLP and GMP Data

Enter this menu to set the Good Laboratory Practices (GLP) data and the Good Manufacturing Practices data.

### 5.7.1 GLP Data Sub-menu

Header

Balance Name

User Name

Project Name

#### 5.7.2 Header

Enables the printing of GLP headings. There are up to 5 headings available.

Alphanumeric settings up to 25 characters are available for each Header setting.

#### 5.7.3 Balance Name

Set the balance identification.

Alphanumeric settings up to 25 characters are available. The default setting is **Adventurer**.

#### 5.7.4 User Name

Set the user identification.

Alphanumeric settings up to 25 characters are available. The default setting is **blank**.

#### 5.7.5 Project Name

Enter this menu to set the Project identification.

Alphanumeric settings up to 25 characters are available. The default setting is **blank**.

## 5.8 Factory Reset

Use this sub-menu to reset the all menu settings to their Factory default settings.

Note: Calibration data is not affected.

- Reset All = resets all menus to their factory default settings.
- Exit = return to application main screen without resetting any menus.

## 5.9 Lockout

Use this sub-menu to lock/unlock certain menus. Once locked, a small lock will appear on the menu icon and the user will not be able to enter that menu.

Click on a menu and choose:

- OFF** = the menu is unlocked
- ON** = the menu is locked

## 6. LEGAL FOR TRADE (LFT)

When the balance is used in trade or a legally controlled application it must be set up, verified and sealed in accordance with local weights and measures regulations. It is the responsibility of the purchaser to ensure that all pertinent legal requirements are met.

### 6.1 Settings

Before verification and sealing, perform the following steps in order:

1. Verify that the menu settings meet the local weights and measures regulations.
2. Verify the units turned **On** meet the local weights and measures regulations.
3. Perform a calibration as explained in Section 5.
4. Set the position of the Security Switch to the locked position.
5. Set Approved Mode to ON in the Balance Setup menu.

**Note:** When Approved Mode is set to ON, external calibration can't be performed.

### 6.2 Verification

A weights and measures official or authorized service agent must perform the verification procedure.

### 6.3 Sealing

After the Balance has been verified, it must be sealed to prevent undetected access to the legally controlled settings. Before sealing the device, ensure that the security switch is in the Locked position and the Approved Mode setting in the Balance Setup menu has been set to ON.

If using a paper seal, place the seal over the security switch and Bottom Housing as shown.

If using a wire seal, pass the sealing wire through the holes in the security switch and Bottom Housing as shown.

Un-Locked

Locked with Paper Seal

Locked with Wire Seal

## 7. PRINTING

### 7.1 Connecting, Configuring and Testing the Printer/Computer Interface

Use the built-in RS-232 Port to connect either to a computer or a printer. If connecting to a computer, use HyperTerminal or similar software like SPDC described below.

(Find HyperTerminal under **Accessories/Communications** in Windows XP.)

Connect to the computer with a standard (straight-through) serial cable.


Choose **New Connection**, “connect using” COM1 (or available COM port).

Select **Baud=9600; Parity=8 None; Stop=1; Handshaking=None**. Click **OK**.

Choose Properties/Settings, then ASCII Setup. Check boxes as illustrated: (**Send line ends...**; **Echo typed characters...**; **Wrap lines...**)

Use RS232 Interface Commands (Section 9.6.1) to control the balance from a PC.


**Note:** When the HyperTerminal configuration is complete, it will automatically print the results of a **Cal Test** operation, and echo print commands sent to the scale.
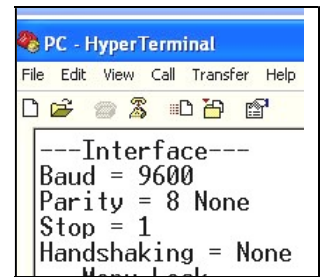
### SPDC Software

The Serial Port Data Collection or SPDC software is provided by Ohaus and can be used on operating systems that do not have the HyperTerminal software mentioned above.

Choose the export file type and export file path and then press Run as shown below.

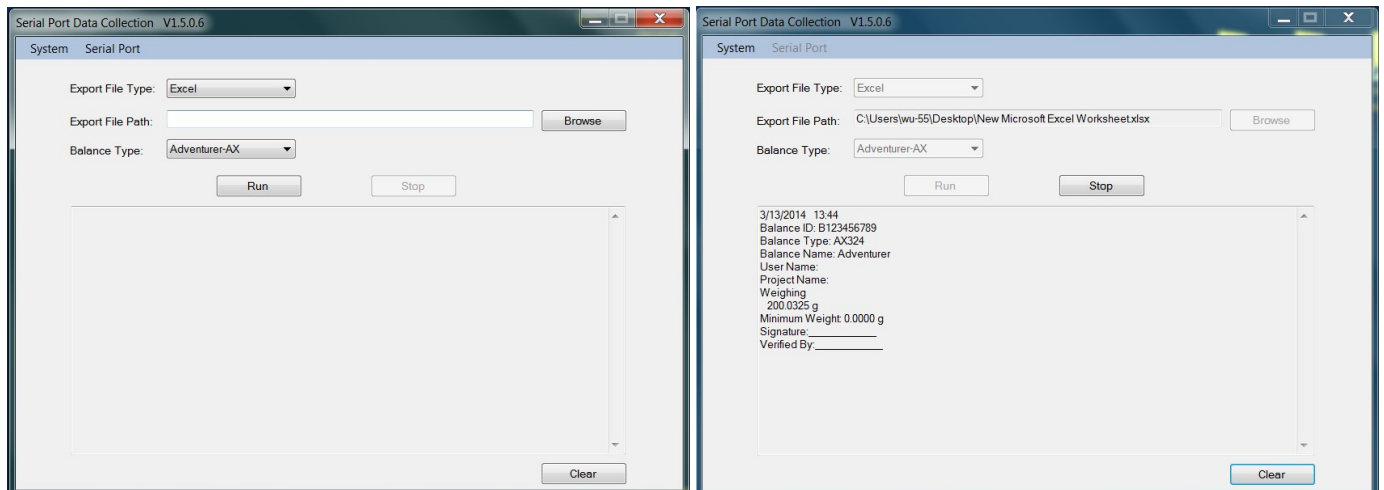**Note:** The SPDC software only supports English language.

### 7.2 Output Format

The Result Data, and G/N/T data, is output in the following format.

| Field:  | Label <sup>1</sup> | Space <sup>2</sup> | Weight <sup>3</sup> | Space <sup>2</sup> | Unit <sup>4</sup> | Space | Stability <sup>5</sup> | Space | G/N <sup>6</sup> | Space | Term. Characters <sup>7</sup> |
|---------|--------------------|--------------------|---------------------|--------------------|-------------------|-------|------------------------|-------|------------------|-------|-------------------------------|
| Length: |                    | 1                  | 11                  | 1                  | 5                 | 1     | ≤ 1                    | ≤ 1   | ≤ 3              | 0     | ≤ 8                           |

1. The length of the label field is not fixed.
2. Each field is followed by a single delimiting space (ASCII 32).
3. The Weight field is 11 right justified characters. If the value is negative, the “-” character is located at the immediate left of the most significant digit.
4. The Unit field contains the unit of measure abbreviation up to 5 characters, right justified.
5. The Stability field contains the “?” character if the weight reading is not stable. The Stability field and the following Space field are omitted if the weight reading is stable.
6. The G/N field contains the net or gross indication. For net weights, the field contains “NET”. For gross weights, the field contains nothing, “G”.
7. The Termination Characters field contains CRLF, Four CRLF or Form Feed (ASCII 12), depending on the LINE FEED menu setting.

7.3 Printout Examples

| Weighing  | Parts Counting  | Percent Weighing   |
|---|---|--|
| <p>Header 1<br/>Header 2<br/>Header 3<br/>Header 4<br/>Header 5<br/>1/15/2014 13:16<br/>Balance ID: B234567890<br/>Balance Type: AX224N<br/>Balance Name: Adventurer<br/>User Name: ohaus<br/>Project Name: ax<br/>Weighing<br/>Sample Name: apple<br/>1.3651 g NET<br/>Gross: 3.9199 g G<br/>Net: 1.3651 g NET<br/>Tare: 2.5548 g T<br/>Minimum Weight: 0.0000 g<br/>Signature: _____<br/>Verified By: _____</p> | <p>Header 1<br/>Header 2<br/>Header 3<br/>Header 4<br/>Header 5<br/>1/15/2014 13:19<br/>Balance ID: B234567890<br/>Balance Type: AX224N<br/>Balance Name: Adventurer<br/>User Name: ohaus<br/>Project Name: ax<br/>Parts Counting<br/>Sample Name: apple<br/>Quantity: 4 PCS NET<br/>Gross: 94.3343 g G<br/>Net: 91.7795 g NET<br/>Tare: 2.5548 g T<br/>APW: 23.09999 g<br/>Sample Size: 23 PCS<br/>Signature: _____<br/>Verified By: _____</p> | <p>Header 1<br/>Header 2<br/>Header 3<br/>Header 4<br/>Header 5<br/>1/15/2014 13:19<br/>Balance ID: B234567890<br/>Balance Type: AX224N<br/>Balance Name: Adventurer<br/>User Name: ohaus<br/>Project Name: ax<br/>Percent Weighing<br/>Sample Name: apple<br/>Percentage: 91.7795 % NET<br/>Gross: 94.3342 g G<br/>Net: 91.7795 g NET<br/>Tare: 2.5548 g T<br/>Reference Weight: 100.0000 g<br/>Signature: _____<br/>Verified By: _____</p> |

| Dynamic Weighing   | Density  | Density   |
|--|--|---|
| <p>Header 1<br/>Header 2<br/>Header 3<br/>Header 4<br/>Header 5<br/>1/15/2014 13:22<br/>Balance ID: B234567890<br/>Balance Type: AX224N<br/>Balance Name: Adventurer<br/>User Name: ohaus<br/>Project Name: ax<br/>Dynamic Weighing<br/>Sample Name: cat<br/>Final wt.: 90.4146 g<br/>Gross: 94.3362 g G<br/>Net: 90.4144 g NET<br/>Tare: 3.9218 g T<br/>Averaging Time: 5 s<br/>Signature: _____<br/>Verified By: _____</p> | <p>Type==Solid<br/>Auxiliary liquid==water<br/>Porous material==off</p> <p>Header 1<br/>Header 2<br/>Header 3<br/>Header 4<br/>Header 5<br/>1/15/2014 13:31<br/>Balance ID: B234567890<br/>Balance Type: AX224N<br/>Balance Name: Adventurer<br/>User Name: ohaus<br/>Project Name: ax<br/>Density Determination<br/>Density Determination: 34.1592 g/cm3<br/>Gross: 97.1644 g G<br/>Net: 93.2426 g NET<br/>Tare: 3.9218 g T<br/>Weight in air: 96.0491 g Header 1<br/>Header 2<br/>Header 3<br/>Header 4<br/>Header 5<br/>1/15/2014 13:37<br/>Balance ID: B234567890<br/>Balance Type: AX224N<br/>Balance Name: Adventurer<br/>User Name: ohaus<br/>Project Name: ax<br/>Density Determination<br/>Density Determination: 13.6849 g/cm3</p> | <p>Type==Solid<br/>Auxiliary liquid==water<br/>Porous material==on)</p> |

|  |  |  |
|--|--|--|
|  | <p>Gross: 95.7991 g G<br/>                 Net: 91.8773 g NET<br/>                 Tare: 3.9218 g T<br/>                 Oiled Weight: 98.8827 g<br/>                 Weight in liquid: 91.8773 g<br/>                 Auxiliary liquid: Water<br/>                 Liquid Density: 0.9982 g/cm<sup>3</sup><br/>                 Water Temp.: 20.0 °C<br/>                 Porous Material: On<br/>                 Oil Density: 0.8000 g/cm<sup>3</sup><br/>                 Dry Weight: 96.0490 g<br/>                 Signature: _____<br/>                 Verified By: _____<br/>                 Weight in liquid: 93.2426 g<br/>                 Auxiliary liquid: Water<br/>                 Liquid Density: 0.9982 g/cm<sup>3</sup><br/>                 Water Temp.: 20.0 °C<br/>                 Porous Material: Off<br/>                 Signature: _____<br/>                 Verified By: _____</p> |  |
|--|--|--|

|   |
|---|
| <p><b>Density</b><br/>                 Type==Solid<br/>                 Auxiliary liquid==other<br/>                 Porous material==on</p>  |
| <p>Header 1<br/>                 Header 2<br/>                 Header 3<br/>                 Header 4<br/>                 Header 5<br/>                 1/15/2014 13:50<br/>                 Balance ID: B234567890<br/>                 Balance Type: AX224N<br/>                 Balance Name: Adventurer<br/>                 User Name: ohaus<br/>                 Project Name: ax<br/>                 Density Determination<br/>                 Density Determination:<br/>                 4.7794 g/cm<sup>3</sup><br/>                 Gross: 93.2556 g G<br/>                 Net: 89.3338 g NET<br/>                 Tare: 3.9218 g T<br/>                 Oiled Weight: 110.5639 g<br/>                 Weight in liquid: 89.3338 g<br/>                 Auxiliary liquid: Other<br/>                 Liquid Density: 1.0000<br/>                 g/cm<sup>3</sup><br/>                 Porous Material: On<br/>                 Oil Density: 0.8000 g/cm<sup>3</sup><br/>                 Dry Weight: 101.7253 g<br/>                 Signature: _____<br/>                 Verified By: _____</p> |

|   |
|---|
| <p><b>Density</b><br/>                 Type==liquid<br/>                 Sinker volume==10ml<br/>                 Liquid temp==26°C</p>   |
| <p>Header 1<br/>                 Header 2<br/>                 Header 3<br/>                 Header 4<br/>                 Header 5<br/>                 1/15/2014 13:56<br/>                 Balance ID: B234567890<br/>                 Balance Type: AX224N<br/>                 Balance Name: Adventurer<br/>                 User Name: ohaus<br/>                 Project Name: ax<br/>                 Density Determination<br/>                 Density Determination: 0.7171<br/>                 g/cm<sup>3</sup><br/>                 Gross: 97.5185 g G<br/>                 Net: 93.5967 g NET<br/>                 Tare: 3.9218 g T<br/>                 Sinker weight in air: 100.7676 g<br/>                 Sinker weight in liquid: 93.5963 g<br/>                 Sinker Volume: 10.0 ml<br/>                 Liquid Temp.: 26.0 °C<br/>                 Signature: _____<br/>                 Verified By: _____</p> |

|   |
|---|
| <p><b>Check Weighing</b><br/>                 Type==liquid<br/>                 Sinker volume==10ml<br/>                 Liquid temp==26°C</p>  |
| <p>Header 1<br/>                 Header 2<br/>                 Header 3<br/>                 Header 4<br/>                 Header 5<br/>                 1/15/2014 13:57<br/>                 Balance ID: B234567890<br/>                 Balance Type: AX224N<br/>                 Balance Name: Adventurer<br/>                 User Name: ohaus<br/>                 Project Name: ax<br/>                 Check Weighing<br/>                 Sample Name: apple<br/>                 93.5966 g NET<br/>                 Result: Accept<br/>                 Gross: 97.5184 g G<br/>                 Net: 93.5966 g NET<br/>                 Tare: 3.9218 g T<br/>                 Over Limit: 4199.9900 g<br/>                 Under Limit: 0.1000 g<br/>                 Signature: _____<br/>                 Verified By: _____<br/>                 Signature: _____<br/>                 Verified By: _____</p> |

| Display Hold             |
|--------------------------|
| Header 1                 |
| Header 2                 |
| Header 3                 |
| Header 4                 |
| Header 5                 |
| 1/15/2014 13:59          |
| Balance ID: B234567890   |
| Balance Type: AX224N     |
| Balance Name: Adventurer |
| User Name: ohaus         |
| Project Name: ax         |
| Display Hold             |
| Sample Name: apple       |
| Hold Weight: 93.5968 g   |
| Gross: 97.5185 g G       |
| Net: 93.5967 g NET       |
| Tare: 3.9218 g T         |
| Mode: Display Hold       |
| Signature: _____         |
| Verified By: _____       |

| Totalization                 |
|------------------------------|
| Header 1                     |
| Header 2                     |
| Header 3                     |
| Header 4                     |
| Header 5                     |
| 1/15/2014 14:11              |
| Balance ID: B234567890       |
| Balance Type: AX224N         |
| Balance Name: Adventurer     |
| User Name: ohaus             |
| Project Name: ax             |
| Totalization                 |
| Total: 734.6187 g Net        |
| Gross: 93.2557 g G           |
| Net: 89.3339 g NET           |
| Tare: 3.9218 g T             |
| Samples: 8                   |
| Average: 91.8273 g           |
| Standard Deviation: 1.9790 g |
| Minimum: 89.3339 g           |
| Maximum: 93.5965 g           |
| Range: 4.2626 g              |
| -----Sample Data (g)----     |
| 1 93.5964                    |
| 2 93.5964                    |
| 3 93.5964                    |
| 4 93.5965                    |
| 5 92.2312                    |
| 6 89.3340                    |
| 7 89.3339                    |
| 8 89.3339                    |
| Signature: _____             |
| Verified By: _____           |

| Formulation              |
|--------------------------|
| Header 1                 |
| Header 2                 |
| Header 3                 |
| Header 4                 |
| Header 5                 |
| 1/15/2014 14:22          |
| Balance ID: B234567890   |
| Balance Type: AX224N     |
| Balance Name: Adventurer |
| User Name: ohaus         |
| Project Name: ax         |
| Formulation              |
| Comp. Total: 11.4528 g   |
| Filler: 2.8063 g         |
| Total: 14.2590 g         |
| Gross: 18.1806 g ? G     |
| Net: 2.8063 g ? NET      |
| Tare: 15.3742 g T        |
| -----Sample Data (g)---- |
| Item 1: 1.7529 g         |
| Item 2: 2.5569 g         |
| Item 3: 1.3553 g         |
| Item 4: 1.3070 g         |
| Item 5: 1.6743 g         |
| Item 6: 2.8062 g         |
| Signature: _____         |
| Verified By: _____       |

## 8. MAINTENANCE

### 8.1 Calibration

Periodically verify calibration by placing an accurate weight on the balance and viewing the result. If calibration is required, refer to section 5.2 for instructions.

### 8.2 Cleaning

**WARNING:** Disconnect the Adventurer Balance from the power supply before cleaning. Make sure that no liquid enters the interior of the balance. Clean the Balance at regular intervals.

Housing surfaces may be cleaned with a lint-free cloth slightly dampened with water or a mild cleaning agent.

Glass surfaces may be cleaned with a commercial glass cleaner. Please follow the steps below on how to remove and install the sliding doors.

**Attention:** Do not use solvents, harsh chemicals, ammonia or abrasive cleaning agents.

#### Removing and reinstalling the glass doors for cleaning:


##### Step 1.

On the back of the balance, press the pin and slide out the door.

##### Step 2.

After cleaning, slide the doors into the slot while pressing the pin mentioned in step 1.
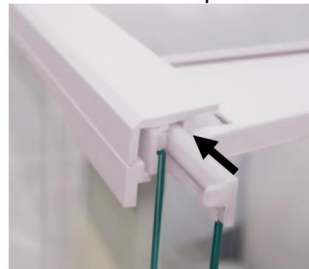
##### Step 3.

Slide the doors into the slot until the back stopper aligns with the other door.

### 8.3 Troubleshooting

TABLE 8-1. TROUBLESHOOTING

| Symptom / Display           | Possible Cause   | Remedy  |
|-----------------------------|--|---|
| Balance will not turn on    | No power to Balance  | Verify connection and voltage   |
| Poor accuracy               | Improper calibration<br>Unstable environment   | Perform calibration<br>Move balance to suitable location  |
| Cannot calibrate            | Calibration Menu locked<br>Approved Mode set to on<br>Unstable environment<br>Incorrect calibration masses | Turn Calibration menu lock off<br>Turn Approved Mode off<br>Move balance to suitable location<br>Use correct calibration masses |
| Cannot change menu settings | Sub-menu locked<br>Approved Mode set to on   | Unlock sub-menu<br>Turn Approved Mode off   |
| Low Reference weight        | Reference weight too small<br>The weight on the pan is too small to define a valid reference weight.       | Increase sample size  |
| Invalid Piece Weight        | Average piece weight is too small  | Increase average piece weight   |
| Operation Timeout           | Weight reading is not stable   | Move balance to suitable location   |
| -----                       | Busy (tare, zero, printing, waiting for a stable weight)   | Wait until completion   |

### 8.4 Service Information

If the troubleshooting section does not resolve your problem, contact an Authorized Ohaus Service Agent. Please visit our website [www.ohaus.com](http://www.ohaus.com) to locate the Ohaus office nearest you. An Ohaus Product Service Specialist will be available to assist you.

## 9. TECHNICAL DATA

### 9.1 Specifications

#### Ambient conditions

- Indoor use only
- Altitude: Up to 2000 m
- Specified Temperature range: 10°C to 30°C
- Humidity: maximum relative humidity 80 % for temperatures up to 30°C decreasing linearly to 50% relative humidity at 40°C
- Mains supply voltage fluctuations: up to  $\pm 10\%$  of the nominal voltage
- Installation category II
- Pollution degree: 2

#### Materials

- Bottom Housing; die-cast Aluminum, Painted
- Top Housing: Plastic (ABS)
- Weighing Platforms: 18/10 stainless steel
- Draft Shield: Glass, plastic (ABS)
- Feet: Plastic (ABS)

TABLE 9-1. SPECIFICATIONS

| Model:                                 | AX124  | AX124/E                  | AX224                      | AX224/E                       | AX324              |
|--|--|--------------------------|----------------------------|-------------------------------|--------------------|
| Capacity                               | 120g   | 120g                     | 220g                       | 220g                          | 320g               |
| Readability d                          | 0.0001g  |                          |                            |                               |                    |
| Repeatability (std. dev.) (g)          | 0.0001g  |                          |                            |                               |                    |
| Linearity (g)                          | ±0.0002g   |                          |                            |                               |                    |
| Span Calibration Points (g)            | 25g, 50g<br>75g, 100g  | 25g, 50g<br>75g,<br>100g | 50g, 100g<br>150g,<br>200g | 50g,<br>100g<br>150g,<br>200g | 100g<br>200g, 300g |
| Calibration                            | Auto-Cal   | External                 | Auto-Cal                   | External                      | Auto-Cal           |
| Weighing units                         | Baht, carat, grain, gram, milligram, mesgal, momme, Newton, ounce, pennyweight, pound, Tael (Hong Kong), Tael, (Singapore), Tael (Taiwan), tical, tola, troy ounce, custom (1) |                          |                            |                               |                    |
| Applications                           | Weighing, Parts Counting, Percent Weighing, Check Weighing, Animal Weighing, Formulation, Density Determination, Totalization, Display Hold                                    |                          |                            |                               |                    |
| Stabilization time (typical)           | ≤ 3 seconds  |                          |                            |                               |                    |
| Sensitivity Temperature Drift (PPM/K)  | 1.5  |                          |                            |                               |                    |
| Typical Min-Weight USP (u=0.10%,k=2)   | 0.20 g   |                          |                            |                               |                    |
| Optimized Min-Weight USP (u=0.10%,k=2) | 0.082 g  |                          |                            |                               |                    |
| Display                                | Full-Color WQVGA Graphic LCD   |                          |                            |                               |                    |
| Display size                           | 4.3 in / 10.9 cm (diagonal)  |                          |                            |                               |                    |
| Backlight                              | White LED  |                          |                            |                               |                    |
| Controls                               | 4-wire resistive touch screen + 6 membrane keys  |                          |                            |                               |                    |
| Communication                          | RS-232, USBx2  |                          |                            |                               |                    |
| Balance power input                    | 12 VDC, 0.5A   |                          |                            |                               |                    |
| Power supply                           | AC Adapter Input: 100-240 VAC 0.3A 50-60 Hz<br>AC Adapter Output: 12 VDC 0.84A   |                          |                            |                               |                    |
| Platform size (diameter)               | 90 mm / 3.5 in   |                          |                            |                               |                    |
| Assembled dimensions (W x D x H)       | 230 x 354 x 340 mm<br>9.1 x 13.9 x 13.4 inch   |                          |                            |                               |                    |
| Shipping dimensions (W x D x H)        | 387 x 507 x 531 mm<br>15.4 x 20.0 x 20.9 inch  |                          |                            |                               |                    |
| Net weight                             | 5.1 Kg / 11.3 lb   |                          |                            |                               |                    |
| Shipping weight                        | 7.8 Kg / 17.2 lb   |                          |                            |                               |                    |

TABLE 9-2. SPECIFICATIONS (continued)

| Model:                                 | AX223  | AX223/E                    | AX423                          | AX423/E                  | AX523                    | AX523/E                  |
|--|--|----------------------------|--------------------------------|--------------------------|--------------------------|--------------------------|
| Capacity                               | 220g   | 220g                       | 420g                           | 420g                     | 520g                     | 520g                     |
| Readability d                          | 0.001g   |                            |                                |                          |                          |                          |
| Repeatability (std. dev.) (g)          | 0.001g   |                            |                                |                          |                          |                          |
| Linearity (g)                          | ±0.002g  |                            |                                |                          |                          |                          |
| Span Calibration Points (g)            | 50g,<br>100g<br>150g,<br>200g  | 50g, 100g<br>150g,<br>200g | 100g,<br>200g<br>300g,<br>400g | 100g, 200g<br>300g, 400g | 200g, 300g<br>400g, 500g | 200g, 300g<br>400g, 500g |
| Calibration                            | Auto-Cal   | External                   | Auto-Cal                       | External                 | Auto-Cal                 | External                 |
| Weighing units                         | Baht, carat, grain, gram, milligram, mesgal, momme, Newton, ounce, pennyweight, pound, Tael (Hong Kong), Tael, (Singapore), Tael (Taiwan), tical, tola, troy ounce, custom (1) |                            |                                |                          |                          |                          |
| Applications                           | Weighing, Parts Counting, Percent Weighing, Check Weighing, Animal Weighing, Formulation, Density Determination, Totalization, Display Hold                                    |                            |                                |                          |                          |                          |
| Stabilization time (typical)           | ≤ 2 seconds  |                            |                                |                          |                          |                          |
| Sensitivity Temperature Drift (PPM/K)  | 3  |                            |                                |                          |                          |                          |
| Typical Min-Weight USP (u=0.10%,k=2)   | 2.0 g  |                            |                                |                          |                          |                          |
| Optimized Min-Weight USP (u=0.10%,k=2) | 0.82 g   |                            |                                |                          |                          |                          |
| Display                                | Full-Color WQVGA Graphic LCD   |                            |                                |                          |                          |                          |
| Display size                           | 4.3 in / 10.9 cm (diagonal)  |                            |                                |                          |                          |                          |
| Backlight                              | White LED  |                            |                                |                          |                          |                          |
| Controls                               | 4-wire resistive touch screen + 6 membrane keys  |                            |                                |                          |                          |                          |
| Communication                          | RS-232, USBx2  |                            |                                |                          |                          |                          |
| Balance power input                    | 12 VDC, 0.5A   |                            |                                |                          |                          |                          |
| Power supply                           | AC Adapter Input: 100-240 VAC 0.3A 50-60 Hz<br>AC Adapter Output: 12 VDC 0.84A   |                            |                                |                          |                          |                          |
| Platform size (diameter)               | 130 mm / 5.1 in  |                            |                                |                          |                          |                          |
| Assembled dimensions (W x D x H)       | 230 x 354 x 340 mm<br>9.1 x 13.9 x 13.4 inch   |                            |                                |                          |                          |                          |
| Shipping dimensions (W x D x H)        | 387 507 x 531 mm<br>15.4 x 20.0 x 20.9 inch  |                            |                                |                          |                          |                          |
| Net weight                             | 5.8Kg/12.8lb   | 5.6Kg/12.4lb               | 5.8Kg/12.8lb                   | 5.6Kg/12.4lb             | 5.8Kg/12.8lb             | 5.6Kg/12.4lb             |
| Shipping weight                        | 8.5Kg/18.8lb   | 8.3Kg/18.3lb               | 8.5Kg/18.8lb                   | 8.3Kg/18.3lb             | 8.5Kg/18.8lb             | 8.3Kg/18.3lb             |

TABLE 9-3. SPECIFICATIONS (continued)

| Model:                                | AX422   | AX422/E                 | AX822                    | AX822/E                  |
|---------------------------------------|---|-------------------------|--------------------------|--------------------------|
| Capacity                              | 420g  | 420g                    | 820g                     | 820g                     |
| Readability d                         | 0.01g   |                         |                          |                          |
| Repeatability (std. dev.) (g)         | 0.01g   |                         |                          |                          |
| Linearity (g)                         | ±0.02g  |                         |                          |                          |
| Span Calibration Points (g)           | 100g, 200g<br>300g,400g   | 100g, 200g<br>300g,400g | 200g, 400g<br>600g, 800g | 200g, 400g<br>600g, 800g |
| Calibration                           | Auto-Cal  | External                | Auto-Cal                 | External                 |
| Weighing units                        | Baht, carat, grain, gram, kilogram, mesgal, momme, Newton, ounce, pennyweight, pound, Tael (Hong Kong), Tael, (Singapore), Tael (Taiwan), tical, tola, troy ounce, custom (1) |                         |                          |                          |
| Applications                          | Weighing, Parts Counting, Percent Weighing, Check Weighing, Animal Weighing, Formulation, Density Determination, Totalization, Display Hold                                   |                         |                          |                          |
| Stabilization time (typical)          | ≤ 1.5 seconds   |                         |                          |                          |
| Sensitivity Temperature Drift (PPM/K) | 3   |                         |                          |                          |
| Typical Min-Weight USP(u=0.10%,k=2)   | 20.0 g  |                         |                          |                          |
| Optimal Min-Weight USP(u=0.10%,k=2)   | 8.2 g   |                         |                          |                          |
| Display                               | Full-Color WQVGA Graphic LCD  |                         |                          |                          |
| Display size                          | 4.3 in / 10.9 cm (diagonal)   |                         |                          |                          |
| Backlight                             | White LED   |                         |                          |                          |
| Controls                              | 4-wire resistive touch screen + 6 membrane keys   |                         |                          |                          |
| Communication                         | RS-232, USBx2   |                         |                          |                          |
| Balance power input                   | 12 VDC, 0.5A  |                         |                          |                          |
| Power supply                          | AC Adapter Input: 100-240 VAC 0.3A 50-60 Hz<br>AC Adapter Output: 12 VDC 0.84A  |                         |                          |                          |
| Platform size                         | 175 x 195 mm / 6.9 x 7.7 in   |                         |                          |                          |
| Assembled dimensions (W x D x H)      | 230 x 354 x100 mm<br>9.1 x 13.9 x 4.0 inch  |                         |                          |                          |
| Shipping dimensions (W x D x H)       | 392 x 557 x 301 mm<br>15.5 x 22.0 x 11.9 inch   |                         |                          |                          |
| Net weight                            | 4.6Kg/10.2lb  | 3.9Kg/8.6lb             | 4.6Kg/10.2lb             | 3.9Kg/8.6lb              |
| Shipping weight                       | 6.5Kg/14.4b   | 5.8Kg/12.8b             | 6.5Kg/14.4b              | 5.8Kg/12.8b              |

TABLE 9-4. SPECIFICATIONS (continued)

| Model:                                | AX1502  | AX1502/E                | AX2202                             | AX2202/E                           | AX4202                              | AX4202/E                            | AX5202                              |
|---------------------------------------|---|-------------------------|------------------------------------|------------------------------------|-------------------------------------|-------------------------------------|-------------------------------------|
| Capacity                              | 1520g   | 1520g                   | 2200g                              | 2200g                              | 4200g                               | 4200g                               | 5200g                               |
| Readability d                         | 0.01g   |                         |                                    |                                    |                                     |                                     |                                     |
| Repeatability (std. dev.) (g)         | 0.01g   |                         |                                    |                                    |                                     |                                     |                                     |
| Linearity (g)                         | ±0.02g  |                         |                                    |                                    |                                     |                                     |                                     |
| Span Calibration Points (g)           | 500g<br>1000g,<br>1500g   | 500g<br>1000g,<br>1500g | 500g,<br>1000g,<br>1500g,<br>2000g | 500g,<br>1000g,<br>1500g,<br>2000g | 1000g,<br>2000g,<br>3000g,<br>4000g | 1000g,<br>2000g,<br>3000g,<br>4000g | 2000g,<br>3000g,<br>4000g,<br>5000g |
| Calibration                           | Auto-Cal  | External                | Auto-Cal                           | External                           | Auto-Cal                            | External                            | Auto-Cal                            |
| Weighing units                        | Baht, carat, grain, gram, kilogram, mesgal, momme, Newton, ounce, pennyweight, pound, Tael (Hong Kong), Tael, (Singapore), Tael (Taiwan), tical, tola, troy ounce, custom (1) |                         |                                    |                                    |                                     |                                     |                                     |
| Applications                          | Weighing, Parts Counting, Percent Weighing, Check Weighing, Animal Weighing, Formulation, Density Determination, Totalization, Display Hold                                   |                         |                                    |                                    |                                     |                                     |                                     |
| Stabilization time (typical)          | ≤ 1.5 seconds   |                         |                                    |                                    |                                     |                                     |                                     |
| Sensitivity Temperature Drift (PPM/K) | 3   |                         |                                    |                                    |                                     |                                     | 1.9                                 |
| Typical Min-Weight USP(u=0.10%,k=2)   | 20.0 g  |                         |                                    |                                    |                                     |                                     |                                     |
| Optimal Min-Weight USP(u=0.10%,k=2)   | 8.2 g   |                         |                                    |                                    |                                     |                                     |                                     |
| Display                               | Full-Color WQVGA Graphic LCD  |                         |                                    |                                    |                                     |                                     |                                     |
| Display size                          | 4.3 in / 10.9 cm (diagonal)   |                         |                                    |                                    |                                     |                                     |                                     |
| Backlight                             | White LED   |                         |                                    |                                    |                                     |                                     |                                     |
| Controls                              | 4-wire resistive touch screen + 6 membrane keys   |                         |                                    |                                    |                                     |                                     |                                     |
| Communication                         | RS-232, USBx2   |                         |                                    |                                    |                                     |                                     |                                     |
| Balance power input                   | 12 VDC, 0.5A  |                         |                                    |                                    |                                     |                                     |                                     |
| Power supply                          | AC Adapter Input: 100-240 VAC 0.3A 50-60 Hz<br>AC Adapter Output: 12 VDC 0.84A  |                         |                                    |                                    |                                     |                                     |                                     |
| Platform size                         | 175 x 195 mm / 6.9 x 7.7 in   |                         |                                    |                                    |                                     |                                     |                                     |
| Assembled dimensions (W x D x H)      | 230 x 354 x100 mm<br>9.1 x 13.9 x 4.0 inch  |                         |                                    |                                    |                                     |                                     |                                     |
| Shipping dimensions (W x D x H)       | 392 x 557 x 301 mm<br>15.5 x 22.0 x 11.9 inch   |                         |                                    |                                    |                                     |                                     |                                     |
| Net weight                            | 4.6Kg/<br>10.2lb  | 3.9Kg/<br>8.6lb         | 4.6Kg/<br>10.2lb                   | 3.9K/<br>8.6lb                     | 4.6Kg/<br>10.2lb                    | 3.9Kg/<br>8.6lb                     | 3.8Kg/<br>8.4lb                     |
| Shipping weight                       | 6.5Kg/<br>14.4b   | 5.8Kg/<br>12.8b         | 6.5Kg/<br>14.4b                    | 5.8Kg/<br>12.8b                    | 6.5Kg/<br>14.4b                     | 5.8Kg/<br>12.8b                     | 5.7Kg/<br>12.6lb                    |

TABLE 9-5. SPECIFICATIONS (continued)

| Model:                                 | AX2201  | AX2201/E                    | AX4201                             | AX4201/E                           | AX8201                             | AX8201/E                           |
|--|---|-----------------------------|------------------------------------|------------------------------------|------------------------------------|------------------------------------|
| Capacity                               | 2200g   | 2200g                       | 4200g                              | 4200g                              | 8200g                              | 8200g                              |
| Readability d                          | 0.1g  |                             |                                    |                                    |                                    |                                    |
| Repeatability (std. dev.) (g)          | 0.1g  |                             |                                    |                                    |                                    |                                    |
| Linearity (g)                          | ±0.2g   |                             |                                    |                                    |                                    |                                    |
| Span Calibration Points (g)            | 500g, 1000g<br>1500g, 2000g   | 500g, 1000g<br>1500g, 2000g | 1000g,<br>2000g<br>3000g,<br>4000g | 1000g,<br>2000g<br>3000g,<br>4000g | 2000g,<br>4000g<br>6000g,<br>8000g | 2000g,<br>4000g<br>6000g,<br>8000g |
| Calibration                            | Auto-Cal  | External                    | Auto-Cal                           | External                           | Auto-Cal                           | External                           |
| Weighing units                         | Baht, carat, grain, gram, kilogram, mesgal, momme, Newton, ounce, pennyweight, pound, Tael (Hong Kong), Tael, (Singapore), Tael (Taiwan), tical, tola, troy ounce, custom (1) |                             |                                    |                                    |                                    |                                    |
| Applications                           | Weighing, Parts Counting, Percent Weighing, Check Weighing, Animal Weighing, Formulation, Density Determination, Totalization, Display Hold                                   |                             |                                    |                                    |                                    |                                    |
| Stabilization time (typical)           | ≤ 1.5 seconds   |                             |                                    |                                    |                                    |                                    |
| Sensitivity Temperature Drift (PPM/K)  | 9   |                             |                                    |                                    |                                    |                                    |
| Typical Min-Weight USP (u=0.10%,k=2)   | 200.0 g   |                             |                                    |                                    |                                    |                                    |
| Optimized Min-Weight USP (u=0.10%,k=2) | 82 g  |                             |                                    |                                    |                                    |                                    |
| Display                                | Full-Color WQVGA Graphic LCD  |                             |                                    |                                    |                                    |                                    |
| Display size                           | 4.3 in / 10.9 cm (diagonal)   |                             |                                    |                                    |                                    |                                    |
| Backlight                              | White LED   |                             |                                    |                                    |                                    |                                    |
| Controls                               | 4-wire resistive touch screen + 6 membrane keys   |                             |                                    |                                    |                                    |                                    |
| Communication                          | RS-232, USBx2   |                             |                                    |                                    |                                    |                                    |
| Balance power input                    | 12 VDC, 0.5A  |                             |                                    |                                    |                                    |                                    |
| Power supply                           | AC Adapter Input: 100-240 VAC 0.3A 50-60 Hz<br>AC Adapter Output: 12 VDC 0.84A  |                             |                                    |                                    |                                    |                                    |
| Platform size                          | 175 x 195 mm / 6.9 x 7.7 in   |                             |                                    |                                    |                                    |                                    |
| Assembled dimensions (W x D x H)       | 230 x 354 x 100 mm<br>9.1 x 13.9 x 4.0 inch   |                             |                                    |                                    |                                    |                                    |
| Shipping dimensions (W x D x H)        | 392 x 557 x 301 mm<br>15.5 x 22.0 x 11.9inch  |                             |                                    |                                    |                                    |                                    |
| Net weight                             | 4.6Kg/10.2lb  | 3.9Kg/8.6lb                 | 4.6Kg/10.2lb                       | 3.9Kg/8.6lb                        | 3.8Kg/8.4lb                        | 3.4Kg/7.5lb                        |
| Shipping weight                        | 6.5Kg/14.4b   | 5.8Kg/12.8b                 | 6.5Kg/14.4b                        | 5.8Kg/12.8b                        | 5.7Kg/12.6lb                       | 5.3Kg/11.6lb                       |

TABLE 9-6. SPECIFICATIONS (continued)

| MODEL                                 | AX124M  | AX224M                        | AX324M                | AX223M  | AX423M                         | AX523M                         | AX1502M                                       | AX2202M                           | AX4202M                            | AX5202M                            | AX8201M                            |
|---------------------------------------|---|-------------------------------|-----------------------|---|--------------------------------|--------------------------------|---|-----------------------------------|------------------------------------|------------------------------------|------------------------------------|
| Max                                   | 120g  | 220g                          | 320g                  | 220g  | 420g                           | 520g                           | 1520g   | 2200g                             | 4200g                              | 5200g                              | 8200g                              |
| Min                                   | 0.01g   | 0.01g                         | 0.01g                 | 0.02g   | 0.02g                          | 0.02g                          | 0.5g  | 0.5g                              | 0.5g                               | 0.5g                               | 5g                                 |
| d=                                    | 0.0001g   |                               |                       | 0.001g  |                                |                                | 0.01g   |                                   |                                    |                                    | 0.1g                               |
| e=                                    | 0.001g  |                               |                       | 0.01g   |                                |                                | 0.1g  |                                   |                                    |                                    | 1g                                 |
| Approval Class                        | I   |                               |                       | II  |                                |                                |   |                                   |                                    |                                    |                                    |
| Repeatability (std. dev.) (g)         | 0.0001g   |                               |                       | 0.001g  |                                |                                | 0.01g   |                                   |                                    |                                    | 0.1g                               |
| Linearity (g)                         | ±0.0002g  |                               |                       | ±0.002g                                       |                                |                                | ±0.02g  |                                   |                                    |                                    | ±0.2g                              |
| Span Calibration Points (g)           | 25g, 50g<br>75g,<br>100g  | 50g,<br>100g<br>150g,<br>200g | 100g<br>200g,<br>300g | 50g,<br>100g<br>150g,<br>200g                 | 100g,<br>200g<br>300g,<br>400g | 200g,<br>300g<br>400g,<br>500g | 500g<br>1000g,<br>1500g                       | 500g,<br>1000g<br>1500g,<br>2000g | 1000g,<br>2000g<br>3000g,<br>4000g | 2000g,<br>3000g<br>4000g,<br>5000g | 2000g,<br>4000g<br>6000g,<br>8000g |
| Calibration                           | Auto-Cal  | Auto-Cal                      | Auto-Cal              | Auto-Cal                                      | Auto-Cal                       | Auto-Cal                       | Auto-Cal                                      | Auto-Cal                          | Auto-Cal                           | Auto-Cal                           | Auto-Cal                           |
| Weighing units                        | milligram, gram, carat  |                               |                       |   |                                |                                | kilogram, gram, carat                         |                                   |                                    |                                    |                                    |
| Applications                          | Weighing, Parts Counting, Percent Weighing, Check Weighing, Animal Weighing, Formulation, Density Determination, Totalization, Display Hold |                               |                       |   |                                |                                |   |                                   |                                    |                                    |                                    |
| Stabilization time (typical)          | ≤ 3 seconds   |                               |                       | ≤2 seconds                                    |                                |                                | ≤1.5 seconds                                  |                                   |                                    |                                    |                                    |
| Sensitivity Temperature Drift (PPM/K) | 1.5   |                               |                       | 3   |                                |                                | 3   |                                   |                                    | 1.9                                |                                    |
| Typical Min-Weight USP (u=0.10%,k=2)  | 0.20 g  | 0.20 g                        | 0.20 g                | 2.0 g   | 2.0 g                          | 2.0 g                          | 20.0 g  | 20.0 g                            | 20.0 g                             | 20.0 g                             | 200.0 g                            |
| Optimal Min-Weight USP (u=0.10%,k=2)  | 0.082 g   | 0.082 g                       | 0.082                 | 0.82 g  | 0.82 g                         | 0.82 g                         | 8.2 g   | 8.2 g                             | 8.2 g                              | 8.2 g                              | 82 g                               |
| Display                               | Full-Color WQVGA Graphic LCD  |                               |                       |   |                                |                                |   |                                   |                                    |                                    |                                    |
| Display size                          | 4.3 in / 10.9 cm (diagonal)   |                               |                       |   |                                |                                |   |                                   |                                    |                                    |                                    |
| Backlight                             | White LED   |                               |                       |   |                                |                                |   |                                   |                                    |                                    |                                    |
| Controls                              | 4-wire resistive touch screen + 6 membrane keys   |                               |                       |   |                                |                                |   |                                   |                                    |                                    |                                    |
| Communication                         | RS-232, USBx2   |                               |                       |   |                                |                                |   |                                   |                                    |                                    |                                    |
| Balance power input                   | 12 VDC, 0.5A  |                               |                       |   |                                |                                |   |                                   |                                    |                                    |                                    |
| Power supply                          | AC Adapter Input: 100-240 VAC 0.3A 50-60 Hz<br>AC Adapter Output: 12 VDC 0.84A  |                               |                       |   |                                |                                |   |                                   |                                    |                                    |                                    |
| Platform size (diameter)              | 90 mm / 3.5 in  |                               |                       | 130 mm / 5.1 in                               |                                |                                | 175x195 mm / 6.9x7.7 in                       |                                   |                                    |                                    |                                    |
| Assembled dimensions (W x D x H)      | 230 x 354 x 340 mm<br>9.1 x 13.9 x 13.4 inch  |                               |                       | 230 x 354 x 340 mm<br>9.1 x 13.9 x 13.4 inch  |                                |                                | 230 x 354 x 100 mm<br>9.1 x 13.9 x 4.0 inch   |                                   |                                    |                                    |                                    |
| Shipping dimensions (W x D x H)       | 387 x 507 x 531 mm<br>15.4 x 20.0 x 20.9 inch   |                               |                       | 387 x 507 x 531 mm<br>15.4 x 20.0 x 20.9 inch |                                |                                | 392 x 557 x 301 mm<br>15.5 x 22.0 x 11.9 inch |                                   |                                    |                                    |                                    |
| Net weight                            | 5.1 Kg / 11.3 lb  |                               |                       | 5.8 Kg / 12.8 lb                              |                                |                                | 4.6 Kg / 10.2 lb                              |                                   |                                    | 3.8 Kg / 8.4 lb                    |                                    |
| Shipping weight                       | 7.8 Kg / 17.2 lb  |                               |                       | 8.5 Kg / 18.8 lb                              |                                |                                | 6.5 Kg / 14.4 lb                              |                                   |                                    | 5.7 Kg / 12.6 lb                   |                                    |

**Note:** M = OIML Approved

TABLE 9-7. SPECIFICATIONS (continued)

| Model                                 | AX224N   | AX223N/E                   | AX423N                         | AX 423N/E                | AX523N/E                 | AX1502N/E  | AX2202N/E                   | AX4202N/E                    | AX8201N/E                    |
|---------------------------------------|--|----------------------------|--------------------------------|--------------------------|--------------------------|--|-----------------------------|------------------------------|------------------------------|
| Max                                   | 220g   | 220g                       | 420g                           | 420g                     | 520g                     | 1520g  | 2200g                       | 4200g                        | 8200g                        |
| Min                                   | 0.01g  | 0.02g                      | 0.02g                          | 0.02g                    | 0.02g                    | 0.5g   | 0.5g                        | 0.5g                         | 5g                           |
| d=                                    | 0.0001g  | 0.001g                     |                                |                          |                          | 0.01g  |                             |                              | 0.1g                         |
| e=                                    | 0.001g   | 0.01g                      |                                |                          |                          | 0.1g   |                             |                              | 1g                           |
| Approval Class                        | I  | II                         |                                |                          |                          |  |                             |                              |                              |
| Repeatability (std. dev.) (g)         | 0.0001g  | 0.001g                     |                                |                          |                          | 0.01g  |                             |                              | 0.1g                         |
| Linearity (g)                         | ±0.0002g   | ±0.002g                    |                                |                          |                          | ±0.02g   |                             |                              | ±0.2g                        |
| Span Calibration Points (g)           | 50g, 100g<br>150g,<br>200g   | 50g, 100g<br>150g,<br>200g | 100g,<br>200g<br>300g,<br>400g | 100g, 200g<br>300g, 400g | 200g, 300g<br>400g, 500g | 500g<br>1000g, 1500g   | 500g, 1000g<br>1500g, 2000g | 1000g, 2000g<br>3000g, 4000g | 2000g, 4000g<br>6000g, 8000g |
| Calibration                           | Auto-Cal   | External                   | Auto-Cal                       | External                 | External                 | External   | External                    | External                     | External                     |
| Weighing units                        | gram, milligram, carat, pennyweight,<br>grain, ounce, troy ounce   |                            |                                |                          |                          | gram, kilogram, carat, pennyweight,<br>grain, pound, ounce, troy ounce |                             |                              |                              |
| Applications                          | Weighing, Parts Counting, Percent Weighing, Check Weighing, Animal Weighing, Formulation, Density Determination,<br>Totalization |                            |                                |                          |                          |  |                             |                              |                              |
| Stabilization time (typical)          | ≤ 3<br>seconds   | ≤2 seconds                 |                                |                          |                          | ≤1.5 seconds   |                             |                              |                              |
| Sensitivity Temperature Drift (PPM/K) | 1.5  | 3                          |                                |                          |                          | 3  |                             |                              | 9                            |
| Typical Min-Weight USP (u=0.10%,k=2)  | 0.20g  | 2.0 g                      | 2.0 g                          | 2.0 g                    | 2.0 g                    | 20.0 g   | 20.0 g                      | 20.0 g                       | 200.0 g                      |
| Optimal Min-Weight USP (u=0.10%,k=2)  | 0.082 g  | 0.82 g                     | 0.82 g                         | 0.82 g                   | 0.82 g                   | 8.2 g  | 8.2 g                       | 8.2 g                        | 82 g                         |
| Display                               | Full-Color WQVGA Graphic LCD   |                            |                                |                          |                          |  |                             |                              |                              |
| Display size                          | 4.3 in / 10.9 cm (diagonal)  |                            |                                |                          |                          |  |                             |                              |                              |
| Backlight                             | White LED  |                            |                                |                          |                          |  |                             |                              |                              |
| Controls                              | 4-wire resistive touch screen + 6 membrane keys  |                            |                                |                          |                          |  |                             |                              |                              |
| Communication                         | RS-232, USBx2  |                            |                                |                          |                          |  |                             |                              |                              |
| Balance power input                   | 12 VDC, 0.5A   |                            |                                |                          |                          |  |                             |                              |                              |
| Power supply                          | AC Adapter Input: 100-240 VAC 0.3A 50-60 Hz<br>AC Adapter Output: 12 VDC 0.84A   |                            |                                |                          |                          |  |                             |                              |                              |
| Platform size (diameter)              | 90 mm /<br>3.5 in  | 130 mm / 5.1 in            |                                |                          |                          | 175x195 mm / 6.9x7.7 in  |                             |                              |                              |
| Assembled dimensions (W x D x H)      | 354x230x340 mm<br>13.9x9.1x13.4 inch   |                            |                                |                          |                          | 354x230x100 mm<br>13.9x9.1x4.0 inch                                    |                             |                              |                              |
| Shipping dimensions (W x D x H)       | 557x392x301 mm<br>22.0x15.5x11.9 inch  |                            |                                |                          |                          |  |                             |                              |                              |
| Net weight                            | 5.1 Kg /<br>11.3 lb  | 5.6 Kg /<br>12.4 lb        | 5.8 Kg /<br>12.8 lb            | 5.6 Kg / 12.4 lb         |                          | 3.9 Kg / 8.6 lb  |                             |                              | 3.4 Kg /<br>7.5 lb           |
| Shipping weight                       | 7.8 Kg /<br>17.2 lb  | 8.3 Kg /<br>18.3 lb        | 8.5 Kg /<br>18.8 lb            | 8.3 Kg / 18.3 lb         |                          | 5.8 Kg / 12.8 lb   |                             |                              | 5.3 Kg /<br>11.6 lb          |

**Note:** N = NTEP Approved

TABLE 9-8. SPECIFICATIONS (continued)

| MODEL                                 | AX124AU   | AX224AU                        | AX324AU               | AX223AU                                       | AX423AU                         | AX523AU                         | AX1502AU                                      | AX2202AU                           | AX4202AU                            | AX5202AU                            | AX8201AU                            |
|---------------------------------------|---|--------------------------------|-----------------------|---|---------------------------------|---------------------------------|---|------------------------------------|-------------------------------------|-------------------------------------|-------------------------------------|
| Max                                   | 120g  | 220g                           | 320g                  | 220g  | 420g                            | 520g                            | 1520g   | 2200g                              | 4200g                               | 5200g                               | 8200g                               |
| Min                                   | 0.01g   | 0.01g                          | 0.01g                 | 0.02g   | 0.02g                           | 0.02g                           | 0.5g  | 0.5g                               | 0.5g                                | 0.5g                                | 5g                                  |
| d=                                    | 0.0001g   |                                |                       | 0.001g  |                                 |                                 | 0.01g   |                                    |                                     | 0.1g                                |                                     |
| e=                                    | 0.001g  |                                |                       | 0.01g   |                                 |                                 | 0.1g  |                                    |                                     | 1g                                  |                                     |
| Approval Class                        | I   |                                |                       | II  |                                 |                                 |   |                                    |                                     |                                     |                                     |
| Repeatability (std. dev.) (g)         | 0.0001g   |                                |                       | 0.001g  |                                 |                                 | 0.01g   |                                    |                                     | 0.1g                                |                                     |
| Linearity (g)                         | ±0.0002g  |                                |                       | ±0.002g                                       |                                 |                                 | ±0.02g  |                                    |                                     | ±0.2g                               |                                     |
| Span Calibration Points (g)           | 25g, 50g<br>75g,<br>100g  | 50g,<br>100g,<br>150g,<br>200g | 100g<br>200g,<br>300g | 50g,<br>100g,<br>150g,<br>200g                | 100g,<br>200g,<br>300g,<br>400g | 200g,<br>300g,<br>400g,<br>500g | 500g<br>1000g,<br>1500g                       | 500g,<br>1000g,<br>1500g,<br>2000g | 1000g,<br>2000g,<br>3000g,<br>4000g | 2000g,<br>3000g,<br>4000g,<br>5000g | 2000g,<br>4000g,<br>6000g,<br>8000g |
| Calibration                           | Auto-Cal  | Auto-Cal                       | Auto-Cal              | Auto-Cal                                      | Auto-Cal                        | Auto-Cal                        | Auto-Cal                                      | Auto-Cal                           | Auto-Cal                            | Auto-Cal                            | Auto-Cal                            |
| Weighing units                        | milligram, gram, carat  |                                |                       |   |                                 |                                 | kilogram, gram, carat                         |                                    |                                     |                                     |                                     |
| Applications                          | Weighing, Parts Counting, Percent Weighing, Check Weighing, Animal Weighing, Formulation, Density Determination, Totalization, Display Hold |                                |                       |   |                                 |                                 |   |                                    |                                     |                                     |                                     |
| Stabilization time (typical)          | ≤ 3 seconds   |                                |                       | ≤2 seconds                                    |                                 |                                 | ≤1.5 seconds                                  |                                    |                                     |                                     |                                     |
| Sensitivity Temperature Drift (PPM/K) | 1.5   |                                |                       | 3   |                                 |                                 | 3   |                                    |                                     | 1.9                                 |                                     |
| Typical Min-Weight USP (u=0.10%,k=2)  | 0.20 g  | 0.20 g                         | 0.20 g                | 2.0 g   | 2.0 g                           | 2.0 g                           | 20.0 g  | 20.0 g                             | 20.0 g                              | 20.0 g                              | 200.0 g                             |
| Optimal Min-Weight USP (u=0.10%,k=2)  | 0.082 g   | 0.082 g                        | 0.082                 | 0.82 g  | 0.82 g                          | 0.82 g                          | 8.2 g   | 8.2 g                              | 8.2 g                               | 8.2 g                               | 82 g                                |
| Display                               | Full-Color WQVGA Graphic LCD  |                                |                       |   |                                 |                                 |   |                                    |                                     |                                     |                                     |
| Display size                          | 4.3 in / 10.9 cm (diagonal)   |                                |                       |   |                                 |                                 |   |                                    |                                     |                                     |                                     |
| Backlight                             | White LED   |                                |                       |   |                                 |                                 |   |                                    |                                     |                                     |                                     |
| Controls                              | 4-wire resistive touch screen   |                                |                       |   |                                 |                                 |   |                                    |                                     |                                     |                                     |
| Communication                         | RS-232, USBx2   |                                |                       |   |                                 |                                 |   |                                    |                                     |                                     |                                     |
| Balance power input                   | 12 VDC, 0.5A  |                                |                       |   |                                 |                                 |   |                                    |                                     |                                     |                                     |
| Power supply                          | AC Adapter Input: 100-240 VAC 0.3A 50-60 Hz<br>AC Adapter Output: 12 VDC 0.84A  |                                |                       |   |                                 |                                 |   |                                    |                                     |                                     |                                     |
| Platform size (diameter)              | 90 mm / 3.5 in  |                                |                       | 130 mm / 5.1 in                               |                                 |                                 | 175x195 mm / 6.9x7.7 in                       |                                    |                                     |                                     |                                     |
| Assembled dimensions (W x D x H)      | 230 x 354 x 340 mm<br>9.1 x 13.9 x 13.4 inch  |                                |                       | 230 x 354 x 340 mm<br>9.1 x 13.9 x 13.4 inch  |                                 |                                 | 230 x 354 x 100 mm<br>9.1 x 13.9 x 4.0 inch   |                                    |                                     |                                     |                                     |
| Shipping dimensions (W x D x H)       | 387 x 507 x 531 mm<br>15.4 x 20.0 x 20.9 inch   |                                |                       | 387 x 507 x 531 mm<br>15.4 x 20.0 x 20.9 inch |                                 |                                 | 392 x 557 x 301 mm<br>15.5 x 22.0 x 11.9 inch |                                    |                                     |                                     |                                     |
| Net weight                            | 5.1 Kg / 11.3 lb  |                                |                       | 5.8 Kg / 12.8 lb                              |                                 |                                 | 4.6 Kg / 10.2 lb                              |                                    |                                     | 3.8 Kg / 8.4 lb                     |                                     |
| Shipping weight                       | 7.8 Kg / 17.2 lb  |                                |                       | 8.5 Kg / 18.8 lb                              |                                 |                                 | 6.5 Kg / 14.4 lb                              |                                    |                                     | 5.7 Kg / 12.6 lb                    |                                     |

9.2 Drawings and Dimensions

Fully assembled dimensions

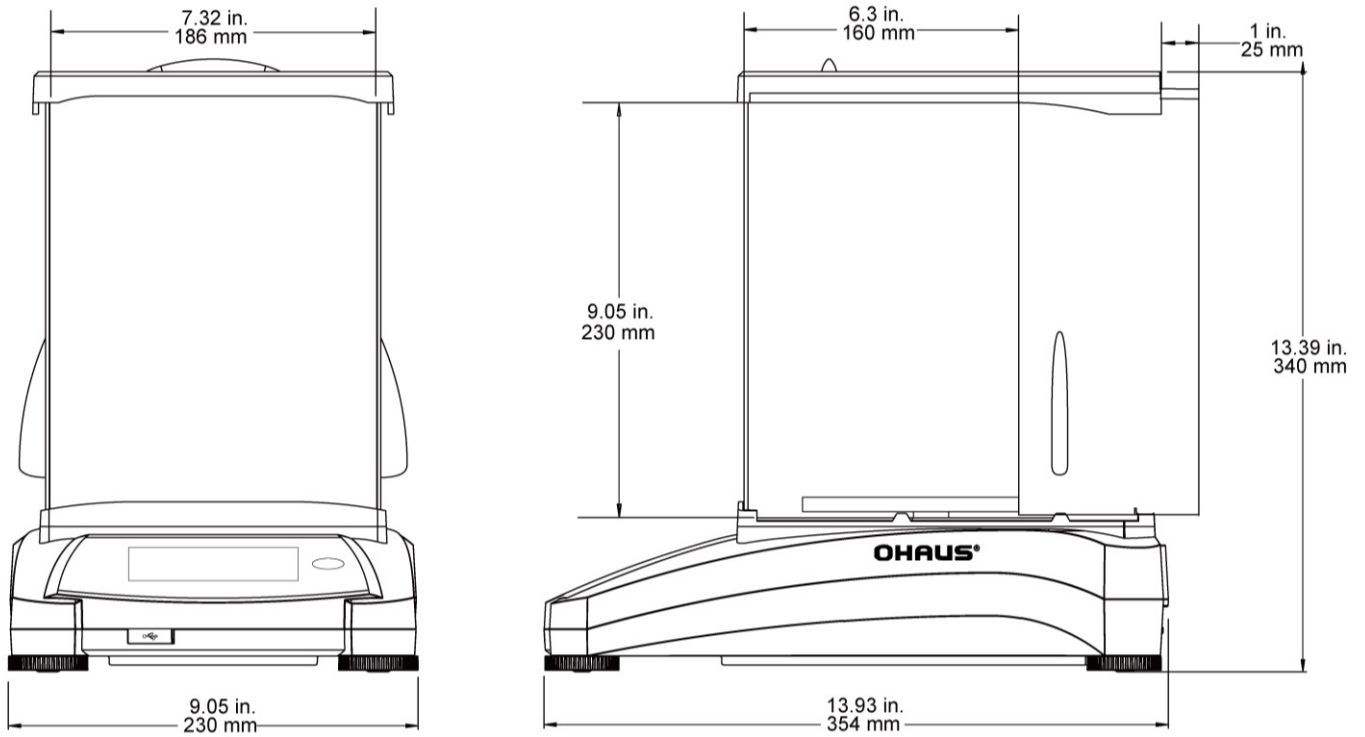
Figure 9-1. Draft Shield models

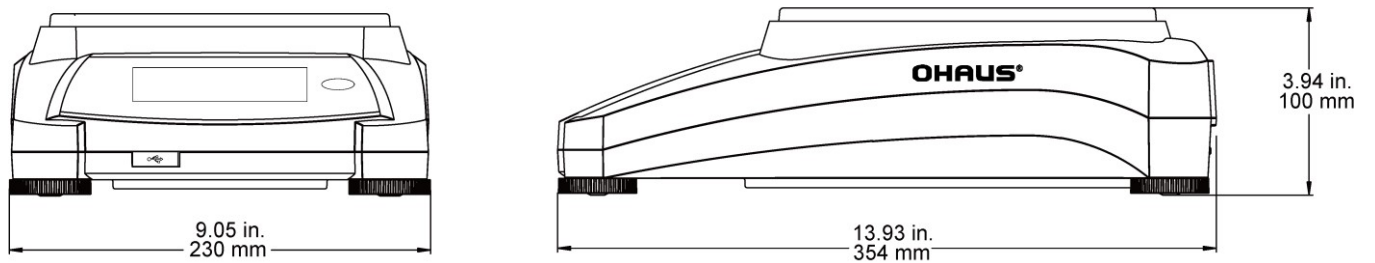
Figure 9-2. Non-Draft Shield models

9.3 Parts and Accessories

TABLE 9-5. ACCESSORIES

| DESCRIPTION                   | PART NUMBER |
|-------------------------------|-------------|
| Auxiliary Display             | 80251396    |
| Density Kit                   | 80253384    |
| Cable, USB Device (Type A-B)  | 83021085    |
| Security Device (Laptop Lock) | 80850043    |
| RS232 Cable, PC 25 Pin        | 80500524    |
| RS232 Cable, PC 9 Pin         | 80500525    |
| Dust Cover                    | 30093334    |

## 9.4 Communication

### 9.4.1 Interface Commands


Commands listed in the following table will be acknowledged by the balance. The balance will return "ES" for invalid commands.

TABLE 9-6. ADVENTURER INTERFACE COMMAND LIST

| Command Characters | Function   |
|--------------------|--|
| IP                 | Immediate Print of displayed weight (stable or unstable).  |
| P                  | Print displayed weight (stable or unstable).   |
| CP                 | Continuous Print.  |
| SP                 | Print on Stability.  |
| SLP                | Auto Print stable non-zero displayed weight.   |
| SLZP               | Auto Print stable non-zero weight and stable zero reading.   |
| xP                 | Interval Print x = Print Interval (1-3600 sec) 0P ends interval Print  |
| 0P                 | See above  |
| H                  | Enter Print Header Lines   |
| Z                  | Same as pressing Zero Key  |
| T                  | Same as pressing Tare Key.   |
| xT                 | Establish a preset Tare value in displayed unit. X = preset tare value. Sending 0T clears tare (if allowed). |
| PT                 | Prints Tare weight stored in memory.   |
| ON                 | Brings out of Standby  |
| OFF                | Goes to Standby.   |
| C                  | Begin Span Calibration   |
| IC                 | Begin internal Calibration, same as trigger from calibration menu.   |
| AC                 | Abort Calibration. Attention: when LFT ON, the operation is not allowed.                                     |
| PSN                | Print Serial Number.   |
| PV                 | Print terminal software version, base software version and LFT ON (if LFT is set ON).                        |
| x#                 | Set Counting APW (x) in grams. (must have APW stored)  |
| P#                 | Print Counting application APW.  |
| x%                 | Set Percent application reference weight (x) in grams. (must have reference weight stored)                   |
| P%                 | Print Percent application reference weight.  |
| xS                 | 0 = print unstable data, same as IP; 1 = print stable only <sup>1)</sup> , same as SP.                       |
| xRL                | 0 = disable response; 1 = enable response. This command only controls the "OK!" response.                    |

### 9.4.2 RS232 (DB9) Pin Connections

- Pin 2: Balance transmit line (TxD)
- Pin 3: Balance receive line (RxD)
- Pin 5: Ground signal (GND)
- Pin 7: Clear to send (hardware handshake) (CTS)
- Pin 8: Request to send (hardware handshake) (RTS)
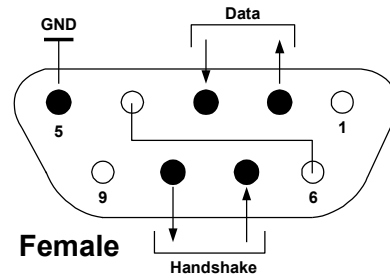
## 9.5 The USB Interface

The Ohaus USB Interface is a unique solution to the problem of connecting a balance to a computer using a Universal Serial Bus (USB). USB devices are categorized into classes such as disk drives, digital cameras, printers, etc. Balances do not have a commonly used class so the Ohaus USB interface uses a generic interface based on the RS232 serial standard.

Data sent from the balance to a computer is in USB format. The USB data is directed to a *virtual port*. This port then appears as an RS232 port to the application program.

When sending a command from a computer to the balance, the application program sends a command to the *virtual port* as if it were an RS232 port. The computer then directs the command from the *virtual port* to the computer's USB connector where the balance is connected. The port receives the USB signal and reacts to the command.

The USB Interface includes a CD with the software drivers to create the required *virtual port* on the computer.

### System Requirements

- PC running Windows 98®, Windows 98SE®, Windows ME®, Windows 2000®, Windows XP®, Windows 7® or Windows 8® (32-bit).
- Available USB port (Type A, 4-pin, female)

### USB Connection

The balance's USB port terminates with a 4-pin, female, USB Type B connector. A USB Cable (type B/male to type A/male) is required (not supplied).

1. Ensure that the balance is powered on and working properly.
2. Power on the computer and verify that its USB port is enabled and working properly.
3. Plug the cable's USB connectors into the computer's USB port and the balance's USB port. Windows® should detect a USB device and the New Hardware Wizard will be initialized.

### Virtual Port Software Installation

1. Insert the supplied CD into the computer's CD drive.  
Different versions of Windows® have slightly different steps to load the driver that is on the CD. In all versions the New Hardware Wizard guides you through the required steps to select the driver that is located on the CD.
2. After clicking Finish, the virtual port should be ready for use.  
Windows® typically adds the virtual port in sequence after the highest number COM port. For example, on PC's equipped with up to 4 COM ports, the virtual port will be COM5.

When using the USB interface with programs that limit the number of COM port designations (e.g. Ohaus MassTracker allows only COM1, 2, 3, & 4), it may be necessary to assign one of these port numbers to the new virtual port.

Example of Windows XP Hardware Wizard

This can be done in the Port Settings of the Device Manager utility, found in the Windows Control Panel.

**USB INPUT**

The balance will respond to various commands sent via the interface adapter. Terminate the following commands when with a [CR] or [CRLF].

**Adventurer Commands**

|             |  |
|-------------|--|
| <b>zC</b>   | perform span calibration                                 |
| <b>0S</b>   | print unstable data                                      |
| <b>1S</b>   | print stable data only                                   |
| <b>P</b>    | same as pressing Print                                   |
| <b>SP</b>   | print stable weight only                                 |
| <b>IP</b>   | immediate print of displayed weight (stable or unstable) |
| <b>CP</b>   | Continuous print of weights                              |
| <b>SLP</b>  | Auto-print stable non-zero weight only                   |
| <b>SLZP</b> | Auto-print stable non-zero weight and zero reading       |
| <b>xP</b>   | Auto-print on 1 to 3600 second intervals (x = 1 to 3600) |
| <b>0P</b>   | Ends interval print                                      |
| <b>T</b>    | same as pressing Tare                                    |
| <b>Z</b>    | same as pressing Zero                                    |
| <b>PV</b>   | print software version                                   |

**Auto-Print Operation**

Once Auto-Print is activated in the menu, the balance will send data as required. If there is data in the print buffer the printer will finish printing this data.

**10. SOFTWARE UPDATES**

Ohaus is continuously improving its balance software. To obtain the latest release, please contact your Authorized Ohaus Dealer or Ohaus Corporation.

## 11. COMPLIANCE

Compliance to the following standards is indicated by the corresponding mark on the product.

| Mark  | Standard  |
|---|---|
|  | This product conforms to the EMC Directive 2004/108/EC, the Low Voltage Directive 2006/95/EC and the Non-automatic Weighing Instruments Directive 2009/23/EC. The complete Declaration of Conformity is available online at <a href="http://europe.ohaus.com/europe/en/home/support/compliance/ce-declaration-of-conformity.aspx">http://europe.ohaus.com/europe/en/home/support/compliance/ce-declaration-of-conformity.aspx</a> . |
|  | AS/NZS CISPR 11, AS/NZS 61000.4.3   |
|  | CAN/CSA-C22.2 No. 61010-1-12<br>UL Std. No. 61010-1 (3rd edition)   |

|   |  |
|---|--|
|  | <p><b>Important notice for verified weighing instruments</b></p> <p>Weighing instruments verified at the place of manufacture bear one of the preceding marks on the packing label and the green 'M' (metrology) sticker on the descriptive plate. They may be put into service immediately.</p>   |
|  | <p>Weighing instruments to be verified in two stages have no green 'M' (metrology) on the descriptive plate and bear one of the preceding identification marks on the packing label.</p> <p>The second stage of the initial verification must be carried out by the approved service organization of the authorized representative within the EC or by the national weights and measures (W+M) authorities.</p> <p>The first stage of the initial verification has been carried out at the manufacturer's work. It comprises all tests according to the adopted European standard EN45501:1992, paragraph 8.2.2.</p> <p>If national regulations limit the validity period of the verification, the user of the weighing instrument must strictly observe the re-verification period and inform the respective W+M authorities.</p> |

### FCC Note


This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.

### Industry Canada Note

This Class A digital apparatus complies with Canadian ICES-001.

**ISO 9001 Registration**

In 1994, OHAUS Corporation, USA, was awarded a certificate of registration to ISO 9001 by Bureau Veritas Quality International (BVQI), confirming that the OHAUS quality management system is compliant with the ISO 9001 standard's requirements. On June 21, 2012, OHAUS Corporation, USA, was re-registered to the ISO 9001:2008 standard. In conformance with the European

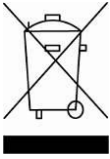**Disposal**

In conformance with the European Directive 2002/96/EC on Waste Electrical and Electronic Equipment (WEEE) this device may not be disposed of in domestic waste. This also applies to countries outside the EU, per their specific requirements.

The Batteries Directive 2006/66/EC introduces new requirements from September 2008 on removability of batteries from waste equipment in EU Member States. To comply with this Directive, this device has been designed for safe removal of the batteries at end-of-life by a waste treatment facility.

Please dispose of this product in accordance with local regulations at the collecting point specified for electrical and electronic equipment.

If you have any questions, please contact the responsible authority or the distributor from which you purchased this device.

Should this device be passed on to other parties (for private or professional use), the content of this regulation must also be related.

For disposal instructions in Europe, refer to <http://europe.ohaus.com/europe/en/home/support/weee.aspx>

Thank you for your contribution to environmental protection.

**LIMITED WARRANTY**

Ohaus products are warranted against defects in materials and workmanship from the date of delivery through the duration of the warranty period. During the warranty period Ohaus will repair, or, at its option, replace any component(s) that proves to be defective at no charge, provided that the product is returned, freight prepaid, to Ohaus.

This warranty does not apply if the product has been damaged by accident or misuse, exposed to radioactive or corrosive materials, has foreign material penetrating to the inside of the product, or as a result of service or modification by other than Ohaus. In lieu of a properly returned warranty registration card, the warranty period shall begin on the date of shipment to the authorized dealer. No other express or implied warranty is given by Ohaus Corporation. Ohaus Corporation shall not be liable for any consequential damages.

As warranty legislation differs from state to state and country to country, please contact Ohaus or your local Ohaus dealer for further details.

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