sartorius stedim
biotech

User Manual | Bedienungsanleitung

Sartocube® | Sartocon® | Sartocon® Slice Sartocon® Slice 200 Cassettes

White Silicone | Weißes Silikon

85032-541-09

Vers. 03 | 2014

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If you require further information on crossflow filtration, other than mentioned in this User Manual, please contact us.

Reproduction, translation and duplication in any form, in part or in full, requires the written permission of the Sartorius Stedim Biotech GmbH.

For technical information write to:

Sartorius Stedim Biotech GmbH
August-Spindler-Strasse 11
D-37079 Goettingen
Internet: www.sartorius-stedim.com

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Contents

1. General	4
1.1 Instructions to this Manual	4
1.2 Storage of unused cassettes	4
1.3 Safety Instructions	5
1.4 Components of a Crossflow System	8
2. Start Up Of Cassettes	9
2.1 Installing Filter Cassettes in a Holding Device	9
2.2 Clamp Force for Sartorius Stedim Biotech Cassettes	10
2.3 Flushing the Cassettes	11
2.4 Integrity Testing	12
2.4.1 Manual Integrity Test	12
2.4.2 Automatic Integrity Test	13
2.5 Determining Clean Water Flux for Filter Cassettes	14
2.6 Thermal Sterilization of Filter Cassettes	15
2.6.1 In Line Steam Sterilization of Sartocon [®] Microfilter Cassettes	16
2.6.2 In Line Steam Sterilization of Sartocon [®] Ultrafilter Cassettes	18
2.6.3 Autoclaving the Filter cassettes	21
2.7 Filtration	22
2.8 Cassette Cleaning	23
2.9 Cassette Disinfection	24
2.10 Cassette Storage	25
2.11 Cassette Return Instructions	26
2.12 Cassette Disposal	26
3. Appendix	27
3.1 Specifications Sartocon [®] Slice 200	27
3.2 Specifications Sartocon [®] Slice	27
3.3 Specifications Sartocon [®]	28
3.4 Specifications Sartocube [®]	28

1. General

1.1 Instructions to this manual

- ▶ Please be sure to read this user manual prior to use of equipment.
- ▶ Please read all safety instructions.
- ▶ This user manual is part of the product. Keep it in a safe place for future reference.
- ▶ Replacement manuals can be requested by telephone or it can be downloaded from our Webpage at: www.sartorius-stedim.com

Symbol

The following symbols are used in this user manual:

This symbol indicates situations in which improper procedure could lead to personal injury.

This symbol indicates that the preservative is not suitable for consumption.

This symbol indicates situations in which improper procedure could lead to equipment damage.

This symbol draws attention to additional information.

1.2 Storage of unused cassettes

The latest date of usage is 3 years from date of manufacture.

The following conditions have to be maintained during storage in order to guarantee full cassette functionality:

Cassettes must be stored in a closed dry area.

Temperature has to be between 5°C to 40°C and absolutely frost free.

No direct solar radiation.

No direct contact with moisture.

Any mechanical influence or damage must be prevented.

Products with damaged packaging should be discarded.

1.3 Safety instructions

The preservative is not suitable for consumption.

During in-line sterilization of the holding device with hot steam, the entire system becomes very hot. Do not touch any part of the system with your bare hands. If you need to handle any part of the system, wear heavy safety gloves to avoid getting severely burned.

During in-line steam sterilization, hot steam exits from the feed inlet and from the permeate outlets. Keep a safe distance to these outlets and do not try to reach inside them; as there is a serious risk of severe burns.

Do not apply more than the torque given in the "**Specifications**". Torquing the filter cassettes too tightly will damage them.

The maximum inlet pressure must not be exceeded (see "**Specifications**"). Over pressurizing will damage the filter cassettes.

Before thermal sterilization, exchange the torque nuts for normal operation to the torque nuts for pressure compensation. These pressure compensation torque nuts keep the torque constant even when the temperature fluctuates. If the pressure compensating torque nuts are not used, the torque will increase beyond the maximum allowable compression during thermal sterilization, causing irreversible damage to the filter cassettes.

Do not allow any type of filter cassette to dry out to avoid loss of filtration properties.

Make sure that the filter cassettes you plan to use have the required chemical compatibility before beginning a filtration run in order to avoid damaging the membranes. If you are not sure about the compatibility of a filter cassette with a certain chemical, please consult your nearest Sartorius Stedim Application specialist.

Do not combine cassettes with white and gray silicone encapsulation in one holding device.

Avoid abrupt changes in pressure (approximately > 1 bar/sec), such as those that can occur when the pump is switched on or when the valves settings are changed.

Keep vibration to the lowest level that is technologically possible.

The effective torque on hydraulic holding devices should be kept constant, especially when the holding device is installed in a large filtration system.

The maximum allowable back pressure on the permeate side must not exceed 0.5 bar for Hydrosart® UF and 1 bar for PESU max and Hydrosart® open channel cassettes.

The outer silicone sealing of cassettes should be washed of with DI water, after submerged storage, before installation into a cassette holding device.

Retentate backpressure is defined as follows:
Permeate pressure is higher than the Retentate pressure

$$P_{\text{perm}} > P_{\text{ret}} ;$$

$$P_{\text{back}} = P_{\text{perm}} - P_{\text{ret}} .$$

Hydrostatic pressure on the permeate side (height and the inner diameter of the filtrate tubing) has an impact on backpressure.
Please make sure that the maximum backpressure pressure is never exceeded.

Note

This has to be distinguished from Transmembrane Pressure.

Transmembrane pressure (TMP):

$$\text{TMP} = [(p_{\text{feed}} + p_{\text{ret}})/2] - p_{\text{perm}}$$

and from Differential Pressure:

$$\Delta p = p_{\text{feed}} - p_{\text{ret}}$$

Feed pressure (p_{feed}) and retentate pressure (p_{ret}), both measurable with a gauge.

Permeate pressure (p_{perm}) is the average pressure on the permeate side, measurable with a gauge.

If the possibility of crystalline particles entering the cassettes cannot be completely ruled out (for example, due to abrasion in the pumps), a prefilter must be used.

Make sure that crystallization caused by changes in temperature or concentration cannot occur during operation or storage.

Do not combine cassettes that have a different feed inlet design. Sartocon cassettes have alternating oval and round ports.

Install the filter cassettes in the appropriate holding device only. You will need either a Sartocon® Slice, Sartocon® 2 plus or a SARTOFLOW® holding device, depending on the type and amount of filter cassette you plan to use. If you wish to use holding devices supplied by other manufacturers, you may need to apply a different torque pressure when installing cassettes. To avoid damaging the cassettes, please contact our application specialist.

Make sure that the correct crossflow conditions are maintained during filtration. Otherwise, the flux rate of the filter cassettes will decrease rapidly during filtration under static conditions (dead-ended filtration) (refer to "Filtration").

The temperature of the retentate must not exceed that given for continuous operation, which is listed in the "Specifications." During crossflow filtration, the temperature will increase over time because the kinetic energy of the pump is transferred to the liquid undergoing filtration. If the temperature increases beyond that indicated in the "Specifications," the pressure resistance of the crossflow system will decrease.

Sartocon® Slice 200

Sartocon® Slice Cassette

Sartocon[®] Cassette

Sartocube[®]

You should use a prefilter. This will enhance the flux rate and increase the service life of the filter cassettes.

You should clean the filter cassettes after every filtration run. During filtration a layer caused by fouling builds up on the membranes. Over time, this layer becomes compacted, making it difficult to remove completely. As a result, the flux rate of your filter cassettes will decrease.

All filter cassettes used in one holding device must have identical membranes.

If you do not need to use filter cassettes for a relatively long period, you should store them in a solution that contains a preservative. This will protect the filter cassettes from microbial contamination and from drying out.


Do not attempt to operate the filtration system with accessories or spare parts other than those supplied by the original manufacturer. In the event that the filtration system is tampered with in any way, we cannot guarantee or take any responsibility for the performance or safety of the system. This may also endanger the safety of the operator.

Use only materials recommended by the manufacturer (e.g. gaskets, tubing or hoses).

1.4 Components of a Crossflow System

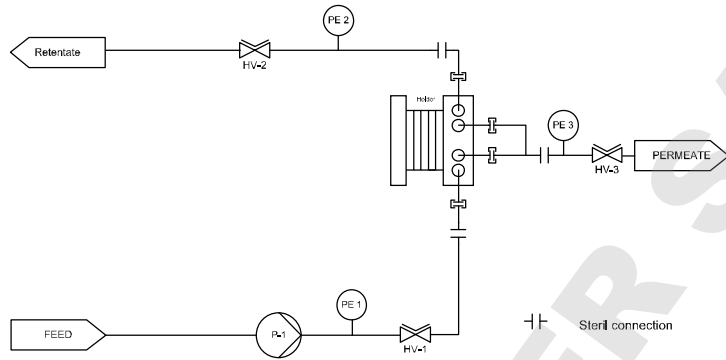Holding Device

Sartocoon Slice 200 (Holder for max. 2 Slice 200 Cassettes)
Mat. No. #17525---01

Sartocoon Slice (for 1-3 Slice Cassette Holder) Mat. No. #17521---001 Need for No. #: 1ZLI-V0001 additional long-spacers to run a Slice 200 cassette

Sartocoon Slice (for 1-5 Slice Cassette Holder [Sartoflow Alpha])
Mat. No. #17521---002

Sartocoon® 2 plus 1-10 Cassettes
[Order Nr.: 17546*]

SARTOFLOW® 10 Holder
1-10 Cassettes

SARTOFLOW® 20 Holder
2-20 Cassettes

2.0 Start up of cassettes

2.1 Installing Filter Cassettes in a Holding Device

Semi-circular cutouts

Sartocor[®] 2 plus
[Order Nr.: 17546*]

Pressure compensating torque nuts

To install the cassettes in a Sartorius Stedim Biotech holding device, proceed as follows:

1. Remove filter cassette from its packaging.
2. Check and record cassette labelling.
3. Place the filter cassette with the semi-circular cutouts facing downward on the two guide rods between the two holding plates.
4. Alternately tighten the compression nuts until the filter cassette(s) is (are) secured in place by the holding plates.
5. Use a torque wrench to tighten the compression nuts alternately to the torque specified. Please see "**Clamp Force for Sartorius Cassettes**".
6. After 15 minutes, check the torque and, if necessary, re-adjust to the torque specified.

Do not torque the holding device beyond the torque specified.

Keep the threads of the holding device clean to ensure adequate torque.

2.2 Clamp Force for Sartorius Stedim Biotech Cassettes

Sartocon® Slice 200 Holding Device, Manual closure

	Required Compression Operating Conditions	Sterilization*	Torque applied on Tie Rods with a Torque Wrench Operating Conditions	Sterilization*
PESU UF MF Sartocon® Slice 200 Cassettes 308146... 308154...	18–21 kN	No thermal sterilization	25 Nm	No thermal sterilization
PESUmax (308146AL...)	25 kN	No thermal sterilization	25 Nm	No thermal sterilization
Hydrosart® UF MF Sartocon® Slice 200 Cassettes 308144... 308186...	14–17 kN	No thermal sterilization	20 Nm	No thermal sterilization

Sartocon® Slice Holding Device, Manual closure

	Required Compression Operating Conditions	Sterilization*	Torque applied on Tie Rods with a Torque Wrench Operating Conditions	Sterilization*
Sartocon® Slice Cassettes MF UF 305146... 305154...	18–21 kN	< 2 kN	25 Nm	first notch
PESUmax 305146AL...	18–21 kN	No thermal sterilization	25 Nm	No thermal sterilization
Hydrosart® MF UF Sartocon® Slice 305144... 305186...	14–17 kN	< 2 kN	20 Nm	first notch

Sartocon® 2 plus Holding Device, Manual closure

	Required Compression Operating Conditions	Sterilization*	Torque applied on Tie Rods with a Torque Wrench Operating Conditions	Sterilization*
Sartocube® and Sartocon® Cassettes 302146... 302154...	18–35 kN	< 10 kN	80 Nm	first notch
PESUmax 302146AL...	21–35 kN	No thermal sterilization	80 Nm	No thermal sterilization
Hydrosart® MF UF Sartocube® and Sartocon® Cassettes 302144... 302186...	18–25 kN	< 10 kN	60 Nm	first notch

Sartoflow® 10 | Sartoflow® 20 and Sartoflow® Systems with Hydraulic Piston, 12.5 cm²

	Required Compression Operating Conditions	Sterilization*	Hydraulic Pressure Operating Conditions	Sterilization*
Sartocube® and Sartocon® Cassettes 302146... 302154...	18–35 kN	< 10 kN	170 bar (2466 PSI) (17 MPa)	80 bar (1160 PSI) (8 MPa)
PESUmax 302146AL...	21–35 kN	No thermal sterilization	170 bar (2466 PSI) (17 MPa)	No thermal sterilization
Hydrosart® MF UF Sartocube® and Sartocon® Cassettes 302144... 302186...	18–25 kN	< 10 kN	170 bar (2466 PSI) (17 MPa)	80 bar (1160 PSI) (8 MPa)

2.3 Flushing the Filter Cassettes

Flush the filter cassettes as indicated below:

- ▶ Before initial use with purified water
- ▶ With buffer or physiological saline before filtration
- ▶ With buffer or physiological saline before cleaning
- ▶ With purified water after cleaning
- ▶ With purified water after disinfection
- ▶ With purified water before checking the functionality of the filter cassettes
- ▶ With purified water before sterilization

Flush Procedure for Cassettes

The installed cassettes must be flushed with purified water or RO water at the recommended approx. ratio of retentate to permeate flow as listed below.

Cassette Type	Ratio of Retentate to Permeate flow
Ultrafilters 1 kD 2 kD	Approx. 5 : 1
Ultrafilters 5 kD	Approx. 3 : 1
Ultrafilters 8 kD up to 300 kD and	
Microfilters 0,1 µm	Approx. 1 : 1
Microfilters 0,2 µm and 0,45 µm and PESU MAX	Approx. 1 : 2

The cassette must be flushed with 10 l purified of water per 0.1 m² filter area.

2.4 Integrity Testing

For multiple cassette integrity testing see "#1 Technical Data Package Integrity testing Sartocor[®] Cassettes".

Ensure the filter cassettes are sufficiently wet before measuring the air diffusion.

2.4.1 Manual Integrity Test

Proceed as follows:

1. Attach a hose to the feed inlet (V_1), and connect the other end of the hose to a compressed air source.
2. Close the retentate outlet (V_2).
3. Adjust the inlet test pressure p_{inlet} to the value specified (see "Specifications" section 2.16), and keep this pressure constant.
4. Close the lower permeate outlet valve V_3 .
5. At a constant inlet pressure p_{feed} , wait four minutes until the system has stabilized.
6. Use a flow meter or an inverted burette filled with water to measure the air diffusion at the upper permeate outlet.

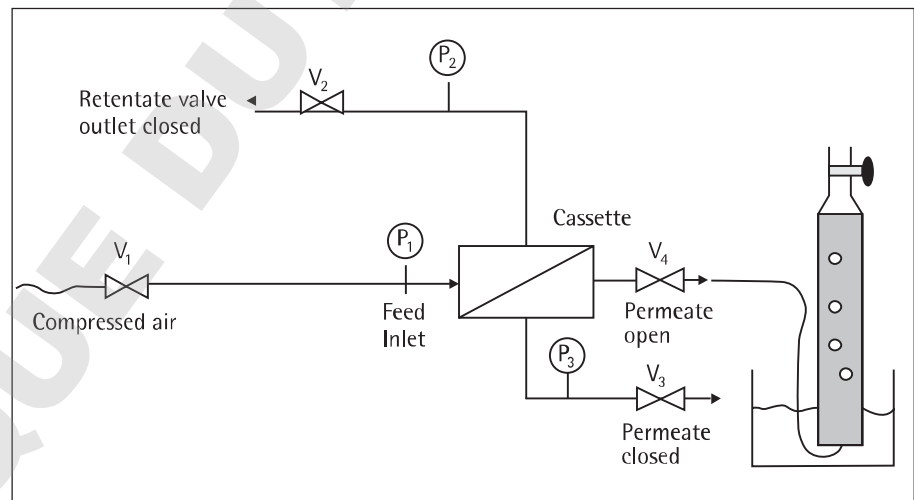
Fig. 1: System set up for manual integrity testing

2.4.2 Automatic Integrity Test (Pressure Decay)

Proceed as follows:

1. Connect the integrity tester [such as Sartochek[®] 3 plus (16290) or Sartochek[®] 4 Plus (26288)] to the power supply and pressure source.
2. Close the retentate outlet (V_2) and open the permeate outlets (V_3, V_4). Open V_5 and close V_1 .
3. Enter the required test parameter according to the Sartochek[®] instruction manual and cassette data.
4. Use the following parameter:
stab. time: 4 min.
test time: 2 min.
5. Enter the maximum allowable diffusion value (see "Specifications").
6. Start the Integrity Test.

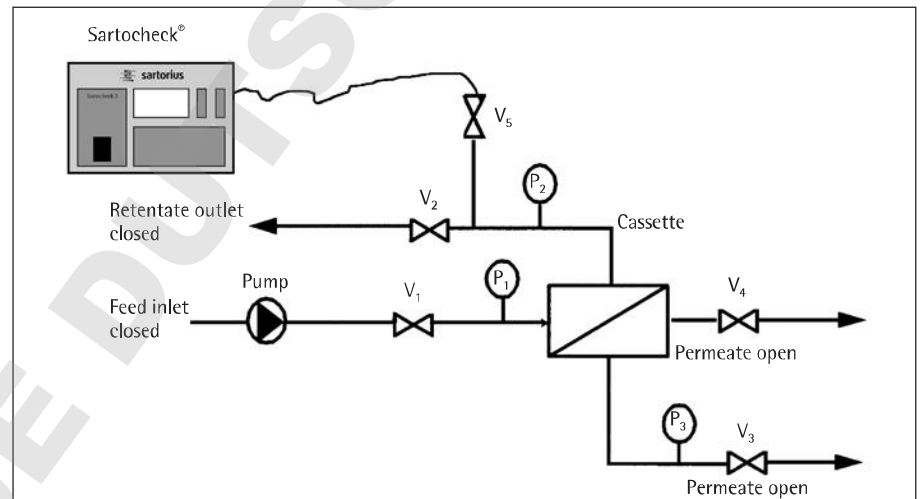
Fig. 2: System set up for automatic integrity testing. V_1 - V_5 : valves; P_1 - P_3 : pressure gauge.

If the hold up volume exceeds 9 l, please contact your Sartorius Stedim application specialist.

2.5 Determining Clean Water Flux for Filter Cassettes (CWF)

The water flux test is an easy and non-destructive test that measures the ability of water to pass through the membrane in a specific system setup. Residuals on the membrane surface will impact water flux rates. Determining the water flux rate enables the user to confirm the ability of a cleaning cycle to restore filtration rate of the system. The water flux should be trended over time and evaluated for potential membrane and system changes.

It cannot be used to assess the removal of preservative, cleaning or process residue. Specific analytical techniques are required to quantify residues.

Determine a clean water flux

Proceed as follows:

1. Rinse new installed cassettes to remove preservatives.
2. Post use sanitization cassettes with 1 N Sodium hydroxide and neutralize.
3. Run the clean water flux with purified water as specified in the table below.
4. Measure retentate rate and permeate rate at constant temperature (20°C).
5. Trend the clean water flux and water flux determined after cleaning to detect membrane and system changes. The clean water flux is used as the base line value.

Any changes in water quality, applied pressure, temperature and viscosity, system design and membrane, will have an effect on the water flux.

In systems containing more than one Sartocube[®] or a Sartocube[®], the interpretation can be simplified by establishing a baseline for a given set of devices in a given system.

Flow path restriction may not allow performing the flux test with the above-mentioned parameters. In that case a different test pressure should be used to perform a water flux test. Once the parameters are defined they have to remain constant during every test to assure comparability.

System Pressure Examples for Water Flux Determination:

Cassette Design	Feed Pressure valve open [psi bar kPa]	Retentate Pressure [psi bar kPa]	Permeate Pressure valve open [psi bar kPa]
Sartocube [®] Slice Cassette and Slice 200	29 2 200	7 0.5 50	0 0 0
Sartocube [®] Cassette	29 2 200	7 0.5 50	0 0 0
Multiple Sartocube [®] Cassettes	29 2 200	10-14 0.7-1 70-100	≤ 7 0.5 50
Sartocube [®] Cassette	29 2 200	10-14 0.7-1 70-100	≤ 7 0.5 50

2.6 Thermal Sterilization of Filter Cassettes

Although the flux rate of used filter cassettes will normally not be fully regained, the flux rate after use and cleaning should approximate the clean water flux.

Trending the water flux helps to evaluate the cleaning status of the system.

Cleaning optimization should be considered:

- ▶ If the permeate flux rate is below 70% of the clean water flux.
- ▶ If the permeate flux rate decreases from run to run.

In case of further questions please contact our Sartorius Stedim application specialist.

Check steamability with "Specifications".

During in-line steam sterilization the entire system, including stainless steel piping and cassette filter holder, will attain $\geq 121^{\circ}\text{C}$. Exercise due caution when working around hot equipment.

Please note that steaming in place (SIP), due to high thermal and mechanical stress, may impact the physical properties of a membrane in terms of flow performance and rejection. Therefore multiple use is limited.

During in-line steam sterilization, hot steam exits from the feed inlet and the permeate outlets. Keep a safe distance away from these outlets.

Before thermal sterilization, exchange the operating torque nuts to the torque nuts for pressure compensation.

This special compensation tool is designed to maintain constant clamp force during thermal sterilization. If the pressure compensation tool is not used, the torque will increase beyond the maximum allowable value as temperature increases. This may cause irreversible damage to the filter cassettes.

Thermal sterilization should only be performed on clean filter cassettes. Contaminants will be irreversibly fixed by thermal sterilization.

The Stabilizing pressure inlet air should be sterile filtered.

Make sure that the steam you will be using is of pharmaceutical quality, free of rust and other particles and does not contain any hydrazine.

Prior to thermal sterilization check the filter cassettes for residual traces of chlorine or disulfide. To denature chlorine, flush with 1% sodium thiosulfate solution followed by flushing with purified water.

Handling safety is increased if you attach thermal condensate run-off valves to the retentate outlet and to the permeate outlets.

2.6.1 In-Line Steam Sterilization of Sartocon® Microfilter Cassettes

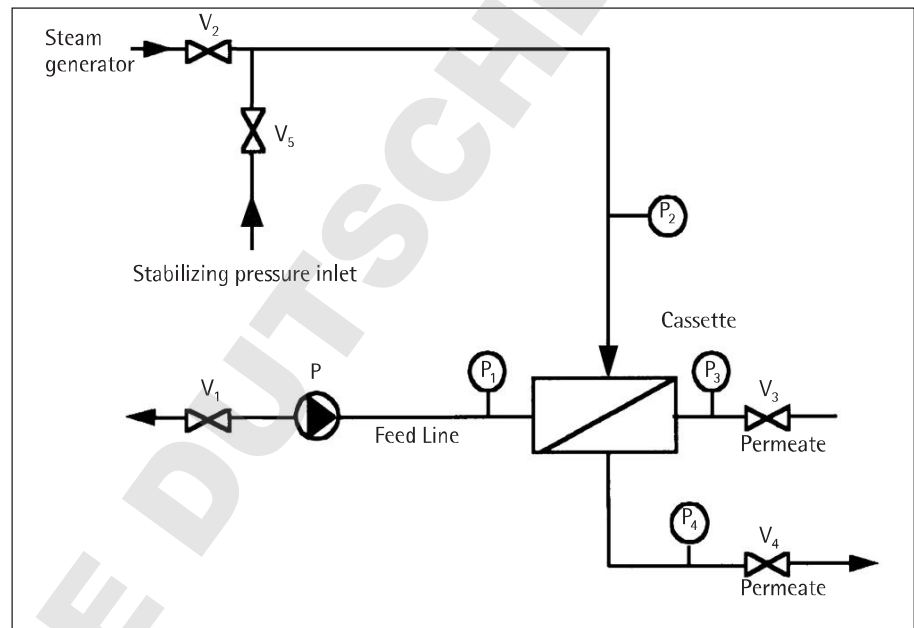
Diagram of the system for in-line steam sterilization of Microfilter cassettes.
V₁-V₅: valves; P₁-P₄: pressure gauge, P: pump

1. Assemble system as shown on the diagram.
2. Open V₂, V₃, V₄ and V₁, close V₅.
3. Install cassettes into filter holder.
4. Attach a RO water source to the feed line.
5. Rinse cassette with RO water for 5 minutes at an inlet pressure of 29 psi|2 bar|200 kPa and retentate pressure of 7 psi|0.5 bar|50 kPa. Adjust valve settings to achieve pressure values if needed.
6. Drain the system.
7. Close valve V₂. Connect V₂ at the retentate line with a steam generator.

8. Exchange the normal operation torque nuts to the pressure compensating torque nuts. For specific information related to the filter holder see user manual of the installed filter holder.
9. Use a torque wrench to tighten the pressure compensation nut to the first notch (see picture).
10. Close valve V_3 and V_4 on the permeate outlets as well as V_5 .
11. Adjust the steam generator pressure to 14.5 psi | 1.0 bar | 100 kPa, and slowly open valve V_2 to initiate steam flow in to the system.
12. We highly recommend the installation of condensate run off valves. The following step has to be incorporated if no condensate run off valves are used. Reduce V_1 valve setting on the feed line to allow a small amount of steam to exit the system. Slowly start to increase pressure build up to 14.5 psi | 1.0 bar | 100 kPa by closing V_1 . A small amount of steam has to remain exiting the system, which means V_1 will not be closed completely.


The retentate pressure must not exceed 14.5 psi | 1.0 bar | 100 kPa at 121°C.

If inline pressure gauges are used, please consider that high temperatures may cause the display to be inaccurate. Usually the actual pressure is lower than indicated. Make sure the pressure is correctly monitored during your validation procedure.

Pressure compensating torque nuts

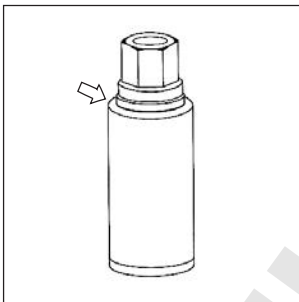first notch (steaming position)

second notch (processing position)

13. We highly recommend the installation of condensate run off valves. The following step has to be incorporated if no condensate run off valves are used. Once the feed pressure P_1 has reached 14.5 psi | 1.0 bar | 100 kPa slowly start to open the permeate valves V_3 | V_4 . A small amount of steam has to escape from each valve.
14. Once the permeate pressure has reached 14.5 psi | 1.0 bar | 100 kPa the actual sterilization begins.
15. After 30 minutes, sterilization is completed. Close valves V_3 | V_4 .
16. Close valve V_1 on the feed line, V_2 on the retentate line. Open V_5 to apply sterile filtered stabilizing pressure of 14.5 psi | 1.0 bar | 100 kPa. At this point, shut off the steam supply.
17. Allow the holding device to cool down for 90 minutes. During cool down V_5 has to remain open.
18. After 90 minutes cool down time close V_5 .
19. Vent system through upper permeate valve V_3 . Do not open V_3 completely. Only vent excess pressure. If sterile air venting is required, please implement a sterilizing air filter.

- Using a torque wrench, alternately retighten the stainless steel pressure compensating nuts to the torque indicated below:

Clamp Force for Sartocor[®] Cassettes

For detailed information see Page 9

Holding Device	Clamp Force
Sartocor [®] Slice	25 Nm Polyethersulfone 20 Nm Hydrosart [®]
Sartocor [®] 2 plus	80 Nm Polyethersulfone 60 Nm Hydrosart [®]
SARTOFLOW [®] 10 and 20	For details see instruction manual.

- After steam sterilization is completed flush filter cassettes with RO water until system has reached room temperature, this will also ensure that the membranes will not dry out.

2.6.2 In-Line Steam Sterilization of Sartocor[®] Ultrafilter Cassettes

Ultrafilter cassettes are sterilized through the retentate and the permeate flowpaths simultaneously to maintain equal pressure within the cassette.

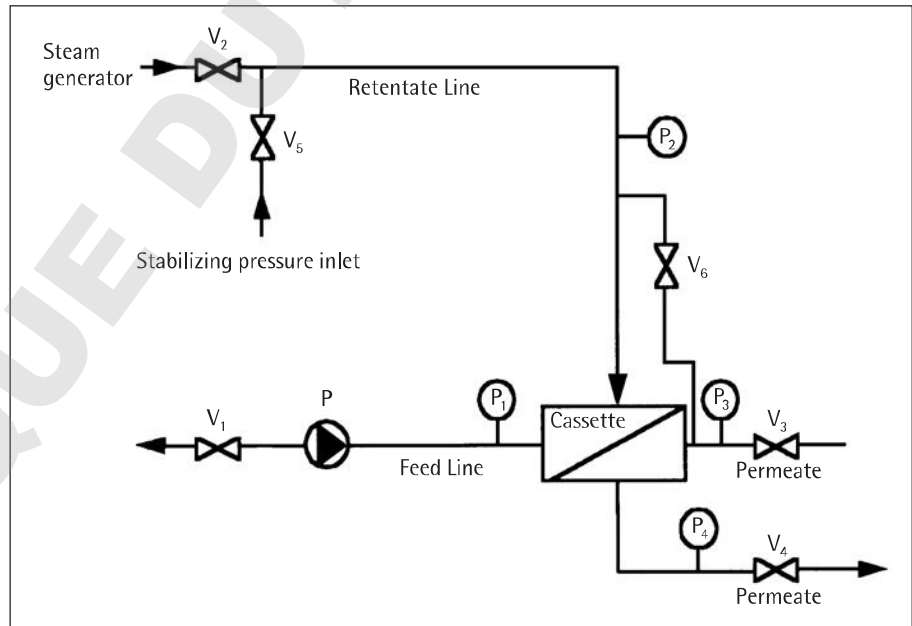
Diagram of the crossflow system for in-line steam sterilization of polyethersulfone Ultrafilter cassettes. V₁-V₆: valves; P₁-P₄: pressure gauge, P: pump

1. Assemble system as shown on the diagramm.
2. Open V_2 , V_3 , V_4 , V_6 and V_1 , close V_5 .
3. Install cassettes into filter holder.
4. Attach a RO water source to the feed line.
5. Rinse cassette with RO water for 5 minutes at an inlet pressure of 29 psi|2 bar|200 kPa and retentate pressure of 7 psi|0,5 bar|50 kPa. Adjust valve settings to achieve pressure values if needed.
6. Drain the system.
7. Close valve V_2 . Connect V_2 at the retentate line with a steam generator.
8. Exchange torque the normal operation nuts to the pressure compensating torque nuts. For specific information related to the filter holder see user manual of the installed filter holder.
9. Use a torque wrench to tighten the compensation nut to the first notch (see picture).
10. Open the connection between the retentate and permeate V_6 .
11. Close Valve V_3 and V_4 on the permeate outlets as well as V_5 .
12. Adjust the steam generator pressure to 14,5 psi|1,0 bar|100 kPa, and slowly open Valve V_2 to initiate steam flow in to the system.
13. We highly recommend the use of condensate run off valves. The following step has to be incorporated in case no condensate run off valves are used. Reduce V_1 valve setting on the feed line to allow a small amount of steam to exit the system. Slowly start to increase pressure build up to 14,5 psi|1,0 bar|100 kPa by partially closing V_1 . A small amount of steam has to remain exiting the system, which means V_1 will not be closed completely.


The retentate pressure must not exceed 14,5 psi|1,0 bar|100 kPa at 121°C.


If inline pressure gauges are used, please consider that high temperatures may cause the display to be inaccurate. Usually the actual pressure is lower than indicated. Make sure the pressure is correctly monitored.

Pressure compensating torque nuts

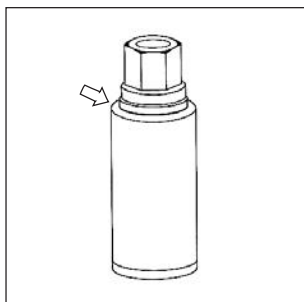First notch (steaming position)

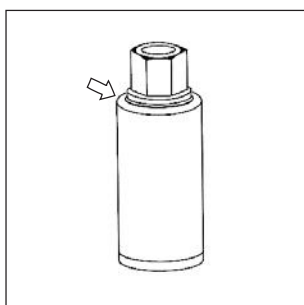Second notch (processing position)

14. We highly recommend the installation of condensate run off valves. The following step has to be incorporated in case that no condensate run off valves are used. Once the feed pressure P_1 has reached 14.5 psi | 1.0 bar | 100 kPa slowly start to open the permeate valves $V_3|V_4$. A small amount of steam has to escape from each valve.
15. Once the permeate pressure has reached 14.5 psi | 1.0 bar | 100 kPa the actual sterilization begins.
16. After 30 minutes, sterilization is completed. Close valves $V_3|V_4$.
17. Close valve V_1 on the feed line, V_2 on the retentate line. Open V_5 to apply sterile filtered stabilizing pressure of 1 bar. At this point, shut off the steam supply.
18. Allow the holding device to cool down for 90 minutes. During cool down V_5 has to remain open.
19. After 90 minutes cool down time close V_5 .
20. Vent system through upper permeate valve V_3 . Do not open V_3 completely only vent excess pressure. If sterile air venting is required, please implement a sterilizing air filter.
21. Using a torque wrench, alternately retighten the stainless steel pressure compensating nuts to the torque indicated below:

Clamp Force for Sartocor[®] Cassettes


Holding Device	Clamp Force
Sartocor [®] Slice and Slice 200	25 Nm Polyethersulfone 20 Nm Hydrosart [®]
Sartocor [®] 2 plus	80 Nm Polyethersulfone 60 Nm Hydrosart [®]
SARTOFLOW [®] 10 and 20	For details see instruction manual.

22. After steam sterilization flush filter cassettes with RO water until system has reached room temperature, this will also ensure that the membranes will not dry out.

2.6.3 Autoclaving the Filter Cassette

1. Tighten the filter cassettes along with the pressure compensation nuts of the holding device.
2. Flush the filter cassettes with purified water as described under "Flushing the Filter Cassettes" section 2.3.
3. Use the torque wrench to tighten the pressure compensation nuts alternately up to the first notch (autoclaving position).
4. Connect a hose to each connection of the holding device.
5. Place the holding device in the autoclave.
6. Fill the hoses with purified water. Immerse all hose ends in a container filled with 5 l of purified water. Place this container in the autoclave as well. **This will prevent the filter cassettes from drying out during autoclaving.**
7. Start the autoclave. To avoid damage of the filter cassette, autoclave the system without final vacuum.
8. Autoclave the holding device and filter cassettes for 30 min. at 121°C.
9. Seal the ends of the hoses immediately after you have opened the autoclave. For this purpose, you can use either hose clamps or valves. **This will prevent contamination and will protect the filter cassettes from drying out.**
10. Allow the holding device to cool off.
11. Use the torque wrench to tighten the pressure compensation nuts alternately up to the second notch:

Pressure compensating torque nuts

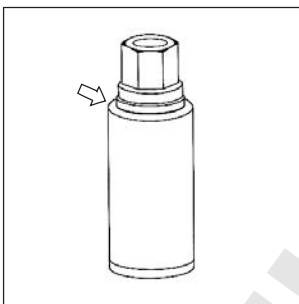First notch (autoclaving position)

Second notch (processing position)

Clamp Force for Sartocor[®] Cassettes

Holding Device	Clamp Force
Sartocor [®] Slice and Slice 200	25 Nm Polyethersulfone 20 Nm Hydrosart [®]
Sartocor [®] 2 plus	80 Nm Polyethersulfone 60 Nm Hydrosart [®]
SARTOFLOW [®] 10 and 20	For details see instruction manual.

2.7 Filtration

A critical time during crossflow filtration is at start up. The process fluid to be filtered initially contacts the clean membrane surface. Until the desired steady state processing conditions are attained, a higher than expected permeate flux rate will prevail. This is especially true for microfiltration.

Please consider:

- ▶ Partly close permeate valve to limit permeate flow.

Starting a Crossflow filtration process:

1. Prime the permeate and retentate areas with one of the liquids listed below:

- ▶ Isotonic physiological saline
- ▶ Processing buffer

To prime the system using one of the liquids listed above, follow the directions given in the section entitled "**Flushing the Filter Cassettes**" section 2.3.

2. Reduce permeate outlet valves setting to allow a small amount of fluid to exit.

3. Generate the desired crossflow rate.

4. Open the permeate outlet valve until the desired TMP is attained.

5. Make sure that the desired crossflow rate is maintained during the complete process.

6. Reset the desired process TMP if required.

The pressure at the feed inlet p_{inlet} must not exceed 58 psi|4 bar.

The maximum allowable backpressure for all cassettes is: ≤ 14 psi|1.0 bar| 100 kPa. Except for Hydrosart® UF.

During filtration, the crossflow rate should be held constant. Viscosity increase during filtration will increase the inlet pressure.

7. At the end of filtration, flush the system (see "**Flushing the Filter Cassettes**").

8. After flushing, clean the system (see "**Cleaning the Filter Cassettes**").

Filter cassettes should be cleaned after every filtration run to prevent a sediment layer build-up on the membrane and in the retentate channels. Over time, this layer becomes compacted, making it difficult to remove completely. As a result, this will decrease the flux rate of your filter cassettes and reduce their lifespan.

We recommend the following cleaning agents and conditions for the various membrane materials. This table gives suggested procedures only; requirements may vary in individual cases:

Cleaning agent	Concentration	pH	Time [min]	Temperature [°C]
Polyethersulfone and Hydrosart®				
Sodium hydroxide	1 N	14	60	50
Phosphoric acid	2% by weight	1.3	30	50

Adjust the following pressures:

$$p_{\text{inlet}} = 29 \text{ psi} | 2.0 \text{ bar} | 200 \text{ kPa}$$

$$p_{\text{ret}} = 0 \text{ psi} | 0 \text{ bar} | 0 \text{ kPa}$$

Microfilter cassettes:

Recirculate cleaning solution with permeate valves closed for 60 minutes. Repeat with fresh cleaning solution and permeate valve open. Repeat as required.

Ultrafilter cassettes:

Recirculate cleaning solution with permeate valves open for 60 minutes. Repeat as required.

Flush the filter cassettes with purified water as described in "**Flushing the Filter Cassettes**". Determine the flux rate of the cleaned cassette as described in "**Determining clean water flux for filter cassettes**". Compare results with the previously determined clean water flux (CWF). It is normal to have a greater drop in clean water flux rate after the first processing and complete cleaning. Good cleaning practice will result in less water flux deviation.

If the flux rate remains below this level after further cleaning attempts, please contact your local Sartorius Stedim application specialist.

Generally, the Sartorius Stedim Biotech product and material limitations, especially regarding pH and temperature stability of the membranes and cassettes, must be adhered to.

The Sartorius INCREASE® Program offers evaluation of product specific cleaning procedures.

Detection Method of NaOH cleaning agent clearance:

- ▶ Neutral pH effluent in both retentate and permeate.
- ▶ Conductivity.

Example of a basic cleaning protocol for Hydrosart® Microfilters:

Step 1:

Flush the filter cassettes with isotonic saline for 5 min. at room temperature.

Conditions:

p_{in} = 36.3 psi|2.5 bar|250 kPa

p_{ret} = 0 psi|0 bar|0 kPa

p_{per} = closed; open after 4 min.

Step 2:

Recirculate 1 N NaOH at room temperature for 10 min; then flush the filter cassettes with the same batch of sodium hydroxide.

Conditions:

p_{in} = 29 psi|2.0 bar|200 kPa

p_{ret} = 0 psi|0 bar|0 kPa

p_{per} = closed; open after 4 min.

Step 3:

Flush with purified water until the pH is neutral.

Conditions:

p_{in} = 29 psi|2.0 bar|200 kPa

p_{ret} = 14.5 psi|1.0 bar|100 kPa

p_{per} = 7.3 psi|0.5 bar|50 kPa

2.9 Cassette Disinfection

Always disinfect clean filter cassettes only. Formaldehyde can denature proteins. As a result the protein might be irreversibly fixed to the membrane, which will result in flux decay.

Disinfectant	Concentration	Time	Temperature
Polyethersulfone and Hydrosart® Formaldehyde	2–3% by weight	30 min	20–30°C
NaOH	1 N	30 min	40°C

Adjust the following pressures:

p_{in} = 29 psi|2.0 bar|200 kPa

p_{ret} = 7.3 psi|0.5 bar|50 kPa

p_{per} = 0 psi|0 bar|0 kPa

Recirculate the disinfectant through the cassette(s).

The disinfectant must flow through both permeate outlets.

Afterwards, flush the filter cassette(s) with purified water as described in "Flushing the Filter Cassettes".

2.10 Cassette Storage

Sartorius Stedim Biotech filter cassettes can be stored after use. To avoid microbial contamination preservatives have to be used.

All membranes will lose their filtration performance if allowed to dry out or freeze. Always store filter cassettes in moist condition.

2.10.1 Short Term Storage

We recommend the following preservatives for short term storage (not exceeding 4 weeks):

Preservative	Concentration
Formaldehyde	2–3% by weight
Ethanol Do not use denatured alcohol	20% by weight*
NaOH	0.1 N

* Hydrosart® can be stored in Etanol 20–40% by weight

For short term storage proceed as follows:

1. Clean the filter cassettes (see "Cleaning the Filter Cassettes").
2. Recirculate the preservative for 5 minutes through the system at the same flow rate used to flush the filter cassettes. **The preservative must flow through both permeate outlets.**
3. Close feed inlet, retentate outlet and both permeate outlets.
4. Do not allow the cassettes to dry out or freeze.

2.10.2 Long Term Storage

We recommend the following preservatives for long term storage:

Preservative	Concentration
Formaldehyde	2–3% by weight
Ethanol Do not use denatured alcohol	20% by weight*
NaOH	0.1 N PESU and Hydrosart up to 12 months

* Hydrosart® can be stored in 20–40% by weight Ethanol

For long term storage proceed as follows:

1. Clean the filter cassettes (see "**Cleaning the Filter Cassettes**").
2. Recirculate the preservative for 5 minutes through the system at the same flow rate used to flush the filter cassettes. **The preservative must flow through both permeate outlets.**
3. Remove the filter cassette from the holding device.
4. Store the filter cassettes in a closed container or bag as follows:
 - ▶ Plastic bag:
Seal the filter cassette together with 50 ml of preservative into the bag.
 - ▶ Container:
Pour preservative over the filter cassette until it is completely immersed, then close the container.

After Storage in NaOH, the storage solution has to be washed from the outer surface with purified water before installation into a holding device.

2.11 Cassette Return Instructions

If cassettes are sent back to Sartorius Stedim for evaluation, please proceed as follows:

1. Inform your local Sartorius Stedim Application Support.
2. Clean the cassettes completely (see "**Cleaning the Filter Cassettes**").
3. Disinfect the cassettes (see current product return form).
4. Include the following information:
 - a. Specify product or media filtered
 - b. Describe failure observed
 - c. Reveal operating parameters such as
 - operating pressure
 - time in use
 - process cycles
 - cleaning conditions
 - process temperature
 - d. Number of cleaning cycles
 - e. Description of disinfectant.

2.12 Disposal (Status: May 2009)

If non-hazardous residual impurities are present, the EWC waste code 150203 (European Waste Catalogue) for absorbents and filter materials, wiping cloths and protective clothing can be applied. If the filters are contaminated with hazardous materials, the EWC waste code 150202* for absorbents and filter materials contaminated with hazardous materials should be used. If you are uncertain, contact your waste disposal authority.

3.0 Appendix

3.1 Specifications Sartocoon® Slice 200

Sartocoon® Slice 200 Ultrafiltration Cassettes

Membrane material	Nominal molecular weight cutoff [NMWCO]	Order numbers	Effective filtration area [m ²]	pH stability	Thermal sterilization	Maximum p _{feed} [bar] at 20°C	Maximum temp. [°C] for continuous operation	Air diffusion value [ml/min] at p _{feed}	Integral Frame	Preservation ¹⁾
Polyether-sulfone (PESU)	1 kD	3081460902E--SW	0.02	1-14	No	4	50	≤ 1.0	PVDF	Glycerin
	5 kD	3081462902E--SW	0.02	1-14	No	4	50	≤ 1.0	PVDF	Glycerin
	8 kD	3081463402E--SW	0.02	1-14	No	4	50	≤ 1.0	PVDF	Ethanol
	10 kD	3081463902E--SW	0.02	1-14	No	4	50	≤ 1.0	PVDF	Ethanol
	Albumin	308146402K--SW	0,02	1-14	No	4	50	≤ 2	PVDF	Ethanol
	30 kD	3081465902E--SW	0.02	1-14	No	4	50	≤ 1.0	PVDF	Ethanol
	50 kD	3081465002E--SW	0.02	1-14	No	4	50	≤ 1.0	PVDF	Ethanol
	100 kD	3081466802E--SW	0.02	1-14	No	4	50	≤ 1.0	PVDF	Glycerin
300 kD	3081467902E--SW	0.02	1-14	No	4	50	≤ 1.0	PVDF	Glycerin	
Hydrosart®	2 kD	3081441902E--SW	0.02	2-14	No	4	50	≤ 1.0	PVDF	Ethanol
	5 kD	3081442902E--SW	0.02	2-14	No	4	50	≤ 1.0	PVDF	Ethanol
	10 kD	3081443902E--SW	0.02	2-14	No	4	50	≤ 1.0	PVDF	Ethanol
	30 kD	3081445902E--SW	0.02	2-14	No	4	50	≤ 1.0	PVDF	Ethanol
	100 kD	3081446802E--SW	0.02	2-14	No	4	50	≤ 1.0	PVDF	Ethanol

¹⁾ The preservation liquid is either Ethanol 20–24% by weight or Glycerin.

3.2 Specifications Sartocoon® Slice

Sartocoon® Slice Ultrafiltration Cassettes

Membrane material	Nominal molecular weight cutoff [NMWCO]	Order numbers	Effective filtration area [m ²]	pH stability	Thermal sterilization	Maximum p _{feed} [bar] at 20°C	Maximum temp. [°C] for continuous operation	Air diffusion value [ml/min] at p _{feed} = 14,5 psi 1 bar	Integral Frame	Preservation ¹⁾
Polyether-sulfone (PESU)	1 kD	3051460901E--SW	0.1	1-14	No	4	50	≤ 15	PVDF	Glycerin
	5 kD	3051462901E--SW	0.1	1-14	No	4	50	≤ 15	PVDF	Glycerin
	8 kD	3051463401E--SW	0.1	1-14	No	4	50	≤ 15	PVDF	Ethanol
	10 kD	3051463901E--SW	0.1	1-14	No	4	50	≤ 15	PVDF	Ethanol
	Albumin	305146401K--SW	0.1	1-14	No	4	50	≤ 5	PVDF	Ethanol
	30 kD	3051465901E--SW	0.1	1-14	Yes	4	50	≤ 15	PVDF	Ethanol
	50 kD	3051465001E--SW	0.1	1-14	No	4	50	≤ 15	PVDF	Ethanol
	100 kD	3051466801E--SW	0.1	1-14	Yes	4	50	≤ 15	PVDF	Glycerin
300 kD	3051467901E--SW	0.1	1-14	Yes	4	50	≤ 15	PVDF	Glycerin	
Hydrosart®	2 kD	3051441901E--SW	0.1	2-14	Yes	4	50	≤ 5	PVDF	Ethanol
	5 kD	3051442901E--SW	0.1	2-14	Yes	4	50	≤ 5	PVDF	Ethanol
	10 kD	3051443901E--SW	0.1	2-14	Yes	4	50	≤ 5	PVDF	Ethanol
	30 kD	3051445901E--SW	0.1	2-14	Yes	4	50	≤ 5	PVDF	Ethanol
	100 kD	3051446801E--SW	0.1	2-14	Yes	4	50	≤ 5	PVDF	Ethanol

¹⁾ The preservation liquid is either Ethanol 20–24% by weight or Glycerin.

Sartocoon® Slice Microfiltration Cassettes

Membrane material	Pore size [µm]	Order numbers	Effective filtration area [m ²]	pH stability	Thermal sterilization*	Maximum p _{feed} [bar] at 20°C	Maximum temp. [°C] for continuous operation	Air diffusion value [ml/min] at p _{feed} = 14,5 psi 1 bar	Integral Frame	Preservation ¹⁾
Polyethersulfone (PESU)	0.1	3051545801W--SW	0.1	1-14	Yes	4	50	≤ 5	PVDF	Ethanol
Hydrosart®	0.2	3051860701W--SW	0.1	2-14	Yes	4	50	≤ 15	PVDF	Ethanol
	0.45	3051860601W--SW	0.1	2-14	Yes	4	50	≤ 15	PVDF	Ethanol

¹⁾ The preservation liquid is either Ethanol 20–24% by weight or Glycerin.

3.3 Specifications Sartococon®

Sartococon® Ultrafiltration Cassettes

Membrane material	Nominal molecular weight cutoff [NMWCO]	Order numbers	Effective filtration area [m ²]	pH stability	Thermal sterilization	Maximum p _{feed} [bar] at 20°C	Maximum temp. [°C] for continuous operation	Air diffusion value [ml/min] at p _{feed} = 14,5 psi 1 bar	Integral Frame	Preservation ¹⁾
Polyethersulfone (PESU)	1 kD	3021460907E-SW	0.7	1-14	No	4	50	≤ 50	PVDF	Glycerin
	5 kD	3021462907E-SW	0.7	1-14	No	4	50	≤ 50	PVDF	Glycerin
	8 kD	3021463407E-SW	0.7	1-14	No	4	50	≤ 50	PVDF	Ethanol
	10 kD	3021463907E-SW	0.7	1-14	No	4	50	≤ 50	PVDF	Ethanol
	Albumin	302146AL07K-SW	0.7	1-14	No	4	50	≤ 20	PVDF	Ethanol
	30 kD	3021465907E-SW	0.7	1-14	Yes	4	50	≤ 50	PVDF	Ethanol
	50 kD	3021465007E-SW	0.7	1-14	No	4	50	≤ 50	PVDF	Ethanol
	100 kD	3021466807E-SW	0.7	1-14	Yes	4	50	≤ 50	PVDF	Glycerin
	300 kD	3021467907E-SW	0.7	1-14	Yes	4	50	≤ 50	PVDF	Glycerin
	Hydrosart ²	2 kD	3021441906E-SW	0.6	2-14	Yes	4	50	≤ 15	PVDF
5 kD		3021442906E-SW	0.6	2-14	Yes	4	50	≤ 15	PVDF	Ethanol
10 kD		3021443906E-SW	0.6	2-14	Yes	4	50	≤ 15	PVDF	Ethanol
30 kD		3021445906E-SW	0.6	2-14	Yes	4	50	≤ 15	PVDF	Ethanol
100 kD		3021446806E-SW	0.6	2-14	Yes	4	50	≤ 15	PVDF	Ethanol

¹⁾ The preservation liquid is either Ethanol 20–24% by weight or Glycerin.

Sartococon® Microfiltration Cassettes

Membrane material	Pore size [µm]	Order numbers	Effective filtration area [m ²]	pH stability	Thermal sterilization*	Maximum p _{feed} [bar] at 20°C	Maximum temp. [°C] for continuous operation	Air diffusion value [ml/min] at p _{feed} = 14,5 psi 1 bar	Integral Frame	Preservation ¹⁾
Polyethersulfone (PESU)	0.1	3021545806W-SW	0.6	1-14	Yes	4	50	≤ 15	PVDF	Ethanol
Hydrosart ²	0.2	3021860706W-SW	0.6	2-14	Yes	4	50	≤ 50	PVDF	Ethanol
	0.45	3021860606W-SW	0.6	2-14	Yes	4	50	≤ 50	PVDF	Ethanol

¹⁾ The preservation liquid is either Ethanol 20–24% by weight or Glycerin.

3.4 Specifications Sartocube®

Sartocube® Ultrafiltration Cassettes

Membrane material	Nominal molecular weight cutoff [NMWCO]	Order numbers	Effective filtration area [m ²]	pH stability	Thermal sterilization	Maximum p _{feed} [bar] at 20°C	Maximum temp. [°C] for continuous operation	Air diffusion value [ml/min] at p _{feed} = 14,5 psi 1 bar	Integral Frame	Preservation
Polyethersulfone (PESU)	1 kD	3021460935E-BSW	3.5	1-14	1-14	4	50	≤ 20	PVDF	Glycerin
	5 kD	3021462935E-BSW	3.5	1-14	1-14	4	50	≤ 20	PVDF	Glycerin
	8 kD	3021463435E-BSW	3.5	1-14	1-14	4	50	≤ 20	PVDF	Ethanol
	Albumin	302146AL35K-BSW	3.5	1-14	1-14	4	50	≤ 20	PVDF	Ethanol
	10 kD	3021463935E-BSW	3.5	1-14	No	4	50	≤ 50	PVDF	Ethanol
	30 kD	3021465935E-BSW	3.5	1-14	Yes	4	50	≤ 50	PVDF	Ethanol
	50 kD	3021465035E-BSW	3.5	1-14	No	4	50	≤ 50	PVDF	Ethanol
	100 kD	3021466835E-BSW	3.5	1-14	Yes	4	50	≤ 50	PVDF	Glycerin
	300 kD	3021467935E-BSW	3.5	1-14	Yes	4	50	≤ 50	PVDF	Glycerin
	Hydrosart ²	2 kD	3021441930E-BSW	3.0	2-14	Yes	4	50	≤ 50	PVDF
5 kD		3021442930E-BSW	3.0	2-14	Yes	4	50	≤ 50	PVDF	Ethanol
10 kD		3021443930E-BSW	3.0	2-14	Yes	4	50	≤ 50	PVDF	Ethanol
30 kD		3021445930E-BSW	3.0	2-14	Yes	4	50	≤ 50	PVDF	Ethanol
100 kD		3021446830E-BSW	3.0	2-14	Yes	4	50	≤ 50	PVDF	Ethanol

¹⁾ The preservation liquid is either Ethanol 20–24% by weight or Glycerin.

Sartocube® Microfiltration Cassettes

Membrane material	Pore size [µm]	Order numbers	Effective filtration area [m ²]	pH stability	Thermal sterilization	Maximum p _{feed} [bar] at 20°C	Maximum temp. [°C] for continuous operation	Air diffusion value [ml/min] at p _{feed} = 14,5 psi 1 bar	Integral Frame	Preservation
Polyethersulfone (PESU)	0.1	3021545830W-BSW	3.0	1-14	Yes	4	50	≤ 50	PVDF	Ethanol
Hydrosart ²	0.2	3021860730W-BSW	3.0	2-14	Yes	4	50	≤ 100	PVDF	Ethanol
	0.45	3021860630W-BSW	3.0	2-14	Yes	4	50	≤ 100	PVDF	Ethanol

¹⁾ The preservation liquid is either Ethanol 20–24% by weight or Glycerin.

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DOMINIQUE DUTSCHER SAS

Sales and Service Contacts

For further contacts, visit www.sartorius-stedim.com

Europe

Germany
Sartorius Stedim Biotech GmbH
August-Spindler-Strasse 11
37079 Goettingen
Phone +49.551.308.0
Fax +49.551.308.3289

Sartorius Stedim Systems GmbH
Robert-Bosch-Strasse 5-7
34302 Guxhagen
Phone +49.5665.407.0
Fax +49.5665.407.2200

France
Sartorius Stedim FMT S.A.S.
ZI Les Paluds
Avenue de Jouques - CS 91051
13781 Aubagne Cedex
Phone +33.442.845600
Fax +33.442.845619

Sartorius Stedim France SAS
ZI Les Paluds
Avenue de Jouques - CS 71058
13781 Aubagne Cedex
Phone +33.442.845600
Fax +33.442.846545

Austria
Sartorius Stedim Austria GmbH
Franzosengraben 12
1030 Vienna
Phone +43.1.7965763.18
Fax +43.1.796576344

Belgium
Sartorius Stedim Belgium N.V.
Leuvensesteenweg, 248/B
1800 Vilvoorde
Phone +32.2.756.06.80
Fax +32.2.756.06.81

Hungary
Sartorius Stedim Hungária Kft.
Kagyló u. 5
2092 Budakeszi
Phone +36.23.457.227
Fax +36.23.457.147

Italy
Sartorius Stedim Italy S.p.A.
Via dell'Antella, 76/A
50012 Antella-Bagno a Ripoli (FI)
Phone +39.055.63.40.41
Fax +39.055.63.40.526

Netherlands
Sartorius Stedim Netherlands B.V.
Phone +31.30.60.25.080
Fax +31.30.60.25.099
filtratie.nederland@sartorius-stedim.com

Poland
Sartorius Stedim Poland Sp. z o.o.
ul. Wrzesinska 70
62-025 Kostrzyn
Phone +48.61.647.38.40
Fax +48.61.879.25.04

Russian Federation
LLC "Sartorius ICR"
Uralskaya str. 4, Lit. B
199155, Saint-Petersburg
Phone +7.812.327.5.327
Fax +7.812.327.5.323

Spain
Sartorius Stedim Spain SA
C/Isabel Colbrand 10,
Oficina 70
Poligono Industrial de Fuencarral
28050 Madrid
Phone +34.90.2110935
Fax +34.91.3589623

Switzerland
Sartorius Stedim Switzerland AG
Ringstrasse 24 a
8317 Tagelswangen
Phone +41.52.354.36.36
Fax +41.52.354.36.46

U.K.
Sartorius Stedim UK Ltd.
Longmead Business Centre
Blenheim Road, Epsom
Surrey KT19 9 QQ
Phone +44.1372.737159
Fax +44.1372.726171

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Asia | Pacific

Australia
Sartorius Stedim Australia Pty. Ltd.
Unit 5, 7-11 Rodeo Drive
Dandenong South Vic 3175
Phone +61.3.8762.1800
Fax +61.3.8762.1828

China
Sartorius Stedim Biotech (Beijing) Co. Ltd.
No. 33 Yu'an Road
Airport Industrial Park Zone B
Shunyi District, Beijing 101300
Phone +86.10.80426516
Fax +86.10.80426580

Sartorius Stedim Biotech (Beijing) Co. Ltd.
Shanghai Branch Office
3rd Floor, North Wing, Tower 1
No. 4560 Jin Ke Road
Pudong District, Shanghai 201210
Phone +86.21.68782300
Fax +86.21.68782332 | 68782882

Sartorius Stedim Biotech (Beijing) Co. Ltd.
Guangzhou Representative Office
Unit K, Building 23
Huihua Commerce Et Trade Building
No. 80 Xianlie Middle Road
Guangzhou 510070
Phone +86.20.37618687 | 37618651
Fax +86.20.37619051

India
Sartorius Stedim India Pvt. Ltd.
#69/2-69/3, NH 48, Jakkasandra
Nelamangala Tq
562 123 Bangalore, India
Phone +91.80.4350.5250
Fax +91.80.4350.5253

Japan
Sartorius Stedim Japan K.K.
4th Fl., Daiwa Shinagawa North Bldg.
8-11, Kita-Shinagawa 1-chome
Shinagawa-ku, Tokyo, 140-0001 Japan
Phone +81.3.4331.4300
Fax +81.3.4331.4301

Malaysia
Sartorius Stedim Malaysia Sdn. Bhd.
Lot L3-E-3B, Enterprise 4
Technology Park Malaysia
Bukit Jalil
57000 Kuala Lumpur, Malaysia
Phone +60.3.8996.0622
Fax +60.3.8996.0755

Singapore
Sartorius Stedim Singapore Pte. Ltd.
1 Science Park Road,
The Capricorn, #05-08A,
Singapore Science Park II
Singapore 117528
Phone +65.6872.3966
Fax +65.6872.2494

South Korea
Sartorius Korea Biotech Co., Ltd.
8th Floor, Solid Space B/D,
PanGyoYeok-Ro 220, Bundang-Gu
SeongNam-Si, GyeongGi-Do, 463-400
Phone +82.31.622.5700
Fax +82.31.622.5799

◀ www.sartorius-stedim.com