

Nested Liquid Handling (LiHa)

Disposable Tips.


BOOST EFFICIENCY, FLEXIBILITY AND WALKAWAY TIME WITH TECAN'S NESTED DISPOSABLE TIPS

The demand for higher throughput triggers many laboratories to look for methods to increase the capacity of their automated systems. Tecan has now addressed this need by creating Nested LiHa Disposable Tips, allowing five trays of tips to be stacked on a single SLAS carrier. This solution allows optimized use of deck space through increasing tip capacity per carrier position and minimizes the need for manual intervention for Freedom EVO® and Fluent® liquid handling workstations. Flexibility of the system is enhanced as the 350 µL tip volume offers greater versatility with the same pipetting precision as the 200 µL LiHa Disposable Tip, whereas the 10 µL Tip allows pipetting in in PCR and 384 well plates.

Main features and benefits

- Five-fold increase in tip capacity on the workdeck – longer walkaway times
- Increased tip volume range & labware compatibility – higher flexibility without compromising on performance
- Quick set-up and selection in the software* – reduced installation time

*From Freedom EVOware® V2.6 and FluentControl® V1.2 onwards

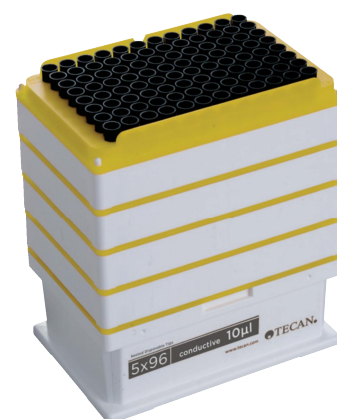

Increased workstation efficiency

The nested design allows tips to be stored in stacks of five racks of 96 tips (480 tips) per SLAS position on the deck of the workstation. This increases tip capacity and/or allows integration of additional analytical modules on the instrument.

Increased tip volume range and improved labware compatibility

A broad range of 3 to 350 μl – with CVs of $\leq 5\%$ and $\leq 0.5\%$ respectively* – makes the new 350 μl nested tip the disposable tip of choice for a wide variety of applications. The 10 μl nested tip is designed in such way, that you can easily reach the bottom of a PCR or 384 well plate.

*Using the Air LiHa Arm with water or serum; methods upon request.

Straightforward implementation

The tip and rack handling for the nested tips is conveniently fully automated and easily implemented through the Freedom EVOware® and FluentControl® software. Empty tip racks can be removed and disposed of using a Robotic Manipulator (RoMa) Arm, Robotic Gripper Arm (RGA) or MCA Gripper option.

Systems currently set up for non-SLAS disposable tips can be easily converted by adding a specific carrier. For further space-saving, carriers are also available with an integrated waste slide for tip and rack disposal.

Quality and purity

Tecan's high quality Nested LiHa Disposable Tips are designed to optimize pipetting performance and instrument reliability. Two purity levels are available to match the requirements of your applications:

Nested LiHa disposable product portfolio

Tecan Pure

Certified to be free from human DNA, DNase, RNase and PCR inhibitors.



Tecan Sterile

Certified to be sterile and free from endotoxins, pyrogens, human DNA, DNase and RNase.



350 μl , Non-filtered

300 83 400

10 μl , Non-filtered

301 04 977*

10 μl , filtered

301 04 978*

350 μl , Non-filtered

300 83 401

10 μl , Non-filtered

301 04 979*

*Coming soon. Contact your local Tecan office for availability.

www.tecan.com



Australia +61 3 9647 4100 **Austria** +43 62 46 89 33 **Belgium** +32 15 42 13 19 **China** +86 21 220 63 206 **Denmark** +45 70 23 44 50 **France** +33 4 72 76 04 80
Germany +49 79 51 94 170 **Italy** +39 02 92 44 790 **Japan** +81 44 556 73 11 **Netherlands** +31 18 34 48 17 4 **Singapore** +65 644 41 886 **Spain** +34 93 490 01 74
Sweden +46 8 750 39 40 **Switzerland** +41 44 922 89 22 **UK** +44 118 9300 300 **USA** +1 919 361 5200 **Other countries** +43 62 46 89 33

© 2016 Tecan trading AG, Switzerland, all rights reserved. For disclaimer and trademarks please visit www.tecan.com

