Corning® Stripettor® Pro Pipet Controller

Quick Start Guide


Catalog Numbers: 4900, 4999

CORNING


To access the Instruction Manual (CLS-AN-817DOC) scan the QR code or visit www.corning.com/resources.

DOMINIQUE DUTSCHER

Battery Charge

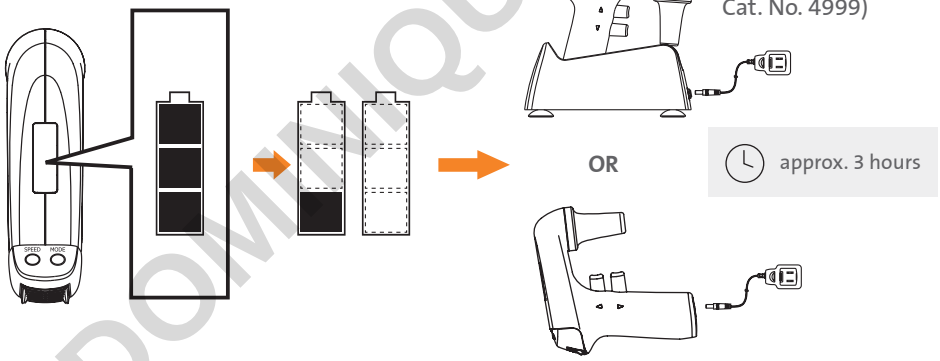

Speed Setting

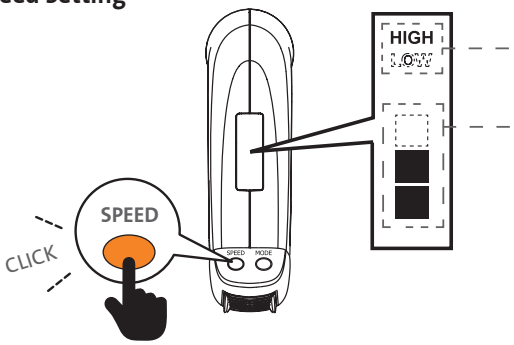


Level (LOW-HIGH)
Click SPEED to change

Speed value (1-3 bars)
· Hold SPEED for 3 sec.
· Press SPEED to change
· Hold SPEED for 3 sec. to accept


The speed can be controlled by pressing the trigger buttons as well.

Mode Control

Click MODE to change between Gravity (drop) and Blow out (arrow) mode

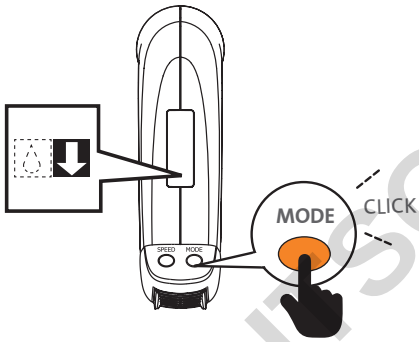


Aspirating and Dispensing

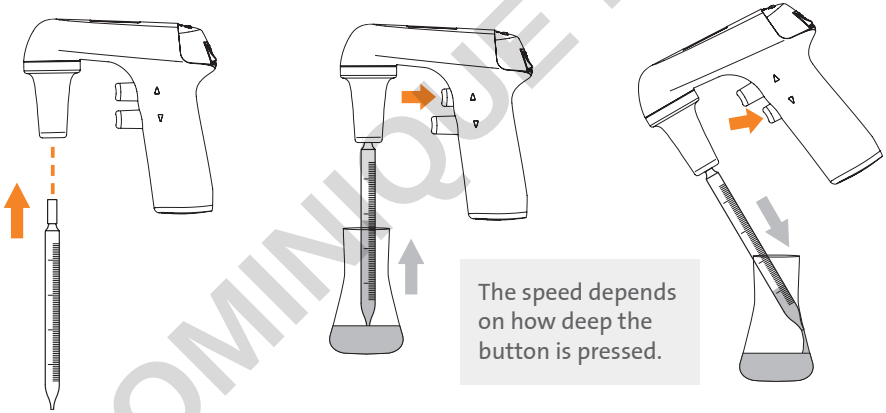


The speed depends on how deep the button is pressed.

Maintenance

Isopropyl alcohol

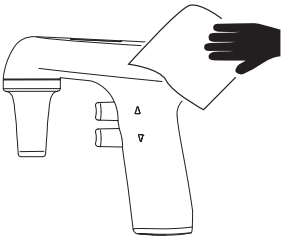

UV sterilization

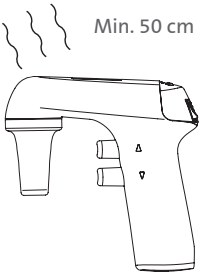

Autoclave

ONLY




121°C for 15 min.




121°C for 20 min.

DOMINIQUE DUTSCHER

Filter Replacement

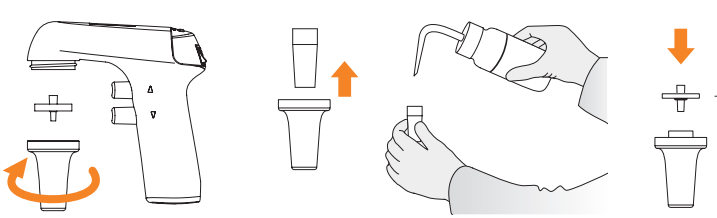

Storage

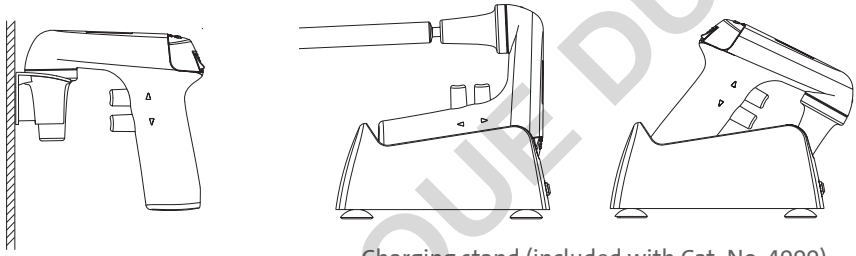

Charging stand (included with Cat. No. 4999)

Safety Information

WARNING: Risk of injury

CAUTION Risk of damage to the device or errors in pipetting of liquids.

Before starting the work with the pipet controller every user should read the operating instructions carefully.

CAUTION:

- Using the device inconsistently with the operating instructions may result in damaging the device.
- The device should be serviced only at an authorized service center, otherwise the manufacturer will be relieved from any liability under the warranty.
- Only original spare parts and accessories, recommended by the manufacturer, shall be used.
- Only the original charger, supplied by the manufacturer, shall be used for charging the battery.
- In case of incorrect functioning of the pipet controller, work shall be stopped. The device shall be cleaned according to the Maintenance Section in the Instruction Manual (CLS-AN-817DOC) and sent for repair to an authorized service center.
- In the case of mechanical damage to the casing, the device shall be immediately sent for repair to an authorized service center.
- The use of excessive force during work shall be avoided.

DOMINIQUE DUTSCHER

WARNING:

- During the work with the pipet controller general safety regulations regarding risks related with laboratory work should be observed. Protective clothing, goggles, and gloves should be worn.
- The pipet controller shall be used only for measuring liquids in conditions specified by the manufacturer, which are limited due to the chemical and mechanical resistance of the device, as well as the user safety.
- The information and instructions provided by the manufacturers of the reagents must be observed.
- Read the Instruction Manual before use. Improper handling of a Li-ion battery may result in the battery catching fire, overheating, or leaking.
- Do not puncture, crush, bend or modify the battery in any way.
- Never attempt to open or modify the battery. This can cause a fire or explosion.
- Do not expose the battery to excessive heat, direct sunlight, or fire.
- Stop charging the battery if charging is not completed within the specified time.
- If there are signs of damage, swelling or leakage, stop use and contact the manufacturer immediately.
- Do not short circuit the (+) and (-) terminals with other metals.

NOTE: The pipet controller is equipped with liquid vapors exhaustion system which protects against corrosion to ensure long instrument life.

Equipment Disposal



According to Directive 2012/19/EU of the European Parliament and of the Council of 4 July 2012 on waste electrical and electronic equipment (WEEE), this product is marked with the crossed-out wheeled bin and must not be disposed of with domestic waste.

Consequently, the buyer shall follow the instructions for reuse and recycling of waste electronic and electrical equipment (WEEE) provided with the products and available at **www.corning.com/weee**.

How to Purchase: For specific availability in your region and purchasing options, terms and conditions of sale, customer/product support, and certificates, visit **www.corning.com/how-to-buy**.

Warranty/Disclaimer: Unless otherwise specified, all products are for research use or general laboratory use only. Not intended for use in diagnostic or therapeutic procedures. Not for use in humans. For a listing of US medical devices, regulatory classifications or specific information on claims, visit www.corning.com/resources.

CORNING

For additional product or technical information, visit **www.corning.com/lifesciences** or call 800.492.1110. Outside the United States, call +1.978.442.2200.

The information contained within is accurate as of the date of publication and subject to change without notice.
For a listing of trademarks, visit www.corning.com/trademarks. All other trademarks are the property of their respective owners.
© 2024-2025 Corning Incorporated. All rights reserved. 9/25 CLS-AN-818DOC REV2

DOMINIQUE DUTSCHER