Lively 300V/ 500V/ 3AP Power Supply

Instruction manual

Catalog No. MP-310
MP-320
MP-510

DOMINIQUE DUTSCHER SAS

www.majorsci.com
service@majorsci.com

Version 08A
Revised on: 2023/02/08

Packing List

MP-310

- 1x MP-310 Power Supply
- 1x Power Cord
- 1x Instruction Manual

MP-320

- 1x MP-320 Power Supply
- 1x Power Cord
- 1x Instruction Manual

MP-510

- 1x MP-510 Power Supply
- 1x Power Cord
- 1x Instruction Manual

Signed by:

Date:

Major Science is liable for all missing or damaged parts / accessories within 7 days after customer received this instrument package. Please contact Major Science immediately regarding this issue. If no response within such time period from consignee party, that will be consignee party's whole responsibility.

DOMINIQUE DUTSCHER SAS

Table of Contents

PACKING LIST **2**

WARNING **4**

SECTION 1 INTRODUCTION **8**

2.1 OVERVIEW 8
2.2 PRODUCT DESCRIPTION & FEATURE 9

SECTION 2 TECHNICAL SPECIFICATION **10**

SECTION 3 INSTALLATION INSTRUCTIONS **11**

SECTION 4 OPERATION INSTRUCTIONS **11**

4.1 CONTROL INTERFACE 11
4.2 GETTING START 12
4.3 GENERAL OPERATING INSTRUCTIONS 13
4.4 CONSTANT MODE 14
4.5 PROGRAM MODE 17
4.6 TYPICAL RUNNING CONDITION 20

SECTION 5 TROUBLE SHOOTING AND MAINTENANCE **21**

REPLACING THE FUSE 24
MAINTENANCE 24

SECTION 6 ORDERING INFORMATION **25**

SECTION 7 WARRANTY **26**

DOMINIQUE DUTSCHER SAS

Warning

Major Science Lively Power Supply has been tested and found to comply with the limits for the CE regulation. Also, it is RoHS compliant to deliver confident product which meets the environmental directive. These limits are designed to provide reasonable protection against harmful interference when the instrument series is operated in a commercial environment. This instrument series used together with power supply unit generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this instrument series in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at their expense. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment. It is strongly recommended for the user to read the following points carefully before operating this equipment.

1. Read and follow the manual instructions carefully.
2. Do not alter the equipment. Failure to follow these directions could result in personal and/or laboratory hazards, as well as invalidate equipment warranty.
3. Use a properly grounded electrical outlet with correct voltage and current handling capacity.
4. Disconnect from power supply before maintenance and servicing. Refer servicing to qualified personnel.
5. Never use this instrument series without having the safety cover correctly in position.
6. Do not use the unit if there is any sign of damage to the external tank or cover. Replace damaged parts.
7. Do not use in the presence of flammable or combustible material; fire or explosion may result. This device contains components which may ignite such materials.
8. Refer maintenance and servicing to qualified personnel.
9. Ensure that the system is connected to electrical service according to local and national electrical codes. Failure to make a proper connection may create fire or shock hazard.
10. Use appropriate materials and operate correctly to avoid possible hazards

DOMINIQUE DUTSCHER SAS

of explosion, implosion or release of toxic or flammable gases arising from overheated materials.

11. The unit shall be operated only by qualified personnel.

Safety Information

Use high level of precaution against any electrical device. Before connecting the electrical supply, check to see if the supply voltage is within the range stated at the rating label, and see to it that the device be seated firmly. Place the unit in a safe and dry location; it must NOT touch the surrounding. Follow the safety precautions for chemicals / dangerous materials. If needed, please contact qualified service representative or service@majorsci.com

Environmental Conditions

Ensure the instrument is installed and operated strictly under the following conditions:

1. Indoor use only
2. ≤95% RH
3. 75 kPa – 106 kPa
4. Altitude must not exceed 2000 meters
5. 4℃~ 40℃ operating temperature
6. Pollution degree: 2
7. Mains supply voltage fluctuations up to ±10% of the normal voltage

Avoiding Electrical Shock

Follow the guidelines below to ensure safe operation of the unit.

The Lively Power Supply has been designed to utilize shielded wires thus minimizing any potential shock hazard to the user. Major Science recommends against the use of unshielded wires.

To avoid electrical shock:

1. In the event of solution spilling on the instrument, it must be dried out for at least 2 hours and restored to NORMAL CONDITION before each operation.
2. Never connect or disconnect wires loading from the power jacks when the red indicator light of power switch is on.
3. WAIT at least 5 seconds after stopping a run before handling output leads or any connected apparatus.
4. ALWAYS make sure that your hands, work area, and instruments are **clean** and **dry** before making any connections or operating the power supply.
5. ONLY connect the power cord to a properly grounded AC outlet.

DOMINIQUE DUTSCHER SAS

Avoiding Damage to the Instrument

1. Do not attempt to operate the device if damage is suspected.
2. Protect this unit from physical damage, corrosive agents and extreme temperatures (direct sunlight, etc.).
3. For proper ventilation and safety concerns, keep at least 10 cm of space behind the instrument, and at least 5 cm of space on each side.
4. Use high level of precaution against the damages on the unit.
5. Do not operate the unit out of environmental conditions addressed above.
6. Do not operate the power supplies in high humidity environments (> 95%), or where condensation may occur.
7. To avoid condensation after operating the power supply in a cold room, wrap the unit in a plastic bag and allow at least 2 hours for the unit to equilibrate to room temperature before removing the bag and operating the unit.
8. Prior to applying any cleaning or decontamination methods other than manufacturer's recommendation, users should check with the manufacturer's instruction to see if the proposed method will damage the equipment.

Equipment Operation

Follow the guidelines below to ensure safe operation of the unit:

1. NEVER access dangerous chemicals or other materials to prevent possible hazard of explosion and damage.
2. Do not operate the unit without lids or covers to prevent possible hazards.
3. A temporary conductivity caused by condensation might occur even though this series is rated Pollution Degree 2 in accordance with IEC 664.

Symbol

Symbols used on the power supply are explained below.



Indicates an area where a potential shock hazard may exist. Consult the manual to avoid possible personal injury or instrument damage.



Indicates disposal instruction.

DO NOT throw this unit into a municipal trash bin when this unit has reached the end of its lifetime. To ensure utmost protection of the global environment and minimize pollution, please recycle this unit.

DOMINIQUE DUTSCHER SAS



Caution/ Warning: Indicates a hazardous situation which, if not avoided, could result in death or serious injury.

MP-310	
	Max. voltage: 300 V Max. current: 700 mA Max. watt: 150 W
MP-320	
	Max. voltage: 300 V Max. current: 3000 mA Max. watt: 300 W
MP-510	
	Max. voltage: 500 V Max. current: 800mA Max. watt: 300W

Potential Risk and Preventive Measure

1. Risk assessment table

Potential Risk / Frequency	Frequent	Likely	Possible	Rare	Unlikely
Bruise			√		
Slash					√
Electrical shock			√		
Power cord plug wrong				√	

2. Preventative measures of risk

Potential Risk	Preventive measures
Bruise	Do not put the machine near the table edge.
Slash	Prevent hard impact on the case.
Electrical Shock	Make sure that your hands, work area, and devices are **clean** and **dry**
Power cord plug wrong	Observe correct adapter plug.

DOMINIQUE DUTSCHER SAS

Section 1 Introduction

2.1 Overview

Major Science Lively Power Supplies are recognized as unique power supplies equipped with very powerful specifications to cover the majority of electrophoresis applications on the market. Sufficient and accurate output voltages, four terminal pairs, and its compacted size can deliver accurate and reliable experimental results from one experiment to another. Lively Power Supplies are RoHS and CE compliant for environmental and safety concerns.



▲Front view



▲Rear view

DOMINIQUE DUTSCHER SAS

2.2 Product Description & Feature

Major Science Lively Power Supply are microprocessor controlled and designed to meet most electrophoresis needs. This manual describes the setup and operation of the Lively Power Supply including important information on safety and maintenance of the unit. The Lively Power Supply are capable of running horizontal & vertical electrophoresis, SDS-PAGE, native PAGE applications, two-dimensional electrophoresis, and electro-blotting. Furthermore, the powerful specifications plus four terminal pairs can be used for multi electrophoresis units simultaneously.

Major Science Lively Power Supply provides Constant Voltage or Constant Current or Constant Power to instruments used in electrophoresis. 4 terminal pairs and the powerful specification equipped enable the maximum capability of Lively Power Supply compared to other existing similar product on the market.

Features of Lively Power Supply:

- Compact size
- Advanced capacity:

MP-310	MP-320	MP-510
150W, 700mA, 300V	300W, 3000mA, 300V	300W, 800mA, 500V

- Microprocessor controller
- Constant voltages, constant currents and constant power
- Four terminal outlets
- LCD display
- Timer with alarm function
- Advanced safety devices
- Stackability
- Wide applications for DNA, RNA and protein electrophoresis

DOMINIQUE DUTSCHER SAS

Section 2 Technical Specification

	MP-310	MP-320	MP-510
Output Voltage / Inc	5 - 300V / 1V	5 – 300 / 1V	5 – 500 / 1V
Output Current / Inc	1 - 700mA / 1mA	10 – 3,000mA / 10mA	1 – 800mA / 10mA
Max. Watt	150W	300W	300W
Rated Voltages	100~ 240V; 47-60Hz, 200W: T2.5A/250V	100-240V: 47-60Hz,420W: T5A/250V	100-240V: 47-60Hz,400W: T5A/250V
Type of Output	1. Voltage or Current with automatic crossover 2. When target constant mode is set, system automatically adjusts the two other parameter to maximum to allow constant run (later could be changed by user)		
Program Storage	30 programmed files		
Program Multi-Step	Up to 6 steps		
Editable Program Function	1.Have typical running conditions program 2.Manual editable program		
Display	2.4” TFT		
Control	Microprocessor controller		
Safety Device	No Load detect		
	Leakage detect		
	Over temperature protection		
	Over load detection		
	Sudden load change detection(could be enabled by proper setting)		
	Shrouded plugs and sockets		
Timer	Constant: 1~9999 mins with alarm, continuous Program: 1~999 mins with alarm, continuous		
Crossover	Yes		
Stackable	Yes		
Automatic Recovery After Power Failure	Yes		
IQ/OQ Protocols	Yes, optional		

DOMINIQUE DUTSCHER SAS

Regulatory	CE, ETL	CE	CE
Operating Temperature	4℃~ 40℃		
Construction Material	Flame retardant ABS faceplate		
Unit Dimension	215 x 335 x 104 mm (W x L x H)		
Weight	approx. 2.1 kg		

Section 3 Installation Instructions

Lively Power Supply is actually a pre-installed instrument. As long as it is placed on a sturdy and level surface in a safe, dry place, and further connects with well-prepared electrophoresis system, it is ready for operation.

Section 4 Operation Instructions


4.1 Control interface



Button Functions

No.	Icon	Function
1	START STOP	Press to activate or stop the unit

DOMINIQUE DUTSCHER SAS

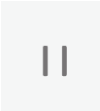2		Press to temporarily interrupt power to an operation in progress; resume power after pausing without resetting the timer
3		Press to go back previous page
4		1. Press to enter the next page 2. Press to enable numeric value set up
5		Press to move cursor left forward between parameters
6		Press to move cursor right forward between parameters
7		1. Press to move cursor up between parameters 2. Press to increase numeric values
8		1. Press to move cursor down between parameters 2. Press to decrease numeric values

4.2 Getting Start

1. Place the unit on a sturdy and level surface in a safe, dry place, away from laboratory traffic.


2. Ensure that the power switch is OFF, and then plug the three-pronged power cord one end into a grounded three-prong AC outlet with appropriate voltage (100V to 240V as indicated on the rating sticker near the AC cord on the back of the unit), and plug the other end into the main power socket.





DOMINIQUE DUTSCHER SAS

3. Connect the DC output jacks from the electrophoresis unit; insert the red lead (+) into the red output jack, and the black lead (-) into the blue output jack.

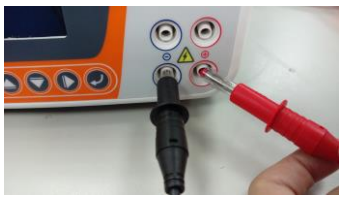

4. Power on the unit by pressing the ON/OFF switch on the back.

5. After a few seconds the screen shows the mode-selecting page.



4.3 General Operating Instructions


There are 3 modes available: **Constant Mode, Program Mode,** and **Typical Running Condition.** Please select according to your needs. Detailed instruction for each mode is explained in the following sections.

DOMINIQUE DUTSCHER SAS


4.4 Constant Mode

Use the **Constant Mode** for applications that require only one specific limit: voltage/ current/ power during the entire operation of electrophoresis.


a. Setting

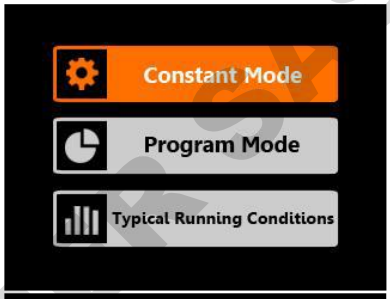

- Use Key and Key to select **Constant Mode**, and then press to enter the next page.

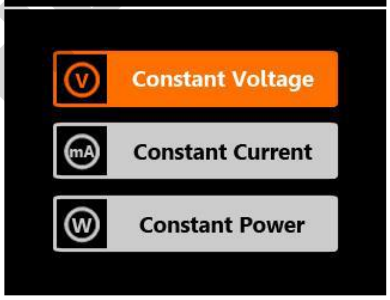

- Use and to select ether Constant Voltage(V) or Constant Current(mA) or Constant Power(W), and then press to enter the next page.

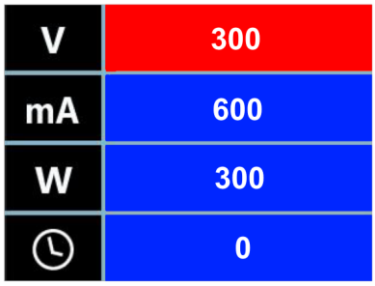

- Use and to move cursor to the parameter such as voltage (V) current (mA) or power (W) or time (min.). Press Key to set the specified parameter.

Note: If the time value is set “0,” it means “infinite.”
The power supply will constantly operate until user manually stops it.


- Use or to set the appropriate value, and then press Key, and move to the next parameter until all the parameters are set in the same operation method.

Note: Under **Constant Mode**, when target value is set, system automatically adjusts the other two parameters to maximum to allow constant run. For example, if constant voltage (V) is set, system will adjust current (mA) and power (W) to the maximum value. User could later lower the other two values by themselves. System will hold the value either at target constant value or the one which has been reached first.

DOMINIQUE DUTSCHER SAS

b. Start an Run

- Press START STOP Key to start operation. The LED is lit, and the screen will show present real-time parameter values.

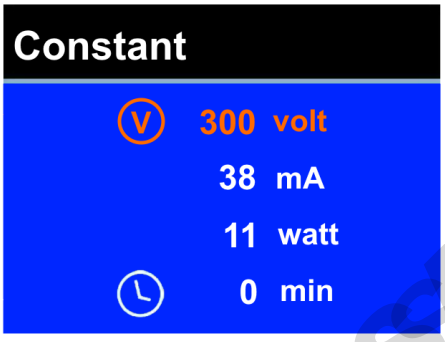

- When parameter reaches the set value, its color would turn orange.

- Press  to switch pages between real-time value display and value set up page.



- When a run is completed, operation stops with alarm and COMPLETE is displayed on the screen. Press START STOP Key to terminate, and turn the power OFF by the switch on the back.

Note: It would take a little time (about 5 seconds) for the unit to power off.

DOMINIQUE DUTSCHER SAS

c. Pause/ Terminate a Run

- Press ⏸ Key to temporarily pause the on-going electrophoresis without terminating the timer. The LED is flashing. Press START STOP Key to resume the run.
- Press the START STOP Key again to stop electrophoresis and to terminate a timed run.
- If you need to modify values of voltage (V)/ current (mA)/ power (W) during operation but do not want to terminate the time run, you may pause electrophoresis by pressing the ⏸ Key. Press  Key to enter the setting page.

  Adjust the values and then press START STOP Key to resume your operation.

- If you want to terminate and start a new run during operation, press START STOP Key instead of ⏸ Key. After adjusting the values press START STOP Key once again to start a new operation run.


Note: If you terminate and start a new run, so does the timer. The timer reset and does not take into account the time before.

DOMINIQUE DUTSCHER SAS

4.5 Program Mode

Instructions for operating Lively Power Supply in the programming operation are provided in this section. The Program Mode allows you to vary levels in voltage (V), current (mA), and power (W) during specified time periods for up to 6 Steps, depending upon your electrophoresis needs. The Lively Power Supply is capable of having 30 different program files storages for user's convenience. Follow the instruction below to set the Program Mode:


a. Setting

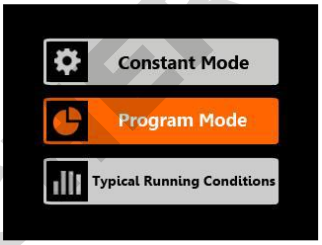


- Use Key and Key to select **Program Mode**, and then press to enter the next page.



- Press or Key to select appropriate file number, and then press Key to enter the next page.

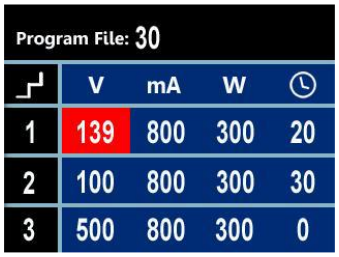


- Use , , or to move cursor to parameters such as voltage (V) current (mA) or power (W) or time (min.).
- Press Key and use and to set the appropriate value, and then press Key. Repeat this step until all the parameters are set properly. There are 6 steps available for use; use to move down for setting Step4-6, and to go back Step1-3.

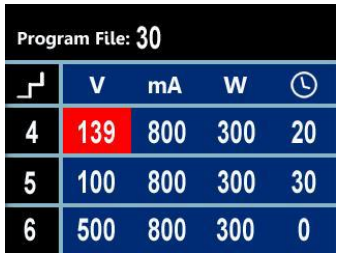

Note: If the time value is set "0" at any stage, it means "stop at previous stage." For example, if Stage 1 is set 5 minutes and Stage 2 is set 0, the system will only complete Stage 1 and stop.

DOMINIQUE DUTSCHER SAS

b. Start a Run

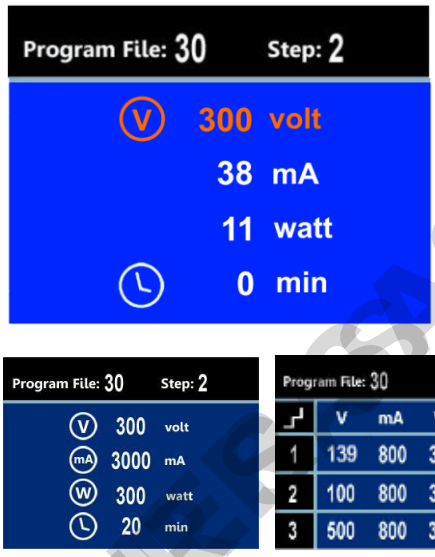

- Press START STOP Key to start electrophoresis. The LED is lit, and the screen will show present real-time parameter values.

- When parameters reach the set value, the color would turn orange.

- Press  to switch pages between real-time value display and value set up page.

- When the run is completed, operation stops with alarm and COMPLETE is displayed on the screen. Press START STOP Key to terminate a timed run, and turn the AC power OFF by the switch on the back.

Note: It would take a little time (about 5 seconds) for the unit to power off.

DOMINIQUE DUTSCHER SAS

c. Pause/ Terminate a Run

- Press ⏸ Key to temporarily interrupt power to ongoing electrophoresis without terminating the operation. The LED is flashing. Press START STOP Key to resume the run.
- Press the START STOP Key again to stop electrophoresis and to terminate a timed run.
- If you need to make changes to the values of parameters during current running but do not want to terminate the time run, you may pause electrophoresis by pressing the ⏸ Key. Press  Key to enter the setting screen.

  Adjust the values and then press START STOP Key to resume your operation.
- If you want to change values of parameters and start a new operation run during current running, press START STOP Key instead of ⏸. After adjusting the values press START STOP Key once again to start a new operation run.


Note: If you terminate and start a new run, so does the timer. The timer resets and does not take into account the time before.

DOMINIQUE DUTSCHER SAS


4.6 Typical Running Condition

Besides **Constant Mode** and **Program Mode**, Lively Power Supply provides a third option, **Typical Running Condition**. It is helpful for those beginners who are not familiar with the parameter setting and is convenient for quick start.


- Turn on the unit; use  key to select and then press  to enter the next page.

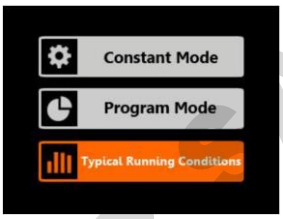

- There are 4 category can be selected. Choose the model according to your experimental need.

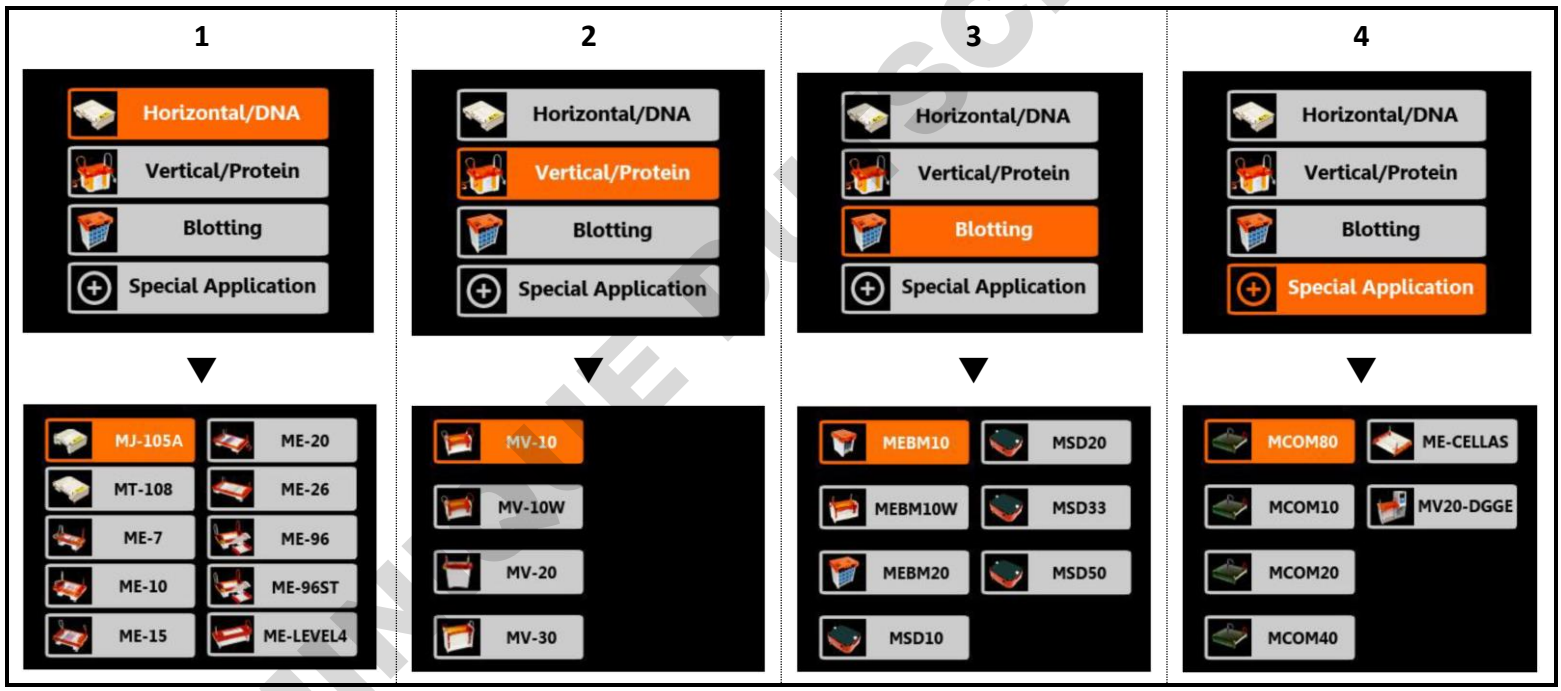

- After selecting model, the pre-set value will display. User could adjust the value or simply just press START STOP to start operation.

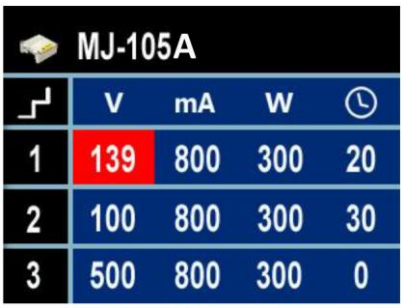

DOMINIQUE DUSCHER SAS

Section 5 Trouble shooting and Maintenance

Many operating problems may be solved by carefully reading and following the instructions in this manual accordingly. Some suggestions for troubleshooting are given below. Should these suggestions not resolve the problem, contact our SERVICE DEPARTMENT or a distributor in your region for assistance. If troubleshooting service is required, please include a full description of the problem.

Problem	Possible Causes	Solution
No Display / lights	No AC power	➢ Check if power supply is plugged, or AC power source has problem.
	AC power cord is not connected	➢ Check AC power cord connections at both ends. Use the correct cords.
	The fuse has blown	➢ Replace the fuse.
Repeated fuse broken	Hardware failure	➢ Contact Major Science service department
Operation stops and the screen displays PLC ERROR	Communication wires on circuit board have loosen or broken.	➢ It is recommended to send the machine back to local distributor or our Service Department for maintenance.
Operation stops with alarm: The screen displays NO LOAD	Electrophoresis leads are not connected to the power supply or to the electrophoresis unit(s), or there is a broken circuit in the electrophoresis cell	➢ Check the connections to the power supply and on your electrophoresis cell to make sure the connection is intact; check condition of wires in electrophoresis unit. Close the circuit by reconnecting the cables. Press **START/STOP** to restart the run.
	High resistance due to tape left on a pre-cast gel, incorrect buffer concentration, or incorrect buffer volumes in the electrophoresis cell	➢ Make sure the tape is removed from the pre-cast gel, buffers are prepared correctly, and the recommended volume of buffer is added to the electrophoresis unit.

DOMINIQUE DUTSCHER SAS

	High voltage application is set to run on a very low current	➢ DISABLE No Load alarm on the Display Screen
Operation stops with alarm. Display shows **LOAD CHANGED**	Bad connections for terminal connectors or damaged wires or damaged platinum wires	➢ Check all the connections to terminators, cables, wires, and gel tanks
Operation stops with alarm: Display shows **OVER VOLTAGE**	Circuit is interrupted	➢ Verify that the running buffer is correct. ➢ Verify the all cables are attached correctly ➢ Turn the Power switch off and on again; restart application. ➢ If you cannot restart the instrument, turn off the power, disconnect the power cord from the outlet, and contact Technical Service.
Operation stops with alarm: Display shows **OVER CURRENT**	Circuit is interrupted	➢ Verify that the running buffer is correct. ➢ Verify the all cables are attached correctly ➢ Turn the Power switch off and on again; restart application. ➢ If you cannot restart the instrument, turn off the power, disconnect the power cord from the outlet, and contact Technical Service.
Operation stops with alarm: Display shows **LEAKAGE**	Ground leak detected during run	➢ Check the electrophoresis system for improper grounding. Restart the power supply by turning the Power switch off and on.

DOMINIQUE DUTSCHER SAS

Operation stops with alarm: Display shows 	Power supply is overheating	➢ Turn off power supply. Check for sufficient airflow around the power supply fan. After cooling down, restart the power supply by turning the Power switch to the on position. ➢ If you cannot restart the instrument, turn off the power, disconnect the power cord from the outlet, and contact Technical Service.
Warning message displays with 5-second beep sound. The screen shows 	The power once been cut and now recover	➢ User does not need to take extra action. The warning sign and beep sound would only last for 5 second; after that, the machine will continue running the unfinished project. ➢ The  sign indicates the machine has been interrupted by sudden power off. Press Enter Key  to clear the sign.

Encountering Problems

1. Check the troubleshooting section.
2. Call Technical Service or e-mail to service@majorsci.com
3. If the unit must be shipped back for repair, contact Major Science or the distributor for a Return Material Authorization number and shipping instructions. The unit will be repaired and returned to you as quickly as possible.

DOMINIQUE DUTSCHER SAS

Replacing the Fuse

1. Turn off the main power switch at the rear of Power Supply and detach the power cord.


2. Open the fuse compartment located at the Main Power Socket by inserting a small flat blade screwdriver into the slot below the ON/OFF switch. Turn the screwdriver to gently pry open the fuse compartment.

   Note: The fuse compartment will not open with the power cord in place.

3. Pull the fuse holder out of the compartment and inspect the fuse. If the fuse is burned or there is a break in the fuse element, replace the fuse with an identical type of fuse as provided in the fuse holder (see figure below).

4. Place the fuse holder back into the compartment.

5. Snap the cover closed.

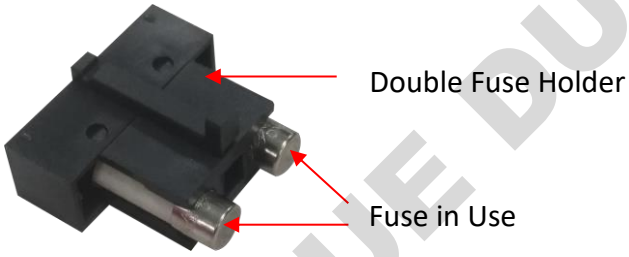

For additional fuses, contact Major Science co. ltd.

Model	Type of Fuse	Ordering Number
MP-310	2.5A/250V	PFV-D-TCF-2.5A-PF-000000
MP-320	5A/250V	PFV-D-000-5A-PF-00000000
MP-510	5A/250V	PFV-D-000-5A-PF-00000000

Maintenance

Lively Power Supply uses all solid-state components and should require no maintenance or recalibration under normal use. If the unit must be returned for repair, contact our **SERVICE DEPARTMENT** or your local distributor for shipping instruction.

DOMINIQUE DUTSCHER SAS

Section 6 Ordering Information

Cat. No.	Description
MP-310	Lively 300V Power Supply
MP-320	Lively 3AP Power Supply
MP-510	Lively 500V Power Supply

DOMINIQUE DUTSCHER SAS

Section 7 Warranty

Major Science warrants apparatus of its manufacture against defects in materials and workmanship, under normal service, for ***one year from the shipping date to purchaser***. This warranty excludes damages resulting from shipping, misuse, carelessness, or neglect. Major Science's liability under the warranty is limited to the receipt of reasonable proof by the customer that the defect is embraced within the terms of the warranty. All claims made under this warranty must be presented to Major Science within one year following the date of delivery of the product to the customer.

Headquarters:
Major Science Co., Ltd.

Contact Information:
Main Office :
No. 156, Sec. 1, Guoji Rd., Taoyuan Dist.,
Taoyuan City 33061,
Taiwan

T/ +886-3-3762878
F/ +886-3-3761310
E-mail：info@majorsci.com

DOMINIQUE DUTSCHER SAS