
SX0004700 Swinnex Filter Holder, 47 mm

- Works only in one direction, as indicated below




Flow

- Adapter, Parts, Assembly

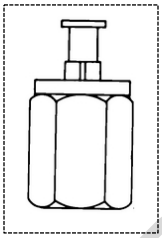

XX3002568 - Adapter, 1/4 in. NPTF to F Luer-Lok

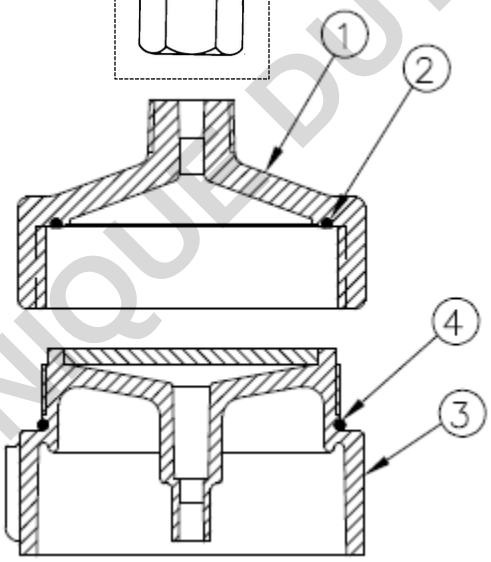





SX0004700 - Swinnex Filter Holder, 47 mm

ITEMS	QTY	DESCRIPTION
4	8	O-RING 2-032
3	8	UNDERDRAIN & BASE ASSY
2	8	O-RING 5-329
1	8	POLYPROPYLENE TOP

DOMINIQUE DUTSCHER SAS

- The adapter is attached by twisting the threaded end into the top (inlet) of the filter holder. The above drawings illustrate the alignment of the adapter with the filter holder.
- Top (or inlet) of holder is item 1. Bottom (outlet) is item 3. Filter fits in between those 2 pieces of the holder.

- To autoclave this assembly (with membrane in place):

Dry Filter Method

1. Holding top upside down, position filter over O-ring seated in top.
2. Replace bottom and tighten by hand. Do not invert the filter holder until unit is tightened to prevent misalignment of the filter.
3. Wrap the unit in autoclavable paper
4. Sterilize at 121 C for 30 minutes
5. Tighten unit with wrenches provided before use

Wet Filter Method

Follow procedure for dry filter method except

1. Wet filter in distilled or DI water before positioning in top half of filter holder
2. Autoclave unit at 121C for 20 minutes

For both of the above methods, MilliporeSigma recommends autoclaving on a “Wet” or “liquid” cycle.

DOMINIQUE DUTSCHER SAS