

OPASTER'ANIOS

BANDELETTES DE CONTROLE - TESTING STRIPS

INDICATIONS

50 bandelettes réactives spécifiques pour le contrôle de la concentration d'ortophtaldéhyde dans la solution d'OPASTER'ANIOS.

INDICATIONS

50 specific reactive strips for testing the concentraction of orthophtaldehyde in the solution OPASTER'ANIOS.

MODE D'EMPLOI

Tremper la bandelette pendant une seconde.

Eliminer l'excès de solution par une légère agitation de la bandelette.

Lire le résultat 60 secondes après trempage. Respecter le délai avant interprétation.

INSTRUCTIONS FOR USE

Soak the strip for one second.

Shake off excess liquid.

Read the result 60 seconds after steeping. Respect this time before reading.

CONDITIONS DE STOCKAGE

Conserver le flacon fermé entre +4°C et +30°C dans un endroit sec pendant 24 mois. Conserver à l'abri de la lumière et de l'humidité.

Si ces conditions sont respectées, les bandelettes sont utilisables jusqu'à 90 jours après ouverture.

STORAGE CONDITIONS

Keep the flask fermly closed at a temperature between +4°C and +30°C, in a dry area for 24 months. Store away from light and moisture.


If these conditions are respected, testing strips are usable until 90 days after opening.

Rose uniforme
CONFORME



Pink uniform
CONFORM

Rose non uniforme
NON CONFORME

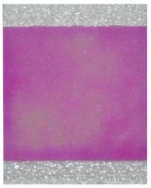

Pink not uniform
NO CONFORM

DOMINIQUE DUTSCHER SAS

laboratoires ANIOS
Le professionnel de la désinfection
Pavé du Moulin
59260 Lille-Hellemmes - France
Tel. +33 3 20 67 67 67 - Fax: + 33 3 20 67 67 68
www.anios.com