misceo by Dynamic | 250F

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Safety instructions

Read the operating instructions in full before starting up and follow the safety instructions.

Keep the operating instructions in a place where they can be accessed by everyone.

Follow the safety instructions, guidelines, occupational health and safety and accident prevention regulations.

Wear your personal protective equipment in accordance with the hazard category of the media to be processed. Otherwise there is a risk from:

- Splashing and evaporation of liquids
- Body parts, hair, clothing and jewellery getting caught.

Check the appliance and accessories beforehand for damage each time you use them. Do not use damaged components.

Only elements approved by MISCEO may be used.

Never operate the appliance without being immersed in liquid.

The tools must never be used outside of a formulation because the sealing system and bearings are destroyed without the tools being cooled by the material.

In the event of unbalance or unusual noises, switch off the appliance immediately. Replace the dispersing element. If there is no difference after the change of the dispersing tool, return it to the dealer or the manufacturer along with a description of the fault.

Warning, this product includes moving parts. Take care when using it and when cleaning it as there is a risk of cuts.

Do not cover the ventilation slots on the block-motor in order to ensure adequate cooling of the motor.

Do not process any flammable or combustible materials.

Do not operate the appliance in explosive atmospheres, with hazardous substances or under water.

Process pathogenic materials only in closed vessels under a suitable extractor hood. Please contact MISCEO if you have any questions.

There may be electrostatic discharges between the medium and the dispersing instrument shaft which could pose a direct danger.

Only process media that will not react dangerously to the extra energy produced through processing. This also applies to any extra energy produced in other ways, e.g. through light irradiation.

Please note that you should avoid standing too close to the device during powder dosing. The powder may be blown away by air turbulence from the drive.

Before handling the device, check that it is unplugged.

The appliance can only be disconnected from the mains supply by pulling out the mains plug or the connector plug.

The voltage stated on the type plate must correspond to the mains voltage.

The socket for the mains cord must be easily accessible.

The appliance does not start up again automatically following a cut in the power supply.

Protect the appliance and accessories from bumps and impacts.

This appliance may be used by children at least 8 years old and by persons with reduced physical, sensory or mental capacities or without experience or knowledge if they are correctly supervised or if they have been given instructions on how to use the appliance safely, and if they correctly understand the risks involved. Children must not play with the appliance. Cleaning and care by the user must not be performed by unsupervised children.

Application

The device offers sufficient power to effectively use the main dispersion and mixing techniques in the laboratory. When combined with a compatible tool, it operates in continuous mode as a disperser or mixer to process liquids of varying viscosities.

ROTOR-STATOR
High shear / low viscosity

This tool is used to :

- Homogenize
- Emulsify
- Disperse powders and pigments in a liquid
- Particle size reduction

MIXING 2 BLADES TOOL

This tool is used to:

- Cut
- Mix
- Dissolve

MIXING 4 BLADES TOOL

This tool is used to:

- Grind
- Mix
- Dissolve

TURBO AGITATOR BLADE
Low shear / low to medium viscosity

This tool is used to:

- Mix
- Stir
- Suspend
- Hydrate

CUTTER BOWL

This tool is used to:

- Grind
- Micronize plants (in the object of extraction)
- Cut

WHISKS TOOL

This tool is used to:

- Mix
- Even distribution
- Stir
- Whip

DISSOLVER DISPERSION TOOL
Medium high speed / low to high viscosity

This tool is used to:

- Fragment
- Disperse particles
- Dissolve powders
- Hydrate

SPIRAL STIRRER TOOL
Low to medium speed / medium to high viscosity

This tool is used to:

- Mix homogeneously
- Distribute evenly

ANCHOR STIRRER
Low to medium speed / medium to high viscosity

This tool is used to:

- Mix
- Stir
- Distribute

Operation

DISMANTLING OF THE TOOL

Once the appliance unplugged, unscrew the tool following the direction of the arrow on the tool (fig. n°1).

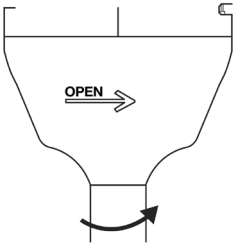
fig. 1) Dismantling of the tool.

MOUNTING OF THE TOOL

Once the appliance unplugged, choose the appropriate tool for the use.

Screw the tool in the reverse direction of the arrow written on the tool (fig. n°2).

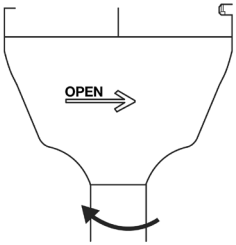
fig. 2) Mounting of the tool.

STARTING-UP

Take into account the ambient conditions listed in the technical data.

Plug the appliance into an electrical wall socket.

When the appliance is plugged in and turned on there is no LED indicating on that the appliance is on.

By pressing the power button (see Fig. 5), the speed setting display indicates '1' and the time display indicates '00: 00'.

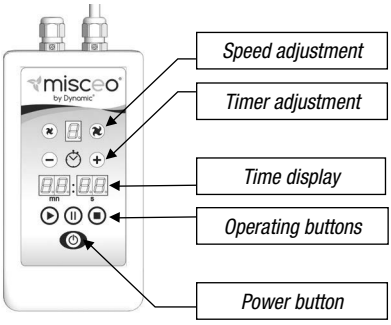
fig. 3) Command unit.

fig. 4) Power button.

fig. 5) Display state at power on.

CHRONOMETER OPERATION

STARTING

The time display scrolls in steps of 1 second.

The appliance automatically switches off after 20mn.

PAUSING

The appliance stops.

The time display is flashing.

The appliance is restarting.

STOPPING THE APPLIANCE

The time display flashes alternately between «STOP» and the value of the operating time.

You cannot restart the appliance any more.

You cannot regulate the speed any more.

The speed setting display indicates '1' and the time display indicates '00: 00'.

TIMER OPERATION

You can use the timer only when the time display indicates «00:00» and when the appliance is stopped.

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One short press: the adjustment step is 1 second.

One long press: the adjustment step is 5 seconds.

One very long press: the adjustment step is 30 seconds.

Maximum time setting: 20mn.

STARTING

The time display scrolls in steps of 1 second.

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The time display is restarting.

STOPPING THE APPLIANCE

The time display flashes alternately between «STOP» and the value of the operating time.

You cannot restart the appliance any more.


You cannot regulate the timer and the speed any more.

The speed setting display indicates '1' and the time display indicates '00: 00' (fig. n°5).

SPEED SETTING

DECREASING SPEED

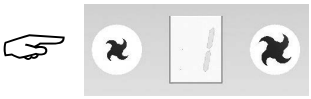You can decrease the speed by pressing the button with the small propeller.

The speed is decreasing by step of 1 for each pressing.

Minimum speed setting : 1.

INCREASING SPEED

You can increase the speed by pressing the button with the big propeller.

The speed is increasing by step of 1 for each pressing.

Maximum speed setting: 8.

OPERATION WITH SUPPORT

CLAMPING OF THE APPLIANCE

1. Spread the fingers of the clamp with the knob.
2. Fit the appliance in the clamp.
3. Tighten the fingers of the clamp.

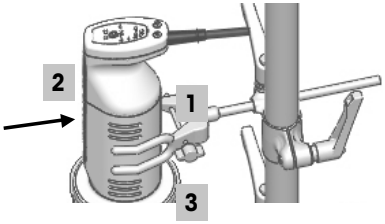
fig. 6) Fitting of the appliance

CLAMPING OF THE CONTAINER / BEAKER

1. Fit the bowl into the V-shaped hold.
2. Fasten the bowl with the belt.

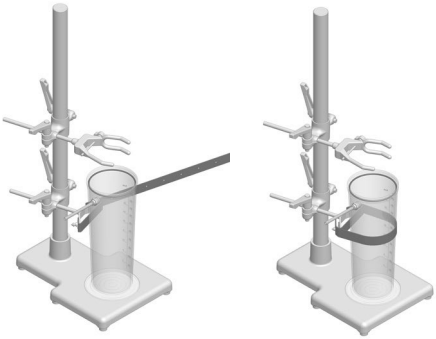
fig. 7) Fitting of the bowl

HEIGHT ADJUSTMENT

1. Release the bosshead with the vertical knob.
2. Adjust to the chosen height.
3. Tighten the bosshead.

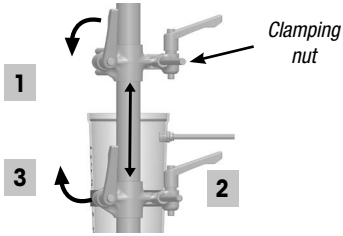
fig. 8) Height adjustment

DEPTH ADJUSTMENT

1. Release the bosshead with the horizontal knob.
2. Adjust to the chosen depth.
3. Tighten the bosshead.

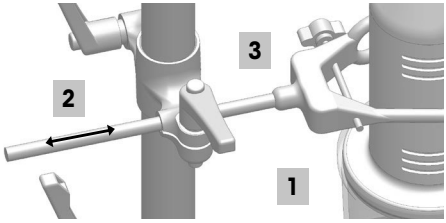
fig. 9) Depth adjustment

PRECAUTIONS OF USE

To avoid unintentional ingress of air into the medium due to excessively turbulent rotation, the whole unit can be placed a little off centre in the stirring vessel.

The safety of the user cannot be guaranteed if the appliance is operated with accessories that are not supplied or recommended by the manufacturer or if the appliance is operated improperly contrary to the manufacturer's specifications or if the appliance or the printed circuit board are modified by third parties.

Maintenance and cleaning

MAINTENANCE

The drive unit is not immune from wear. The motor brushes and the coupling can wear out over time.

It is preferable to regularly check:

- That the power cable is in good condition (the copper wires in the cable should not be visible).
- That the tube is watertight.
- That the blades are sharp.

If the power cable is damaged, it must be replaced by the manufacturer, its after-sales service or people of similar qualification to avoid danger.

CLEANING

Wear protective gloves when cleaning the appliances.

For cleaning disconnect the main plug.

The appliance can be cleaned using a sponge or a slightly damp cloth.

Do not allow moisture to get into the appliance when cleaning.

Please consult with MISCEO before using any cleaning or decontamination methods, other than those recommended here.

SPARE PARTS ORDER

When ordering spare parts, please give:

- Machine type
- Manufacturing number, see type plate
- Item and designation of the spare part

To view our spare parts catalogue, visit our website : https://www.misceo-cosmetics.com/fr/ Then go to the 'Our catalogue' page and click on 'Download'.


Electric requirements

The appliance of a class 2 category (double insulation) powered by a single phase alternating current and does not require a "u" ground plug.

The installation must be protected by a ground fault circuit breaker and a 10 A caliber fuse.

Use a standard 10/16 A single phased socket.

Check that the electricity supply voltage and the value on the appliance specifications plate are compatible.

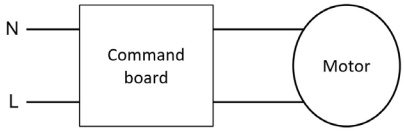Warranty / Declaration of conformity

DECLARATION OF CONFORMITY

We declare under our sole responsibility that this product corresponds to the regulations 2004/108/EC, 2006/42/EC and 2011/65/EU and conforms with the standard IEC 60235-1.

WARRANTY

The warranty does not cover worn out parts, nor does it apply to faults resulting from improper use, insufficient care or maintenance not carried out in accordance with the instructions in this operating manual. In accordance with MISCEO warranty conditions, the warranty period is 12 months. For claims under the warranty please contact MISCEO. You may also send the machine direct to our factory, enclosing the delivery invoice and giving reasons for the claim. You will be liable for freight costs.

misceo by DYNAMIC FRANCE S.A.S.
518 rue Léo Baekeland
85290 Mortagne-sur-Sèvre Cedex - FRANCE
Tel. : +33 (0)2 51 63 02 72
contact@misceo-cosmetics.com
www.misceo-cosmetics.com

REPAIR

The appliance may only be opened by experts.

In case of repair, the appliance has to be cleaned and free from any materials which may constitute a health hazard.

If you require servicing, return the appliance in its original packaging. Use a suitable transport packaging.

Send the appliance to the address below:

misceo by DYNAMIC FRANCE S.A.S.
518 rue Léo Baekeland
85290 Mortagne-sur-Sèvre Cedex - FRANCE
Tel. : +33 (0)2 51 63 02 72
contact@misceo-cosmetics.com
www.misceo-cosmetics.com

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