

# AGANI™

*Hypodermic needle*



## Product Solutions

**Your needles of choice**  
for a wide range of applications



 **TERUMO**

The safety and  
comfort of your patients  
is in your hands.



**Injection & Access**



# Smooth and easy puncture increases patient comfort.

Make the administration of medication as painless as possible with the special geometry of the sharp lancet point and silicone-coated cannula that allows a smooth penetration and prevents any knocking effect.<sup>1</sup>



Precisely shaped needle point  
with silicone coating.

Agani hypodermic needles have a sharp lancet tip and an outside layer of silicone lubricant.

**The sharp point**

- Facilitates the penetration of the skin
- Helps to minimize the trauma

**The silicone coating**

- Provides good gliding properties
- Allows a smooth insertion

# Secure and easy handling

increases safety and user comfort.

The steeper point angle of the short bevel<sup>2</sup> shortens the length of the needle bevel, thereby reducing the risk of vascular perforation.

The outside diameter of the needle can be recognized by the ISO colour-coded hub.

- Easy identification of the required size
- Easy recognition of flashback through the semi-transparent hub

All needles come sterile, individually packed in color-coded blister packaging.

- Supplied in boxes of 100 needles
- Comprehensive information on the packaging



# Choose from a wide range from 18G up to 30G.



## 18G Cannula (ø 1.2 mm)

Length 38 mm  
Colour Pink



## 19G Cannula (ø 1.1 mm)

Length 25, 38, 50 mm  
Colour Cream



## 20G Cannula (ø 0.9 mm)

Length 25, 38 mm  
Colour Yellow



## 21G Cannula (ø 0.8 mm)

Length 16, 25, 38, 50 mm  
Colour Green



## 22G Cannula (ø 0.7mm)

Length 32, 38 mm  
Colour Black



## 23G Cannula (ø 0.6 mm)

Length 16, 25, 32 mm  
Colour Blue



## 24G Cannula (ø 0.55 mm)

Length 25 mm  
Colour Purple



## 25G Cannula (ø 0.5 mm)

Length 16, 25 mm  
Colour Orange



## 26G Cannula (ø 0.45 mm)

Length 13, 23 mm  
Colour Brown



## 27G Cannula (ø 0.4 mm)

Length 16, 19 mm  
Colour Grey



## 30G Cannula (ø 0.3 mm)

Length 13 mm  
Colour Yellow



Does not contain  
components made of PVC  
or natural rubber latex.



For further ordering information, please refer to the AGANI hypodermic needle product information card.

# AGANI™

Hypodermic needle



## Life is easier with needles from Terumo.

Agani hypodermic needles provide you and your patients with the best possible levels of care and comfort.

Sharp lancet point – easy puncture

Silicone coating – smooth needle movement

Thin walled needles – increased flow rate

Short bevel options – reduced risk of vascular perforation

Large portfolio – suitable for a wide range of applications

### Peak performance even in little things – made easy!

Terumo Corporation  
+81 3 3374 8111

Terumo Europe NV  
+32 16 38 12 11

Your Terumo partner

[www.terumo-europe.com](http://www.terumo-europe.com)

© Registered Trademark  
Published by Terumo Europe NV

MP210GB-0919FK-110-19E



## Product Information



### Sterile hypodermic needle for single use

Make the administration of medication as painless as possible with the special geometry of the sharp lancet point and silicone-coated cannula that allows a smooth penetration and prevents any knocking effect<sup>1</sup>.

Description	Raw materials
1 Hub	Polypropylene + master batch
2 Cannula	Stainless steel
Adhesive	Epoxy glue
Lubricant	Silicone oil

Does not contain components made of PVC or natural rubber latex.



#### Sterilization

With ethylene oxide to a Sterility Assurance Level of at least 10<sup>-6</sup>.



#### Storage

- Shelf life of 5 years
- The batch number and expiry date are printed on the packaging
- Do not store at extreme temperature and humidity



#### Precautions

- Check integrity of the package before use
- Do not reuse
- Do not use if unit package is damaged
- Dispose of safely after single use to avoid risk of infection

<sup>1</sup> Data on file

## Your needles of choice for a wide range of applications

### Product range

Product Code	Colour	Cannula ø	Cannula Length	Bevel
AN*1838R1	Pink	18G 1.2 mm	1 ½ inch 38 mm	Regular 11°
AN*1838S1	Pink	18G 1.2 mm	1 ½ inch 38 mm	Short 17°
AN*1925R1	Cream	19G 1.1 mm	1 inch 25 mm	Regular 11°
AN*1938R1	Cream	19G 1.1 mm	1 ½ inch 38 mm	Regular 11°
AN*1938S1	Cream	19G 1.1 mm	1 ½ inch 38 mm	Short 17°
AN*1950R1	Cream	19G 1.1 mm	2 inch 50 mm	Regular 11°
AN*2025R1	Yellow	20G 0.9 mm	1 inch 25 mm	Regular 11°
AN*2038R1	Yellow	20G 0.9 mm	1 ½ inch 38 mm	Regular 11°
AN*2116R1	Green	21G 0.8 mm	⅝ inch 16 mm	Regular 11°
AN*2125R1	Green	21G 0.8 mm	1 inch 25 mm	Regular 11°
AN*2138R1	Green	21G 0.8 mm	1 ½ inch 38 mm	Regular 11°
AN*2150R1	Green	21G 0.8 mm	2 inch 50 mm	Regular 11°
AN*2232R1	Black	22G 0.7 mm	1 ¼ inch 32 mm	Regular 11°
AN*2238R1	Black	22G 0.7 mm	1 ½ inch 38 mm	Regular 11°
AN*2316R1	Blue	23G 0.6 mm	⅝ inch 16 mm	Regular 11°
AN*2325R1	Blue	23G 0.6 mm	1 inch 25 mm	Regular 11°
AN*2332R1	Blue	23G 0.6 mm	1 ¼ inch 32 mm	Regular 11°
AN*2425R1	Purple	24G 0.55 mm	1 inch 25 mm	Regular 11°
AN*2516R1	Orange	25G 0.5 mm	⅝ inch 16 mm	Regular 11°
AN*2525R1	Orange	25G 0.5 mm	1 inch 25 mm	Regular 11°
AN*2613R1	Brown	26G 0.45 mm	½ inch 13 mm	Regular 11°
AN*2623R1	Brown	26G 0.45 mm	⅞ inch 23 mm	Regular 11°
AN*2716R1	Grey	27G 0.4 mm	⅝ inch 16 mm	Regular 11°
AN*2719R1	Grey	27G 0.4 mm	¾ inch 19 mm	Regular 11°
AN*3013R1	Yellow	30G 0.3 mm	½ inch 13 mm	Regular 11°

### Packaging

**Single unit** Blister pack (film/paper)

**Unit box** 100 pieces

Dimensions 105 x 84 x 86 mm

For item numbers AN\*1950R1 and AN\*2150R1 dimensions 105 x 84 x 110 mm

**Shipping carton** 10.000 pieces

Dimensions 450 x 440 x 450 mm

For item numbers AN\*1950R1 and AN\*2150R1 5.000 pieces, dimensions 545 x 240 x 450 mm

### Regulatory Status

Class IIa, sterile product following MDD 93/42/EEC

**Notified body** TÜV Product Service GmbH, TÜV SUD Group 0123


**Manufacturer** Zhejiang Kindly Medical Devices Co. Ltd. (China)

**EC representative** Shanghai International Corp. GmbH (Europe) - Hamburg, Germany

**Distributor** TERUMO EUROPE N.V. Interleuvenlaan 40, 3001 Leuven, Belgium

Choose from a wide range  
from 18G up to 30G

The outside diameter of the needle can be recognized by the ISO colour-coded hub.

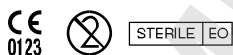

Terumo Corporation  
+81 3 3374 8111

Your Terumo partner

Terumo Europe NV  
+32 16 38 12 11