

Simply better membrane technology, simply better Stericup® filters.

Don't take risks with your media. If your sterile filtration membranes are clogging and/or binding serum proteins in your media, your cells risk:

- Cellular damage
- Exposure to toxins
- Suboptimal cell growth

Trust Stericup® filters with NEW improved Millipore Express® PLUS membrane. Our most popular Stericup® membrane, Millipore Express® PLUS polyethersulfone (PES), now delivers even faster filtration rates while reducing overall protein and reagent loss.


Learn more at: www.merckmillipore.com/sterilefiltration




Membrane design makes a difference.

Membrane	Pore Structure	Membrane Top	Membrane Bottom	Main Features
Millipore Express® PLUS PES (polyethersulfone)	 Asymmetric			<ul style="list-style-type: none"> • High flow rate • Low protein binding • Low volume retention

Protein Binding by Membrane Type

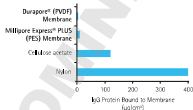

Low Protein Binding: Minimizes Protein Loss

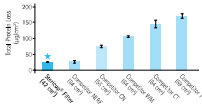


13 mm disks were cut out from different filtration devices (n=3). Protein binding was quantified by incubating disks in 1 mL tracer (1 mg/mL goat IgG in PBS containing ¹²⁵I-goat IgG at 0.1 µCi/mL). After 2 h of shaking at RT, the tracer solution was aspirated. Disks were washed 3X in PBS then assayed for bound radioactivity.


Merck Millipore is a division of MERCK

Simply better Stericup® filters.

Faster filtration.

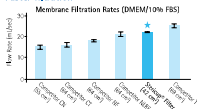

Lower retention volume.

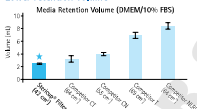

500 mL of DMEM+10%FBS was filtered to completion through various 0.22 µm hydrophilic PES membrane filtration devices. Filtration rates and retention volumes were calculated (n=3).

Stericup® and Steritop® Filter Units

Description	Membrane	Pore Size (µm)	Funnel Capacity (mL)	Receiver Bottle (mL)	Qty/Pk	Catalogue No.
Stericup®-BP Filter Unit	Millipore Express® PLUS (PES)	0.22	150	150	12	SCGFU01RE
		0.22	250	250	12	SCGFU02RE
		0.22	500	500	12	SCGFU06RE
		0.22	1000	1000	12	SCGFU11RE
		0.22	150	150	12	SCGF501RE
Steritop®-BP Filter Unit	Millipore Express® PLUS (PES)	0.22	150	150	12	SCGPT01RE
		0.22	250	250	12	SCGPT02RE
		0.22	250	250	12	SCGPT02RE
		0.22	250	250	12	SCGPT02RE
		0.22	500	500	12	SCGPT06RE
		0.22	500	500	12	SCGPT06RE
		0.22	1000	1000	12	SCGPT10RE

Steriflip® Filter Units

Description	Membrane	Volume (mL)	Pore Size (µm)	Qty/Pk	Catalogue No.
Steriflip®-BP Filter Unit	Millipore Express® PLUS (PES)	50	0.22	25	SCBP00525
Steriflip®-BV Filter Unit	Durapore® (PVDF)	50	0.22	25	SE1M179M6
Steriflip®-BV Filter Unit	Durapore® (PVDF)	50	0.45	25	SE1M003M00

Sterile 33 mm Millex® Syringe Filters

Description	Membrane	Pore Size (µm)	Qty/Pk	Catalogue No.
Millex®-GP Filter Unit	Millipore Express® PLUS (PES)	0.22	1000	SLGP033RK
		0.22	250	SLGP033RB
		0.22	50	SLGP033RS
		0.45	50	SLHP033RB
Millex®-GV Filter Unit	Durapore® (PVDF)	0.22	50	SLGV033RS
		0.45	50	SLHV033RS

To Place an Order or Receive Technical Assistance

In Europe, please call Customer Service:
 France: 0825 045 645
 Germany: 01 805 045 645
 Italy: 848 845 645
 Spain: 901 516 645 Option 1
 Switzerland: 0848 645 645
 United Kingdom: 0870 900 4645
 For other countries across Europe,
 please call: +44 (0) 115 943 0840

Or visit: www.merckmillipore.com/offices

For Technical Service visit:
www.merckmillipore.com/techservice

Get Connected!

Join Merck Millipore Bioscience on your favorite social media outlet for the latest updates, news, products, innovations, and contests!

facebook.com/MerckMilliporeBioscience

twitter.com/Merck4Bio



www.merckmillipore.com/offices

Merck, Millipore and the M logo are trademarks and Stericup®, Millipore Express®, Durapore®, Millex®, Steriflip®, and Steriflip® are registered trademarks of Merck KGaA, Darmstadt, Germany.
 LB No. FR2244/EN/LJ 05-08EN-11-06974 03/2013
 © 2013 EMD Millipore Corporation, Billerica, MA USA. All rights reserved.