




INNOVATIVE SURGICAL GLOVES

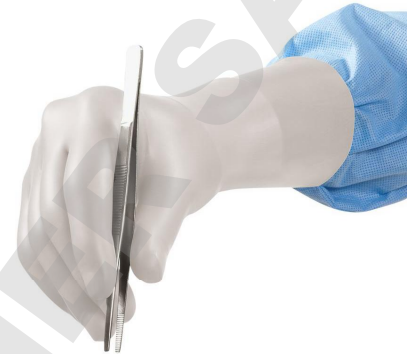

GAMMEX® LATEX

ADVANCED POWDER-FREE LATEX GLOVES FOR GENERAL SURGERY.

Ansell's GAMMEX® Latex glove sets the top surgical standards with their design and latex formulation. An optimal combination of safety and comfort based on Ansell's proprietary technology, with unparalleled levels of ease of use, this glove provides a secure barrier without compromising sensitivity.

Donning is made easy, even with damp hands, thanks to the water-based synthetic inner coating. The adhesive band at the top of the straight cuff allows for a secure hold on the surgical gown: this efficient design relieves pressure on the forearm while operating.

GAMMEX® Latex gloves provide a secure grip on instruments for easy, sure handling in both wet and dry conditions. Grip is supported by the micro-rough outer surface, manufactured using our proprietary textured former combined with closely controlled dosage in the final chlorination stage.

The glove's very low latex protein and chemical accelerator residue content is achieved through lengthy leaching processes and post-production washing cycles, while formulation is safe on skin, as it is free from thiurams and MBT (mercaptobenzothiazole).

Overall, GAMMEX® Latex is the answer to the most demanding surgical requirements, as proven by its success in operating theatres across Europe.

PRODUCT DESCRIPTION

Material	Natural Rubber Latex
Colour	White
Shape	Anatomic
Cuff	Straight cuff with adhesive band
External Surface	Micro-textured, chlorinated & siliconized
Internal Surface	Polyurethane coated & siliconized

SIZE / PRODUCT CODE

5.5	330048055	6	330048060
6.5	330048065	7	330048070
7.5	330048075	8	330048080
8.5	330048085	9	330048090
9.5	330048095		

GAMMEX®, the most innovative and technologically advanced line of surgical gloves available, offers surgeons, OR nurses, surgical technicians and patients the active protection and ultimate performance they need. Through the broad range of available options, GAMMEX® surgical gloves afford the highest level of customised protection against external aggression, allergic skin reactions and bacterial and viral contamination, regardless of the surgical environment.

Natural Latex Powder-Free Sterile Surgical Glove

MD Class IIa & PPE Cat. III
Registering Authority: British Standards Institution (0086): MD, British Standards Institution (0086): PPE

If relevant, read carefully the instructions for use flyer accompanying the product. Single use only. Sterile unless this package is opened or damaged. This product contains natural rubber latex which may cause allergic reactions, including anaphylactic responses.

Ref-XXXX-20180110

Ansell is a world leader in providing superior health and safety protection solutions that enhance human well being. With operations in North America, Latin America, EMEA and Asia, Ansell employs more than 13,000 people worldwide and holds leading positions in the industrial and medical gloves market, as well as in the sexual health and wellbeing category worldwide. Ansell operates in four main business segments: Medical Solutions, Industrial Solutions, Specialty Markets and Sexual Wellness. Information on Ansell and its products can be found at www.ansell.eu.

Ansell Healthcare Europe N.V. (European Head Office)

Riverside Business Park, Block J
Boulevard International 55, 1070 Brussels, Belgium
Tel. +32 (0) 2 528 74 00 • Fax +32 (0) 2 528 74 01 • Fax Customer Service +32 (0) 2 528 74 03
<http://www.ansell.eu> • E-mail info@ansell.eu

