


CHEMICAL RESISTANCE DATA



SHIELDskin XTREME™

SHIELDskin Xtreme Eco Nitrile 300 DI+

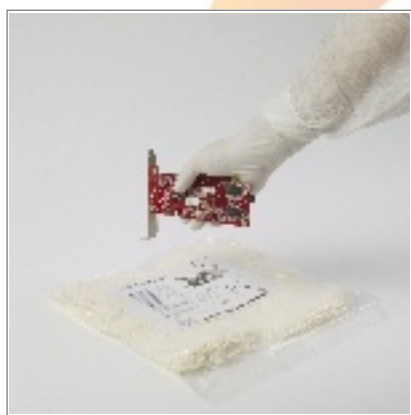

- Cleanroom glove - PPE Complex Design (Category III)
- 100% Nitrile glove - Powder Free
- Unique Skin Nitrile Technology
- Ambidextrous - Textured finger tips
- AQL 1.5 (EN374-2:2003 Level 2)
- Length 300 mm
- Reduced allergy risk (Type IV) - Latex Free glove
- Virus & Microorganisms resistant (ASTM F1671-97b and EN374-2:2003)
- Typical particles level: less than 950 per cm² more or equal 0.5µm

1239-45-8	Ethidium bromide 5%	Breakthrough time : 480 minutes	Class 6
1310-58-3	Potassium Hydroxyde 40%	Breakthrough time : 480 minutes	Class 6
50-00-0	Formaldehyde, 37%	Breakthrough time : 480 minutes	Class 6
75-59-2	Tetramethylammonium hydroxide 2.5%	Breakthrough time : 480 minutes	Class 6
7664-93-9	Sulphuric acid, 50%	Breakthrough time : 480 minutes	Class 6
7647-01-0	Hydrochloric Acid, conc.	Breakthrough time : 79 minutes	Class 3
67-63-0	Isopropanol 70%	Breakthrough time : 45 minutes	Class 2
67-63-0	Isopropanol 100%	Breakthrough time : 35 minutes	Class 2
64-17-5	Ethanol, 70%	Breakthrough time : 29 minutes	Class 1
1336-21-6	Ammonium Hydroxide, 25%	Breakthrough time : 16 minutes	Class 1
64-17-5	Ethanol 100%	Breakthrough time : 11 minutes	Class 1
7664-39-3	Hydrofluoric Acid, 40%	Breakthrough time : 8 minutes	Class 0
67-56-1	Methanol 99,9%	Breakthrough time : 0 minutes	Class 0
67-64-1	Acetone 100%	Breakthrough time : 0 minutes	Class 0