



trackman™


Guiding You to Success!



VERSATILE FUNCTIONALITY

SMART AND SIMPLE

PORTABLE AND RELIABLE

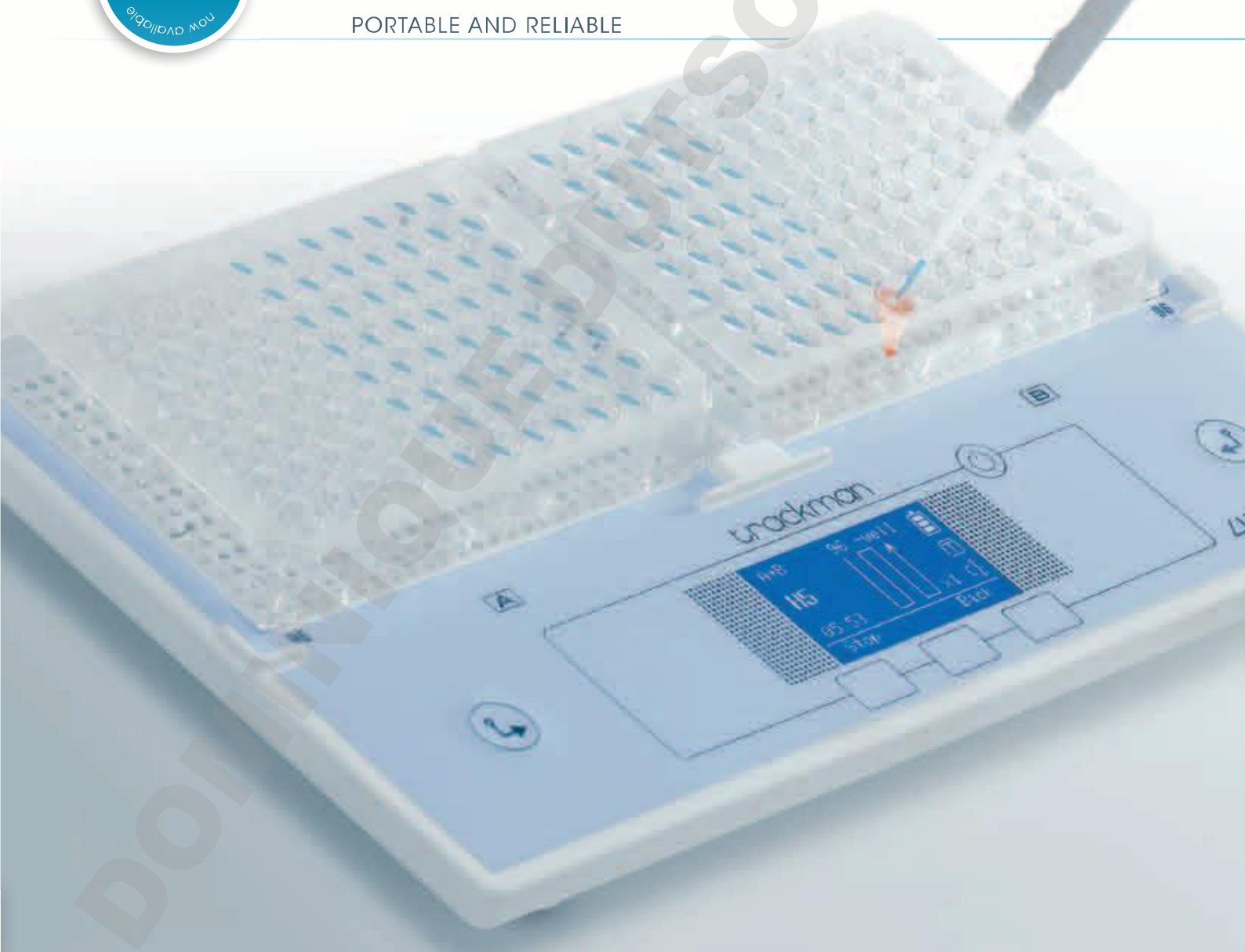

PERSONAL LABORATORY ASSISTANT

VERSATILE FUNCTIONALITY

SMART AND SIMPLE

PORTABLE AND RELIABLE



trackman™

Guiding You to Success!

TRACKMAN is your new personal laboratory assistant. Designed specifically to keep you on track as you pipette from one microplate to another. TRACKMAN simplifies your work, improves your productivity and helps avoid pipetting errors and cross-contamination.




Versatile functionality

For a wide range of applications!

SWITCH FROM USING ONE SIZE PLATE TO ANOTHER EASILY



Thanks to its two plate locations and its plate position indicators, TRACKMAN can be used with both 96 and 384 well microplates, without the need for additional accessories.

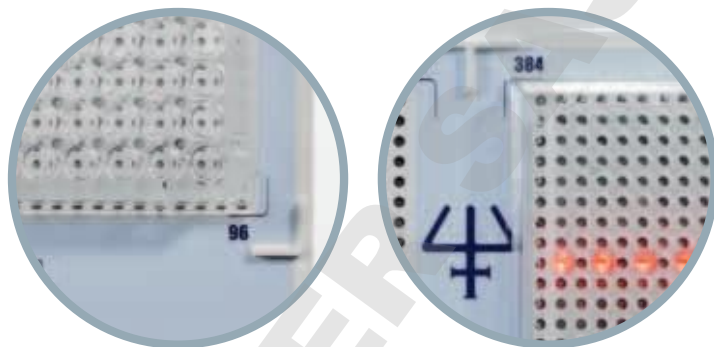

ADAPTS TO YOUR LAB EQUIPMENT


No need to purchase new lab equipment because TRACKMAN works with most microplate types* and pipette types on the market.

- Use with single and multichannel pipettes
- Use with 96 and 384 well microplates
- Use with 2.0, 1.5 and 0.5 mL tubes thanks to the optional smart microtube-holder (accessory)

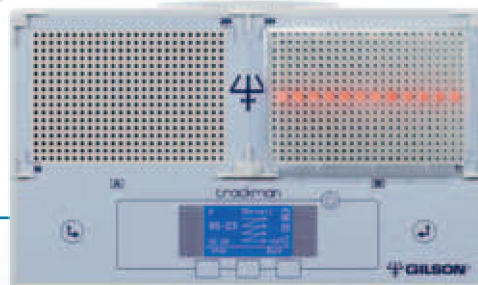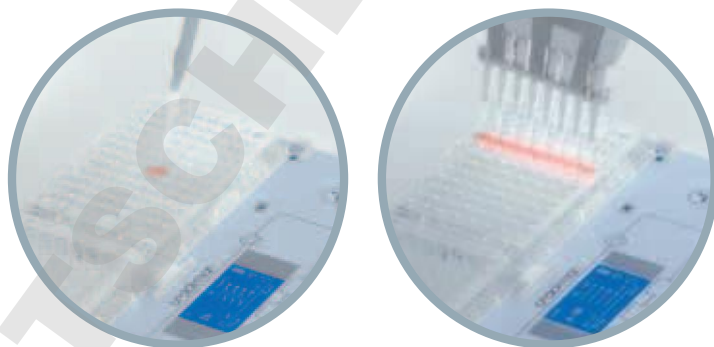

16 programming possibilities for a variety of applications

different pipetting patterns, (individually through rows, columns or other configurations.), different pipette types and vessel types!

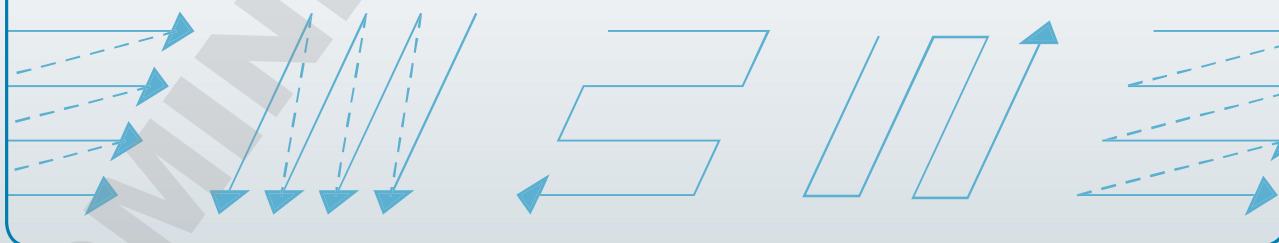

*transparent, translucent, ISO/ANSI SBS 1-2004, 2-2004, 3-2004 and 4-2004 are accepted

Smart and simple

Makes your daily work easier!

SIMPLIFY EVERYDAY PIPETTING!

TRACKMAN, the light box that keeps track of your pipetting, has intuitive software that requires no training. Using TRACKMAN helps avoid pipetting mistakes and eliminates cross-contamination.


SECURE YOUR WORK RIGHT FROM THE START

TRACKMAN's LED lights help you to place your 96 and 384 well microplates in the correct position on the light panel, while fixed notches secure your plates firmly in place.

ADAPTS TO YOUR NEEDS

TRACKMAN's timer function enables hands-free operation, allowing you to stay focused on your work.

For optimal positioning on the bench, TRACKMAN can be used in a flat position or adjusted to a slight incline using its retractable stand.

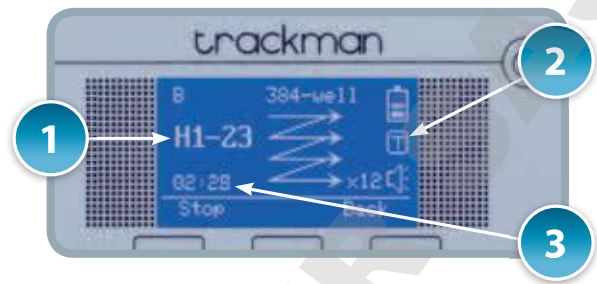

INTUITIVE SOFTWARE INTERFACE

- 1 Indicator of well or tube position
- 2 Timer mode
- 3 Timer



Portable and reliable

Take your TRACKMAN with you!

ALWAYS AT YOUR SIDE



TRACKMAN is a small and lightweight device. Because of its compact size, TRACKMAN can easily find a place on any bench, in any laboratory. Its long-lasting battery allows cordless operation for up to 8 hours.

SECURE TRANSPORT

When travelling with your TRACKMAN, be sure to use the safe transport function. This feature prevents damage to the touch keys, keypad and battery, allowing you to enjoy your TRACKMAN for many years to come!

Complete your TRACKMAN

Expand the functionality of your TRACKMAN with the versatile microtube-holder to transfer liquids from one microtube to another (0.5 mL, 1.5 mL and 2 mL).

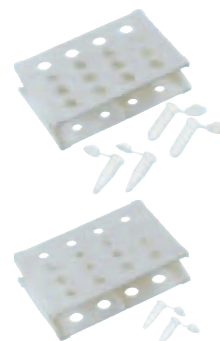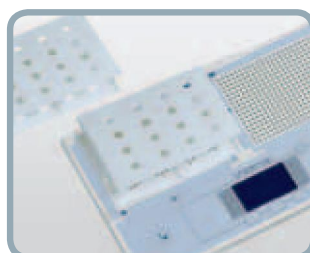

Enjoy hands-free operating and individual speed control, when you equip your TRACKMAN with Gilson's optional foot pedal.



***TRACKMAN simplifies your work,
improves your productivity
and will guide you to success!***



Ordering information

Model	Reference number
TRACKMAN	F70301
Power supply with AC adapter and cable	F30705
3 reagent reservoirs (25 mL)	F267660

Optional Accessories	Reference number
Pedal for TRACKMAN	F70261
0.5 mL, 1.5 mL and 2 mL microtube-holder (vessel capacity: 16 microtubes)	F70105

www.gilson.com | sales@gilson.com | service@gilson.com | training@gilson.com

Gilson, Inc. World Headquarters

3000 Parmenter Street | P.O. Box 620027 | Middleton, WI 53562-0027 USA
Tel: (1) 800-445-7661 or (1) 608-836-1551 | Fax: (1) 608-831-4451

Gilson S.A.S.

19, avenue des Entrepreneurs | BP 145, F-95400 Villiers-le-Bel, FRANCE
Tel: (33-1) 34 29 50 00 | Fax: (33-1) 34 29 50 20



ISO9001 Certified

LT800500. Copyright © 2009. Gilson S.A.S. All rights reserved.
Specifications subject to change without notice. All trademarks
are the property of Gilson Inc. and its subsidiaries.