
OZ Biosciences' Magnetic Rack


following the magnetic field…

How to use our Magnetic Rack


Magnetic Separation Rack Instructions
Read before use

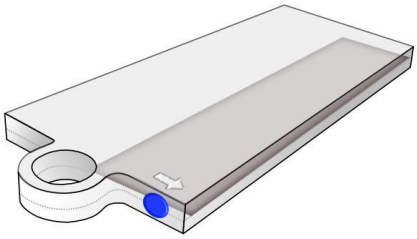


Magnetic insert
The colored dot and the arrow on the magnetic insert indicates the location of the magnetic field.

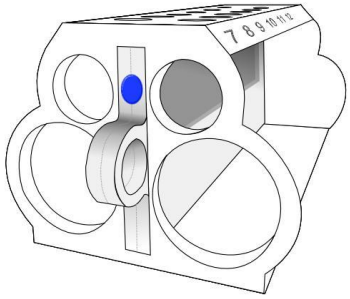

POSITION OF THE MAGNETIC FIELD FOR 1.5mL or 15mL* TUBES

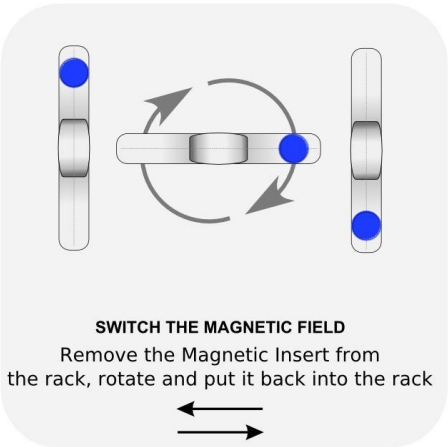

SWITCH THE MAGNETIC FIELD
Remove the Magnetic Insert from the rack, rotate and put it back into the rack

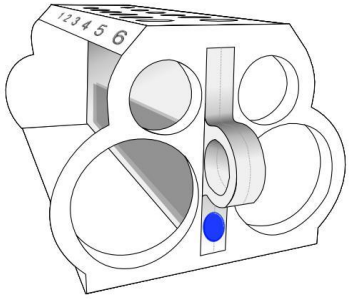

POSITION OF THE MAGNETIC FIELD FOR 50mL* TUBES

* use it in vertical position.

Our Magnetic rack should be stored in a cool, dry environment and should be cleaned with a damp cloth and mild detergent or 70% ethanol. Do not autoclave. This magnetic rack should be handled with care, avoid direct contact with magnet material; pacemakers, implants, watches, electronic devices, all ferrous material and magnetic media have to be keep away.

Contact information FRANCE & Europe	Contact information USA
▪ **Technical service:** +33(0) 486 948 514 tech@ozbiosciences.com	▪ **Technical service:** +1(858) 246 7840 techUSA@ozbiosciences.com
▪ **Customer service:** +33(0) 486 948 516 contact@ozbiosciences.com	▪ **Customer service:** +1(855) 246 7840 contactUSA@ozbiosciences.com

We wish you the best for your experiments with OZ Biosciences' reagent.
OZ Biosciences' team