
Operation

The vacuum is controlled by the pressure sensitive turquoise control button on the top of the hand piece.
By means of the plastic screw opposite the control button, the residual suction level can be set, which will prevent dripping of waste liquids from the tip or pipette. The aspirator is factory pre-set to allow a low level of constant suction; by turning the screw one full turn clockwise, the valve is shut, and residual vacuum is stopped.

Cleaning and servicing

The Hand Operator may be decontaminated with H_2O_2 gas (max. 35 %) for 60 min and is entirely autoclavable. When working with biohazardous material, disinfect the VACUBOY with a 70 % solution of isopropanol or ethanol or a suitable disinfectant before disassembling for autoclaving. Please refer to the explosion drawing.

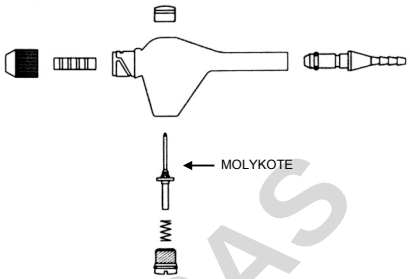To clean, immerse the pieces in mild laboratory detergent or 70 % ethanol, but do not use abrasives, as a roughened surface may affect the smooth operation of the control valve. Apply a thin coat of grease MOLYKOTE as needed to the control valve to keep the action smooth (arrow).

When working with salt-rich solutions, rinse out the aspirator hand piece with distilled water after use to avoid build-up of dried deposits.

Chemical compatibility


Avoid prolonged contact (days or weeks) with:
- concentrated inorganic acids
- sodium hypochlorite
- methylene chloride
- concentrated aldehydes
- other halogenated hydrocarbons


Resistant to:
- DMSO
- ethanol/methanol
- dimethylformamide
- acetone

Note:
- Flushing the aspirator with ethanol immediately after contact with organic solvents is the best safeguard in all cases.
- For O-ring replacement, use INTEGRA O-rings only

Accessories and consumables

VACUBOY hand operator set 155500

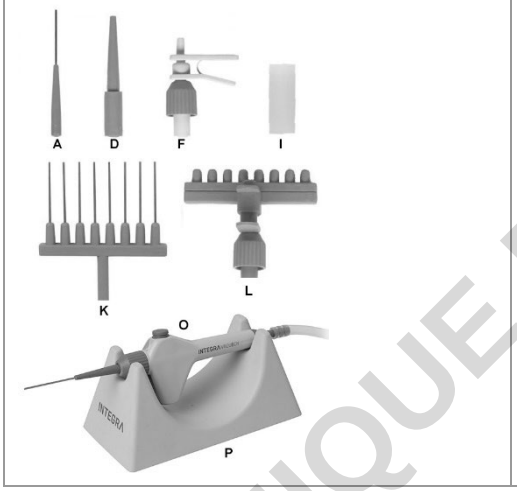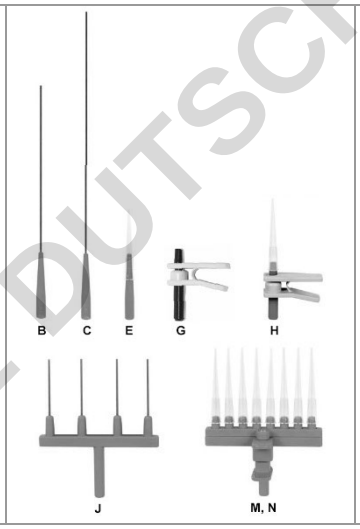	VACUBOY hand operator set small: hand operator (O), two aspiration adapters (A and F), rubber adapter for Pasteur pipettes (I), grease (Dow Corning MOLYKOTE HP-870)	155510
	VACUBOY hand operator set large: hand operator (O), aspiration adapter D, A, K, L, F, rubber adapter for Pasteur pipettes (I), stand for hand operator (P), silicone tubing 1.8 m, grease	155500
A	Adapter 1 channel, 40 mm, stainless steel	155502
B	Adapter 1 channel, 150 mm, stainless steel	155522
C	Adapter 1 channel, 280 mm, stainless steel	155525
D	Adapter 1 channel, with locking for tips/GRIPTIPS, plastic	155504
E	Adapter 1 channel, for tips/GRIPTIPS, plastic, pack of 5	159023
F	Adapter 1 channel, with ejector and locking for tips, plastic	155526
G	Adapter 1 channel, with ejector for tips, plastic	159026
H	Adapter 1 channel, with ejector for GRIPTIPS, plastic	159027
I	Adapter for Pasteur pipettes, rubber	155505
J	Adapter 4 channel, 40 mm, stainless steel	155524
K	Adapter 8 channel, 40 mm, stainless steel	155503
L	Adapter 8 channel, with ejector and locking for tips, plastic	155520
M	Adapter 8 channel, with ejector for tips, plastic	159024
N	Adapter 8 channel, with ejector for GRIPTIPS, plastic	159025
P	Stand for VACUBOY hand operator	155501
	Silicone tubing, Ø 8 mm, 1.8 m, VACUBOY to bottle	158332
Q	Stand for PIPETBOY (*acu* or *pro*) and VACUBOY	155065

FAC SIMILE TECHNIQUE DUTSCHER SAS

464083_V09_OI_VACUBOY

Manufacturer: INTEGRA Biosciences AG
CH-7205 Zizers, Switzerland
www.integra-biosciences.com
info@integra-biosciences.com

Customer service: Please contact your local INTEGRA Biosciences representative. You can find the name and address on www.integra-biosciences.com.

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FRANCIQUE DUTSCHER SAS
Fabricante: INTEGRA Biosciences AG