
Quick Start Guide For MagnetaPure 16

- **The drawer should be released before power on, otherwise, the alarm will appear and the instrument may be damaged.**
- **This Guide can not instead of the User Manual, read the User Manual before operation is necessary.**
- **Check the voltage before using.**

1 Preparation

Take out the instrument from carton, remove the tape that holds cabin door and display screen, unscrew the screw that fixed the lifting part at the bottom of the instrument, release the lifting part. Then install the screw into the placement hole (to prevent the screw loss, for next use).

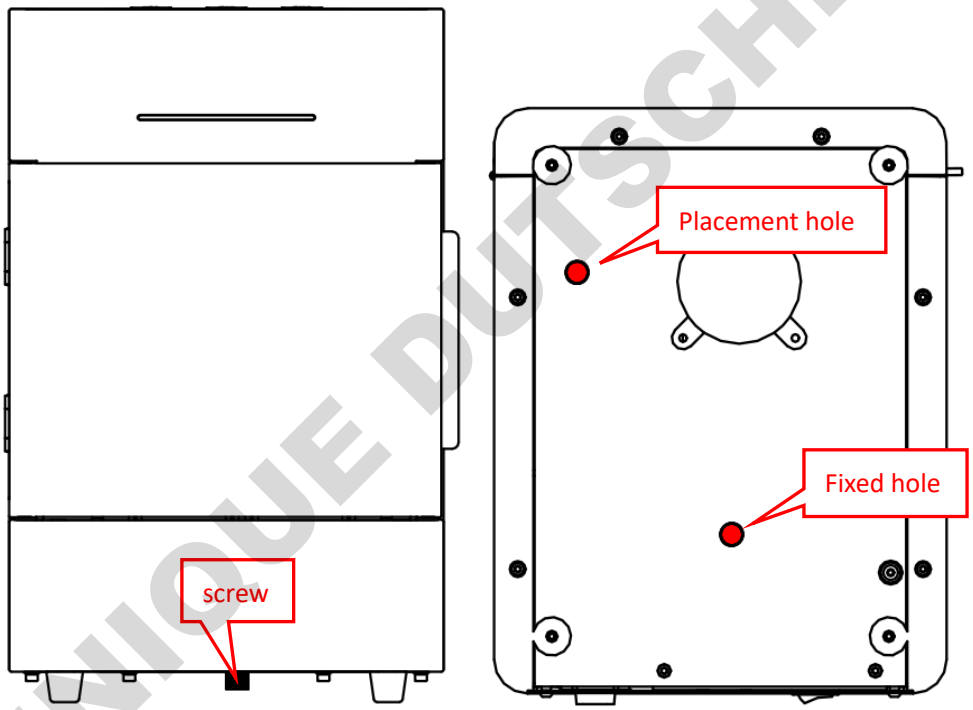

Release the lifting part

Note: Attention to the magnetic rod in case any damage during demolition!

2 Connect the power


Power the instrument on(AC 100～240V).

3 Magnetic rods and kits installation

3.1 Reagent kits installation

Open the transparent lid, place kits to test area, pay attention to direction of kit (For installation, put kits in locating slot according to arrows①, then press it in according to arrows②. Note: attention to the direction of 96-deep well plate).


DOMINIQUE DUTSCHER SAS

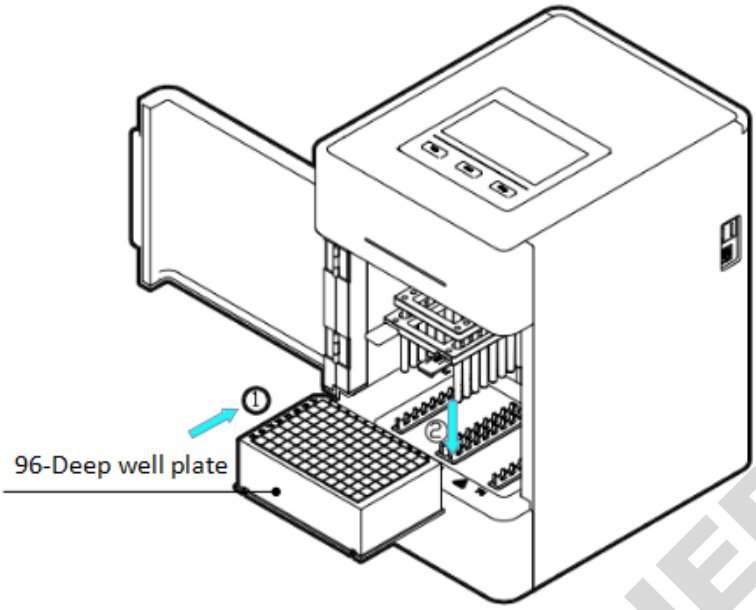

Reagent kits installation

3.2 Magnetic rods installation

Insert magnetic rod's tip into fixed slot of the magnetic rod's tip holder according to the arrow direction.

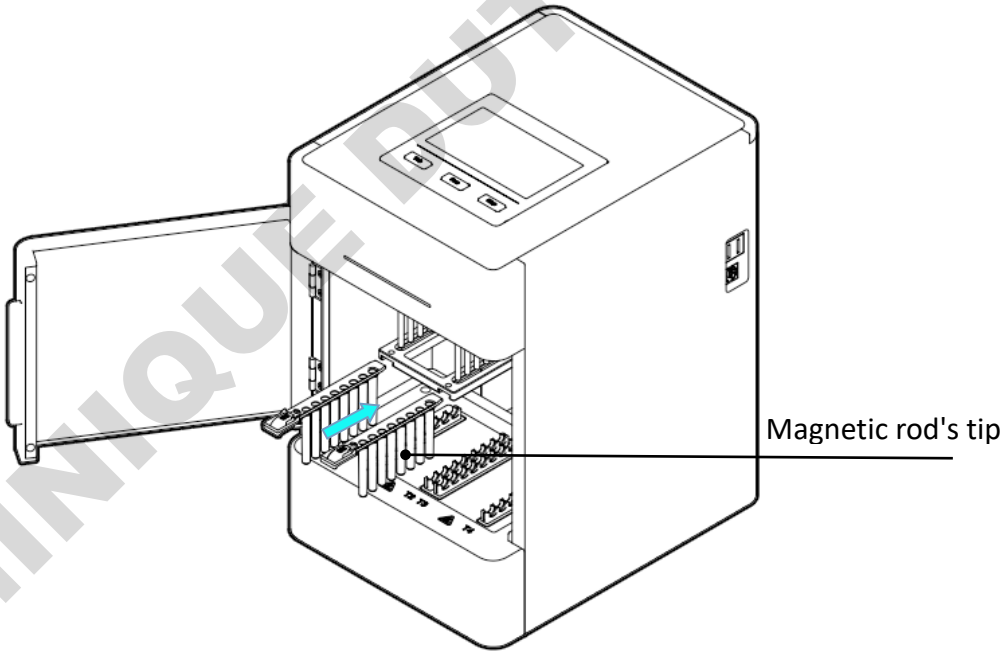

Magnetic rod's tip installation

4 APP Installation

Please copy the "APK" file to mobile phone memory in the root directory(The Kernel version of a Mobile should be over 6.0), then go to "file management" of the mobile, next, move to internal storage to find the APK file, install the software please.

5 Instrument Operation

For details of setting and operations, please refer to Operation Manual.

DOMINIQUE DUTSCHER SAS