

#### **CLEARLINE LADDER**

Range: 100 bp - 1000 bp

**Cat. No.: 257738** 100 lanes of 5 μl

MADE IN **DENMARK** 

-	ClearLine Ladder
ID No.	CL0.050-0042
Cap colour	Black
Content	1 x 0.5 ml

## **Features and General Description**

ClearLine DNA Ladder is a broad range dsDNA ladder with bands from 100 bp to 1000 bp (figure 1).

The ladder is supplied in loading buffer, ready-to-use on agarose and polyacrylamide gels. It is suitable with TBE, TAE, SB and LB electrophoresis systems.

The 500 bp band of ClearLine DNA Ladder comes with a higher intensity and can be used as a reference point.

### Suggestions for use of ClearLine PCR Ladder

- Important: Mix ladder briefly before use. Do not heat the ladder.
- Load 5  $\mu$ l of ladder per lane (for gel wells of ~5 10 mm width).
- If band intensity is too high, then load a smaller volume of the ladder per lane.
- Agarose gel electrophoresis: Prepare 1-1.5 % gel. The dye should migrate 60 70 % the length of the gel.

- Polyacrylamide gel electrophoresis: Prepare 8 % gel. The dye should migrate approx. 90 % the length of the gel.
- Ethidium bromide (0.5  $\mu$ g/ml) is the recommended gel stain.

#### Storage, Stability and Shipment

ClearLine DNA ladder may be kept safely at room temperature for at least 6 months from date of shipment. ClearLine DNA Ladders are guaranteed for 12 months when stored at 4°C. Aliquot product if necessary, to avoid repeated freezing and thawing cycles.

Shipped at ambient temperature. For long term storage, keep the ClearLine DNA Ladder at -20 °C for up to 3 years.

#### **Quality control**

Agarose gel analysis shows that all bands are present at the expected location and band intensity.

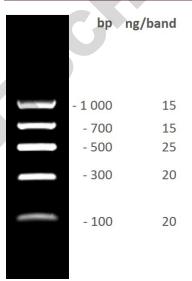


Figure 1: ClearLine DNA Ladder. 5  $\mu$ l ClearLine PCR Ladder was loaded on a 1.5 % agarose in 1x TBE and stained with ethidium bromide.

Other product sizes, combinations and customized solutions are available. Please look at www.dutscher.com or ask for our complete product list for PCR Enzymes. For customized solutions please contact us.

Made in Denmark

Issued 08/2021



# **ClearLine DNA Ladders**

