

IncuSafe

Multigas Incubators

50 L



MCO-50M-PE



Personal-size compact Multigas Incubator

Patient-specific, compact design, triple stackable **IncuSafe** Multigas Incubator with variable CO₂ and O₂ control to simulate *in vivo* conditions. During cell culturing, the inCu-saFe germicidal interior and SafeCell UV lamp work continuously to prevent contamination.

Reproduction of *in vivo* conditions


With a unique, solid-state zirconia sensor for precise oxygen control (1-18%; 22-80%), the MCO-50M Multigas Incubator is able to reproduce low oxygen concentrations found in many tissues and organs.

Time-Saving Decontamination

The high-speed decontamination system uses vaporised hydrogen peroxide and UV light. It cleans the chamber of the incubator safely in less than three hours, achieving a minimal 6 log reduction of major contaminants.

Precise Control & Intelligent Monitoring

An OLED alphanumeric keypad allows convenient but secure user control. It can display internal conditions, such as CO₂ level and temperature. Transfer of data is easy via a USB port.

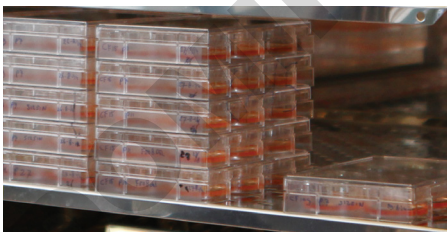

Sensitive Cell Culturing

Culturing cells at physiological oxygen levels allows them to grow faster and live longer, and reduces the frequency of mutations.



Individual Cell Culturing

Compact and stackable these incubators are ideal for individual cell cultures from patient samples or small scale research projects.



Easy to Use

Adjustable audible and visual alarms are standard, along with integrated system diagnostics and predictive performance supervision. The password-protected control panel provides security and minimizes risk of accidental changes in setpoint.

IncuSafe Multigas Incubators

Zirconia O₂ Sensor

The incubator's unique, solid-state, Zirconia O₂ sensor delivers precise control of physiological oxygen levels to simulate *in vivo* conditions.



Dual IR CO₂ Sensor



The incubator's DUAL IR sensor and P.I.D. control enables ultra-fast CO₂ recovery without overshoot, even following multiple door-openings.

Active Background Decontamination

The exclusive inCu-saFe copper-enriched stainless steel alloy interior offers the germicidal properties of copper and the durability of stainless steel. The optional, isolated, SafeCell UV lamp decontaminates circulating air and water in the humidifying pan, without harming cultured cells.

Condensation Management

With a unique antibacterial coating, the 'dew stick' – controlled by Peltier technology – condenses water on its surface, which then drips into the humidifying pan, preventing unwanted condensation in the chamber and possible contamination.

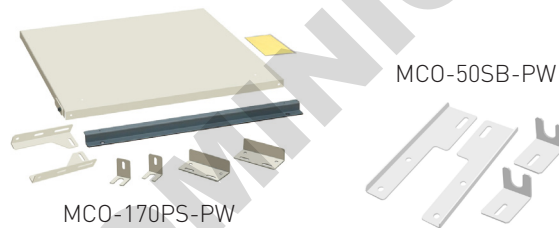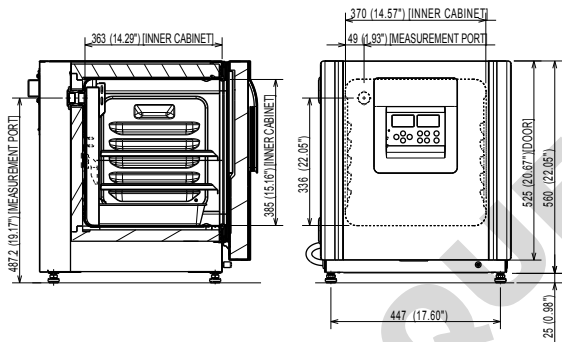

EEA, Switzerland and Turkey only



For medical use

The MCO-50M-PE is certified as a Class IIa Medical Device (93/42/EEC and 2007/47/EC). Applicable countries: Austria, Belgium, Cyprus, Denmark, Finland, France, Germany, Ireland, Italy, Liechtenstein, Luxembourg, Malta, the Netherlands, Spain, Switzerland and the United Kingdom only



For laboratory use

Applicable countries: EEA countries, Switzerland and Turkey

| Model Number | | MCO-50M-PE |
|---|--------|---|
| External Dimensions (W x D x H) ¹⁾ | mm | 480 x 550 x 585 |
| Internal Dimensions (W x D x H) | mm | 370 x 363 x 385 |
| Volume | liters | 50 |
| Net Weight | kg | 46 |
| Performance | | |
| Temperature Control Range & Fluctuation | °C | AT +5 - +50, ±0.1 |
| Temperature Uniformity ²⁾ | °C | ±0.25 |
| CO ₂ Control Range & Fluctuation ³⁾ | % | 0 - 20, ±0.15 |
| O ₂ control range & Fluctuation ⁴⁾ | % | 1 - 18 and 22 - 80, ±0.2 |
| Humidity Level & Fluctuation | %RH | 95, ±5 |
| Control | | |
| Temperature Sensor | | Thermistor |
| CO ₂ Sensor | | Dual IR |
| O ₂ Sensor | | Stabilized Zirconia Sensor |
| Display | | Digital (White graphic OLED) readable to 0.1 increments |
| Construction | | |
| Exterior Material | | Painted Steel |
| Interior Material | | Stainless Steel Copper-Enriched Alloy |
| Insulation Material | | Styrene AcryloNitrile copolymer |
| Heating Method | | Direct Heat & Air Jacket System |
| Outer Door | qty | 1 |
| Outer Door Lock | | N/A |
| Field Reversible Door | | Included |
| Inner Door | qty | 1 |
| Shelves | qty | 2 x Stainless Steel Copper-enriched Alloy |
| Shelf Dimensions (W x D x H) | mm | 353 x 308 x 12 |
| Max. Load per Shelf | kg | 7 |
| Max. Shelf Capacity | qty | 6 |
| Access Port | qty | 1 |
| Access Port Position | | Rear upper left |
| Access Port Diameter | ∅ mm | 30 |
| Alarms | | (R = Remote Alarm, V = Visual Alarm, B = Buzzer Alarm) |
| Power Failure | | R |
| Out of Temperature Setting | | V-B-R |
| High Temperature | | V-B-R |
| Out of CO ₂ Setting | | V-B-R |
| Out of O ₂ setting | | V-B-R |
| Door open | | V-B |
| Electrical and Noise Level | | |
| Power Supply | V | 220-240 |
| Frequency | Hz | 50 |
| Noise Level ⁵⁾ | dB | 29 |
| Options | | |
| SafeCell UV [®] System | | MCO-170UVS-PE |
| H ₂ O ₂ Decontamination Board | | MCO-50HB-PW |
| Electric door lock with password | | MCO-170EL-PW |
| H ₂ O ₂ Vapor Generator | | MCO-50HP-PW |
| H ₂ O ₂ Reagent, pack of 6 bottles | | MCO-50H202-PE |
| CO ₂ Gas Pressure Regulator | | MCO-010R-PW |
| N ₂ Gas Pressure Regulator | | MCO-010R-PW |
| Automatic CO ₂ Cylinder Changeover System | | MCO-5GC-PW |
| InCu-saFe [®] Shelf | | MCO-50ST-PW |
| Double stacking bracket | | MCO-170PS-PW |
| Stacking plate | | MCO-50SB-PW |
| Roller Base | | MCO-50RB-PW |

Appearance and specifications are subject to change without notice.

¹⁾ Exterior dimensions of main cabinet only, excluding handle and other external projections

²⁾³⁾⁴⁾ Ambient temperature 25°C, setting 37°C, CO₂ 5%, O₂ 5%, no load

⁵⁾ Nominal value

* MCO-50M-PE requires MCO-50HB-PE, MCO-170EL-PW, MCO-50HP-PW and SafeCell UV option for H₂O₂ decontamination



PHC Europe B.V.

Nijverheidsweg 120 | 4879 AZ Etten-Leur | Netherlands
T: +31 (0) 76 543 3839 | F: +31 (0) 76 541 3732

www.phchd.com/eu/biomedical