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embecta Medical II LLC
300 Kimball Drive
Parsippany, NJ 07054 USA


embecta Netherlands B.V.
c/o Regus
Weerdestein 97
1083 GG Amsterdam

EC REP
embecta Penel Limited
Pottery Road
Dun Laoghaire
Co. Dublin, IE
ec-rep@embecta.com

CH REP
embecta Switzerland Sàrl
Route de Crassier 11
Terre-Bonne Park Z1
1262 Eysins,
Switzerland

UK REP
embecta U.K. Limited
Winnersh Triangle
Building 220, Wharfedale Road
Wokingham
United Kingdom, RG41 5TP

EN Defined per ISO 15223-1
SV Definierad enligt ISO 15223-1
DA Definition nach ISO 15223-1
NO Gedefinieerad volgens ISO 15223-1
FI Definitsione coma da norma ISO 15223-1
EL Definición según ISO 15223-1
ES Definido de acordo com a norma ISO 15223-1
PT Ορίζεται σύμφωνα με το πρότυπο ISO 15223-1
IT Definiert gemäß ISO 15223-1
DE Definert iht. ISO 15223-1

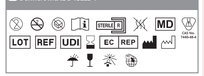CH REP
- **EN** Swiss Authorized Representative
- **SV** Representant agréé en Suisse
- **DA** Bevollmächtigter Vertreter in der Schweiz
- **NO** Erkende vertegenwoordiger in Zwitserland
- **FI** Rappresentante autorizzato per la Svizzera
- **EL** Representante autorizado na Suíça
- **ES** Εξουσιοδοτημένος αντιπρόσωπος για την Ελβετία
- **PT** Auktoriserad representant i Schweiz
- **IT** Autorisert representant i Sveits

UK REP
- **EN** UK Responsible Person
- **SV** Ansvarig för Storbritannien
- **DA** Verantwortliche Person im Vereinigten Königreich
- **NO** Verantwoordelijke persoon VK
- **FI** Responsabile per il Regno Unito
- **EL** Responsable en Reino Unido
- **ES** Pessoa responsável no Reino Unido
- **PT** Υπεύθυνος στο Ηνωμένο Βασίλειο
- **IT** Ansvarig person i Storbritannien
- **DE** Ansvarlig person i Storbritannia

UK CA 0050
- **EN** UK Conformity Assessed
- **SV** Evaluation de la conformité au Royaume-Uni
- **DA** UK-Konformität bewertet
- **NO** Conformiteit VK beoordeeld
- **FI** Conformità valutata per il Regno Unito
- **EL** Evaluación de conformidad del Reino Unido
- **ES** Conformidade avaliada no Reino Unido
- **PT** Αξιολογήθηκε η συμμόρφωση στο Ηνωμένο Βασίλειο
- **IT** Overensstemmelse for Storbritannien bedömd
- **DE** UKCA-merket

- **EN** Open Here
- **SV** Ouvrir ici
- **DA** Hier öffnen
- **NO** Hier openen
- **FI** Aprire qui
- **EL** Abrir aquí
- **ES** Abrir aqui
- **PT** Ανοίξτε εδώ
- **IT** Öppna här
- **DE** Åpnes her

CE 2797 PAH

embecta, the embecta logo, and Micro-Fine are trademarks of Embecta Corp. and registered or used in the United States or other countries. © 2025 Embecta Corp. All rights reserved.

embecta.com 20250315 PC00000

BACK PAGE

FRONT PAGE

Micro-Fine™

embecta™

- **EN** Insulin syringes with sterile interior **U-100** INSULIN
- **FR** Seringues à insuline à intérieur stérile INSULINE **U-100**
- **DE** Insulinspritzen mit sterilem Innenraum **U-100** INSULIN
- **FI** Insulinespuiten met steriele binnenkant **U-100** INSULINE
- **SV** Sirinåe per insulina con l'interno sterile INSULINA **U-100**
- **NO** Jeringas de insulina con interior estéril INSULINA **U-100**
- **PT** Seringas de insulina com interior estéril INSULINA **U-100**
- **EL** Σύριγγες ινσουλίνης με αποστειρωμένο εσωτερικό ΙΝΣΟΥΛΙΝΗ **U-100**
- **ES** Insulinsprutor med steril insida **U-100**-INSULIN
- **IT** Insulinsprøyter med steril innside **U-100** INSULIN

EN Intended Purpose:
Insulin Syringes are sterile, single-use devices designed for subcutaneous administration of insulin. These syringes have only been tested and validated for the intended use.

Indications:
Insulin Syringes are intended to be used for subcutaneous administration of insulin at a concentration of 100U/ml (U-100) for the treatment of diabetes mellitus.

Contraindications:
None known.

Patient target group(s):
The injections of a desired dose of insulin into the subcutaneous layer is intended for patients with diabetes mellitus.

Intended user(s):
Laypersons, caregivers or health care professionals.

Clinical benefits:
The general clinical benefit of insulin syringe products is the ability to administer insulin into the subcutaneous tissue. The administration of insulin is critical for patients with diabetes mellitus.

Precautions:
- Consult with your healthcare professional if you need any advice on recommended needle length or injection technique, have any difficulties using the device, or have potentially experienced an injury or cross-contamination.
- Do not use if the unit package or unit is damaged. Before use, verify the integrity of each unit.
- The following residual risks are related to insulin syringes: pain at injection site, inflammation, bleeding, bruising, induration, local infection, allergic reaction, hypoglycemia, hyperglycemia, lipodystrophy, and medical intervention (broken needle in skin). Needle stick injury can result in cross-contamination (infection) with blood pathogen. A pyrogenic (fever producing) needle may cause fever.
- Once the syringe is prepared, the medication should be administered as soon as possible.
- Never re-shield the syringe after injection.
- Do not share syringes with others.
- Keep out of reach of children.
- Do not re-use device. Re-use can be linked to hardening of the fat under the skin. This may change insulin absorption, making blood glucose control more unpredictable.[1,2]
- EU Only: users should report any serious incident related to the device to the Manufacturer and National Competent Authority.
- Used Insulin Syringes should be disposed of in a special "sharps" container, or a puncture-proof, clearly marked container.

Instructions For Use (Read these 12 steps first):
1. The sterility barrier is composed of the white cap and coloured shield. Before use, verify the integrity of each unit. Do not use if the product sterility barrier is compromised.
2. Wash hands and gather supplies. Twist white cap then pull straight off and discard.
3. Wipe top of insulin bottle with an alcohol swab. If you are taking cloudy insulin, roll the bottle between your hands until it is uniformly cloudy. To avoid the formation of air bubbles do not shake the bottle of insulin.
4. To expose the needle, twist the coloured shield, pull straight off, being careful not to bend the needle or let the needle touch anything. Discard coloured shield.
5. Pull the syringe plunger down to the desired number of units. You need air in the syringe equal to the amount of insulin you will take.
6. Push the needle through the center of the rubber top of the insulin bottle and push plunger down completely.
7. Leave the needle in the insulin bottle. Carefully turn the bottle and syringe upside down, so the bottle is on top.
8. Pull the plunger down slowly, aligning the thin black line of the plunger with the desired number of units on the syringe.
9. If air bubbles appear in syringe, push the plunger up, injecting insulin back into vial and redraw insulin to desired number of units. Confirm the dose is correct and then pull the syringe out of the bottle.
10. Confirm the dose is correct and then clean a small area of skin. Ensure the skin surface is completely dry before injecting.
11. Hold the syringe like a pencil. Pinch up your skin and push the needle quickly through the skin at 90° (straight in) to the skin surface. Push the insulin in with the plunger. Pull the needle out of your skin. Release the skin pinch.
12. Dispose of the product according to local regulations and medical facility biohazardous waste protocols.

Notes:

Insert Style: 4 page Accordion Fold
Unfolded: 22" x 8.5" - Folded: 5.5" x 4.25" (2 fold)

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