

DIRECTION FOR USE

READY-TO-USE

exeol GTA 2%

High level disinfectant

INSTRUMENTS ENDOSCOPY


Ready-to-use glutaraldehyde-based solution for final disinfection by immersion of pre-cleaned semi-critical and critical reusable invasive and non-invasive heat-sensitive medical devices.




Before use, check the expiry date on the label. After opening, the product can be used until the expiry date stated on the label.
exeol GTA 2% can be used for up to 30 days maximum.
Regularly check the validity of the bath using the test strips **exeol strips GTA 2%.**
Before using **exeol GTA 2%**, carry out pre-cleaning or complete mechanical cleaning and intermediate rinsing of the medical devices to be treated.

READY TO USE | MANUAL SOAKING | TEST STRIPS


5L container

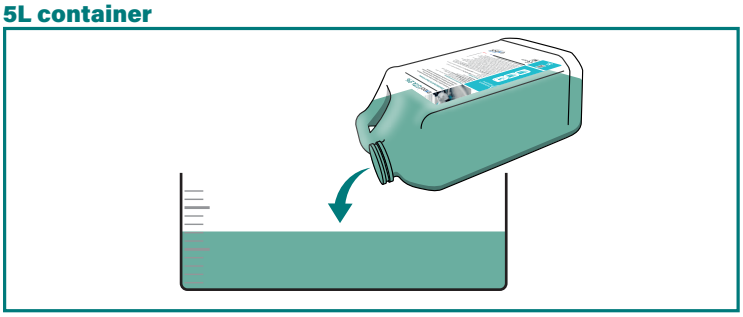


1 When using for the first time, remove the tamper-proof ring before unscrewing the cap. Unscrew the container cap and transfer the contents into the soaking tank. Close the container after use.

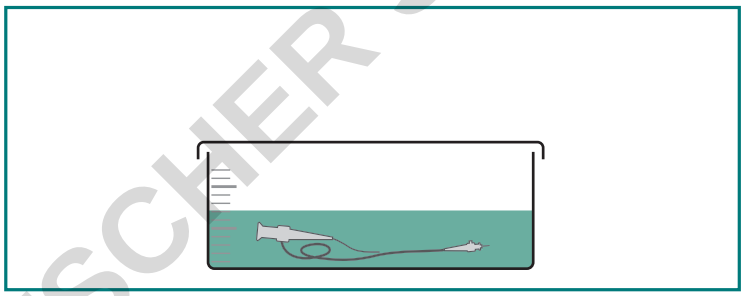


2 Open and disassemble the medical devices and place them in the solution, ensuring that they are completely immersed. Ensure that the channels are irrigated, if these are present. Cover the soaking tank.

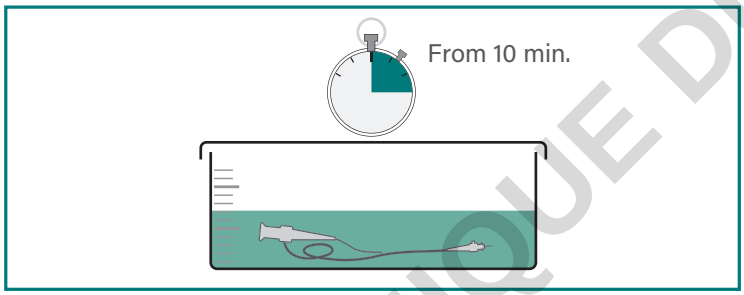


3 Contact time: from 10 min. depending on the required activity spectrum.

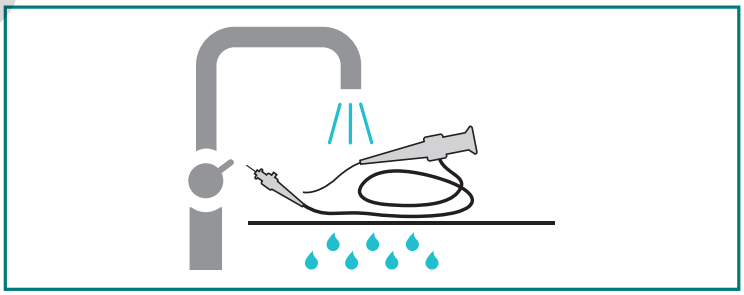

4 Rinse thoroughly using sterile or filtered water (0.2µm) before reuse.

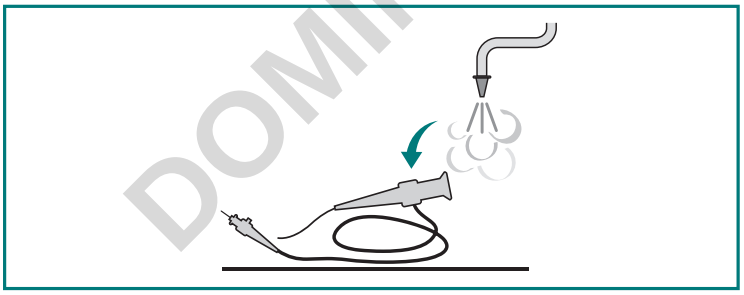


5 Before putting in storage, dry the device using medical air on a dry surface cleaned in advance.

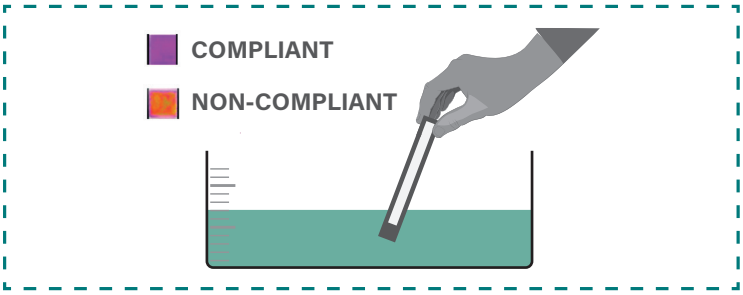

1. Immerse the control strip **exeol strips GTA 2%** for 1 second in **exeol GTA 2%.**
2. Immediately shake the strip lightly to remove excess liquid.
3. Exactly 15 seconds after removal of **exeol GTA 2%**, observe the color of the test area and compare it to the reference visuals above.

The **exeol strips GTA 2%** test strip is 100% sensitive and specific when instructions for use are followed.

Always read the label and product information before use. PPE :

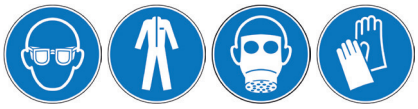

PROFESSIONAL USE ONLY

Product file available on www.exeol.fr

Sodel
190 rue René Barthélemy
14100 Lisieux, France.
TEL: +33 (0)2 31 31 10 50
www.exeol.fr

CE marking obtained in 2023
MD class IIb
CE 0459
MD

MADE IN FRANCE

SANTÉ exeol
A SODEL DIVISION

PRT EXEOL GTA2%-EN-02_2026-01

DOMINIQUE DUTSCHER SAS