TECHNICAL DATA SHEET
SCHEDA TECNICA

FORENSIC SWABS
TAMPONI FORENSI

DOMINIQUE DUTSCHER SAS

DESCRIPTION

Forensic swabs available with wooden shaft and cotton tip or with plastic shaft and viscose tip. They are contained in a labelled cylindrical test tube Ø12x150 mm, made of unbreakable medical polypropylene in order to preserve the sample from possible contamination or degradation, for an easier use and to guarantee the safety of the operators. The pre-applied label on the test tube and on the cap serves both for the identification of the sample and as a guarantee seal that the product has not been previously used. The swabs are sterilized with ethylene oxide in order to guarantee the complete elimination of the nucleic acids present on the swab. Each batch of forensic swabs is analyzed by an accredited laboratory that verifies the absence of human DNA, DNAse, RNAse and PCR inhibitors.

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INTENDED PURPOSE

Forensic swabs for use in the genetic and forensic field for tests such as legal paternity/maternity, determination of the genetic profile, etc. The swabs are sterile and ready to use systems, designed for the collection of biological samples and compatible with the main DNA extraction systems. Not suitable for diagnostic purposes. Microbiotech is not responsible for any unauthorized use. The device has no CE marking (it is not a medical device) as it is not intended for diagnostic purposes.
For professional use only.

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APPLICATOR

TAMPONE

The applicator shaft of the swab is mounted on the tube plug which serves as an optimal grip for an efficient swab handling. Available in different versions:

- **Wooden shaft and cotton tip:** made in 100% cotton tip (European Pharmacopea) and atoxic wooden stick.
- **Plastic shaft and viscose tip:** made in anti-shock polystyrene, very flexible but easily breakable as needed, with atoxic viscose tip.

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TEST TUBE

PROVETTA

Swab is contained in a cylindrical test tube Ø 12 x 150 mm, in transparent, unbreakable medical polypropylene (PP), to preserve sample and user's safety. With pressure cap in polyethylene (HDPE) with specific colour according to stick type. Pre-labelled on test tube to allow a correct sample identification. Item code, lot number and expiry date printed on each label for a full traceability.

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
WOODEN STICK

ASTA LEGNO

Sterile swabs with wooden shaft and cotton tip, in polypropylene Ø 12 x 150 mm test tube with label.

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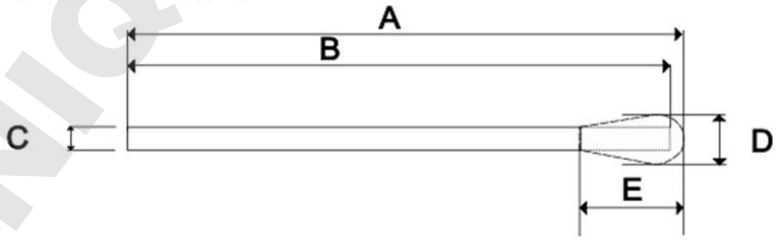100% cotton tip (European Pharmacopea) and handle in biodegrable atoxic white birch wood, completely round, straight and polished, resistant to sudden changes of temperature. The tip is bonded to a handle, atoxic glue to fix cotton tip. Very absorbent, soft, and safe, cotton-tipped applicators perform consistently and reliably.

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
PLASTIC STICK

Sterile swabs with plastic shaft and viscose tip, in polypropylene Ø 12 x 150 mm test tube with label.

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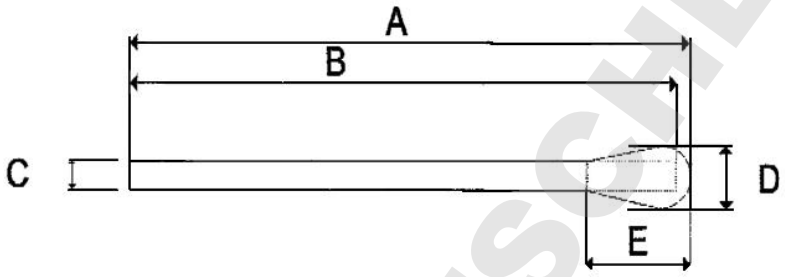This specialty swab is made with viscose fibers. It is bonded securely to a rigid white birch polystyrene handle by an aqueous based adhesive.

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PRECAUTIONAL MEASURES

- The swab must be suitable only for the defined sample and the defined sampling site.
- Do not deviate from the intended use.
- Do not use the product if it is expired or the package is opened/damaged. Sterility guaranteed if unopened.
- Use the swab following aseptic procedures.
- Single-use device; do not reuse. Reusing the device could contaminate the sample.
- The fiber is only guaranteed to adhere to the shaft for instant sampling.
- The shaft is breakable, exert moderate pressure during sample collection.
- Store in a cool, dry place at a temperature between +2°C (35.6°F) and +30°C (86°F). Do not freeze.
- After use, the device may contain infectious microorganisms. Use appropriate PPE and dispose of the

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DOMINIQUE DUTSCHER SAS

test tube and swab according to current regulations for medical waste.
- Keep out of flame or heat sources which might damage the product.
- Keep out of reach of children.
- Strictly follow the user's instructions.
- The device is intended only for the collection and transport of samples for nucleic acid analysis, not intended for diagnostic tests.

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QUALITY ASSURANCE

All raw materials, components and every lot of finished product undergo strict quality inspections carried out by N.A. in order to check conformance to specifications.

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PRODUCT DETERIORATION

The content of the unopened and undamaged units is guaranteed to be sterile. Do not use the units in case of damage, dehydration or contamination. The Devices must be stored at +2°C (35.6°F) / +30°C (86°F) until used. Improper storage will result in a loss of efficacy. Do not use after expiration date indicated on label.

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HOW TO USE

1. Remove the applicator from the test tube.
2. Collect the sample from the concerned site applying moderate pressure during sample collection, the shaft is breakable.
3. Put the applicator into the test tube, seal tightly with the cap and write the sample's data on the label.

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