Rev 2026/05

Instant FLYFood - Ready-to-Use Vials

Pre-dosed Instant Medium fo *Drosophila melanogaster* - Biosigma

NEW

Description:


Instant FlyFood - Vials is the first instant powdered medium, pre-dosed directly into vials, designed for the rapid and standardized preparation of nutritive medium for *Drosophila melanogaster* rearing.
The mixture, pre-dosed and packaged directly in the final rearing vial, only requires the addition of room-temperature or slightly warm water to form a stable and homogeneous gel.

The specially stabilized formulation ensure:
- Immediate preparation compared to traditional methods, with no weighing, transferring, or boiling required
- Reduced preparation time and operational costs
- Compatibility with thermosensitive compounds, which can be directly added to the hydration water
- Protection against microbial contamination through integrated propionic acid
- High reproducibility and standardization between vials and between batches

The product is available in the following formats:
- Ø 25 mm – Code BFF025I
- Ø 28 mm – Code BFF028I

Instructions for Use:

1. Gently tap the vial to allow the powder to settle completely at the bottom, then remove the cap.
2. Place the vial on a vortex mixer and start agitation. The powder must be well mixed and free of agglomerates.
3. While **maintaining agitation**, add:
 - 5 mL of water for Ø 25 mm vials
 - 8 mL of water for Ø 28 mm vials
4. Vortex vigorously until the powder is completely and uniformly hydrated.
5. Allow to stand until full water absorption and uniform gel formation are achieved (about 60 minutes).
 ***Note**: Gel formation time may vary depending on environmental conditions. Ensure the gel is compact and dry before adding flies, to avoid trapping.*
6. The product is ready for Drosophila rearing: apply the suitable cap (not supplied) for immediate use or store the vials in a refrigerator until use.

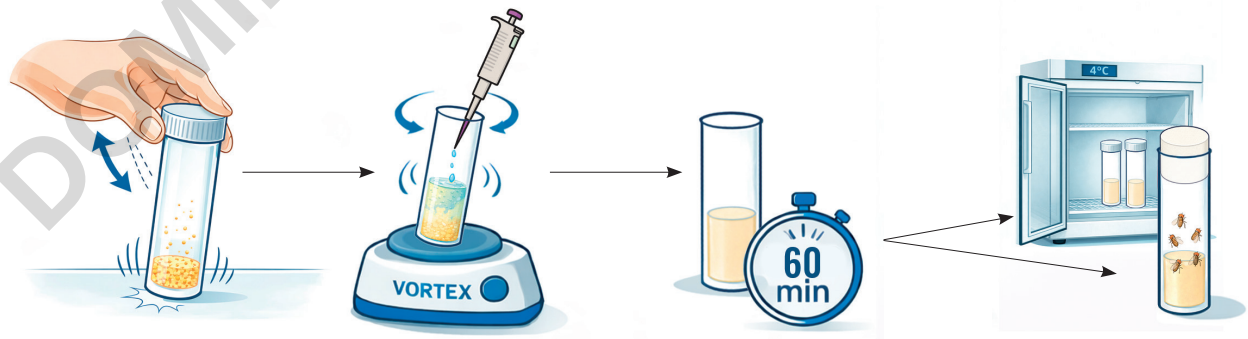Storage and Stability:

- **Unhydrated product**: Store in a dry place at 15–25 °C, protected from light and heat sources.
- **Hydrated product**: Store at 4 °C until use.
- **Shelf-life**: 24 months for the unhydrated product when stored under the indicated conditions.

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DOMINIQUE DUTSCHER SAS
Biosigma S.p.A Dominique Dutscher Group