

Low Temp PCR[®] Rack 96 Well





Slide, Set, Work & Store PCR[®] Samples

The PCR[®] Storage Rack, by Heathrow Scientific, offers a sleek and modern design that's packed full of features. The patent pending design offers a unique lid the rack slides out from, and then rotates 180°, to become a convenient stand the rack sits on. Ergonomic textured finger grips provide a secure hold, reducing the chance of dropping the PCR[®] rack, and makes handling and transportation easy, especially when using gloves. The lid has a side opening that enables faster cooling of samples to ensure integrity and is made from translucent polycarbonate for easy sample and well identification.

The PCR[®] rack base is made from robust ABS construction that helps keep samples safe so you don't have to spend extra money and time on expensive retesting. A location ID of numbers and letters allows for organized sample placement which results in easy and quick sample retrieval. The PCR[®] rack has a closed base design that helps samples integrity by staying cold longer.


The PCR[®] storage rack is versatile; it can be used on a bench top or for freezer storage. The compact footprint is stackable, which helps save valuable space in the lab. The PCR[®] storage rack is offered in multiple colors for coding options.

- The lid is ambidextrous; it can be switched to accommodate left or right hand orientation
- The lid is detachable
- Holds 0.2mL tubes, individual or strip
- 96 well 8x12 configuration
- Can be used for freezer storage down to -80°C (-112°F)

Part No.	Tubes	Wells	Well Dia.	L x W x H	Pkg.	Color
HS120538	0.2mL	96	5.2mm	12.5 cm x 9.2 cm x 3.7 cm (4.9 x 3.6 x 1.5 in)	5/pk	Assorted Colors 
HS120539	0.2mL	96	5.2mm	12.5 cm x 9.2 cm x 3.7 cm (4.9 x 3.6 x 1.5 in)	5/pk	Blue 
HS120540	0.2mL	96	5.2mm	12.5 cm x 9.2 cm x 3.7 cm (4.9 x 3.6 x 1.5 in)	5/pk	Green 
HS120541	0.2mL	96	5.2mm	12.5 cm x 9.2 cm x 3.7 cm (4.9 x 3.6 x 1.5 in)	5/pk	Purple 

* Patent Pending

Heathrow Scientific[®]

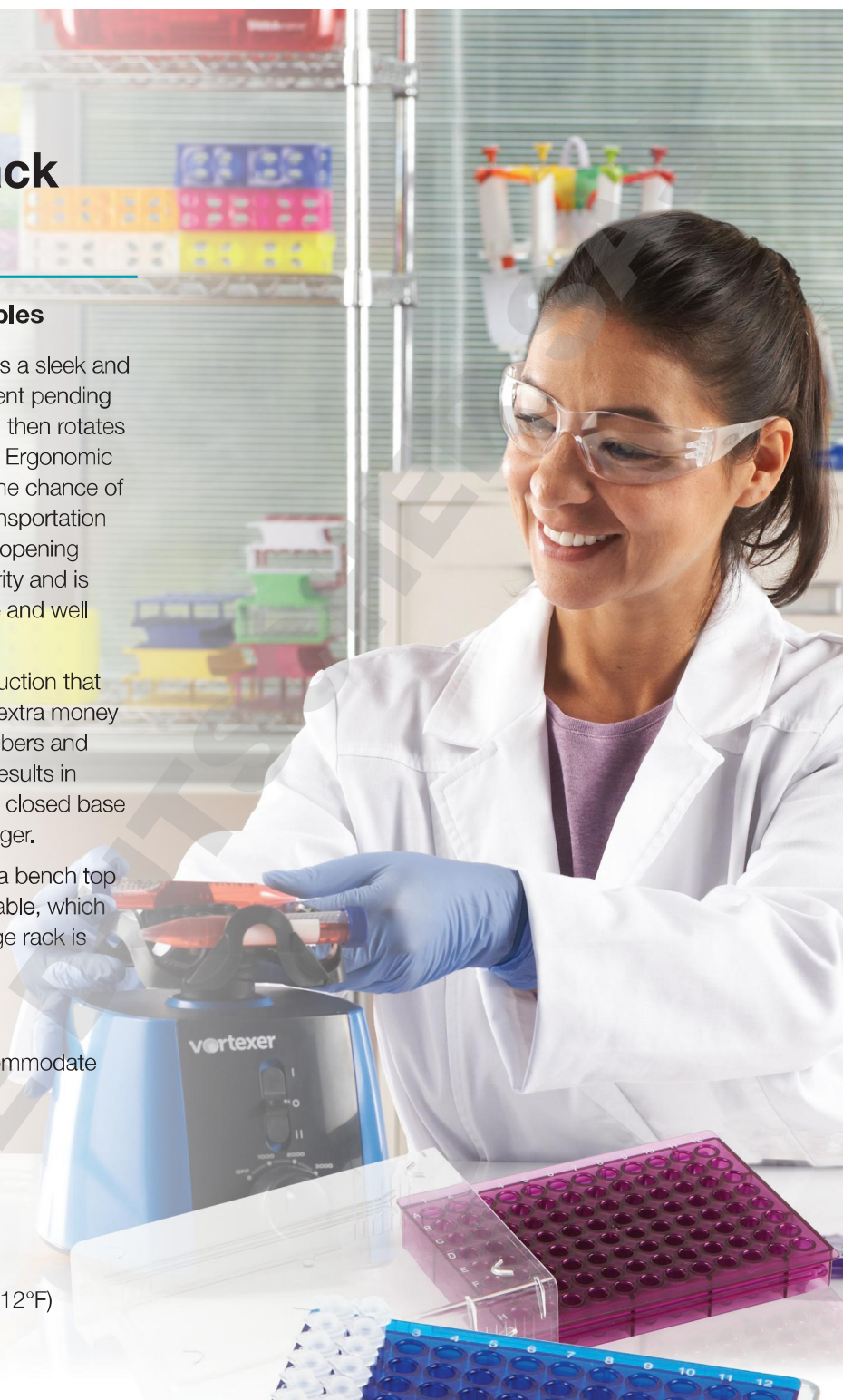



ID of numbers and letters allows for organized sample placement

Side opening enables faster cooling of samples to ensure integrity



A closed base design that helps samples integrity by staying cold longer

Ergonomic textured finger grips



Slide, Set, Work & Store

The patent pending design offers a unique lid the rack slides out from, and then **rotates 180°**, to become a convenient stand the rack sits on.

