

DATA SHEET

202840en	Date: 16/04/15	
Review	Page	

Codes: 202840, 202841, 202842, 202843, 202848, 202844 Security screw cap containers

Catalogue page 2014: 160

Code	Ø mouth mm	Height mm	Volume ml	Case quant.	Case weight kg	Case volume m3
202840	34	43	20	1000	7.20	0.102
202841	45	43	40	600	8.75	0.097
202842	45	60	60	500	8.20	0.097
202843	45	75	90	400	7.80	0.097
202848	53	68	90	300	6.35	0.097
202844	53	78	120	300	7.70	0.100

Directive 98/79/EC. "In vitro" diagnostic medical devices



INTENDED USE

Hospital use, as a containment container with safety closure for transporting biological samples of human origin (urine, aspirates, sputum, biopsies, ...) between hospital departments (for example: from the operating room to pathology, from microbiology to biochemistry, from emergency to different laboratories, etc.).

Intended for the transport of both solid and liquid samples, as well as anatomical specimens with Formaldehyde (Formalin).

Elaborated	Elaborated Verified	
Genna	Chuallh	

www.deltalab.es



DATA SHEET


202840en	Date: 16/04/15	
Review	Page	
,	2 do 2	


DESCRIPTION

Manufactured in translucent polypropylene with yellow polyethylene cap. Container neck with continuous thread. Body graduated in ml. They are provided with **double internal security closure**, so the leakproof closure makes them able to transport liquid samples without risks. Supplied uncapped with the caps and containers in the same box. Ideal for the pneumatic transport of liquid samples, including those that have Formaldehyde (Formalin).



The base and the cap are ribbed so the handling with gloves is comfortable and reliable. There is a ring in the upper part of the body, in order to avoid the dripping in case of liquid decanting. For other colors (blue, orange, magenta), ask for the minimum orders quantities.

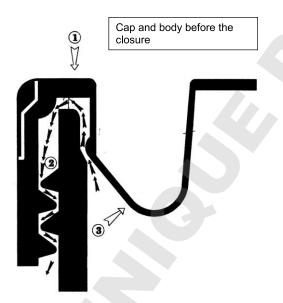

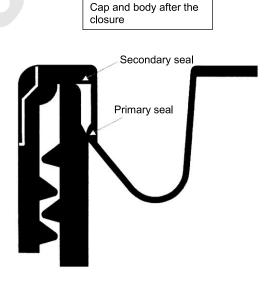

Elaborated	Verified		
Genna	Chwallh	San Land	

www.deltalab.es