
Oosafe® 35 mm Dish, Non-Treated Surface

Description (D, H) (mm)	Prod. No.	Qty/pc	Qty/cs
40.0, 13.1	OOPW-TF06	10	500

Intended Use:

General Intended Use: For the examination of specimens, including blood and tissue donations, derived from the human body.

Specific Intended Use: Aid in small volume cells culturing and observation of cells growth under microscope.

Product compartments: 1

Dish Volumes:

Well compartment 1:9 ml

Safety

Observe the expiry date and the storage condition on the label.
Do only use optically perfect and undamaged plasticware.


Handling

Prior to use, open the sealing and leave the dishes off-gas before use.

Oosafe® 35 mm Wall Dish, High Wall, Label Area Grip packages may be off-gassed by partially opening the sleeve or taking the dishes entirely out of the packaging. In either case, the opened sleeves or dishes must be kept in a sterile environment at room temperature and humidity. The dishes should remain covered with the lids until they are loaded with the sample and medium.

Place small droplets on the bottom of the dish. The recommended volume per drop is 20-50 µl. Overlay the drop with oil. The recommended volume is 5ml. It is also possible to use the dish entirely instead of only making drops. Use the 2 finger grips to hold the dish securely.

When the dish is covered by the lid, the air can still go in and out of the dish.

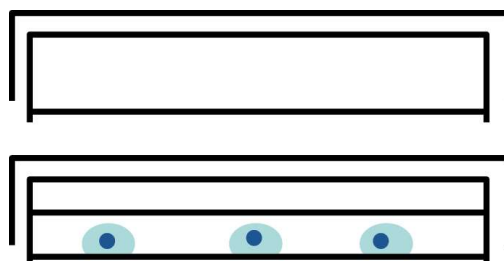

Discard after use.

Storage: Store at room temperature.