
Oosafe® Andrology Tube, Non-Treated Surface

Description (volume) (ml)	Prod. No.	Qty/pc	Qty/cs
5ml	OOPW-AT10	10	500

Intended use


General Intended Use: For the examination of specimens, including blood and tissue donations, derived from the human body.

Specific Intended Use: Aid in andrology use, examination and diagnosis of infertility with or without gassing.

Product compartments: 1

Tube Volumes:

Tube compartment 1: 5 ml

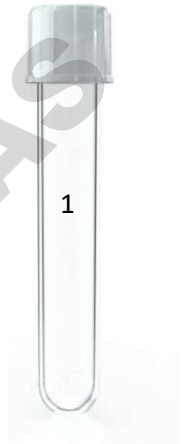

Safety


Observe the expiry date and the storage condition on the label.
Do only use optically perfect and undamaged plasticware.

Handling

The package of Oosafe® Andrology Tube is breathable material with no VOC release. After opening the sealing, the tube can be directly put into use.

Although it can be filled up to 5 ml, the recommended volume for use is under 4 ml. Follow the scale on the tube.

The cap can be loosely put over the tube for gas changing. The cap can also be pushed down thoroughly for tight sealing.

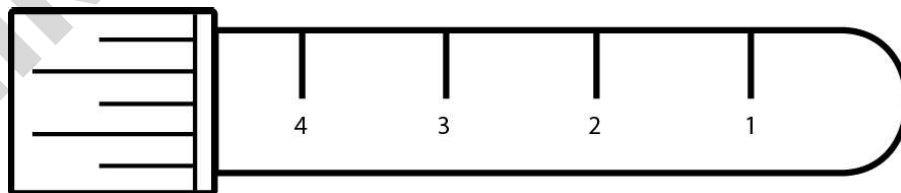

Discard after use.

Storage: Store at room temperature.