User Guide

Amicon® Ultra-4 10K Centrifugal Filter Devices

for volumes up to 4 mL

UFC801008D
UFC801024D
UFC801096D

IVD For *in vitro* diagnostic use

Languages available

English 2
Français 8
Italiano 14
Deutsch 20
Español 26
Português 32
Ελληνικά 38
Nederlands 44
Dansk 50
Svenska 56
Latviski 62
Lietuviškai 68
Magyar 74
Čeština 80
Polski 86
日本語 92
中文 98
한국어 104
Norsk 110
Slovenčina 116
Türkçe 122

DOMINIQUE DUTSCHER SAS

SigmaAldrich.com

Millipore®

English

Introduction

Amicon® Ultra-4 10K centrifugal filter devices provide fast ultrafiltration, with the capability for high concentration factors and easy concentrate recovery from dilute and complex sample matrices. The vertical design and available membrane surface area provide fast sample processing, high sample recovery (typically greater than 90% of dilute starting solution), and the capability for 80-fold concentration. Typical processing time is 10 to 20 minutes. Solute polarization and subsequent fouling of the membrane are minimized by the vertical design, and a physical deadstop in the filter device prevents spinning to dryness and potential sample loss. The concentrate is collected from the filter device sample reservoir using a pipettor, while the ultrafiltrate is collected in the provided centrifuge tube. The device can be spun in a swinging-bucket (for optimal performance) or fixed-angle rotor. Amicon® Ultra-4 10K devices are supplied non-sterile and are for single use only.

Intended Use


Amicon® Ultra-4 centrifugal filter devices are disposable, single-use ultrafiltration filter devices for *in vitro* diagnostic use for processing biological solutions in the 1-4 mL volume range. They are intended for concentrating biological samples such as serum, urine, cerebrospinal fluid and other body fluids prior to *in vitro* diagnostic analysis. Device intended to be single-use and used by a laboratory professional.

Applications

- Concentration of biological samples containing antigens, antibodies, enzymes, nucleic acids (DNA/RNA samples, either single- or double-stranded), microorganisms, column eluates, and purified samples
- Purification of macromolecular components found in tissue culture extracts and cell lysates, removal of primer, linkers, or molecular labels from a reaction mix, and protein removal prior to HPLC
- Desalting, buffer exchange, or diafiltration

Materials Supplied

The Amicon® Ultra-4 10K device is supplied with a cap, a filter device, and a centrifuge tube.

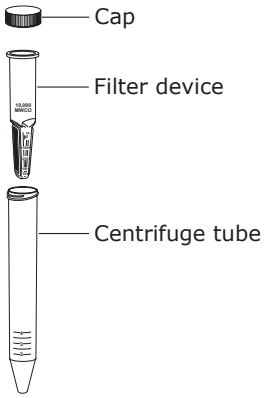Required Equipment

- Centrifuge with swinging-bucket (preferred) or fixed-angle rotor with wells/carriers that can accommodate 17 mm × 124 mm 15 mL conical-bottomed tubes

 CAUTION: To avoid damage to the device during centrifugation, check clearance before spinning.
- Pipettor with 200 microliter (µL) tip for concentrate recovery

Suitability

Preliminary recovery and retention studies are suggested to ensure suitability for intended use. See the "How to Quantify Recoveries" section.

Storage and Stability

Refer to product label for storage conditions and shelf life.

Rinsing Before Use

The ultrafiltration membranes in Amicon® Ultra-4 10K devices contain trace amounts of glycerine. If this material interferes with analysis, rinse the device with buffer or Milli-Q® water before use. If interference continues, rinse with 0.1 N NaOH followed by a second spin of buffer or Milli-Q® water.

CAUTION: Do not allow the membrane in Amicon® Ultra filter devices to dry out once wet. If you are not using the device immediately after rinsing, leave fluid on the membrane until the device is used.

DOMINIQUE DUTSCHER SAS

How to Use Amicon® Ultra-4 Centrifugal Filter Devices

1. Add up to 4 mL of sample (3.5 mL if using a 23° fixed-angle rotor) to the Amicon® Ultra filter device.
2. Place capped filter device into centrifuge rotor (swinging-bucket preferred); counterbalance with a similar device. Check Centrifuge clearance before use.
3. **When using a swinging-bucket rotor,** spin the device at 4,000 × g maximum for approximately 10–20 minutes.
4. **When using a fixed-angle rotor,** orient the device with the membrane panel facing up and spin at 7,500 × g maximum for approximately 10–20 minutes.


 NOTE: Refer to Figure 1 and Table 1 for typical spin times.
5. To recover the concentrated solute, insert a pipettor into the bottom of the filter device and withdraw the sample using a side-to-side sweeping motion to ensure total recovery. The ultrafiltrate can be stored in the centrifuge tube.

 NOTE: For optimal recovery, remove concentrated sample immediately after centrifugation.

 NOTE: Filter may not function correctly if allowed to dry out after wetting.

Desalting or Diafiltration

Desalting, buffer exchange, or diafiltration are important methods for removing salts or solvents in solutions containing biomolecules. The removal of salts or the exchange of buffers can be accomplished in the Amicon® Ultra-4 device by concentrating the sample, then reconstituting the concentrate to the original sample volume with any desired solvent. The process of "washing out" can be repeated until the concentration of the contaminating microsolute has been sufficiently reduced. See example below.

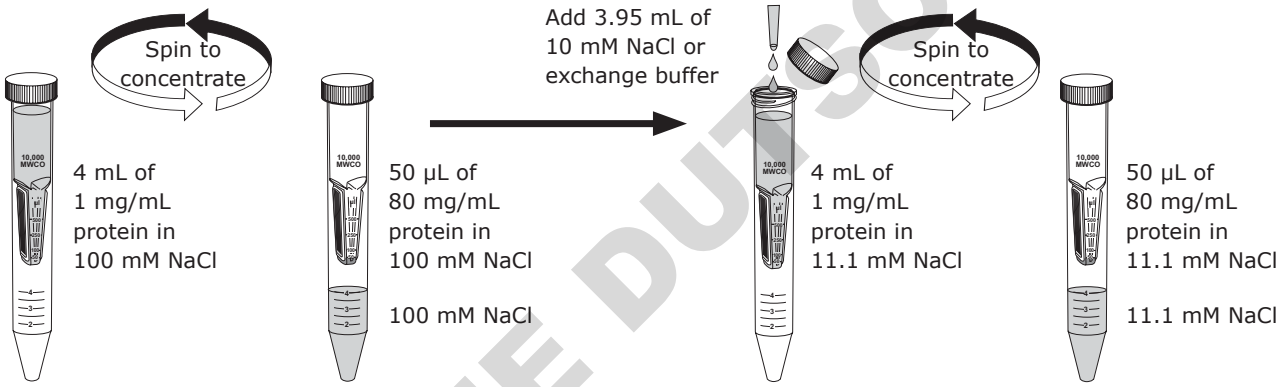
DOMINIQUE DUTSCHER SAS

Performance

Flow Rate

Factors affecting flow rate include sample concentration, starting volume, chemical nature of solute, relative centrifugal force, centrifuge rotor angle, membrane type, and temperature. Figure 1 and Table 1 can be used to estimate the time required to achieve a given volume of filtrate or concentrate. A typical spin time for a 4 mL sample is approximately 10 to 20 minutes. While most of the sample is filtered in the first 10 to 15 minutes of centrifugation, the lowest concentrate volume (30–75 µL) is reached after spinning for 15 to 20 minutes.

Figure 1. Typical Filtrate Volume vs. Spin Time

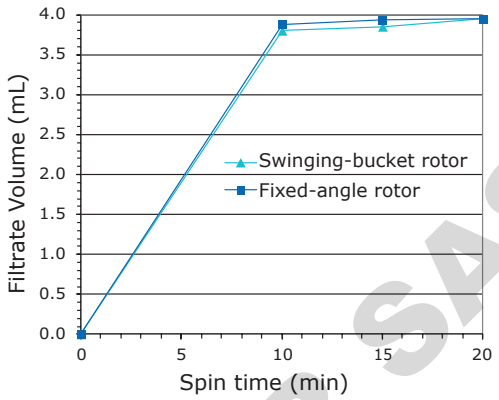Spin Conditions: Swinging-bucket rotor, 4,000 × g, or fixed-angle rotor, 7,500 × g, room temperature, 4 mL starting volume.
Protein marker used: Cytochrome c, n=6.

Table 1. Typical Concentrate Volume vs. Spin Time

	Concentrate volume (µL)	
Spin time (min)	**Swinging-bucket rotor 4,000 × g**	**35° Fixed-angle rotor 7,500 × g**
10	176	97
15	76	54
20	58	35

Spin Conditions: Room temperature, 4 mL starting volume.
Protein marker used: Cytochrome c, n=6 (mean value of 3 device lots).
Shaded volumes were used for the calculation of protein recovery in Table 3.

Protein Retention and Concentrate Recovery

The membranes used in Amicon® Ultra devices are characterized by a molecular weight cutoff (MWCO); that is, their ability to retain molecules above a specified molecular weight. Solutes with molecular weights close to the MWCO may be only partially retained. Membrane retention depends on the solute's molecular size and shape. For most applications, molecular weight is a convenient parameter to use in assessing retention characteristics. For best results, use a membrane with a MWCO at least two times smaller than the molecular weight of the protein solute that one intends to concentrate. Refer to Table 2.

Table 2. Typical Retention of Protein Markers

Marker/Concentration	Molecular Weight	Device MWCO	% Retention Swinging-bucket	% Retention Fixed-angle	Spin Time (min)
α-Chymotrypsinogen (1 mg/mL)	25,000	10K	>95	>95	15
Cytochrome c (0.25 mg/mL)	12,400		>95	>95	15
Vitamin B-12 (0.2 mg/mL)	1,350		<15	<15	15

Spin Conditions: Swinging-bucket rotor, 4,000 × g, or 35° fixed-angle rotor, 7,500 × g, 4 mL starting volume, room temperature, n=6 (mean value of 3 device lots).

Factors that determine sample recovery include the nature of the protein solute relative to the device MWCO chosen, starting concentration, and concentration factor. Table 3 provides typical recoveries for the Amicon® Ultra-4 10K device.

Table 3. Typical Concentrate Recovery

			Concentrate Volume (µL)		Concentration Factor (x)		Concentrate Recovery (%)	
Marker/Concentration	**Device MWCO**	**Spin Time (min)**	**Swinging-bucket**	**Fixed-angle**	**Swinging-bucket**	**Fixed-angle**	**Swinging-bucket**	**Fixed-angle**
Cytochrome c (0.25 mg/mL)	10K	15	76	54	52.3	76.6	97.3	98.5

Spin Conditions: Swinging-bucket rotor, 4,000 × g, or 35° fixed-angle rotor, 7,500 × g, 4 mL starting volume, room temperature, n=6 (mean value of 3 device lots). The shaded volumes were taken from Table 1.

DOMINIQUE DUTSCHER SAS

Maximizing Sample Recovery

Low sample recovery in the concentrate may be due to adsorptive losses, over-concentration, or passage of sample through the membrane.

- Adsorptive losses depend upon solute concentration, its hydrophobic nature, temperature and time of contact with filter device surfaces, sample composition, and pH. To minimize losses, remove concentrated samples immediately after centrifugal spin.
- If the starting sample concentration is high, monitor the centrifugation process in order to avoid over-concentration of the sample. Over-concentration can lead to precipitation and potential sample loss.
- If the sample appears to be passing through the membrane, choose a lower MWCO Amicon® Ultra-4 device.

How to Quantify Recoveries

Calculate total recovery, percent concentrate recovery, and percent filtrate recovery using the method below. The procedure provides a close approximation of recoveries for solutions having concentrations up to roughly 20 mg/mL.

NOTE: Appropriate assay techniques include absorption spectrophotometry, radioimmunoassay, refractive index, and conductivity.

Direct Weighing Procedure

The density of most dilute proteins is nearly equal to the density of water (i.e., 1 g/mL). Using this property, the concentrate and filtrate volumes can be quantified by weighing them and converting the units from grams to milliliters. This technique is valid only for solutions with concentrations of approximately 20 mg/mL or less.

- Before use, separately weigh the empty filter device, the centrifuge tube, and an empty tube for concentrate collection.
- Fill filter device with solution and reweigh.
- Assemble device and centrifuge per instructions.
- Collect the concentrate with a pipettor and dispense it into the preweighed concentrate collection tube.
- Remove the device from the centrifuge tube and weigh the centrifuge tube and concentrate collection tube.
- Subtract weight of empty device/tubes to calculate weights of starting material, filtrate, and concentrate.
- Assay the starting material, filtrate, and concentrate to determine solute concentration.
- Calculate recoveries using the weight/volume data and the measured concentrations as follows:

$$\%\ \text{concentrate recovery} = 100 \times \frac{W_c \times C_c}{W_o \times C_o}$$

$$\%\ \text{filtrate recovery} = 100 \times \frac{W_f \times C_f}{W_o \times C_o}$$

$$\%\ \text{total recovery} = \%\ \text{concentrate recovery} + \%\ \text{filtrate recovery}$$

W_c = total weight of concentrate before assay
W_o = weight of original starting material
W_f = weight of filtrate

C_c = concentrate concentration
C_o = original starting material concentration
C_f = filtrate concentration

Specifications

Maximum initial sample volume	
Swinging-bucket and fixed-angle rotors (45° and 35°)	4.0 mL
Fixed-angle rotor (23°)	3.5 mL
Typical final concentrate volume	50–100 µL
Maximum relative centrifugal force	
Swinging-bucket rotor	4,000 × g
Fixed-angle rotor	7,500 × g
Active membrane area	3 cm^2

Dimensions	
Filter device in tube (capped)	
Length: 124 mm (4.9 in.)	Diameter: 17.3 mm (0.7 in.)
Filter device	
Length: 73.4 mm (2.9 in.)	Diameter: 17.2 mm (0.7 in.)
Materials of construction	
Filter device	Copolymer styrene/butadiene
Membrane	Ultracel® low binding regenerated cellulose
Filtrate tube	Polypropylene
Filtrate cap	Polypropylene
Filtrate liner	Polyethylene
Filter	MWCO 10K

DOMINIQUE DUTSCHER SAS

Chemical Compatibility

Amicon® Ultra centrifugal devices are intended for use with biological fluids and aqueous solutions. Before use, check the sample for chemical compatibility with the device. Go to SigmaAldrich.com/FilterChemicalCompatibility for more information.

Symbol Definitions

Symbol	Definition	Symbol	Definition
IVD	*In vitro* diagnostic medical device		Do not use if package is damaged and consult instructions for use
REF	Catalogue number		Date of manufacture
LOT	Batch code		Manufacturer
eIFU	Consult electronic instructions for use		Importer
WWW	Download product documentation online	CE	CE conformity marking
NON STERILE	Non-sterile	UK CA	UK Conformity Assessed
	Do not reuse	CH REP	Swiss authorized representative
	Use-by date	UDI	Unique device identifier
	Temperature limit		

Product Ordering

Purchase products online at SigmaAldrich.com.

Description	Qty/Pk	Catalogue Number
Amicon® Ultra-4	8	UFC801008D
	24	UFC801024D
	96	UFC801096D

Notice

We provide information and advice to our customers on application technologies and regulatory matters to the best of our knowledge and ability, but without obligation or liability. Existing laws and regulations are to be observed in all cases by our customers. This also applies in respect to any rights of third parties. Our information and advice do not relieve our customers of their own responsibility for checking the suitability of our products for the envisaged purpose.

Collection and Disposal

All samples must be clearly labelled. Suitable instruments must be used for obtaining and preparing samples.

NOTE: Follow precautions for disposal of items contaminated with potentially infectious or hazardous biological material according to all applicable international, federal, state, and local regulations.

Technical Assistance

Visit the tech service page on our web site at SigmaAldrich.com/techservice.

Any serious incident of this device should be reported to manufacturer and competent authority of country where user is established.

Standard Warranty

The applicable warranty for the products listed in this publication may be found at SigmaAldrich.com/terms.

DOMINIQUE DUTSCHER SAS

Revision History

2021-OCT	• IFU PR05779 Date of Issue OCT 2021 - Replaced PR05139. • IFU, Packaging Damage and Temperature Limit symbols added. • Chemical Compatibility and Ordering Information linked to website. • Disposal and Complaints information added. • UK Responsible Person and UKCA symbol information added
2024-MAR	• Figure 1 corrected • "Materials of construction": Filtrate cap corrected text to polypropelene
2024-OCT	• Added to Symbol Definitions table: CH-REP, Importer, UDI • Added UKCA symbol on title page

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Made in Ireland

Merck Millipore Ltd.
Tullagreen, Carrigtwohill,
Co. Cork, Ireland
T45 KD29

an affiliate of Merck KGaA,
Darmstadt, Germany

UK Responsible Person:

Millipore (UK) Ltd.
Kirkton Campus
Fleming Road,
Livingston, UK
EH54 7BN

Merck, Millipore, Amicon, Milli-Q, Ultracel, and Sigma-Aldrich are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.
© 2001-2025 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved.

The life science business of Merck operates as MilliporeSigma in the U.S. and Canada.