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Medicom SAS
Boulevard de la chanterie | 49124 Saint Barthélemy d'Anjou - France
Tel +33 0 241 963 434 | sales.support.pro.eu@medicom.com | www.medicom-eu.com

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SafeMask® Fit+Pro

Instructions for use of respiratory masks with and without expiratory valve FFP2 NR and FFP3 NR for single use.

Description: SafeMask® Fit+Pro meets the requirement EN 149:2001+A1:2009 on filtering half masks to protect against particles. They provide effective respiratory protection for users in industries where workers will be exposed to solid (dust) particles and / or non volatile liquid particles. The product also meets the requirements of Personal Protective Equipment Regulation (EU) 2016/425.

Taking into account the normative indications of conformity and assessment: - FFP2 NR and FFP3 NR disposable masks are category III personal protective equipment (PPE) (death or irreversible health damage), which can be used by industry professionals.

FFP Type	With or without Valve	Reference
FFP2 NR	Without valve	M62110
	With valve	M62110V
FFP3 NR	Without valve	M63110
	With valve	M63110V

Type-examination (Module B) by: INSPEC International B.V. (NB 2849) Beechavenue 54-62, 1119 PW, Schiphol-Rijk, Netherlands
Conformity to type (Module D) by: BSI Group The Netherlands B.V. (NB 2797) Say Building, John M. Keynesplein 9, 1066 EP, Amsterdam, Netherlands

EU Declaration of conformity : The internet address where the EU declaration can be accessed is: **www.medicom-eu.com.**

Warning and Use Limitations:

Always be sure that the product is:
- Suitable for the application;
- Fitted correctly;
- Worn during all period of exposure;
- Replaced when necessary.

Proper selection, training, and use are essential in order for the product to help protect the wearer from certain airborne contaminants. Failure to follow all instructions for use and/or failure to properly wear the face mask during all periods of exposure may adversely affect the wearer's health.

This mask may only be used by people who are trained to use it.

This mask does not provide any protection against gases, vapours or fumes.
This half mask does not supply any oxygen. Do not use in atmospheres which do not contain at least 17% oxygen. Do not use in an explosive atmosphere.
This mask does not require any maintenance. It has been designed for one single continuous work session, i.e. for 8 hours as a maximum (one shift).

The mask must be used on intact skin only.
It is unlikely that requirements for leakage will be achieved if facial hair passes under the face seal;
Leave the premises if:
Dizziness or other distress occurs
Breathing becomes difficult

Dispose of and replace the mask if it becomes damaged, breathing resistance becomes excessive or at the end of the shift.


This device is for single use only (**see Fig. 7**). The reuse of a used device generates risks of contamination of the wearer or the surrounding people due to the degradation of the performance of the device.

The manufacturer declines all liability for damage resulting from misuse of this product or use not in accordance with the following instructions.

.

Donning Instructions
Before fitting the device, ensure hands are clean. All respirators components should be inspected for damage prior to use.

1. Carefully stretch out the headband using both hands
2. Unfold the mask (**see Fig 1**)
3. Preform the nasal bar by slightly folding it at the center
4. Position the lower band around your neck
5. Place the mask below your chin with the nose clip pointing up (**See Fig 2**)
6. Hold the upper band with your hand and place over the back of your head above your ears (**See Fig 3**).
7. Secure the headband and adjust if necessary with the toggle (**See Fig 4**).
8. Adjust the nose clip and make it fit your nose as much as possible (**See Fig 5**).

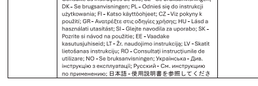Cover the front of the respirator with both hands being careful not to disturb the fit of the respirator.

Unvalved respirator – Exhale sharply;
Valved respirator – Inhale sharply;

If air leaks around the nose, re-adjust the noseclip to eliminate leakage. Repeat the above fit test.

If air leaks at the respirator edges, work the straps back along the sides of the head to eliminate leakage and tighten elastic with the toggle. Repeat the above fit check.

Proper selection, training, and use are essential in order for the product to help protect the wearer from certain airborne contaminants. Failure to follow all instructions for use and/or failure to properly wear the face mask during all periods of exposure may adversely affect the wearer's health.

This mask may only be used by people who are trained to use it.

This mask does not provide any protection against gases, vapours or fumes.
This half mask does not supply any oxygen. Do not use in atmospheres which do not contain at least 17% oxygen. Do not use in an explosive atmosphere.
This mask does not require any maintenance. It has been designed for one single continuous work session, i.e. for 8 hours as a maximum (one shift).

The mask must be used on intact skin only.
It is unlikely that requirements for leakage will be achieved if facial hair passes under the face seal;
Leave the premises if:
Dizziness or other distress occurs
Breathing becomes difficult

Dispose of and replace the mask if it becomes damaged, breathing resistance becomes excessive or at the end of the shift.

This device is for single use only (**see Fig. 7**). The reuse of a used device generates risks of contamination of the wearer or the surrounding people due to the degradation of the performance of the device.

The manufacturer declines all liability for damage resulting from misuse of this product or use not in accordance with the following instructions.

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