ClearLine® PapPen

Pen for Immunostaining

The PapPen forms a water-repellent barrier around the specimen and enables simple, clean and economical immunostaining



Water-repellent
Creates a clear water-repellent barrier

Precise application
Easy and targeted marking around the specimen


Saves reagents
Reduces the consumption of valuable antisera

Time-saving
Clean handling without wiping steps

Versatile
Suitable for immunostaining, ASD staining and enzyme staining

Application (Paraffin Sections)

1. Deparaffinize the section with xylene, then wash with ethanol. Remove excess liquid with tissue paper
2. Draw a circle around the specimen on the slide with the PapPen and allow to dry for 1–2 minutes at room temperature
3. Then incubate in phosphate buffer solution (PBS or TBS) for about 5 minutes

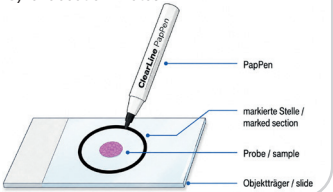

DOMINIQUE D