Betriebsanleitung
Operating instructions
Notice d'instructions

Hei-SHAKE Salto Core

heidolph
research made easy

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Schnittstellenbefehle RS232 (ASCII)

Generic commands

Send	Receive	Description
"IDENT\r\n"	"< Name > < MaterialNr > < SerialNr > < FertigungsAuftrag > \r\n"	Read identification numbers
"SW_VERS\r\n"	"MC X1.Y1.Z1:X2.Y2.Z2, HMI A1.B1.C1:A2. B2.C2\r\n"	Read all software versions
"STATUS\r\n"	"STATUS X\r\n"	Read actual status
"IN_SP_ALL\r\n"	"SP_3 = X, SP_10 = Y, SP_11 = Z\r\n"	Read all setpoints
"IN_PV_ALL\r\n"	"PV_5 = X, PV_12 = Y S\r\n"	Read all actual values
"CC_ON\r\n"	"CC_ON\r\n"	Set connection control
"CC_OFF\r\n"	"CC_OFF\r\n"	Reset connection control
"RESET\r\n"	"RESET X\r\n"	Reset communication

Speed control commands

Send	Receive	Description
"OUT_SP_3 X\r\n"	"OUT_SP_3 X\r\n"	Set target speed (rpm)
"OUT_SP_12 X\r\n"	"OUT_SP_12 X\r\n"	Set maximum motor moment
"IN_PV_5\r\n"	"IN_PV_5 X\r\n"	Get actual speed (rpm)
"START_2\r\n"	"START_2\r\n"	Start movement
"STOP_2\r\n"	"STOP_2\r\n"	Stop movement
"OUT_MODE_10 0\r\n"	"OUT_MODE_10 0\r\n"	Select speed controller mode

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Contents

Introduction

About this document ... 32
- Typographic conventions ... 32
- Copyright protection ... 32
- Safety instructions in the official languages of the European Union ... 32

General product information

Directives applied, product certification ... 33
Copyright notice ... 33
California Residents ... 33
Residual risk ... 33
Intended use ... 33
Compliant use ... 33
Reasonably foreseeable misuse ... 34
Transportation ... 34
Storage ... 34
Acclimatization ... 34
Permissible ambient conditions ... 34
Energy saving potential and cost efficiency ... 35

Safety

General safety instructions ... 35
Electrical safety ... 35
Interfaces ... 35
Data security ... 36
Operational safety ... 36
Occupational safety ... 36
Electrostatic charging of liquid or powdered media ... 36
Personal protective equipment (PPE) ... 37
Environmental protection ... 37
Biohazard ... 37
Other regulations ... 37

Device description

Mechanical design– control panel– terminals ... 38
- Emergency off button ... 39

Positioning the device ... 39
Power supply ... 40
Mounting laboratory bottles ... 40

DOMINIQUE DUTSCHE

Operation 42
Switching the device on/off 42
Remote control status indicator 42
Menu navigation 42
Settings 43
Starting/deceleration ramp 43
Automatic restart 43
Rotation mode 44
Operating modes 45
Continuous operation 45
Timer mode 46
Remote mode 47
Device activation via PC 47
Device activation via server-based software 47
Operation behavior in remote operation 47
Troubleshooting 48
Warning codes 48
Error codes 49

Appendix

Technical specifications 50
RS232 (ASCII) interface commands 51
Scope of delivery 52
Accessories 52
Device service 52
General cleaning instructions 52
Repairs - Return of equipment 53
Maintenance 53
Disposal 53
Warranty statement 54
Contact details 54
Certificate of decontamination 55

DOMINIQUE DUTSCHE

About this document

This operating manual describes the features and operation of overhead shakers, Hei-SHAKE Salto Core type. The operating instructions are an integral part of the described device!

Typographic conventions

This document uses the following symbols, signal words, and highlights:

	Description
	Warning symbols in combination with a signal word indicate dangers: **DANGER** Indicates an immediate dangerous situation. Failure to respect the indications will result in death or serious injury. **WARNING** Indicates a potential danger. Failure to respect the indications will result in serious injuries. **CAUTION** Indicates a potential hazard. If not avoided, damage to property and minor to moderate injuries can occur.
	Mandatory signs are used to indicate important and useful information on handling a product. This information is used to ensure operational safety and to maintain the value of the product.
[…]	Square brackets identify labels of control elements on the device as well as labels and entries in software masks and user interfaces.
→	The arrow indicates specific instructions to be followed to ensure operational safety when handling the product.

Copyright protection

This document is protected by copyright and is intended for use by the purchaser of the product only.

No transfer to third parties, reproduction in any form, including excerpts, and by any means, as well as utilization and/or disclosure of the contents is permitted without the prior written consent of Heidolph Scientific Products GmbH. Any violation is subject to compensation for damage.

Safety instructions in the official languages of the European Union

A summary of all safety instructions in the official languages of the European Union can be found in our Safety Guide for the product group Hei-SHAKE (shakers and mixers, Ref. 01-005-006-99). This document is available for download on our homepage in the most up-to-date version.

DOWNLOAD DEUTSCHE

Directives applied, product certification

CE marking

The device meets all requirements of the following directives:

- Machinery Directive 2006/42/EC
- EMC Directive 2004/108/EC
- RoHS 2 directive 2011/65/EU, including amendment 2015/863

Copyright notice

The software implemented in this product is protected by copyright laws. The rights holder is Heidolph Scientific Products GmbH, Walpersdorfer Straße 12, 91126 Schwabach/Germany. Any open source software components in this software are excluded from our copyright. Further information is available in the service area on our website www.heidolph.com.

Violations of copyright (e.g. unauthorized use or modification of the software) may result in civil claims (e.g. omission, compensation for damages) and/or entail criminal penalties within the meaning of the legislation of the destination country.

California Residents

Important information for California residents regarding Prop 65. Please visit www.P65Warnings.ca.gov for more information.

Residual risk

The device has been designed and manufactured in accordance with the state-of-the-art standards at the time of development and the recognized safety regulations. During mounting and use, as well as during maintenance, repair and cleaning work, there are nevertheless certain residual risks associated with the described device.

These are identified and described at the appropriate points in this document.

Intended use

The described overhead shaker unit has been specially developed for the following tasks:

- Mixing
- Emulsifying
- Suspending
- Separating

The areas of application of the described device include chemical, biological and environmental-analytical laboratory and research applications, applications of basic research and comparable facilities.

Due to its design, the device in its delivery condition may only be used in analytical processes or in laboratory-like conditions in the food, cosmetics, and pharmaceutical industries as well as other comparable industries that manufacture products intended for consumption by humans or animals, or for use on humans or animals.

Any other use of this device is not considered as intended!

Compliant use

The compliance of each individual application must always be evaluated by the user. Any additional measures necessary to ensure compliance are always the responsibility of the user.

DOMINIQUE DUTSCHER

Reasonably foreseeable misuse

Additional measures may be necessary for use under conditions or for purposes deviating from the intended use, and/or specific guidelines and safety regulations must be observed. Corresponding requirements must be evaluated and implemented by the operator in each individual case.

Compliance with and implementation of all relevant directives and safety measures for the respective field of application is the sole responsibility of the operator.

All risks resulting from improper use are borne solely by the operator.

The device may only be operated by authorized and instructed personnel. Training and qualification of the operating personnel as well as ensuring that the device is handled responsibly are the sole responsibility of the operator!

Transportation

During transportation, avoid severe shocks and mechanical stresses that can cause damage to the device. Keep the original packaging in a dry and protected place for later use.

For safety reasons, the device must always be transported by two people (exception: variant Salto 4).

If possible, use the two guide rods on the device as holding points. When setting down the device, make sure to hold both thumbs above the guide rod as shown in the picture to avoid crushing.

Storage

Always store the product in its original packaging. To protect against damage and excessive material aging, store the device in an environment that is as dry, temperature-stable and dust-free as possible.

Acclimatization

After each transportation and after storage under critical climatic conditions (e.g., high temperature difference between inside and outside), allow the device to acclimatize at room temperature for at least two hours to prevent possible damage due to condensation before commissioning it in the place of use. If necessary, extend the acclimatization phase if the temperature differences are very high.

Make all supply connections only after the device has been acclimatized!

Permissible ambient conditions

The device is designed for indoor use only. The device is **NOT** suitable for outdoor use! The device is **NOT** suitable for use in potentially explosive areas!

When used in corrosive atmospheres, the service life of the device may be reduced depending on the concentration, duration and frequency of exposure.

DOMINIQUE DUTSCHER

Energy saving potential and cost efficiency

Even in standby mode, small amounts of energy are used to supply individual device functions. However, especially over longer periods of nonuse, this resource consumption adds up and causes avoidable costs.

Therefore, always switch off the device when not in use, especially for long periods of time (devices with main switch: Main switch in position [O]; devices without main switch: Unplug) to effectively avoid unnecessary energy consumption and the resulting costs.

General safety instructions

Before commissioning and using the device, familiarize yourself with all the safety regulations and occupational safety guidelines applicable at the place of use and observe them at all times.

Only operate the device if it is in faultless technical condition. In particular, ensure that there is no visible damage on the device itself and, where applicable, on connected devices or the supply connections.

If there is missing or misleading information on the device or regarding occupational safety, contact the responsible safety specialist or our technical service.

Only use the device in accordance with the regulations for intended use.

Electrical safety

Before connecting the device to the power supply, ensure that the voltage indicated on the rating plate matches the specifications of the local power utility company.

Ensure that the power supply circuit provided is protected by means of a residual-current device (RCD).

Always use the power supply cord provided with the device.

Prior to use, check that the device and the power supply cord are free of visible damage.

Make sure that the power plug can be reached directly at any time, in order to disconnect the device from the power supply at any time without delay.

Have repairs and/or maintenance work on the device carried out exclusively by an authorized electrician or by the technical service department of Heidolph Scientific Products GmbH.

Always disconnect the device from the power supply system before carrying out any maintenance, cleaning or repair work.

Interfaces

Safely separate extra low voltage inputs and outputs rated 25 V AC or 60 V DC or above according to EN 61140, or use double or reinforced insulation according to EN 60730-1 or DIN 60950-1.

Use only shielded connection cables. Connect the shield to the connector housing. Always cover unused interfaces (front and back) to prevent liquids from entering the device.

DOMINIQUE DUTSCHER

Data security

The user is responsible for ensuring data security when transferring data between the described device and other devices.

Only use secure networks for the data transfer and avoid use of critical infrastructure.

Only use high-quality shielded data cables for the data transfer.

For data transfer via USB, an industrial standard USB hub should be preferably used to ensure the most stable connection possible.

Operational safety

Operate the device under a closed ventilated fume hood when working with potentially hazardous substances (see EN 14175 and DIN 12924).

Do not make any unauthorized changes or modifications to the device!

Only use genuine spare parts and accessories, or those expressly approved by the manufacturer!

Rectify malfunctions or faults on the device immediately.

Switch off and disconnect the device from the power supply, preventing reconnection, if it is not possible to eliminate the malfunction or rectify the fault immediately.

Observe all other applicable regulations such as laboratory and workplace guidelines, recognized safety technology rules and special local regulations.

Occupational safety

Always use the prescribed personal protective equipment (PPE) such as protective clothing, safety goggles, protective gloves, safety shoes, etc.

Do not operate any other devices in the immediate vicinity of the device …

- which can generate electromagnetic fields in the frequency range between 9×10^3 Hz to 3×10^{11} Hz,
- which generate emission or radiation sources in the frequency range 3×10^{11} Hz to 3×10^{15} Hz (in the optical spectral range wavelengths from 1,000 µm to 0.1 µm),
- which generate ultrasonic or ionizing waves.

Do not process any substances that could release energy in an uncontrolled manner (e.g. self-ignition).

Do not process substances in which the energy input through mixing poses a danger.

Wipe off any liquid that may have spilled on the device immediately.

Electrostatic charging of liquid or powdered media

Liquid or powdered media as well as the laboratory vessels used can become dangerously electrostatically charged due to the shaking movement or when emptying or cleaning the vessel. The intensity of the charge generated depends on the properties of the medium, the flow rate, the process parameters as well as the size, geometry and material of the laboratory vessels (source: TRGS 727, Section 4).

The evaluation of a potential explosion hazard due to electrostatic charging and the definition of suitable safety measures are always the responsibility of the operator and are always based on the substances to be processed and the individual application.

DOMINIQUE DUTSCHER

Personal protective equipment (PPE)

The operating company must determine and provide the necessary PPE, depending on the respective application and the media and chemicals used.

The required instruction of the personnel is solely within the operating company's responsibility.

Environmental protection

When processing environmentally hazardous substances, take appropriate measures to avoid hazards to the environment.

The evaluation of appropriate measures such as the marking of a hazardous area, their implementation, and the training of the relevant personnel is the sole responsibility of the operator!

Biohazard

When processing biohazardous substances, take appropriate measures to prevent hazards to persons and the environment, including:

- Instruction of the personnel regarding the necessary safety measures.
- Provision of personal protective equipment (PPE) and instruction of the personnel in its use.
- Marking of the device with the biohazard warning symbol.

The evaluation of appropriate measures such as the marking of a hazardous area, their implementation, and the training of the relevant personnel is the sole responsibility of the operator!

Other regulations

In addition to the notes and instructions in this document, observe all other applicable regulations such as laboratory and workplace guidelines, hazardous substances ordinances, recognized rules of safety engineering and occupational medicine as well as particular local regulations!

Noncompliance will invalidate any warranty claims against Heidolph Scientific Products GmbH.

The operating company is solely liable for all damage resulting from unauthorized changes or modifications to the device, from the use of unapproved or non-genuine spare parts and accessories, or from disregarding the safety instructions and hazard warnings or the manufacturer's instructions!

DOMINIQUE DUTSCHER

Mechanical design– control panel– terminals

1 Adapter with clamping units for four, eight or twelve laboratory bottles.

2 Emergency off button

3 [Standby] button
Used to activate or switching the device into standby mode:
The LED assigned to the button lights up green briefly when the operating voltage is applied and changes then to permanent white: status indicator [device in standby mode].
Press the button once to activate the device: the LED changes to permanent green, status indicator [device in normal operation].
Press the button again to switch the device back to standby mode: LED switches back to permanent white.


4 Display of the operating parameters and messages.

5 Push-and-rotary control
Rotary function: select menu items, adjust operating parameters.
Push function: confirm entries.

6 [Menu] button: Use this button to switch between the Start screen and the Home screen with the various sub-menus.

7 [Timer] button: Pushbutton to activate the [program timer] function: The setpoint display changes to orange. In this operating state, the device timer can be programmed with the push-and-rotary control in a range from 1 min – 999 h 59 min. When the value is confirmed by pressing the button again, the assigned LED lights up.

8 [Start/Stop] button: Pushbutton to activate/deactivate the rotation function. The button LED lights up when the rotation function is activated.

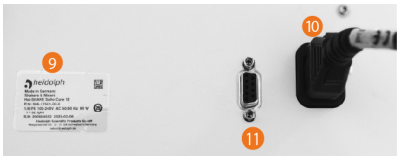9 Rating plate: Device and manufacturer data, serial number (example).

10 Appliance inlet: Power supply connection

11 RS232 interface (Sub-D, 9-pin)
The RS232 interface of the device can be used for maintenance and service purposes (e.g. firmware update) as well as for controlling the device and for reading out device data. Observe the instructions in section "Remote mode" on page 47 and "RS232 (ASCII) interface commands" on page 51.

DOMINIQUE ITSCHE

Emergency off button

By pressing the emergency off button, the rotation movement is stopped immediately, the device generates an acoustic warning signal and the LED ring of the push-and-turn control flashes red. At the same time, the error message [Error: 800A] is displayed.

To resume operation, the emergency off button must first be unlocked (turn the button slightly to the left). Then press the [STANDBY] button to briefly put the device into standby mode.

Press [Standby] again to put the device into the previously set operating mode.

WARNING

Before switching the device back on after an emergency off, it is to correct the cause(s) of the disturbance. See also section "Troubleshooting" on page 48.

Positioning the device

CAUTION

Improper positioning of the device may result in direct and indirect damage to property due to falling and/or spilling of fluids!

The surface must be even, non-slip, proper and sufficiently stable.

The device must be freely accessible at all times.

Make sure that a sufficient safety distance to moving components can be maintained on all sides.

Observe the generally valid rules for the safe positioning of work equipment.

DANGER

When using glass vessels, there is an increased risk of injury and contamination due to cracking or glass breakage (drop, falling containers during assembly, etc.).

Take care when loading the appliance with glass vessels!

The professional mounting and positioning of the device including the supplied accessories is within the sole responsibility of the operator!

Heidolph Scientific Products GmbH shall under no circumstances be liable for any personal injury or damage to property resulting from failure to observe the instructions for the correct installation and positioning of the device.

DOMINIQUE DUTSCHER

Power supply

DANGER

Observe the notes and instructions in the section on "Electrical safety" on page 35.

Power supply cord with locking system

Only use the supplied power supply cord with locking system!

Connect the cable coupling of the cord [2] to the IEC appliance inlet [1] on the back of the device. Make sure that the locking system engages [3] correctly.

Connect the power supply cord to a properly secured mains socket-outlet.

Disconnecting the power supply cord

Disconnect the power supply cord from the mains socket-outlet.

Press in the locking lever on the cable coupling and disconnect the power supply cord from the IEC appliance plug on the back of the device.

When the mains voltage is applied, the device is in standby mode (the [STANDBY] button LED lights up white).

Mounting laboratory bottles

In the delivery state, the device is equipped with attachments for 2,000 ml bottles. To mount the bottles, the adapter must be moved to an ergonomic position:

- → Make sure that the device is 'on' (press Standby, the button LED lights up green).
- → Use the push-and-rotary control to position the adapter appropriately at minimum rotation speed: Turn the push-and-rotary control to the left (negative speed range) to move the adapter counter-clockwise (looking from the left at the motor/control housing); turn the push-knob to the right (positive speed range) to move the adapter clockwise (looking from the left at the motor/control housing).
- → Turn the device off (press Standby again, button LED lights up white) to avoid an undesired start when mounting the bottles.
- → Loosen the clamping lever [1] (rotation 90 ° counterclockwise).


- → Pull the knurled handle of the retaining bolt upwards on the released clamping lever [2] and insert the bottle.
 - Make sure that the pressure plate is positioned in the middle of the bottle neck.

DOMINIQUE DUTSCHER

- → Press the retaining bolt downwards to secure the bottle with the pressure plate.
- → Lock the released clamping lever again (rotation 90 ° clockwise).
- → If necessary, the position of the clamping lever can be adjusted: To do this, pull it outwards, turn the clamping lever to a suitable position and release the clamping lever again.
 - Repeat the procedure if necessary until the bottle can be clamped correctly.

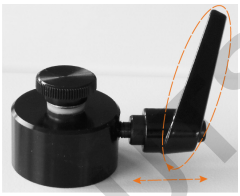Special extensions are required for smaller bottles. These are available as optional accessories:

- Mounting kit M, P/N: 549-26000-00, see www.heidolph.com

To convert the system, proceed as follows:

- → Use the supplied hex wrench to remove the mounted attachments.

- → Screw the extension for the 500 ml or 1000 ml bottle onto the retaining bolt.

- → Using the supplied hexagon wrench, screw the pressure plate to the extension on the retaining bolt.

During initial installation, the threads of the mounts for the pressure plates are glued into the threaded holes of the housing with Loctite 2400 to ensure a permanently secure hold.

We recommend securing both the mounts and the extensions in a similar manner during each installation as part of a retrofit.

DOMINIQUE DUTSCHER


Switching the device on/off

When the mains voltage is applied, the device is in standby mode (the [STANDBY] button LED lights up white).

Briefly press [Standby] to put the device on. The button LED changes to permanent green.

The manufacturer logo appears briefly on the display, followed by the start screen with the following values:

- Setpoint and actual speed [1]
- Timer value [2]
- On time [3] in [hhh:mm:ss]

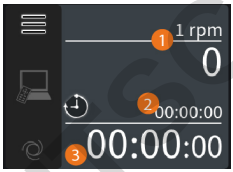
Press [Standby] again to put the device off. "Shutting down" is briefly displayed on the display, and the LED of the button changes back to white.

Unplug the device to de-energize it! When de-energized, the [STANDBY] button LED goes out.

Remote control status indicator

When the device is connected to a remote controller via the RS 232 interface, the remote status is indicated on the home screen as follows:

- Remote icon gray: device is controlled locally, no input signal.
- Remote icon magenta-coloured: device operates in remote mode and is controlled by input signals. In this operating state, the LED ring of the push-and-rotary control also lights up magenta-coloured.

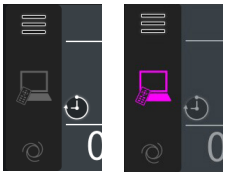
Menu navigation

Press the [Menu] button [4] on the control panel below the display to display and to hide the menu bar.

The [Menu] key returns you from any view to the Home screen.

The individual menu items are selected by turning the push-and-turn control. The selected menu item is highlighted. Press the push-and-rotary control to open the selected menu.

- [5]: [Home] button, opens the home screen from each view.
- [6]: Menu item [Device settings], for more information see section "Settings" on page 43.
- [7]: Menu item [Error list], opens the list of all recorded errors; for more information, see section "Troubleshooting" on page 48.
- [8]: Menu item [Info], touch this function button to open the info menu with the following options: Serial number, HMI version, micro-controller version, QR code for further information and technical documents.

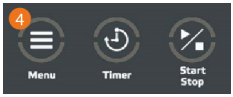
DOMINIQUE DUTSCHER

Settings

Expand the menu bar as described in the section "Menu navigation" on page 42 and open the [Device Settings] menu, [6].

In the device settings, you can define the starting and braking behavior of the device [9] and activate/deactivate the [Automatic restart] function [10], see the following sections.

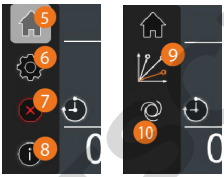Starting/deceleration ramp

This parameter defines the starting and deceleration behavior of the device when the rotation movement is switched on or off.

→ In the [Device Settings] menu, use the push-and-rotary control to highlight [Starting/Deceleration Ramp], [1] and press the push-turn knob.

The following options are available:

- [Fast]: The speed is increased or decreased as quickly as possible to the set value.
- [Medium]: The speed is increased or decreased at normal speed to the set value.
- [Slow]: The speed is gradually increased or decreased to the set value.

→ Use the push-and-rotary control to highlight the desired mode and confirm the selection by pressing the push-and-rotary control. The icon associated with the selected mode is displayed in the menu bar area.

The setting remains active even if the power supply is interrupted!

In general, the actual start-up behavior depends on other parameters such as the total load and the load distribution. The ideal setting must be determined and defined by the user in individual cases. It is not possible to give a general manufacturer recommendation for the different scenarios due to the application-specific parameters.

Automatic restart

The [Automatic restart] function can be activated and deactivated via this menu item.

When the [Automatic restart] function is activated, the device switches on automatically upon return of the power supply after a power interruption and the speed is augmented to the last set value(exception: emergency off). If the function remains deactivated, the device must be activated by the operator from standby mode when the mains voltage returns.

WARNING

Before activating the function, take appropriate measures to make sure that, in case of a possible unmonitored restart, no hazard for persons or risk of damage to property can occur.

→ In the [Device Settings] menu, use the push-and-rotary control to highlight [Automatic restart].

→ Enable or disable the function by pressing the push-and-rotary control.

DOMINOUE DUTSCHER

Status display in the menu bar area

- Icon [Rotation] white: Function disabled. Press the push-and-rotary control to enable the function. The icon changes to orange.
- Icon [Rotation] orange: Function enabled. Press the push-and-rotary control to disable the function. The icon changes to white.

In the as-delivered condition, the [Automatic restart] function is disabled and has to be enabled manually.

The function can be activated and deactivated at any time, i.e. without interrupting a running process. The evaluation of potential risks by a possible unmonitored automatic restart as well as of appropriate security measures is the sole responsibility of the user!

Operating behavior when the [Automatic restart] function is activated

When the power is restored after a power interruption, a warning signal is emitted, simultaneously the LED ring of the push-and-rotary control flashes yellow and the warning symbol for the automatic restart appears for five seconds on the display (see figure on the right). After this warning time has expired, the motor speed is increased again to the previously set value! The timer is automatically restarted after each mains power interruption.

Rotation mode

WARNING

- Always place individual bottles in the center. Always ensure that the bottles are evenly distributed.
- Use suitable attachments for the used bottles.
- Before starting the rotation movement, make sure that all bottles are properly closed and securely clamped.
- Lower the rotation speed when the base vibrates during operation.
- Always use the necessary and appropriate personal protective equipment!
- Never touch rotating parts with your hands. Do not bend over rotating parts.
- Wear tight-fitting work clothes, remove jewelry before starting work, tie long hair under a hair net or under a suitable head covering

DOMINIQUE DUTSCHER

→ Use the push-and-rotary control to set the required rotation speed (range: -32 to +32 rpm).
 - In the range -32 to -01 rpm, the adapter rotates clockwise (looking at the motor/control housing from the left)
 - In the range +01 to +32 rpm, the adapter turns counter-clockwise (looking from the left at the engine/control housing)

→ Press [Start/Stop] to turn on the rotation movement.
 - The rotation movement starts at the lowest rotation speed and is increased to the set speed according to the set starting ramp (see section "Starting/deceleration ramp" on page 43).
 - During startup, ensure the stability of the device until the set rotation speed is reached.
 - Rotational speed and direction can be continuously adjusted at any time during operation.

Operating modes

The device can be used in continuous operation or with a timer (timer mode). Regardless of the selected operating mode, observe the following specific safety instructions:


WARNING

- Gradually increase the rotation speed and observe the device behavior. Lower the rotation speed or reduce the total load when the device starts to vibrate during operation.
- Always ensure that the bottles are evenly distributed.
- Pay particular attention to even distribution of the weight load when bottles of different size and/or differently filled bottles are mounted at the same time!
- Before switching on the device, make sure that all bottles are properly fixed.
- Always use the necessary and appropriate personal protective equipment!

Continuous operation

In continuous operation, the rotational movement is started and stopped manually:

→ To start or stop the rotational movement, press the [Start/Stop] button on the control panel of the device.

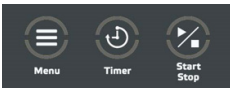
In continuous mode, the LED of the [Start/Stop] button and the LED ring of the push-and-rotary control light up white.

In continuous operation without a timer, the elapsed time since the start of the process is shown on the display. If the maximum run time of 999 hours, 59 minutes and 59 seconds is exceeded, the value [>999] is displayed permanently until the process is switched off.

DOMINIQUE DUTSCHE

Timer mode

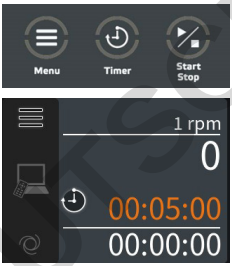In timer mode, a specific period of time can be set for a continuous rotational movement.

- → To activate the timer function, briefly press the [Timer] button. The [Timer] setpoint changes to orange. In this state a timer value in the range of 1 min – 999 h 59 min can be defined by means of the push-and-rotary control. To do this, turn the push-and-rotary control until the desired value is displayed (in the example, five minutes.
- → Briefly press the push-and-rotary control or the [Timer] button to accept the input value. The new value must be confirmed within three seconds, otherwise the entry will be discarded!

After the input has been accepted, the set value display of the timer is reduced again and displayed in white font.

Status indicator

As soon as the timer is programmed and set, the timer icon will appear in orange on the display.

Starting rotational mode with timer

To start or stop the rotational movement, press the [Start/Stop] button on the control panel of the device.

The timer starts automatically, the LED of the [Start/Stop] and [Timer] buttons and the LED ring of the push-and-rotary control light up permanently white.

The remaining time until the expiration of the programmed time period is shown as [hh:mm:ss].

When the rotational mode stops, the timer is also stopped and reset.

When the timer expires, the LED ring of the push-and-rotary control flashes white and the device generates an acoustic signal. Press the push-and-rotary control to reset. In this state, a new process can be started with the last set value in timer mode.

End timer mode

To end the timer function, press and hold the [Timer] button for two seconds. The LED is switched off.

DOMINIQUE DUTSCHER

Remote mode

WARNING

Secure the device in remote mode with a clearly visible warning sign and, if necessary, take further ambient-specific protective measures that protect against damage to property and injuries in the event of unexpected/unmonitored starting of the device.

To activate the remote access of the device, a message of type “SetTargetSpeed“ or “SetMotorState“ has to be sent. Afterwards, the remote symbol appears on the display and remote commands are executed normally.

The transmission of the messages is made via ASCII strings. The protocol is active after the start. By a defined ADIP “Connect“ command, you can switch to the ADIP mode.

RS232 interface parameters

- 9600 Baud
- Parity: none
- Data: 8 Bit
- Stop: 1 Bit

Command termination

All ASCII strings have to be finished with \r\n (0x0D 0x0A):

- …`\r\n`

For further information on interface handling, please contact our technical service, see "Contact details" on page 54.

Device activation via PC

The device can be addressed via the integrated RS232 interface via an external control. The data interface enables activation of the device and recording of the process data in remote mode. A prerequisite for this is a suitable terminal program or software solution such as the Hei-PROCESSING Solutions from Heidolph Scientific Products GmbH

Device activation via server-based software

The device can be addressed via the integrated RS232 interface via a server-based software. If required, contact an authorized dealer or our technical service, see section "Contact details" on page 54.

Operation behavior in remote operation

In remote mode, the individual menus of the device control can be opened as described and the parameters/values can be read out, but not changed.

It is also not possible to increase or decrease the rotational speed intensity on the device (push-and-rotary control without function).

If necessary, remote operation can be terminated by the operator by pressing the [Standby] button.

DOMINIQUE DUTSCHER

Troubleshooting

Failure	Cause	Possible remedy
Display remains dark after switching on	No mains voltage	Check the power supply cord for damages/correct connection
	Power supply cord defective	Check building fuses
	Display defective	Contact Heidolph Sales or Service
LED of the [Start/Stop] button does not light up when the function is activated	Button LED defective	Contact Heidolph Sales or Service
Rotational movement does not stop/start	Motor overheat protection has tripped	Wait approx. 20 minutes and reduce the load
	Electrical fault (no engine noise)	Contact Heidolph Sales or Service
	Mechanical defect (engine noise audible)	Contact Heidolph Sales or Service
	Emergency off button pressed	Unlock when the cause of the fault is eliminated
	Timer activated	Check timer settings

In case of recurring errors please contact the responsible sales department or our technical service (see section "Contact details" on page 54).

Warning codes

Code	Description	Device behavior, remedy
W100B	Permissible PCB temperature exceeded	The device changes to a safe operating state - let it cool down, restart the process.
W100A	Permissible motor temperature exceeded	
W100C W100D W100E W100F	Overvoltage	Device changes to a safe operating state - restart the process.
W1010 W1011 W1012 W1013 W1014	Undervoltage	Device changes to a safe operating state - restart the process.
W1015	Engine blockage	Device changes to a safe operating state - restart the process.
W1017	Overload	The device changes to a safe operating state - mount bottles evenly, restart the process.

DOMINIQUE DUTSCHER

Error codes

Code	Description	Device behavior, remedy
E8001	Internal memory error	Device changes to safe operating state - restart the device.
E8006 E8007 E8008	Motor error	
E8009	Undervoltage	
----	HMI communication error	Unit switches off, contact technical service.
E800A	Emergency stop has been triggered	Unit changes to safe operating state - unlock emergency off button, switch the unit off and on again with the [Standby] button, restart the process.

- Pending error messages are displayed in orange font.
- Acknowledged error messages are displayed in white font.

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Technical specifications

Hei-SHAKE Salto Core	Variant 4	Variant 8	Variant 12
Dimensions (W × H × D, in mm)	536 × 520 × 458	816 × 520 × 458	1,095 × 520 × 458
Weight	approx. 20 kg	approx. 26 kg	approx. 31 kg
Interface(s)	1 × RS232 (Sub-D, 9 pole)		
Permissible load	max. 30 kg		
Capacity	4 bottles	8 bottles	12 bottles
Movement type/device type	Overhead shaker		
Speed range	1 – 32 rpm (clockwise/counter-clockwise rotation)		
Drive	Gear stepper motor		
Degree of protection	IP21		
Protection class	I ⏚		
Overvoltage category	II		
Degree of pollution	2		
Acoustic pressure	< 60 dB (A), in accordance with IEC 61010		
Rated voltage	1/N/PE 100-240V AC 50/60 Hz		
Power input	max. 90 W		
Operating and storage temperature	5 °C – 31 °C at up to 80 % rel. humidity 32 °C – 40 °C at up to 50 % rel. humidity (decreasing linearly)		
Installation altitude	up to 2,000 m asl		

DOMINIQUE DUTSCHE

RS232 (ASCII) interface commands

Generic commands

Send	Receive	Description
"IDENT\r\n"	"< Name > < Material No. > < Serial No. > < Production Order > \r\n"	Read identification numbers
"SW_VERS\r\n"	"MC X1.Y1.Z1:X2.Y2.Z2, HMI A1.B1.C1:A2.B2.C2\R\N"	Read all software versions
"STATUS\r\n"	"STATUS X\r\n"	Read actual status
"IN_SP_ALL\r\n"	"SP_3 = X, SP_10 = Y, SP_11 = Z\R\N"	Read all setpoints
"IN_PV_ALL\r\n"	"PV_5 = X, PV_12 = Y S\r\n"	Read all actual values
"CC_ON\r\n"	"CC_ON\r\n"	Set connection control
"CC_OFF\r\n"	"CC_OFF\r\n"	Reset connection control
"RESET\r\n"	"RESET X\r\n"	Reset Communication

Speed control commands

Send	Receive	Description
"OUT_SP_3 X\r\n"	"OUT_SP_3 X\r\n"	Set target speed (rpm)
"OUT_SP_12 X\r\n"	"OUT_SP_12 X\r\n"	Set maximum motor torque
"IN_PV_5\r\n"	"IN_PV_5 X\r\n"	Get actual speed (rpm)
"START_2\r\n"	"START_2\r\n"	Start movement
"STOP_2\r\n"	"STOP_2\r\n"	Stop movement
"OUT_MODE_10 0\r\n"	"OUT_MODE_10 0\r\n"	Select Speed Controller Mode

DOMINIQUE DUTSCHE

Scope of delivery

Item	Quantity	Product no.
Hei-SHAKE Salto Core 4, or	1	546-17301-00
Hei-SHAKE Salto Core 8, or	1	546-17401-00
Hei-SHAKE Salto Core 12	1	546-17501-00
Power supply cord	1	country specific
Warranty registration	1	01-006-002-78

Accessories

information on available accessories can be found on our website at www.heidolph.com.

Device service

When carrying out service work on the device (cleaning, maintenance, repair), observe the general instructions and safety information described in this section.

DANGER

Switch the device's main switch off and disconnect it from the power supply before carrying out maintenance work, cleaning, or repairs.

When cleaning, avoid the penetration of liquids.

Before replacing the fuses, switch off the device and disconnect the power supply cord.

Always replace the two device fuses in pairs with original manufacturer fuses. Further information on available accessories can be found on our website at www.heidolph.com.

After fuse replacement, check the device for a safe condition according to IEC 61010-1.

General cleaning instructions

Wipe all surfaces and the control panel with a damp cloth if necessary. Persistent contamination can be removed with mild soapy water.

CAUTION

Clean the device's surfaces with a soft, lint-free and only slightly moistened cloth.

Never use any aggressive or abrasive cleaning agents or aids.

DOMINIQUE DUTSCHE

Repairs - Return of equipment

Repairs to the device may only be carried out by authorized experts! Unauthorized repairs during the warranty period will result in the loss of the warranty claim. Regardless of the warranty claim, the owner is solely liable for damage caused by unauthorized repairs.

→ In case of repair and before returning your device, contact our technical service at the following e-mail address:
 - service@heidolph.de.

→ In your message, please provide us with the following information in addition to an error description:
 - Item number
 - Serial no.

The required data can be found on the rating plate of the device.

A service representative will contact you as soon as possible to agree on the next steps.

Include the completed certificate of decontamination with every device return, see "Certificate of decontamination" on page 55.

Maintenance

The device contains no user-serviceable components. If necessary (in the event of abnormal operating behavior such as excessive noise or heat generation, for example) contact our technical service.

Disposal

- When disposing of the device, observe the provisions of the WEEE Directive 2012/19/EU and its transposition into national law in the country of use.
- When disposing of portable batteries, observe the provisions of the European Battery Directive 2013/56/EU and their implementation in national law in the country of use.
- Check the device and all components for residues of substances that are hazardous to health, the environment and biohazardous before disposal.
- Properly remove and dispose of residues of substances that are hazardous to health, the environment and biohazardous!

Registration in accordance with the Electrical and Electronic Equipment Act

Heidolph Scientific Products GmbH, headquartered in 91126 Schwabach/Germany, address Walpersdorfer Str. 12, is registered under the number DE 50705753 in the electronic waste equipment register of the foundation ear (www.ear-system.de).

DOMINIQUE DUTSCHE

Warranty statement

Heidolph Scientific Products GmbH provides a three-year warranty against material and manufacturing defects.

Glass and wear parts, transportation damage, and damage resulting from improper handling or non-intended use of the product are excluded from the warranty.

The warranty period for registered products begins on the date of purchase. Register the product with the enclosed warranty card or on our homepage www.heidolph.com.

For non-registered products, the warranty period begins with the date of the serial production (to be determined by the serial number).

In the event of material or manufacturing defects, the product will either be repaired or replaced free of charge within the warranty period.

Contact details

Heidolph Scientific Products GmbH

Technical Service
Walpersdorfer Str. 12
D-91126 Schwabach/Germany

Email: service@heidolph.de

Representations

To find your local Heidolph distributor please visit www.heidolph.com

DOMINIQUE DUTSCHE

Certificate of decontamination

Enclose the certificate of decontamination, duly completed, with your device return. Submissions without a certificate of decontamination cannot be processed!

CERTIFICATE OF DECONTAMINATION

IN CASE OF RETURNS

heidolph
research made easy

Please fill in the required fields.
Note: The sender must package the goods properly and appropriately for transport.

Heidolph Scientific Products GmbH
Walpersdorfer Straße 12
91126 Schwabach
Phone: +49 (0) 9122 9920-380
Fax: +49 (0) 9122 9920-19
E-Mail: service@heidolph.de

SENDER

Name _____ First name _____

Company/institution _____ Department _____

_____ Workgroup _____

Address _____

ZC/City _____

Country _____ Phone _____

Email _____

DEVICE DETAILS

Article number _____ **Serial no.** _____

Ticket number _____

Reason for sending in _____

Has the device been cleaned, decontaminated/disinfected? **Yes** **No** (Please mark as applicable)

If yes, which measures were carried out?

Does this device pose a risk to people and/or the environment due to the processing of substances that are hazardous to health, the environment and/or are biohazardous? **Yes** **No** (Please mark as applicable)

If yes, with which substances did the device come into contact?

LEGALLY BINDING DECLARATION

The principal/consignor is aware that they are liable to the agent/consignee for losses or damage incurred due to incomplete and incorrect information.

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hereby declare, that the product (from serial number on 200066855) designated below is in compliance with the basic requirements of all applicable EU-directives stated below with regard to design, type of model sold and manufactured by us. This certificate will be invalid if the product is modified without the prior written consent and agreement of the manufacturer.

Reax 2	541-21001-XX
Reax top	541-10000-XX
Reax control	541-11000-XX
Reax 20/4	541-20004-XX
Reax 20/8	541-20008-XX
Reax 20/12	541-20012-XX
Titramax 100	544-11100-XX, 544-11200-XX
Duomax 1030	543-32105-XX, 543-32205-XX, 543-32210-XX
Polymax 1040	543-42105-XX, 543-42110-XX, 543-42205-XX, 543-42210-XX
Polymax 2040	542-40005-XX, 542-40010-XX
Multi Reax	545-10000-XX
Vibramax 100	544-21200-XX
Rotamax 120	544-41200-XX
Titramax 1000	544-12200-XX
Unimax 1010	543-12310-XX
Unimax 2010	542-10020-XX
Hei SHAKE Orbital Core	546-11XXX-XX
Hei SHAKE Salto Core	546-17XXX-XX
Titramax 101	544-11300-XX
Vibramax 110	544-31200-XX
Inkubator 1000	549-90010-XX
Promax 1020	543-22332-XX
Promax 2020	542-20020-XX

Maschinenrichtlinie / Machinery Directive 2006/42/EG
EMV-Richtlinie / Electromagnetic Compatibility Directive 2014/30/EU
Delegierte (RoHS-) Richtlinie / Delegated (RoHS) Directive 2015/863/EU
Angewandte (harmonisierte) Normen / (Harmonized) Standards applied:
EN ISO 12100:2010, EN 61326-1:2013, EN 61010-1:2010, EN IEC 63000:2018

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UKCA Declaration of Conformity

Declaration of Conformity

In accordance with UK Government guidance

This declaration is issued under the sole responsibility of the manufacturer,
Heidolph Scientific Products GmbH
Walpersdorfer Straße 12
91126 Schwabach / Germany

Product: Shaker & Mixer
Model:

Reax top /control	541-1XXXX-XX
Reax 20/4…12	541-200XX-XX
Reax 2	541-21001-XX
Duomax 1030	543-32XXX-XX
Unimax X010	54X-1XXX0-XX
Promax X020	54X-2XXXX-XX
Polymax X040	54X-4XXXX-XX
Titramax 1XXX	544-1XXXX-XX
Vibramax 1XX	544-X1200-XX
Rotamax 120	544-41200-XX
Multi Reax	545-10000-XX
Inkubator 1000	549-90010-XX
Hei SHAKE Orbital Core	546-11XXX-XX
Hei SHAKE Salto Core	546-17XXX-XX

Description:
Shaker & Mixer with different working directions

The object of the declaration described above is in conformity with the relevant UK Statutory Instruments (and their amendments):

2008 No. 1597	*The Supply of Machinery (Safety) Regulations*
2008 2016 No. 1091	*The Electromagnetic Compatibility Regulations 2016*
2012 No. 3032	*The Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012*

and complies with the following technical standards :
EN ISO 12100:2010, EN 61326-1:2013, EN 61010-1:2010, EN IEC 63000:2018

UK Authorised Representative (for authorities only)
ProductIP (UK) Ltd.
8. Northumberland Av.
London WC2N 5BY

Signed for and on behalf of Heidolph Scientific Products GmbH
Walpersdorfer Straße 12, 91126 Schwabach / Germany

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With this declaration, we confirm (according to current knowledge) that all sold laboratory devices and componentes by Heidolph Scientific Products GmbH fulfill the requirements of the EU directive 2011/65/EU (RoHS) and the extension 2015/863.
All devices are compatible with the requirements of the RoHS Directive for the following materials:
Max. 0,01% of weight in homogenes material for cadmium and max. 0,1% of the weight in homogeneous material for Lead, Mercury, Hexavalent chromium, Polybrominated biphenyls (PBB), Polybrominated diphenyl ethers (PBDE), Bis(2-ethylhexyl) phthalate (DEHP), Butyl benzyl phthalate (BBP), Dibutyl phthalate (DBP), Diisobutyl phthalate (DIBP).
In the case of individual assemblies, maximum concentrations maybe exceeded within the permissible exceptions of the Directive.

Schwabach, 04.September.2024

Andreas Hahn
Managing Director

Jörg Ziel
Quality Manager

COMMON DEUTSCHE

heidolph
research made easy

China RoHS DECLARATION OF CONFORMITY

Heidolph Scientific Products GmbH has made reasonable efforts to ensure that hazardous materials and substances may not be used in its products.

In order to determine the concentration of hazardous substances in all homogeneous materials of the subassemblies, a “Product Conformity Assessment” (PCA) procedure was performed. As defined in GB/T 26572 the “Maximum Concentration Value” limits (MCV) apply to these restricted substances:

- Lead (Pb): 0.1%
- Mercury (Hg): 0.1%
- Cadmium (Cd): 0.01%
- Hexavalent chromium (Cr(VI)): 0.1%
- Polybrominated biphenlys (PBB): 0.1%
- Polybrominated diphenyl ether (PBDE): 0.1%

Environmental Friendly Use Period (EFUP)

EFUP defines the period in years during which the hazardous substances contained in electrical and electronic products will not leak or mutate under normal operating conditions. During normal use by the user such electrical and electronic products will not result in serious environmental pollution, cause serious bodily injury or damage to the user’s assets.

The Environmental Friendly Use Period for Heidolph Scientific Products GmbH products is 25 years.

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This table is created according to SJ/T 11364-2014

MATERIAL CONTENT DECLARATION FOR Heidolph Scientific Products GmbH PRODUCTS

部件名称 Part name	有毒有害物质或元素 Hazardous substances: 铅 Pb	汞 Hg	镉 Cd	六价铬 Cr(VI)	多溴联苯 PBB	多溴二苯醚 PBDE	环保期限标识 EFUP
包装 Packaging	O	O	O	O	O	O	
塑料外壳/组件 Plastic housing / parts	O	O	O	O	O	O	
电池 Battery	O	O	O	O	O	O	
玻璃 Glass	O	O	O	O	O	O	
电子电气组件 Electrical and electronic parts	X	X	X	O	O	O	
控制器/测量设备 Controller / measuring device	X	O	X	O	O	O	
金属外壳/组件 Metal housing /parts	X	O	O	O	O	O	
电机 Motor	X	O	O	O	O	O	
配件 Accessories	X	O	O	O	O	O	

DOMINIQUE DUTSCHE

China RoHS Certification

heidolph
research made easy

注释：此表格适用于所有产品。以上列出的原件或组件不一定都属于所附产品的组成。
Note: Table applies to all products. Some of the components or parts listed above may not be part of the enclosed product.

O: 表示该有毒有害物质在该部件所有均质材料中的含量均在GB/T 26572规定的限量要求以下。
O: Indicates that the above mentioned hazardous substance contained in all homogeneous materials of the part is below the required limit as defined in GB/T 26572.

X: 表示该有毒有害物质至少在该部件某一均质材料中的含量超出GB/T 26572规定的限量要求。
X: Indicates that the above mentioned hazardous substance contained in at least one of the homogeneous materials of this part is above the required limit as defined in GB/T 26572.

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Apart from the disclosures in the above table, the subassemblies are not intentionally manufactured or formulated with lead (Pb), mercury (Hg), cadmium (Cd), hexavalent chromium (CrVI), polybrominated biphenyls (PBB), and polybrominated diphenyl ethers (PBDE).

Products manufactured by Heidolph Scientific Products GmbH may enter into further devices or can be used together with other appliances .
With these products and appliances in particular, Heidolph Scientific Products GmbH will not take responsibility for the EFUP of those products and appliances.

Schwabach, 04.September.2024

Andreas Hahn
Managing Director

Jörg Ziel
Quality Manager

DOMINIQUE DUTSCHE

© Heidolph Scientific Products GmbH

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Subject to change without notice. The printed version of this document is not regularly updated. The latest issue of this document can be found by visiting our homepage.

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