Vari-Plates[™]

Vari-Strip™ PCR Tube Strips in Frames

Vari-Strip^{TM} low profile PCR tube strips are moulded from virgin polypropylene under clean-room conditions. They can be used as a stand-alone product, see page 44, or with a frame for ease of handling.

Vari-StripTM is available in clear polypropylene for standard PCR techniques. For fluorescent detection, like qPCR, Vari-StripTM is available with white well tubes which give the highest sensitivity and the highest consistency as all the fluorescence signal is reflected back to the detector.

The VaristripTM can be inserted into a frame and there are three frames available, one for the Roche Lightcycler® 480, one recommended for use with ABi Fast Block cyclers and a standard 96 well skirted frame, to create Vari-PlatesTM.

Vari-Plates[™] are also available pre-loaded for maximum convenience.



Vari-Plate™



4ti-0950-F Vari-Strip™ 480/96 frame for Roche LightCycler 480




4ti-0910-F Vari-Strip™ FastBlock frame for ABi Fast Block



4ti-0757-F Universal fully skirted plate

Code	Description	Quantity
4ti-0753	Vari-Strip™ low profile 8 strip tubes, white	120 tube strips
4ti-0753/C	Vari-Strip™ low profile 8 strip tubes, clear	120 tube strips
4ti-0754	Vari-Strip™ low profile 8 strip tubes, white,	120 tube strips
	with flat optical caps	+ 120 cap strips
4ti-0754/C	Vari-Strip™ low profile 8 strip tubes, clear,	120 tube strips
	with flat optical caps	+ 120 cap strips
4ti-0950-F	Vari-Strip™ 480/96 frame, clear	5 frames
4ti-0910-F	Vari-Strip™ FastBlock frame, clear	5 frames
4ti-0757-F	Vari-Strip™ Skirted frame, white	5 frames
4ti-0753/757	Vari-Strip™ , white, pre-loaded in fully skirted frame	50 plates
4ti-0753/C/757	Vari-Strip™, clear, pre-loaded in fully skirted frame	
4ti-0753/950	Vari-Strip [™] , white, pre-loaded in semi-skirted LC480 fran	ne 50 plates
4ti-0753/C/910	Vari-Strip™ , clear, pre-loaded in semi-skirted	
	ABi Fast Block frame	50 plates

For thermal cycler compatibility, see 4ti-0950, 4ti-0910 and 4ti-0960 in the plate compatibility table

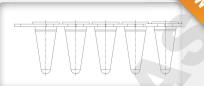





Colour options



Other options

- Available barcoded
- Available with flat optical caps



MAX = 200μl WORKING = 150μl

- Use as individual strips or in a frame, for easier handling
- · Low profile
- · Barcoding of frames possible
- Polypropylene thin-walled tubes for optimum PCR results
- · Suitable for cap sealing
- RNase, DNase, Human Genomic DNA Free

Also see

PLATE COMPATIBILITY TABLE (PAGE 46-48)
BARCODING (PAGE 54)
STRIP CAPS (PAGE 45)



SCAN QR CODE TO CHECK PLATE COMPATIBILITY FOR YOUR INSTRUMENT OR SEE PAGE 46-48

www.4ti.co.uk/plate-compatibility

+44 (0)1306 884885 www.4ti.co.uk | 43