

Restriction Enzyme BsaW I



Cat.# FG-BsaWI Size 250 units Conc. 10 units/µl

Store at -20℃


Supplied with: 10X FastGene® Buffer IV (FG-REB4)

10X FastGene® FastCut Buffer (FG-REBHF)

6X DNA Loading Buffer

Sterile water

Recognition site

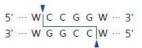

For Research Use Only. Not for use in diagnostic procedures.

ISO9001

Dilution buffer:

FastGene® Diluent A

Heat Inactivation

BsaW I can be inactivated at 80°C for 20 min.

Methylation sensitivity

dam methylation: Not sensitive dcm methylation: Not sensitive CpG methylation: Not sensitive

Prolonged incubation

A minimum amount of enzyme required to digest 1 μg substrate DNA for 16 hr; 0.13 U.

Relative activity in FastGene® Buffers

 FastGene® Buffer I:
 50%

 FastGene® Buffer II:
 100%

 FastGene® Buffer III:
 100%

 FastGene® Buffer IV:
 100%

 FastGene® FastCut Buffer:
 100%

Note

It is not affected by *dam, dcm,* or mammalian CpG methylation. It is an isoschizomer of Bet I.

Source: Bacillus stearothermophilus W1718

Reaction conditions

1X FastGene® Buffer IV 60℃ 1X FastGene® FastCut Buffer, 60℃

FastGene® FastCut Buffer

FastGene® restriction enzyme can cut substrate DNA in 5-15 with FastGene® FastCut Buffer.

1X FastGene® Buffer IV

20 mM Tris-acetate (pH 7.9 at 25°C) 50 mM potassium acetate 10 mM magnesium acetate 100 μg/ml BSA

Unit definition

One unit is defined as the amount of enzyme required for complete digestion of 1 μ g bacteriophage λ at 37°C for 1 hr in 50 μ l reaction mixtures.

Quality control

- Unit definition assay
- Overdigestion assay
- Endonuclease assay
- Extreme pure assay

Standard reaction condition

- Normal protocol

Component	Final Conc.	Volume
Substrate DNA	1 μg	Χ μΙ
10X FastGene® Buffer IV	1 X	5 μΙ
BsaW I	10 unit	1 μΙ
Sterile water		up to 50 μl

→ Incubate at 60°C for 1 hr

Fast protocol

Component	Final Conc.	Volume
Substrate DNA	1 μg	Χ μΙ
10X FastGene® FastCut Buffer	1 X	5 μΙ
BsaW I	10 unit	1 μΙ
Sterile water		up to 50 μl

→ Incubate at 60°C for 15 min

 $\ensuremath{\mathbb{X}}$ We recommend 5-10 units of enzyme per μg DNA and 10-20 units for genomic DNA in a 1 h digest..