



# AUTOMATION MADE EASY

STERILE

FREE OF  
detectable  
DNase

FREE OF  
detectable  
human DNA

FREE OF  
detectable  
RNase

 non-  
cytotoxic

 non-  
pyrogenic

## THINCERT®

## 96 WELL HTS INSERT

### For high-throughput transport studies and co-cultures


The new ThinCert® 96 well HTS insert is the optimal tool for scientists who want to generate tissue models (e.g., endothelia and epithelia) for transport, uptake, or co-culture studies in submerged as well as air-lift culture with an automation-friendly format.

By combining high membrane permeability and best possible transparency, Greiner Bio-One is offering a unique 96 well insert version, exclusively designed for applications that require both high diffusion rates as well as microscopic readouts. This is assured by a special pore arrangement and a high pore density of  $2 \times 10^7/\text{cm}^2$ .

In addition, a product version featuring a translucent membrane and a pore density of  $1 \times 10^8/\text{cm}^2$  complements the portfolio for experiments that need highest diffusion rates while imaging is not required.


#### KEY FACTS

- / Optimal for transport and permeability studies, air-lift cultures, co-cultures, and toxicity screening
- / 96 well system for high-throughput applications
- / Automation-friendly design
- / Transparent polycarbonate membrane with  $0.4 \mu\text{m}$  pore size
- / Reduced wicking effect
- / High membrane flatness for reproducible cell culture conditions

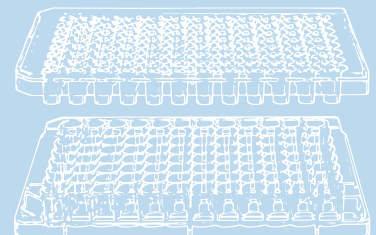

## PRODUCT OVERVIEW

FOR FURTHER INFORMATION AND/OR SAMPLE ORDERING  
PLEASE VISIT OUR WEBSITE OR CONTACT US.

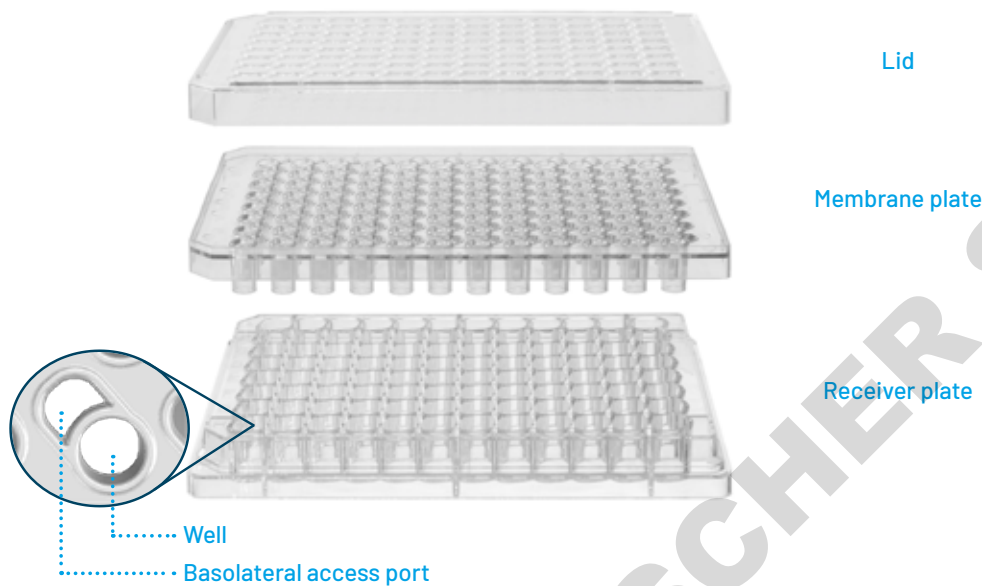

The two-component system comprises a 96 well porous membrane plate manufactured entirely from polycarbonate (PC) and a corresponding polystyrene (PS) receiver plate. The tissue culture treated PC membrane allows for an ideal exchange of nutrients and substances, creating in vivo-like cultivation conditions to facilitate optimal cell growth, monolayer formation, and tissue differentiation.

We have optimised the ports to ensure a convenient accessibility of the lower (basolateral) compartment with pipettes, automated liquid handling robots as well as electrodes for transepithelial electric resistance (TEER) measurements. The precise, well-centred fit of the membrane plate prevents wicking\* and provides stable and reproducible assay conditions.

\* Wicking: Undesired formation of a liquid bridge between upper and lower compartment due to capillary forces within the narrow space between membrane and receiver plate.

### ThinCert® 96 Well HTS Insert (Membrane plates and receiver plates)

Growth area: 14 mm<sup>2</sup>, Working volume (well of membrane plate): 15 - 160 µl,

Working volume (well of receiver plate): 120 - 300 µl, Lid: yes, condensation rings

Item no.	Pore density	Ø Pores	Optical features of membrane	Surface treatment	Sterile	Qty. inner / outer
655640	1 x 10 <sup>8</sup> /cm <sup>2</sup>	0.4 µm	translucent	TC	+	1 / 5
655641	2 x 10 <sup>7</sup> /cm <sup>2</sup>	0.4 µm	optimised transparency	TC	+	1 / 5

### Receiver plates for ThinCert® 96 Well HTS Insert

Working volume (well of receiver plate): 120 - 300 µl, Lid: yes, condensation rings

Item no.	Growth area	Surface treatment	Sterile	Qty. inner / outer
655169	-	non-treated	+	8 / 32
655167	53 mm <sup>2</sup>	TC	+	8 / 32