
NEUTRAX PF

ACID RINSE AGENT FOR MECHANICAL TREATMENT

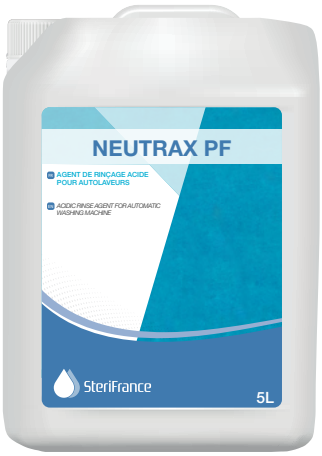

PRODUCT BENEFITS

- Dual descaling action: dissolves limescale deposits and prevents the formation of new deposits
- Neutralising action of detergents used in the wash phase
- Rinsing action: breaks down residual foam

www.sterifrance.com

SteriFrance

CE

NEUTRAX PF

ACID RINSE AGENT FOR MECHANICAL TREATMENT

INSTRUCTIONS FOR USE

ACID NEUTRALISATION IN LAVEUR INSTRUMENT DISINFECTOR

Setting the intermediate rinse phase of the LDI cycle :

- **Recommended concentration:** .... 0.3 to 0.6% depending on detergent used and water hardness
- **Recommended temperature:** ..... 40 to 50°C
- **Recommended contact time:** ..... 4 to 5 minutes

FEATURES

- Aspect : **Liquid**
- Colour : **Colorless**
- Pure pH : **2**

PRECAUTIONS OF USE

Refer to the Safety Data Sheet

COMPOSITION

- Organic acids
- Excipients

PACKAGING


❶ Box of 3 low x 5L.....réf. **1032015**

❷ Box of 4 x 5L.....réf. **1032020**

❸ Box of 2 x 10L.....réf. **1032022**





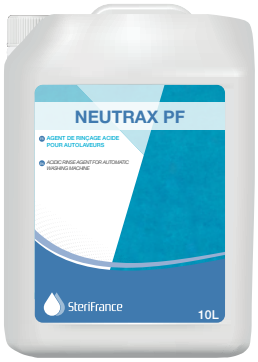

Product for professional use - please observe the precautions for use. Read the label and product information before use. Class 1 medical device.

Non-contractual document - Franklab - Printed in France

SteriFrance

FRANKLAB SAS – Zone Industrielle du Chemin Vert,
40 Rue du Chemin Vert 78610 Le Perray-en-Yvelines ■ France.
Tél. : **+ 33 (0)1 39 44 93 40** ■ Fax : **+ 33 (0)1 39 44 93 41**
contact@sterifrance.com ■ **www.sterifrance.com**

FT 260
v6_03.2026