Dutscher
TOGETHER, HELPING SCIENCE ADVANCE

Technical Data Sheet




Description

- Mikrozid® universal liquid, disinfecting detergent

Part number

- 727414

Pieces by sales : 1 L

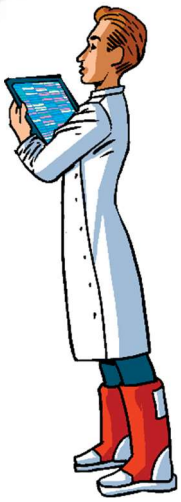

Description

- Uniquely versatile, low-alcohol formulation for the rapid cleaning and disinfection of surfaces & non-invasive medical devices.
- Alcohol-based rapid cleaning and disinfection of medical devices and surfaces when short contact times are required and where a material friendly disinfectant is needed.
- Apply mikrozid® universal liquid undiluted to objects and surfaces and ensure complete wetting and allow to take effect. Make sure that all visible soiling is removed before disinfection. Wipe up any excess with a non-dyed disposable cloth. schülke guarantees a shelf life of 1 month after the package has been opened.
- 100 g solution contains the following active substances: 17,4g Propan-2-ol, 12,6g Ethanol (94% w/w)
- Labelling according to Regulation (EC) No. 648/2004: <5% anionic surfactants
- Color: colourless
- Form: liquid
- Flash point:26 °C / Method : DIN 51755 Part 1
- Viscosity, dynamic: No data available
- Density:ca. 0.95 g/cm3
- pH:3 - 3.6 / 100 %
- The product is highly compatible with metals and synthetics. Please follow carefully any national fire and explosion protection guidelines when using alcohol-based disinfectants. Only treat small surfaces. The surface coverage depends on the room temperature and on the structure of the surface to be disinfected. More information avaible on request.