FR

Sonde 52 31

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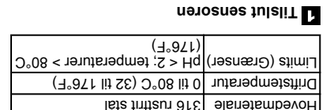
www.hach.com

© Hach Company/Hach Lange GmbH, 2016, 2017, 2021. All rights reserved.

Printed in France 02/2021, Edition 4

HACH

52 31

DOC022.98.90545

EN

52 31 probe

Additional information is available on the manufacturer's website.

Caution

Personal injury hazard. Glass components can break. Handle with care to prevent cuts.

Chemical exposure hazard. Refer to the current material safety data sheets (MSDS) for safety protocols.

Introduction

This puncture gel-filled electrode is used for pH measurements in semisolid samples (e.g., cheeses, meats, fish, fruits, bread doughs).

Technical information

Range	2 to 14 pH
Junction	Open, Ceramic
Reference	Ag/AgCl
Minimum depth	15 mm (0.6 in.)
Body material	316 stainless steel
Operating temp	0 to 80 °C (32 to 176 °F)
Limits	pH < 2; temperatures > 80 °C (176 °F)

1 Connect the probe

Do not turn or twist the probe connector.

2 Calibration

Refer to the meter documentation for the calibration procedure.

3 Measurement

Caution

The internal body is made from glass and can break. Do not bend the electrode. Do not put the electrode directly in thick samples or samples containing solid materials (e.g., bones). If necessary, first open a path in the sample with a stainless steel perforator.

Refer to the meter documentation for the mesaurement procedure.

Maintenance and storage

Cleaning: soak the probe in a cleaning solution. Rinse with deionized water.

Contamination	**Cleaning agent**
Proteins	Pepsin cleaning solution
Grease, oils, fats	Electrode cleaning solution
Limescale	0.1 N HCl solution

Hydration: soak the dry probe in a dilute HCl solution for several hours.

Storage: put the probe in the protector with KCl 3 M and glycerol solution.

Troubleshooting

Be sure that the correct probe type is used. Do not hold the tubes by hand to prevent temperature changes.

Warranty

6 months for manufacturing faults only. Damage from use is not covered.

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DOMINIODE DUTSCHER SAS