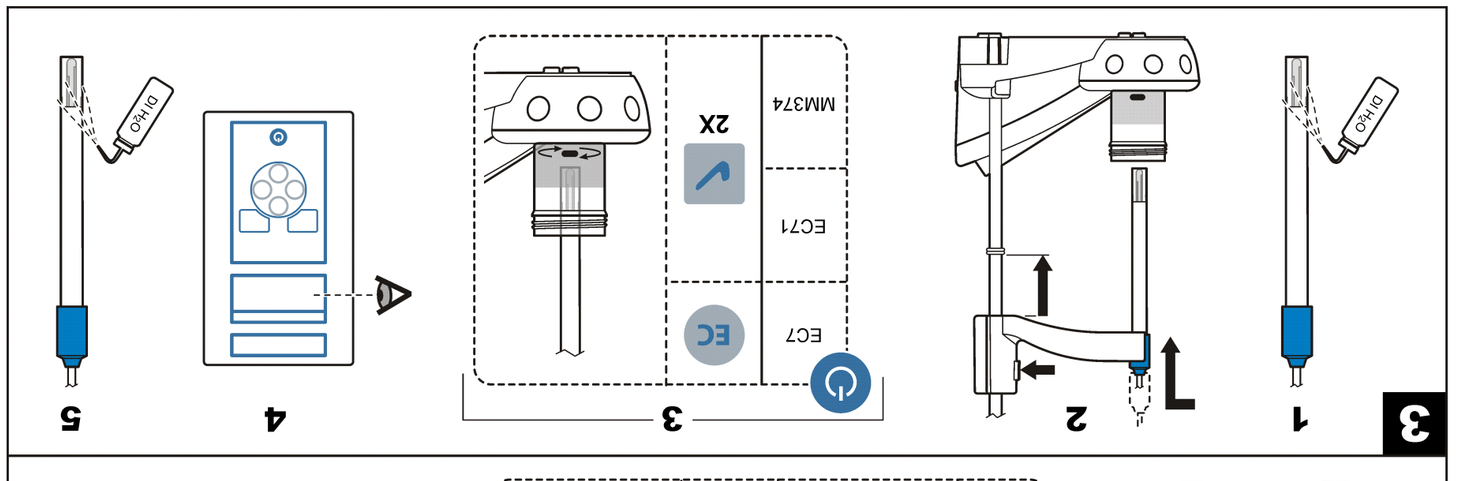
EN

50 70 probe

Refer to the meter manual for more information.

Caution

⚠ *Personal injury hazard. Glass components can break. Handle with care to prevent cuts.*

⚠ *Chemical exposure hazard. Refer to the current material safety data sheets (MSDS) for safety protocols.*

Introduction

This platinum cell with temperature sensor is used for conductivity measurements in general water samples.

Technical information

Range	0.2 µS/cm to 200 mS/cm
Constant	0.7 cm^{-1}
Electrode	Platinum (3x)
Temperature	Pt 1000
Minimum depth	25 mm (0.98 in.)
Operating temp	–30 to 80 °C (–22 to 176 °F)
Limits	Conductivity over 200 mS/cm; samples that contain HF

Preparation

Before first use or after long-term storage: put the electrode in ethanol for 15 s, then rinse with distilled water.

1 Connect the probe

2 Calibration

⚠ *Air bubbles under the probe tip when submerged can cause slow stabilization or error in measurement. If bubbles are present, gently shake the probe until bubbles are removed.*

3 Measurement

Maintenance and storage

To clean inorganic contamination, soak the cell in a dilute HCl solution. Rinse with deionized water.

To clean organic contamination, soak the cell in solution. Put in ethanol for 15 s, then rinse with deionized water.

Contamination	Cleaning agent
Grease, oils, fats	Electrode cleaning solution
Limescale	0.1 N HCl solution

Storage: use dry storage for the cell.

Troubleshooting

Be sure that the samples are suitable for the application.

Warranty

1 year for manufacturing faults only. Damage from use is not covered.

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© Hach Company/Hach Lange GmbH, 2011-2012. All rights reserved.