Chlorine, Total, Ultra Low Range

DOC316.53.01545

Fluorescence Method **Method 10316**

3 to 100 µg/L Cl_2 **Reagent**

Scope and application: For testing total chlorine in clean water applications. Accurate to less than 20 ppb to prevent reverse osmosis (RO) membrane damage.

Test preparation

Before starting

Analyze the samples immediately. The samples cannot be preserved for later analysis.
Do a reagent blank calibration before a set of new reagent bottles is used. Refer to Reagent blank calibration on page 4.
Always do tests in 16-mm sample vials. Do not put the instrument in the sample or pour the sample into the sample port.
If the sample temperature is higher or lower than 20 ± 3 °C (68 ± 5.5 °F), refer to Temperature correction on page 4 to calculate the actual reading. ***Note:*** *As an alternative, keep the prepared sample at room temperature. Measure the sample when it is 20 ± 3 °C (68 ± 5.5 °F).*
Pretreat glass sample containers, sample vials and caps to remove chlorine demand. Refer to Sample collection on page 2.
Rinse the sample vial and cap with the sample three times before the vial is filled.
Make sure that the vials are clean and there are no scratches where the light passes through them.
Make sure that there are no fingerprints or liquid on the external surface of the vials. Wipe with a no-lint cloth before measurement.
Cold waters can cause condensation on the vial or bubbles in the vial. Examine the vial for condensation or bubbles. Remove condensation with a no-lint cloth. Invert the vial to remove bubbles.
Keep the cap on the vials when not in use to prevent ammonia contamination.
Install the vial cover on the vial before READ is pushed.
Do not use the same vials for free and total chlorine. If trace iodide from the total chlorine reagent is carried over into the free chlorine determination, monochloramine will interfere. It is best to use separate, dedicated vials for free and total chlorine measurements.
If the concentration of the sample is below or above the specified range of the method, the instrument will show a warning on the display.
If the test result is over-range, dilute the sample with a known volume of high quality, chlorine demand-free water and do the test again. Some loss of chlorine may occur due to the dilution. Multiply the result by the dilution factor.
After the test, immediately empty and rinse the sample cell. Rinse the sample cell and cap three times with deionized water.
Review the Safety Data Sheets (MSDS/SDS) for the chemicals that are used. Use the recommended personal protective equipment.
Dispose of reacted solutions according to local, state and federal regulations. Refer to the Safety Data Sheets for disposal information for unused reagents. Refer to the environmental, health and safety staff for your facility and/or local regulatory agencies for further disposal information.

Items to collect

Description	Quantity
Fluorescence Total Chlorine Reagent 1	5 drops
Fluorescence Total Chlorine Reagent 2	5 drops
Glass vial with cap, 16 mm	1

DOMINIQUE DUTSCHER SAS

Items to collect (continued)


Description	Quantity
Glass beaker, 100 mL (optional)	1
Sample vial holder	1

Refer to Consumables and replacement items on page 6 for order information.


Sample collection

- Analyze the samples immediately. The samples cannot be preserved for later analysis.
- Chlorine is a strong oxidizing agent and is unstable in natural waters. Chlorine reacts quickly with various inorganic compounds and more slowly with organic compounds. Many factors, including reactant concentrations, sunlight, pH, temperature and salinity influence the decomposition of chlorine in water.
- Collect samples in clean glass bottles. Do not use plastic containers because these can have a large chlorine demand.
- Pretreat glass sample containers, sample vials, caps and other items that the sample may touch to remove chlorine demand. Soak the glass sample containers, sample vials, caps and other items that the sample may touch in a weak bleach solution (1 mL commercial bleach to 1 liter of deionized water) for at least 1 hour. Rinse fully with deionized or distilled water. If sample containers, sample vials and caps are rinsed fully with deionized or distilled water after use, only occasional pretreatment is necessary.
- Make sure to get a representative sample. If the sample is taken from a spigot or faucet, let the water flow for at least 5 minutes. Let the container overflow with the sample several times and then put the cap on the sample container so that there is no headspace (air) above the sample.
- If the sample is collected in a sample vial from a spigot or faucet, adjust the water flow to add the sample carefully. As an alternative, collect the sample in a chlorine pretreated beaker. Rinse the beaker 3 times with the sample before the sample is collected.


Test procedure

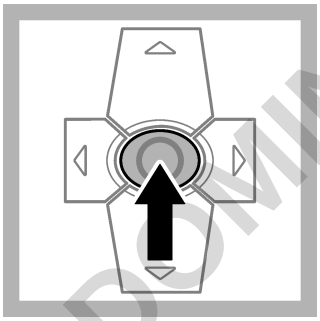

1. Set the instrument to on.

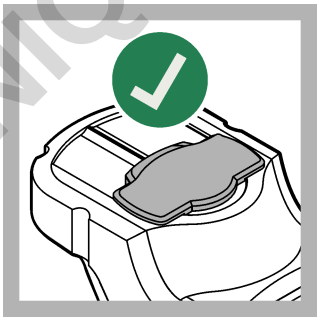

2. Keep the sample port cover closed during startup.

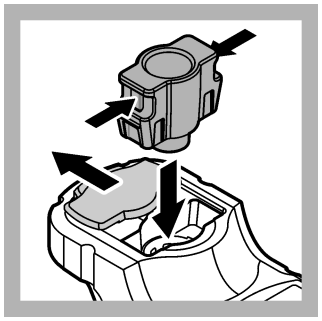

3. After startup, install the vial adapter. Push in the tabs on the sides of the vial adapter as shown.




4. Select **TESTS**.


DOMINIQUE DUTSCHER SAS




5. Select **TOTAL CHLORINE** .

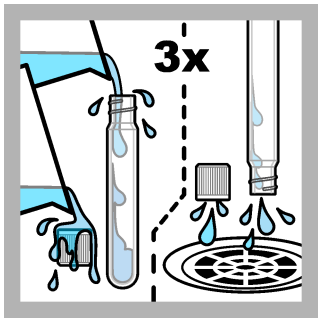


6. Rinse the sample vial and cap three times with the sample.

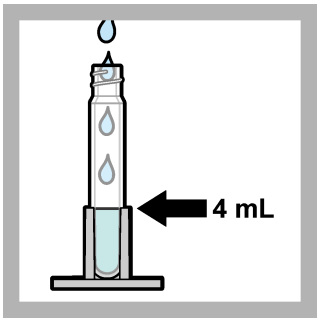


7. Put the sample vial in the sample vial holder. Add the sample to the top of the sample vial holder (4 mL).

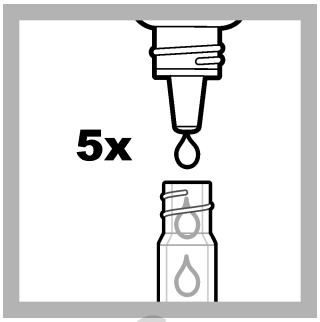


8. Immediately add **5** drops of Fluorescence Total Chlorine **Reagent 1**. Hold the bottle vertically.

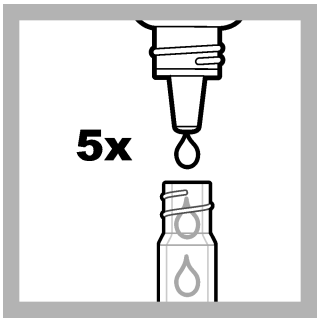

9. Immediately add **5** drops of Fluorescence Total Chlorine **Reagent 2**. Hold the bottle vertically.

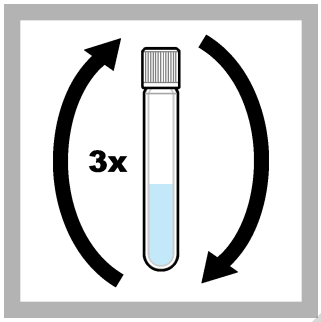

10. Put the cap on the vial. Invert the vial 3 times to mix. Put the vial in the sample vial holder.

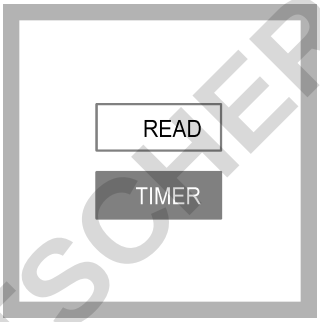


11. Select **TIMER**.




12. Select **START**.



13. When the timer ends, clean the external surface of the sample vial with a no-lint cloth.

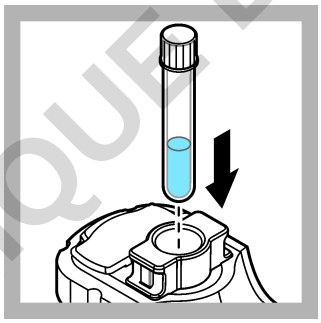

14. Put the sample vial in the vial adapter.

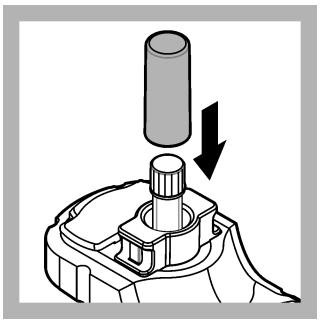

15. Put the vial cover on the sample vial.



16. Select **BACK**.

DOMINIQUE DUTSCHER SAS



17. Select **READ**. The results show in ppb (µg/L).

If the sample temperature is higher or lower than 20 ± 3 °C (68 ± 5 °F), refer to Temperature correction to calculate the actual reading.

Temperature correction

- If the sample temperature is **more** than 20 °C (68 °F):
  Actual reading = reading - [(0.015 × reading) × (sample temperature in °C - 20)]
- If the sample temperature is **less** than 20 °C (68 °F):
  Actual reading = reading + [(0.015 × reading) × (20 - sample temperature in °C)]

Interferences

Interfering substance	Interference level
Alkalinity	More than 500 mg/L $CaCO_3$. Adjust to pH 6–7 with 1 N Sulfuric Acid. Measure the amount to add on a separate sample aliquot, then add the same amount to the sample that is tested. Correct the test result for the dilution from the volume addition.
Bromine, Br_2	Positive interference at all levels
Chlorine Dioxide, ClO_2	No interference until 0.15 ppm
Hardness	No effect at less than 1000 mg/L as $CaCO_3$. For higher hardness samples, the sample may become cloudy after the reagents are added. However, a cloudy sample will have only a very small effect on the result.
Manganese, Oxidized (Mn^{4+}, Mn^{7+}) or Chromium, Oxidized (Cr^{6+})	Mn^{7+}—Interferes at all levels. 10 ppb positive interference for 50 ppb Mn^{7+} Mn^{4+}—No interference until 1 ppm Cr^{6+}—No interference until 0.2 ppm
Ozone	Interferes at all levels
Highly buffered samples or extreme sample pH	Can prevent the correct pH adjustment (of the sample) by the reagents. Sample pretreatment may be necessary. Adjust to pH 6–7 with acid (Sulfuric Acid, 1 N) or base (Sodium Hydroxide, 1 N). Correct the test result for the dilution caused by the volume additions.
Chloride (Cl^-)	Above 1200 ppm

Reagent blank calibration

Do a reagent blank calibration before a set of new reagent bottles is used.

1. Set the instrument to on. Keep the sample port cover closed during startup.
2. After startup, make sure that the sample port is clean. Use a soft cloth or no-lint paper tissue (or air duster) to clean the sample port as necessary.
3. Install the vial adapter in the instrument. Push in the tabs on the sides of the vial adapter.
4. Select **CALIBRATION** > **TOTAL CHLORINE** .

DOMINIQUE DUTSCHER SAS

5. Select **CAL-RB**.
6. Prepare a reagent blank. Replace the sample with deionized water (or distilled water or organic-free water) in the test procedure. Do not use the timer.
7. Clean the external surface of the sample vial with a no-lint cloth.
8. Install the sample vial in the vial adapter.
9. Put the vial cover on the sample vial.
10. Push the **Enter** key to measure the reagent blank.
11. Push the **Enter** key to save the results.
12. Immediately empty and rinse the sample vial. Rinse the sample vial and cap three times with deionized water (or distilled water).

Slope calibration check

Do a calibration check

Use the DR1300 FL Fluorescence Standards Kit to do a calibration check at regular intervals. Refer to Consumables and replacement items on page 6 for ordering information.

1. Make sure that the DR1300 FL Fluorescence Standards Kit is within the expiration date.
2. Get the Certificate of Analysis for the DR1300 FL Fluorescence Standards Kit.

   The Certificate of Analysis is available at https://app.hach.com/coaweb/customer_coa_request.asp.
3. Push **Enter** to set the instrument to on. Keep the sample port cover closed during startup.
4. After startup, make sure that the sample port is clean. Use a soft cloth or no-lint paper tissue (or air duster) to clean the sample port as necessary.
5. Install the vial adapter in the instrument. Push in the tabs on the sides of the vial adapter.
6. Select **CALIBRATION** > **TOTAL CHLORINE** .
7. Select **CALIBRATION CHECK**.
8. Measure deionized water as follows:
   a. Fill a sample vial with deionized water.
   b. Clean the external surface of the sample vial with a no-lint cloth.
   c. Install the sample vial in the vial adapter.
   d. Put the vial cover on the sample vial.
   e. Push **Enter** to select **ZERO**.
9. Measure the Standard 1 vial as follows:
   a. Clean the Standard 1 vial with a no-lint cloth.
   b. Install the Standard 1 vial in the vial adapter.
   c. Put the vial cover on the sample vial.
   d. Push **Enter** to select **READ**.
10. Measure the Standard 2 vial as follows:
    a. Clean the Standard 2 vial with a no-lint cloth.
    b. Install the Standard 2 vial in the vial adapter.
    c. Put the vial cover on the sample vial.
    d. Push **Enter** to select **READ**.
11. If the Standard 1 reading is not within ±2 ppb of the Total Chlorine value from the Certificate of Analysis and the Standard 2 reading is not within ±5 ppb of the Total Chlorine value from the Certificate of Analysis, do a slope calibration. Refer to the DL1300 FL User Manual for instructions.

DOMINIQUE DUTSCHER SAS

Method performance

The method performance data that follows was derived from laboratory tests that were measured on a DR1300 FL fluorometer during ideal test conditions. Users can get different results under different test conditions

Precision (95% confidence interval)
50.0 ± 1.0 µg/L Cl_2

Summary of method

Chlorine can be in water as free chlorine and as combined chlorine. Both forms can be in the same solution and can be determined together as total chlorine. Free chlorine is in a solution as hypochlorous acid or hypochlorite ion. Combined chlorine represents a combination of inorganic and organic chloramines.

The buffer is added to the sample to adjust the pH. The non-fluorescent indicator reacts with the chlorine in the sample to make fluorescent molecules, which increases the fluorescence intensity. The meter calculates the chlorine concentration based on the fluorescence intensity change.

To determine the concentration of combined chlorine, complete a free chlorine test and a total chlorine test. Subtract the results of the free chlorine test from the total chlorine test to get the combined chlorine concentration.

Consumables and replacement items

Required reagents

Description	Quantity	Item no.
Reagent kit, ULR Fluorescence Total Chlorine, includes: Sample vial holder, two 16-mm vials with caps and two reagents in dropper bottles	1	34252000
Refill kit, ULR Fluorescence Total Chlorine Reagents, 100 tests, includes: Two reagents in dropper bottles	1	34252001

Required apparatus

Description	Quantity/test	Unit	Item no
Fluorescence sample vials with caps, 10-mL round, 16 mm x 100 mm	1	6/pkg	100866

Optional apparatus

Description	Unit	Item no.
Beaker, glass, 100 mL	each	50042H
Sodium Hydroxide Standard Solution, 1.0 N	100 mL MDB	104532
Sulfuric Acid Standard Solution, 1 N	100 mL MDB	127032
Pipet, TenSette, 1.0–10.0 mL	each	1970010
Pipet tips for TenSette Pipet, 1.0–10.0 mL	50/pkg	2199796
Paper, pH, 0–14 pH range	100/pkg	2601300
Sample vial holder	1	3563500
Water, organic-free	500 mL	2641549
Deionized water	500 mg	27248

DOMINIQUE DUTSCHER SAS

FOR TECHNICAL ASSISTANCE, PRICE INFORMATION AND ORDERING:
In the U.S.A. – **Call toll-free 800-227-4224**
Outside the U.S.A. – **Contact the HACH office or distributor serving you.**
On the Worldwide Web – **www.hach.com**; E-mail – **techhelp@hach.com**

HACH COMPANY
WORLD HEADQUARTERS
Telephone: (970) 669-3050
FAX: (970) 669-2932

© Hach Company/Hach Lange GmbH, 2022, 2024. All rights reserved.
05/2024, Edition 3