



# Foot Valve Complete System & System Network

## Foot Valve Complete System

Cat # 59-121C Includes: Flypad; Flypad Frame; Blowgun; Foot Pedal; "T" Fitting; 5 feet 1/8" I.D. tubing. All components also available separately. For details, visit [www.flystuff.com](http://www.flystuff.com).

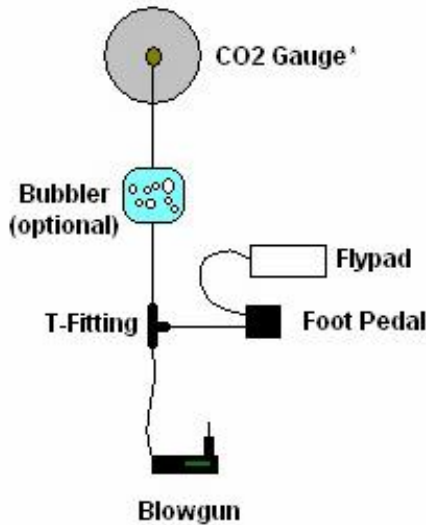

## Foot Valve System Network

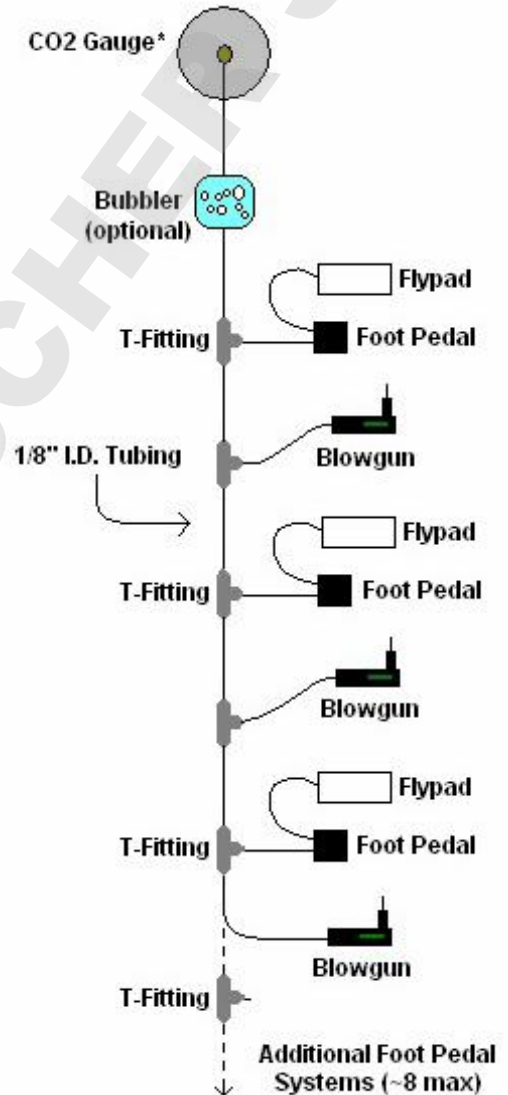

### Notes:

- 1.) Foot valve systems (cat # 59-121C) offer the best efficiency of CO<sub>2</sub> usage since introduction of CO<sub>2</sub> into the flypad (cat # 59-114) is user regulated.
- 2.) The flypad frame (cat# 59-118) will further increase the efficiency of CO<sub>2</sub> usage, since CO<sub>2</sub> (heavier than air) will be better retained in the flypad's staging area.
- 3.) A tube reducer (cat# 59-152, 59-153) may be necessary to connect 1/8" I.D. tubing to the CO<sub>2</sub> tank gauge, as most gauge outputs are not configured for 1/8" I.D. tubing (cat# 59-124).