

BradyJet J2000 Colour Label Printer

Identify samples by colour to
increase productivity



www.bradyeurope.com/J2000



The more samples laboratories need to process, the more time consuming it becomes to retrieve one specific sample. Brady Corporation proposes sample colour coding to solve this and offers the BradyJet J2000 Colour Label Printer to do it, to stay productive and to avoid spending too much time on sample hunting.

Sort samples by colour

The BradyJet J2000 Colour Label Printer can print any colour on reliable laboratory sample labels. By setting up a colour code, laboratory specialists will be able to recognise sample categories at a glance. With over 16 million colour variations, users may switch colour every day, or use a different colour per application, storage or urgency type depending on the needs of their laboratory. Next to colour, the BradyJet J2000 also offers 4800 dpi photo quality print to add highly legible data or barcodes a decent scanner could never misread.



Print reliable labels

The BradyJet J2000 prints on Brady's reliable laboratory sample labels that come in matching sizes and shapes for most sample containers including tubes, vials, cryo containers, dishes, plates and general glassware. These labels are designed to stay attached and to remain legible when printed with a Brady printer, even if they are stored in liquid nitrogen and in freezers, or when used in hot water baths or autoclaves. Brady's laboratory sample labels can also resist common laboratory chemicals like acetone, ethanol, toluene, xylene and IPA.



Easy to use

With 26.4 by 38.9 cm dimensions, the BradyJet J2000 has a small footprint and is very practical and easy to use. Spring loaded Brady label rolls can be switched in less than 20 seconds and with automatic label set up, the printer is ready to print in no time. And with up to 101.6 mm per second print speed, even larger colour coded sample labels can quickly be printed on demand.



For more information or a
live demo, visit
www.bradyeurope.com/J2000

Printer specifications

PRINTER	
Specifications	BradyJet J2000 Colour Label Printer
Includes	This kit contains a J2000 printer, a power cord and AC adapter, a USB cable, a Quick Start Guide, an ink cartridge, a material roll for printhead alignment and a Brady Workstation Suite
Print Speed (mm/sec)	101
Print Technology	Inkjet
Print Resolution (dpi)	Up to 4800
Stand-alone Printing (Detached from PC)	Peripheral
Colour Capability	Full Colour
Weight (kg)	3.20
Size W x H x D (mm)	264 x 180 x 389
Cutter	Auto Cutter
Warranty	1 year
MEDIA	
Continuous vs Die-cut Labels	Continuous Tape & Die Cut Labels
Custom Labels	Yes
Materials Supported	B-2595, B-2569, B-2551, B-7425J
Min. Label Length (mm)	6.35
Min. Label Width (mm)	19
Max. Label Width (mm)	101.60
Max. Liner Width (mm)	108.00
Max. Print Width (mm)	101.60
Recommended Usage per Day	Around 3000 labels per day, depending on specific label size, text printed and other factors
SOFTWARE & CONNECTIVITY	
Connectivity Options	USB
Software Compatibility	Brady Workstation, Windows based driver for 3rd party software use
COMPLIANCE	
Approval/Compliance	USA: UL, FCC Class A; Canada: cUL, IC/ISED; European Union : CE, RoHS 2011/65/EU, 2014/35/EU, 2015/863, LVD, 2014/30/EU, RED 2014/53/EU; Australia/New Zealand : RCM

Materials

Diagram	Order Reference	Art. No.	Colour	Label Width A (mm)	Label Height B (mm)	Diameter C (mm)	Total Qty/ Pack
Inkjet Receptive Polypropylene							
▲	J20-257-7425J	150007	White	20.96	6.35	—	1500
▲	J20-151-7425J	150011	White	50.80	6.35	—	2000
▲	J20-152-7425J	150005	White	25.40	9.53	—	2000
⊕	J20-249-7425J	150023	White	25.40	9.53	9.53	1500
▲	J20-59-7425J	150006	White	25.40	12.70	—	2000
⊕	J20-251-7425J	150033	White	25.40	12.70	11.18	1500
▲	J20-37-7425J	150020	White	38.10	12.70	—	1500
▲	J20-53-7425J	150024	White	50.80	12.70	—	2000
▲	J20-256-7425J	150016	White	33.02	15.24	—	1500
▲	J20-88-7425J	150015	White	25.40	19.05	—	2000
▲	J20-136-7425J	150019	White	38.10	19.05	—	1500
▲	J20-179-7425J	150025	White	25.40	25.40	—	1500
▲	J20-235-7425J	150022	White	38.10	25.40	—	1250
▲	J20-137-7425J	150018	White	50.80	25.40	—	1500
▲	J20-258-7425J	150028	White	38.10	44.45	—	900

Printers

Article No.	Order Reference	Description
150159	J2000-EU	BradyJet J2000 Colour Label Printer - EU
199966	J2000-EU-LABS	BradyJet J2000 Colour Label Printer - EU with Brady Workstation LAB Suite
150157	J2000-UK	BradyJet J2000 Colour Label Printer - UK
199967	J2000-UK-LABS	BradyJet J2000 Colour Label Printer - UK with Brady Workstation LAB Suite
150162	J2000	BradyJet J2000 Colour Label Printer - US
199971	J2000-US-LABS	BradyJet J2000 Colour Label Printer - US with Brady Workstation LAB Suite



The Basic Suite is included with every printer, all you need to do is download it.
Workstation.BradyID.com/Free




Easily identify any sample with the most wanted laboratory sample label design apps in our Brady Workstation Lab suite.

Accessories

Article No.	Order Ref.	Description
150154	J20-CMY	CMY pigment based ink cartridge for J2000 Printer
150160	J20-CUTTER	Replacement Cutter Blade for J2000
150163	J20-PAD	Replacement Ink Pads for J2000



For more information on supported LAB materials, visit www.bradyeurope.com/J2000

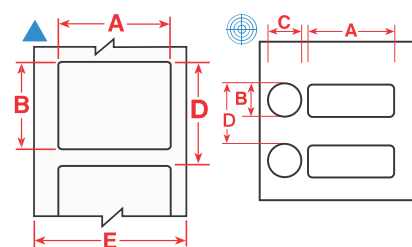

CMY pigment based ink cartridge for J2000 Printer

Polypropylene Labels for J2000 printer