

5 REASONS for choosing ROCKER vacuum pumps



Product Protection:



1 Prote™ overflow protection

Rocker series pumps are equipped with Prote™ overflow protection which can prevent a large amount of liquid from being sucked into the pump.



2 Moisture trap with filter cartridge

Rocker series pumps are supplied with a filter cartridge at the air inlet to absorb particle and moisture which can prolong the life of pump.

Environmentally Friendly Production:



3 RoHS Compliance

ROCKER is committed to manufacturing environment-friendly products, so our pumps are already RoHS compliant as manufactured.



4 No air pollution, maintenance free

Rocker series pumps are driven by piston, oil free, no regular oil changes and no oil pollution to the Lab.




International safety certification :

5

- CE certification
- CSA certification
- IP 30 protection class

How does the Prote™ overflow protection work?

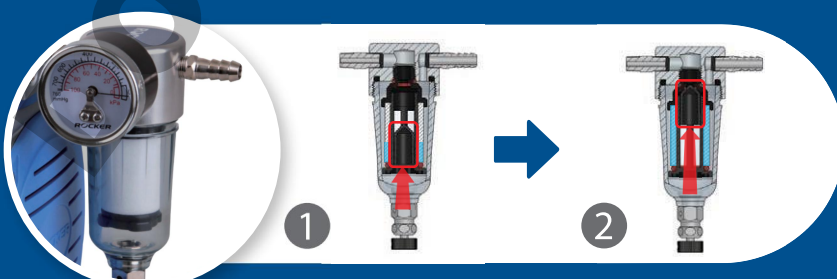

1 Once the liquid is sucked into the moisture trap accidentally, the buoy of Prote™ overflow protection will float up.

2 Then block the air inlet to prevent the liquid from flowing into the pump.

ROCKER

Rocker Scientific Co., Ltd.

6F., No.6, Zhulin Rd., Linkou Dist.,
New Taipei City 244, Taiwan (R.O.C.)

T/ 886-2-26033311 F/ 886-2-26036622

www.rocker.com.tw