


Simport



CryoSette™
Frozen Tissue Storage Containers

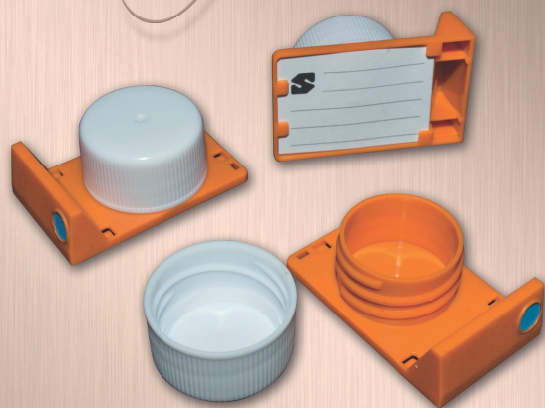

2588 Bernard-Pilon
Beloeil, Quebec
J3G 4S5
Canada

T - 450 464-1723
F - 450 464-3394
info@simport.com
www.simport.com

Distributed by:

MKT-M956LL1P4 — 06/07/2018

Printed in Canada



Simport

The CryoSette™

Frozen Tissue Storage Containers



www.simport.com

M956 CryoSette™

Frozen Tissue Storage Containers

M9560

Base made of acetal
Screw closure made of high-density polyethylene

- Offers sample specimen storage in upright chest freezers and vapor-phase Dewar flasks at -196 °C*
- 2D Barcoding available for efficient sample tracking
- Maintains the integrity of tissue samples
- Choice of two storage boxes
- 95 kPa tested

The Simport CryoSette™ is designed for frozen tissue collection, transport and storage. It can also be used for storing a multitude of other specimens in all laboratories. The robust design of the cryosette ensures that tissue morphology will be well preserved.

This 2.5 ml container features a wide-mouth opening and a high integrity screw closure. It is designed with a flat bottom and straight sides for easy tissue removal. The specific shape of the Cryosette allows its manipulation to be optimized, whether you are using gloves or not. The CryoSette™ efficiently stores tissue samples in -86 °C ULT freezers or in the gas phase of liquid nitrogen in Dewar flasks at -196 °C but not in LN₂ liquid phase.

The CryoSette™ screw closure and base are designed so they will not seize at low temperatures. Rapid visual 1 1/4 turn, even during repeated freeze / thaw cycles. Color-coded inserts are available in 5 colors.

The CryoSette™ offers two surfaces for handwriting or applying barcode labels: an anterior surface and a larger space underneath the base. A thin cardboard I.D. label and a plastic transparent cover are available for the latter. A 2D barcode insert can also be placed on the anterior writing surface for rapid visual sample identification. The CryoSette™ is compatible with common tracking methods. The anterior writing surface is improved for ink adherence.

The large size of the CryoSette™ chamber is suitable for most storage applications.
Chamber dimensions: 21 mm O.D x 18 mm I.D. x 10.5 mm H.
Total CryoSette™ dimensions : 24 mm x 38 mm x 11 mm H.



M956-21 and M956-40


CryoSette™ Frozen Tissue Storage Boxes

M956-21 made of polypropylene

- Compatible with leading frozen sample storage systems
- Thick walls minimizing deformation due to temperature changes
- Can be handled by robotic equipment
- Rack orientation marks on sides
- Latches for locking and unlocking rack lid
- Autoclavable

The M956-21 storage rack is designed to contain 21 frozen tissue CryoSette™ containers. The box is compact and has an SBS standard footprint to make suitable to several robotic storage systems. These innovative boxes can be stored in low-temperature freezers.


Temperature range -196 °C to 121 °C. Suitable for cryogenic storage, but only in gas phase of liquid nitrogen.
Dimensions : 127,76 mm x 85,48 mm x 41,6 mm H.

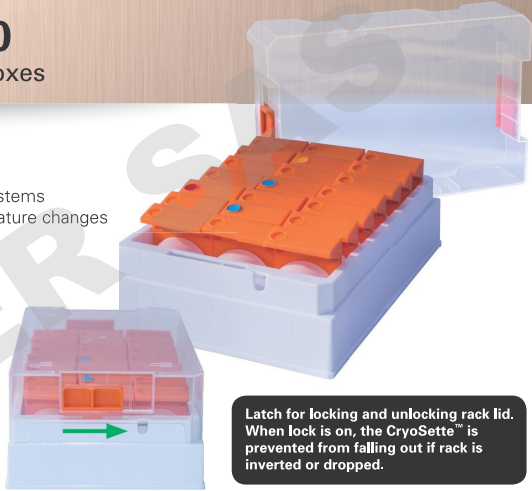

M956-40 made of polycarbonate

The M956-40 CryoSette™ storage box can accommodate up to 40 Tissue Storage Containers in a space as small as 133 x 133 x 52 mm H. Made of extra strong polycarbonate, this durable cryogenic storage box is designed to be used at temperatures between -196 °C and +121 °C and are autoclavable at 120 °C, 15 psig (1 bar) for 20 minutes. Only in gas phase of LN₂.

The frosted cover allows the user to write directly on top of the box, and is keyed to the base to prevent misalignment.

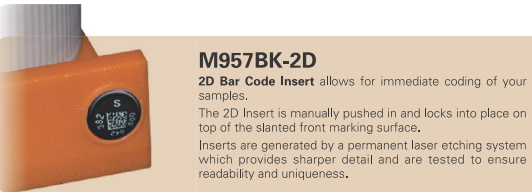

Cat. #	Description	Qty/Pk	Qty/Cs
M9560	CryoSette™ Tissue Storage Container	50	250
M956-21W	CryoSette™ 21 Place Storage Box	—	10
M956-40B	CryoSette™ 40 Place Storage Box (Separator included)	5	10
M956SEP	Cardboard Separator kit for M956-40B	—	6
M956LAB	Cardboard I.D. Label for M956	100	—
M956COV	Transparent Plastic I.D. Label Cover for M956	100	—
M957BK-2D	CryoSette™ 2D Barcode Insert	100	500
M957B	CryoSette™ Color Insert – Blue	100	500
M957G	CryoSette™ Color Insert – Green	100	500
M957L	CryoSette™ Color Insert – Lilac	100	500
M957R	CryoSette™ Color Insert – Red	100	500
M957Y	CryoSette™ Color Insert – Yellow	100	500

Available options: Sterile, pre-capped with screw caps