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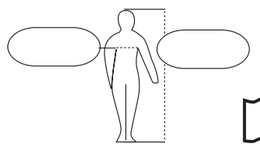
ENGLISH

1. Trademark
2. Model identification
3. CE marking : PPE category III
4. European Protection Standards
5. Size and body measurements
6. Refer to instructions for use
7. Fabric care symbols
8. Single use

PHYSICAL PERFORMANCE OF WEEJET®

Fabric Physical Properties	*Test Method*	*Class*
Abrasion resistance	EN 530	6
Flex cracking	ISO 7854 B	6
Trapezoidal Tear MD / XD	ISO 9073-4	3
Tensile Strength MD / XD	ISO 13934-1	2
Puncture resistance	EN 863	2
Seam Strength	ISO 13935-2	4

RESISTANCE TO PENETRATION BY LIQUIDS EN ISO 6530

Fabric Chemical Properties	*Penetration*	*Repellency*
Sulphuric Acid 30%	CLASS 3	CLASS 3
Sodium Hydroxide 10%	CLASS 3	CLASS 3
Butan-1-ol	CLASS 3	CLASS 3
o-xylene	CLASS 3	CLASS 3

FABRIC AND TAPED SEAMS RESISTANCE TO PERMEATION BY LIQUIDS EN ISO 6529

Fabric Chemical Properties	*Breakthrough time (min)*	*Class*
Sulphuric Acid 50 %	480 min	6

FABRIC RESISTANCE TO PENETRATION OF INFECTIVE AGENTS EN 14126

Test Method	*Class*
Resistance to penetration by blood-borne pathogen ISO 16604	6/6
Resistance to penetration by contaminated liquids EN ISO 22610	6/6
Resistance to penetration by biologically contaminated aerosols EN ISO 22611	3/3
Resistance to penetration by contaminated solid particles ISO 22612	3/3

WHOLE SUIT TEST PERFORMANCE

Liquid tight suits and Spray Tight Suits **Type 3B** EN 14605	PASS
High level spray test **Type 4B** EN 14605	PASS
Dry particles **Type 5 B** Method as defined by EN ISO 13982-1:2004 + A1:2010	PASS
Light splashes **Type 6 B** Method as defined by EN 13034 : 2005 +A1: 2009	PASS
Nominal protection factor according to **EN 1073-2** : 2002	CLASS 1
Surface resistivity **EN1149-5** : 2018	PASS

Model and typical area of use

Yellow coverall with elasticated hood, wrists, ankles and back, adhesive flap over the zipper thumb loop. Coverall is designed to protect against chemical risks (jets – splashes – sprays – airborne solid particulates) and against infectious agents and against contamination by radioactive particles and electrostatics (dissipative properties).

How to wear protective clothing

Remove the coverall from its packaging, open the center zipper and put on the garment, taking care not to damage it. Close the zipper completely and cover it with the adhesive flap, without creasing. If there is a risk of airborne solid particles, it is advisable to tape additionally protective gloves with ends of the sleeves.

Limitations of use

Do not use the coverall if any defects are noticed (e.g.: holes, damaged seams, ….). The coverall meets the requirements: Ljmn 82/90 ≤30% - Ls 8/10 ≤ 15%. Type 6 chemical protective suits have been tested to the whole suit test. Wear for long periods of time can cause heat stress. Heat stress and discomfort can be reduced or eliminated by using appropriate undergarments or suitable ventilation equipment. Coverall is for single use only and must be disposed after any job. The user is the sole judge of the adequacy between his risk analysis and the type of protection required, including the additional equipment necessary for a complete protection of the body. If necessary, use additional PPE with the same level of protection (such as gloves, breathing apparatus, etc.) in order to provide full body protection. The manufacturer cannot accept responsibility for any improper use or disposal of garments produced by them. This garment is flammable, keep away from fire. The person wearing the electrostatic dissipative protective clothing shall be properly earthed. The resistance between the person skin and earth shall be less than 10^8 Ω e.g. by wearing adequate footwear on dissipative or conductive floors. Electrostatic dissipative protective clothing shall not be open or removed whilst in presence of flammable or explosive atmospheres or white handling flammable or explosive substances. Electrostatic dissipative protective clothing is intended to be worn in Zones 1, 2, 20, 21 and 22 (see EN 60079-10-1 and EN 60079-10-2) in which the minimum ignition energy of any explosive atmosphere is not less than 0,016 mJ. Electrostatic dissipative protective clothing shall not be used in oxygen enriched atmospheres, or in Zone O (see EN 60079-10- 1) without prior approval of the responsible safety engineer. The electrostatic dissipative performance of the electrostatic dissipative protective clothing can be affected by wear and tear, laundering and possible contamination electrostatic dissipative protective clothing shall be worn in such a way that it permanently covers all non-complying materials during normal use. Use the product within 5 years of the date of manufacturing.

Storage and disposal

The coveralls must be stored in their original packaging in a dry place, away from sunlight. They must be disposed of in compliance with current regulations.

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