




**THE No. 1 System  
for Solvent Safety**

# The next Level! Safety Caps 2.0

Safe Solvent Supply for  
any HPLC / UHPLC System.

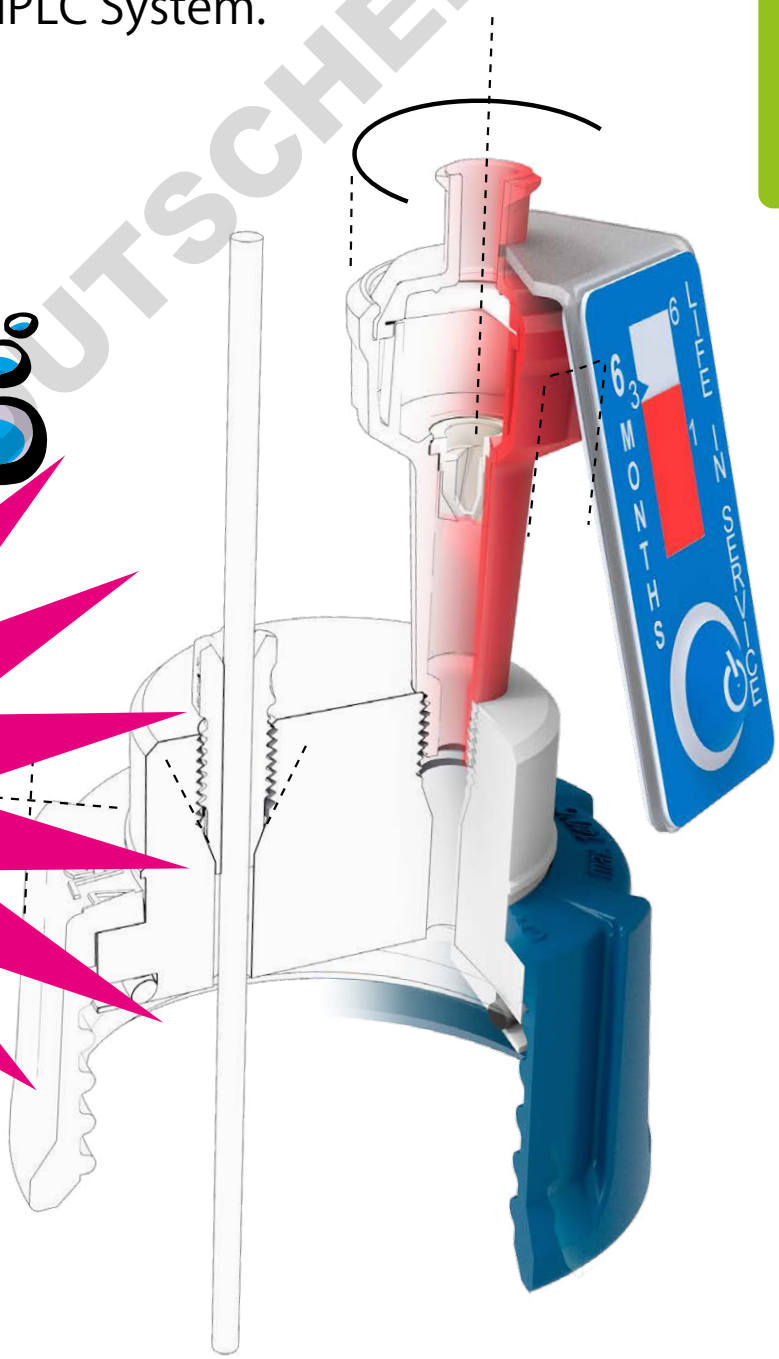

[www.scatt-europe.com](http://www.scatt-europe.com)

# The next Level!

Safety Caps 2.0

- ✓ **Luer Lock Adapter**  
Useful for many applications
- ✓ **Air valve**  
for analytical and preparative applications
- ✓ **PFA Fittings**  
Improved design for easy tubing connection. Excellent chemical resistance. Flammability classification V-0 (UL-94).
- ✓ **Medical Grade PTFE**  
Excellent resistance. Food safe. FDA conform. Flammability classification V-0 (UL-94).
- ✓ **Insertion Label**  
Mark your bottles with individual descriptions
- ✓ **PPS Screw Cap**  
Autoclavable / sterilisable up to 200 °C. Flammability classification V-0 (UL-94).
- ✓ **PTFE Filter**  
For clean solvents
- ✓ **Exchange Indicator**  
Starts working when the button is pressed: never miss an exchange again!
- ✓ **Labelling Field**  
Simply turn the clip around and label individually!
- ✓ **Valve membrane**  
Blocks solvent vapours, opens at underpressure
- ✓ **PTFE-Body**  
360° freely rotatable. Provides easy container change without distortion of tubing.
- ✓ **PTFE Body with integrated FKM Seal**
- ✓ **Ergonomic Design**  
Enhanced stability and construction. Improved handling. Easy installation. Smart exchange.

\* Air valve and Fittings for tubing connection are included in delivery.



Part No.	Description	Unit	Thread size	Connections for Ø 3.2 mm (1/8") OD tubing	Blind plugs for tubing connections
307 019	Safety Cap I (V 2.0) *	1	GL 45	1x	-
307 909	Safety Cap II (V 2.0) *	1	GL 45	2x	1x
307 910	Safety Cap III (V 2.0) *	1	GL 45	3x	2x
307 410	Safety Cap IV (V 2.0) *	1	GL 45	4x	2x
307 520	Safety Cap VI (V 2.0) *	1	GL 45	6x	2x
317 010	Air valve for Safety Caps	1	UNF 1/4" 28G	-	-
397 010	Air valve for Safety Caps	10	UNF 1/4" 28G	-	-
107 063	Fitting for Ø 3.2 mm OD tubing, blue	5	UNF 1/4" 28G	-	-
107 060	Fitting for Ø 3.2 mm OD tubing, black	5	UNF 1/4" 28G	-	-
107 065	Fitting for Ø 3.2 mm OD tubing, yellow	5	UNF 1/4" 28G	-	-
160 501	Blind plug for tubing connections	10	UNF 1/4" 28G	-	-

# HPLC Starter Kits

- Price advantage compared to single order
- Suitable for all common HPLC systems (independent of manufacturer)

- Use blind plugs (included) to lock unused connections



Fig.	Part No.	Description	Content	Unit	Thread size	Connections for $\varnothing$ 3.2 mm OD tubing
<b>A</b>	399 200	<b>Starter-Kit-Supply-S1</b>	Safety Cap I (V 2.0)	3x	GL 45	1x
			Safety Cap II (V 2.0)	1x	GL 45	2x
			Blind plug for tubing conn.	1x	UNF 1/4" 28G	-
<b>B</b>	399 201	<b>Starter-Kit-Supply-S2</b>	Safety Cap II (V 2.0)	4x	GL 45	2x
			Blind plug for tubing conn.	4x	UNF 1/4" 28G	-



## SCAT - Europe GmbH

Waldecker Straße 7  
64546 Mörfelden-Walldorf  
Germany

 **+49 (0) 61 05 / 30 55 86 - 0**

 **+49 (0) 61 05 / 30 55 86 - 99**

 **info@scat-europe.com**

 **www.scat-europe.com**



## SCAT - Americas, Inc.

1704 East Blvd., Suite 101  
Charlotte, NC 28203  
USA

 **+1.704.438.9477**

 **+1.704.373.2603**

 **info@scat-americas.com**

 **www.scat-americas.com**

